



# NOELLE<sup>®</sup> S554.100 User Guide

---

## Maternal and Neonatal Birthing Simulator

NOELLE is an interactive educational system developed to assist a certified instructor. It is not a substitute for a comprehensive understanding of the subject matter and not intended for clinical decision making.

### **User Guide 12.3.1**

© Gaumard Scientific Company, 2012

All Rights Reserved

[www.Gaumard.com](http://www.Gaumard.com)



# Contents

<b>CONTENTS.....</b>	<b>3</b>
<b>END USER LICENSE AGREEMENT .....</b>	<b>5</b>
<b>CARE AND CAUTIONS .....</b>	<b>7</b>
OVERALL WARNINGS .....	8
<b>GETTING STARTED .....</b>	<b>9</b>
OVERVIEW .....	10
TERMINOLOGY .....	10
<b>EQUIPMENT SET-UP.....</b>	<b>11</b>
<i>Leg Assembly .....</i>	<i>12</i>
<i>Power Supply.....</i>	<i>12</i>
<i>Control Computer .....</i>	<i>12</i>
<i>Virtual Monitor.....</i>	<i>13</i>
<b>WORKING WITH GAUMARDUI.....</b>	<b>15</b>
INITIALIZING THE SIMULATOR .....	16
<i>User Profile Menu .....</i>	<i>16</i>
THE ENVIRONMENT .....	17
<i>Details .....</i>	<i>19</i>
<i>Palette.....</i>	<i>23</i>
<i>Lab .....</i>	<i>24</i>
<i>Scenarios.....</i>	<i>30</i>
<i>Labor.....</i>	<i>46</i>
<i>Model (Newborn) .....</i>	<i>67</i>
<i>Speech.....</i>	<i>68</i>
<i>Log.....</i>	<i>70</i>
<i>Evaluation Form.....</i>	<i>75</i>
MENUS .....	81
<i>File .....</i>	<i>81</i>
<i>Setup .....</i>	<i>83</i>
<i>Monitors .....</i>	<i>89</i>
<i>CPR .....</i>	<i>92</i>
<i>A/V .....</i>	<i>94</i>
<i>Help .....</i>	<i>95</i>
<b>WORKING WITH NOELLE.....</b>	<b>96</b>
AIRWAY .....	97
<i>Nasal and Oral Intubation.....</i>	<i>97</i>
BREATHING.....	97
<i>Pulmonary Ventilation.....</i>	<i>97</i>
<i>CPR .....</i>	<i>97</i>
CIRCULATION .....	97
<i>Bilateral IV arms .....</i>	<i>97</i>
OBSTETRICS .....	99
<i>Normal Labor and Delivery.....</i>	<i>99</i>

<i>Vacuum-Assisted Delivery</i> .....	102
<i>Leopold Maneuver</i> .....	103
<i>Shoulder Dystocia</i> .....	104
<i>Cesarean Delivery</i> .....	104
<i>Prolapse of the Umbilical Cord</i> .....	105
<i>Placenta Previa</i> .....	105
<i>External Version</i> .....	105
<i>Breech Birth</i> .....	106
<i>Placenta Delivery</i> .....	107
<i>Postpartum Activity</i> .....	108
<i>Post-Partum Hemorrhage</i> .....	108
<i>Episiotomy Repair</i> .....	109
SYSTEMIC .....	109
OTHER .....	110
<i>Vital Signs Monitor</i> .....	110
<b>WORKING WITH NEWBORN</b> .....	<b>111</b>
AIRWAY .....	112
BREATHING .....	112
CARDIAC .....	112
CIRCULATION .....	112
CEPHALIC .....	112
SYSTEMIC .....	112
OTHER .....	113
<b>APPENDIX</b> .....	<b>114</b>
FACTORY PRESET LABOR SCENARIOS .....	115
<i>Newborn Flowcharts</i> .....	167
TIPS ON CREATING SCENARIOS .....	179
TROUBLESHOOTING .....	181
<i>General Troubleshooting Guide</i> .....	181
<i>Microphone Boost (Windows® 7)</i> .....	183
<i>Connecting to the Gaumard Monitors</i> .....	184
<i>Diagnostics</i> .....	188
SELECTED CONSUMABLES AND REPLACEMENTS PARTS .....	189
<i>Selected Parts List</i> .....	189
<i>Replacing Common Consumable and Replacement Parts</i> .....	191
WARRANTY .....	193
CONTACT US .....	194

# End User License Agreement

This is a legal agreement between you, the end user, and Gaumard® Scientific Company, Inc. ("Gaumard"). This software is protected by copyright laws and remains the sole property of Gaumard. By installing the GaumardUI simulation (the "Software") media, you agree to be bound by the terms of this agreement. If you do not agree to the terms of this agreement, promptly return the uninstalled media and accompanying items to Gaumard at the address indicated below.

**1. Grant of License.** Gaumard hereby grants to you (an individual or institution) the right to install and activate the Software on one computer for use with one Code Blue III Interactive patient simulator system. The software may also be installed on any number of other computers at the same institution so that students may access the learning resources. One copy of the software may be made for backup purposes. You may not network this Software, or allow multiple users unless you purchased a multi-user workstation license. Sharing this Software with other individuals or allowing other individuals to view the contents of this Software is in violation of this license.

**2. Copyright.** The Software is owned by Gaumard and protected by United States copyright laws and international treaty provisions. Therefore, you must treat this Software like any other copyrighted material. You may not make this Software or copies thereof available in any manner or form or use, copy or transfer the Software, in whole or in part, except as provided herein.

**3. Other Restrictions.** You may not rent or lease this Software to any other party. You may not alter, merge, modify, adapt, reverse engineer, decompile or disassemble the software, or disclose the contents of this Software to any other party.

**4. Electronic Transmission of Software.** If you received the Software by electronic transmission or by Internet delivery, by installation of the Software, you acknowledge that you have read and understand this license agreement and agree to be bound by its terms and conditions.

**5. Term of Agreement.** The term of this Agreement and the license granted to you pursuant hereto shall commence upon installation of this Software. This Agreement and the license granted herein may otherwise be terminated by Gaumard in the event that you are in breach of any provision of this Agreement. In the event of termination, you agree to immediately return this Software, accompanying items, and any copies thereof to Gaumard.

## **6. LIMITED WARRANTY**

**(A)** THE CD-ROM MEDIA (THE "MEDIA") WHICH CONTAINS THIS SOFTWARE IS WARRANTED, FOR A PERIOD OF 30 DAYS FROM THE DATE OF PURCHASE, TO BE FREE FROM DEFECTS IN MATERIAL AND WORKMANSHIP. ELECTRONIC TRANSMISSION IS WARRANTED TO BE FREE FROM DEFECTS AT THE MOMENT OF TRANSMISSION. YOUR SOLE AND EXCLUSIVE REMEDY, AND GAUMARD'S SOLE LIABILITY, IS TO REPLACE THE DEFECTIVE MEDIA OR TO REPEAT THE ELECTRONIC TRANSMISSION PROVIDED THAT YOU NOTIFY GAUMARD IN WRITING OF SUCH DEFECT OR DEFECTIVE TRANSMISSION AND RETURN THE DEFECTIVE MEDIA, IF ANY, DURING THE 30-DAY WARRANTY PERIOD.

**(B)** EXCEPT AND TO THE EXTENT EXPRESSLY PROVIDED IN PARAGRAPH (A), THE SOFTWARE AND ACCOMPANYING WRITTEN MATERIALS ARE PROVIDED ON AN "AS IS" BASIS, WITHOUT ANY WARRANTIES OF ANY KIND, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE. NO ORAL OR WRITTEN INFORMATION OR ADVICE GIVEN BY GAUMARD, ITS DEALERS, DISTRIBUTORS, AGENTS OR EMPLOYEES SHALL CREATE A WARRANTY OR IN ANY WAY INCREASE THE SCOPE OF THIS WARRANTY, AND YOU MAY NOT RELY ON ANY SUCH INFORMATION OR ADVICE. GAUMARD DOES NOT WARRANT, GUARANTEE, OR MAKE ANY REPRESENTATIONS REGARDING THE USE OR THE RESULTS OF USE, OF THE SOFTWARE OR WRITTEN MATERIALS IN TERMS OF CORRECTNESS, ACCURACY, RELIABILITY, CURRENTNESS, OR OTHERWISE, AND THE ENTIRE RISK AS TO THE RESULTS AND PERFORMANCE OF THE SOFTWARE IS ASSUMED BY YOU. IF THE SOFTWARE OR WRITTEN MATERIALS ARE DEFECTIVE, YOU AND NOT GAUMARD OR ITS DEALERS, DISTRIBUTORS, AGENTS, OR EMPLOYEES, ASSUME THE ENTIRE COST OF ALL NECESSARY SERVICING, REPAIR OR CORRECTION OTHER THAN EXPRESSLY DESCRIBED ABOVE.

**(C)** NEITHER GAUMARD NOR ANYONE ELSE WHO HAS BEEN INVOLVED IN THE CREATION, PRODUCTION OR DELIVERY OF THIS PRODUCT SHALL BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL OR INCIDENTAL DAMAGES (INCLUDING DAMAGES FOR LOSS OF BUSINESS PROFITS, BUSINESS INTERRUPTION, LOSS OF BUSINESS INFORMATION, AND THE LIKE) ARISING OUT OF THE USE OR INABILITY TO USE SUCH PRODUCT OR RELATED TO THIS AGREEMENT EVEN IF GAUMARD HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. GAUMARD SHALL NOT BE LIABLE TO YOU FOR ANY INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES OR LOST PROFITS ARISING OUT OF OR RELATED TO THIS AGREEMENT OR YOUR USE OF THE SOFTWARE AND/OR THE RELATED DOCUMENTATION, EVEN IF GAUMARD HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN NO EVENT SHALL GAUMARD'S LIABILITY HERE UNDER, IF ANY, EXCEED THE PURCHASE PRICE PAID BY YOU FOR THE SOFTWARE.

ALL RIGHTS NOT EXPRESSLY GRANTED IN THIS LICENSE AGREEMENT ARE RESERVED BY GAUMARD.

#### **ACKNOWLEDGMENT**

BY INSTALLATION OF THIS SOFTWARE, YOU ACKNOWLEDGE THAT YOU HAVE READ AND UNDERSTAND THE FORE GOING AND THAT YOU AGREE TO BE BOUND BY ITS TERMS AND CONDITIONS. YOU ALSO AGREE THAT THIS AGREEMENT IS THE COMPLETE AND EXCLUSIVE STATEMENT OF AGREEMENT BETWEEN THE PARTIES AND SUPERSEDES ALL PROPOSED OR PRIOR AGREEMENTS, ORAL OR WRITTEN, AND ANY OTHER COMMUNICATIONS BETWEEN THE PARTIES RELATING TO THE LICENSE DESCRIBED HEREIN.

## Care and Cautions

---

## Overall Warnings

Remember that damage caused by misuse is not covered by your warranty. It is critical to understand and comply with the following guidelines:

### Procedures

Do not attempt to intubate without lubricating the airway adjunct with silicone lubricant (provided). Failure to do so will make intubation very difficult and is likely to result in damage.

When simulating drug administration via endotracheal tube, providers must use an empty syringe. Passing liquids into the trachea or esophagus may cause internal damage.

Mouth to mouth resuscitation without a barrier device is not recommended, as it will contaminate the airway. Treat the simulator with the same precautions that would be used with a real patient.

Always adjust the motor arm prior to every delivery.

Always keep clear of the birthing mechanism while the system is on.

Never operate the birthing mechanism without the tummy cover in place.

### NOELLE IV arm

Only use Gaumard's provided simulated blood. Any other simulated blood containing sugar or any additive may cause blockage and/or interruption of the vasculature system.

The use of needles larger than 22 gauge will reduce the lifetime of the lower arms' skin and veins.

When the arm veins require replacement, contact Gaumard to arrange for a lower arm exchange. For a small fee, we will deliver reconditioned and warranted lower arm assemblies to your facility. After receiving the replacement arms, use the same box and the enclosed shipping label to return the old arms to Gaumard. For international and express service, additional fees may be charged. Refer to the Consumables and Replacement Parts section of this guide, and contact customer service for more information.

### Storage

Store NOELLE in a cool, dry place. Extended storage above 85 degrees Fahrenheit (29 Celsius) will cause the simulator to soften and slowly warp. It is acceptable to *operate* NOELLE at an ambient temperature of 95 degrees Fahrenheit (35 Celsius).

Do not store the simulator with a discharged battery. It is good practice to re-charge the battery at the end of every simulation session. In addition, make sure the battery is re-charged at least once every 6 months even if the simulator is not being used; otherwise permanent loss of capacity might occur because of self-discharge.

### Cleaning

The simulator should be cleaned with a cloth dampened with diluted liquid dishwashing soap. If medical adhesives remain on the skin, clean with alcohol wipes. **DO NOT USE "GOO GONE"** as the citric acid in the formula will cause pitting of the various materials comprising your simulator.

NOELLE is "splash-proof" but not water-proof. Do not submerge or allow water to enter the interior of the simulator. Do not expose the tablet computer to water or excessive dust.

### Set Up

NOELLE will only power on when connected to the power supply.

NEVER disconnect the communications module while the GaumardUI software is running. The software will halt, and the module may be damaged.

### Birth Canal Maintenance

Ball point pens, ink and markers permanently stain the birth canal insert.

Do not wrap this or any other Gaumard product in newsprint.

The birth canal insert can be cleaned by wiping with a mild solution of soap and water. After cleaning, dust with talcum powder.

Store the simulator in a cool, dry place.

After exercise is completed, **DO NOT** leave birthing baby in contact with the birth canal.

### Post-partum hemorrhage check list and Warnings

Use only Gaumard's provided simulated blood. Any other simulated blood brand containing sugar or any additive may cause blockage and/or interruption of the vasculature system.

At the end of every simulation, always flush the system with distilled water to prevent clogging.

Always position the simulator so post-partum hemorrhage fluid flows away from the birth canal and the simulator itself.

Do not allow PPH fluid to puddle beneath the simulator or reach the lower back.

To prevent staining or molding, always clean NOELLE using diluted soap and water. Remove the birth canal and clean thoroughly.

### Labor

Always lubricate the fetus and the birth canal before every delivery. Failure to do so will result in damage to the birthing mechanism and the birth canal.

Do not pull the baby upward in contrast to the birthing mechanism's linear trajectory. Doing so can bend the motor arm and cause damage to the birthing mechanism.



## Getting Started

---

## Overview

The NOELLE S554.100 is a comprehensive package of simulator technology, scenario-based training, and performance measurement and debriefing tools designed to build competencies needed to help manage OB emergencies.

### Highlights

- Prepare for OB emergencies; evaluate and report on training and clinical outcomes
- Full size NOELLE maternal and neonatal birthing simulator with eclampsia and hemorrhage capabilities
- Mobility allows training to be in L+D and postpartum units
- Build team and technical competencies
- Maternal vital signs monitor
- Fetal heart tones and neonatal vital signs monitor
- Set up and run OB emergency simulations for shoulder dystocia, postpartum hemorrhage, eclampsia, umbilical cord prolapse, breech vaginal delivery, operative vaginal delivery, neonatal resuscitation

### Genuine NOELLE Simulator

- IV arms for meds/fluids
- Intubatable airway with chest rise
- Removable stomach cover
- Programmable eclampsia
- Advanced birthing mechanism
- Programmable postpartum hemorrhage
- Birthing fetus with placentas and umbilical cords

### Newborn

- Full term intubatable newborn with cyanosis and umbilical pulse
- Chest compressions and ventilations are measured and logged
- Realistic heart and lung sounds
- Realistic crying

### Simulated vital signs monitor

- Single large 23 inch touchscreen monitor with desktop controller and wireless communication to laptop
- Display up to 8 numeric values including HR, ABP, RR, CO2, SpO2, temperature, NIBP, and time
- Select up to 5 dynamic waveforms including ECG II, ABP, respiration, CO2, and pulse oximetry
- Fetal Heart Rate monitor

### Other

- Simulation transport case for electronic products
- FCC, IC, CE Certifications
- One year warranty; extend to three years
- Installation and training available

## Terminology

**Apply** - In the context of a simulation, to apply settings is to send details of the patient's condition to the simulator itself. When settings are successfully applied, NOELLE's condition should match that shown on the Status panel.

**GaumardUI** - the Gaumard User Interface software application, used to control the simulator and evaluate care providers.

**Facilitator** - the person conducting the simulation; an instructor or lab staff member.

**Palette Item** - Any full or partial set of physiological parameters saved together under a single name.

**Profile** - a unique GaumardUI configuration, including custom palettes, scenarios, and options. Each profile acts as a separate program whereby changes made to one profile have no effect on the others.

**Provider** - a person participating in the simulation as a healthcare provider.

**Scenario** - a saved sequence of physiological states, which flow like a "play list." Scenarios provide a level of automation that unburdens the facilitator and allows standardized presentation of symptoms.

**Scenario Item** - a Palette Item that is part of a scenario. Scenario Items may also represent a fixed delay period ("Wait") or a pause ("Wait Indefinitely").

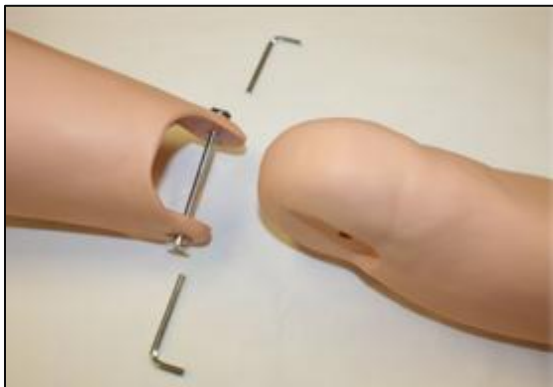
## Equipment Set-up

---

## Leg Assembly

Follow the steps below to install the lower legs. Always, remove the lower legs when transporting NOELLE inside the protective case.

1. Remove the fixed bolts from the knee joints using the hexagonal wrench included.




2. Position the lower legs and insert the bolts. Use the two provided hexagonal wrenches to secure the knee bolts. Do not over tighten.



## Power Supply

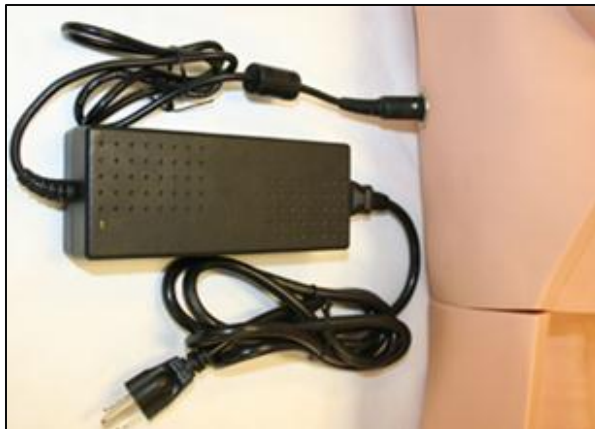Connect the AC adapter labeled NOELLE to the power input located on NOELLE's right side.

**Always operate NOELLE with the power supply connected.**



The NOELLE is equipped with an internal backup battery for transporting during scenarios. Always reconnect NOELLE to the charger after arriving to the designated simulation area.

The backup battery level is displayed on the GaumardUI status panel. For more information about the battery indicator, refer to page 18.

Newborn is not equipped with an internal battery.

**WARNING:** Do not store the simulator with a discharged battery. Re-charge the system at least once every 3 months if the simulator is not being used; otherwise permanent loss of capacity might occur.

## Control Computer

Control commands are sent to the simulator via a wired USB communication module. Follow the steps below to connect the wired communication module to the simulator and the control laptop computer.

GaumardUI controls only one simulator at a time. To switch between NOELLE and Newborn, first exit the GaumardUI software and reconnect the communications cable to the simulator you wish to control.

**WARNING:** Do not connect the simulator to Ethernet cards, LAN networks or unauthorized diagnostic equipment. Doing so may cause damage to the system.

1. Connect the Ethernet cable to the USB communication module.



2. Connect the USB communication module to an available USB port and power on the Micro+ computer.



3. Connect the Ethernet cable to the communications port located on the simulator's right side.



## Virtual Monitor

The mobile virtual monitor system works with GaumardUI to display the vital signs of NOELLE, birthing baby and Newborn.



## Equipment installation

Refer to the manufacturer's documentation included with the virtual monitor system components for important safety, installation and start-up information before turning on the computer.

### Virtual monitor wireless Connectivity

The Micro+ laptop and the all-in-one virtual monitor computer **automatically** establish a Wi-Fi connection at **startup**. To establish the connection, first power on the Micro+ laptop followed by the virtual monitor computer.

The simulator's vital signs are generated by the GaumardUI software and sent to the Gaumard Monitors software via Wi-Fi. To verify the connection between the computers, click on the **wireless icon** located on the task tray. The wireless network name is randomly generated at the factory. To troubleshoot connection issues, please navigate to page 181.



## Gaumard Monitors

After the Wi-Fi connection is established, double click or tap the **Gaumard Monitors** icon located on the virtual monitor's home screen.

The Gaumard Monitors software is now ready to receive vital signs information once GaumardUI is initialized.



Continue to the next section to begin working with the simulator and the GaumardUI control software.

## Working with GaumardUI

---

# Initializing the simulator

After reading the care and cautions information included with the system components, power on control computer.

Following a brief loading screen, the simulator selection screen is shown. Select the simulator connected to the communication module and click **start**. The simulator will power on in approximately 30 seconds. To troubleshoot startup issues, go to page xxx.

**GaumardUI controls one simulator at a time. To switch between NOELLE and Newborn, first exit the GaumardUI software and reconnect the Ethernet cable to the simulator you wish to control.**

## User Profile Menu

Each profile stores an independent library of customizable palettes, scenarios, and options settings. Changes made to items in one profile do not affect items stored in others.

The **Quick Start** profile was created in conjunction with experienced healthcare instructors and working medical professionals. The Quick Start NOELLE profile contains prebuilt labor scenarios for birth simulation and post-partum complications. It serves a convenient starting point and it can be customized to fit most simulation objectives.

Select a profile and click **Load** to continue.

- **Quick Start NOELLE** – Contains a total of twelve labor scenarios.
- **NOELLE Advanced** – Contains fourteen linear scenarios, twenty-two labor scenarios and one branched scenario.
- **Quick Start Newborn** – contains a total of eleven scenarios.

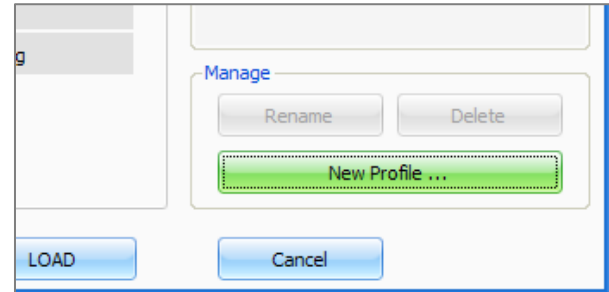
For more information on the items included on the Quick Start profiles, navigate to Page 115. To change profiles from inside the GaumardUI environment, select Profiles from the File menu.

Profiles are used to organize and protect software settings. As you begin to customize NOELLE, it will become clear how profiles can best serve your needs. For example:

- It may be appropriate to assign one profile to each user of your NOELLE system.
- Others may choose to create a profile dedicated to a specific academic course, which might be taught by multiple instructors.

- For the most detailed exercises, it is sometimes useful to devote an entire profile to one particular subject area, or even one particular scenario.

To create a user made profile, click on the **New Profile** button on the right panel.



Enter a name for the new profile followed by a description.

To include scenarios and palettes from other profiles, click the applicable check box. For security, enable PIN protection, which will require a user to enter a four-digit key before loading the protected profile.

Finally, click **create** to save the new profile. To import and export profiles, navigate to Page 82.



# The Environment

The GaumardUI environment is used by the facilitator as the simulation command center. In addition to managing the simulator's functionality and vital signs, the facilitator can evaluate student performance and act as the simulator's voice. This powerful combination of tools provides the facilitator with an indispensable tool for simulation, training and testing.

## The Status Panel

The Status panel is visible along the left edge of the GUI window at all times. Backup battery level information, signal strength, volume levels and vitals are displayed in real time. Vital sign parameters such as heart rate, blood pressure and respiratory rate are highlighted in yellow while undergoing change.

Click the **arrow** to expand the viewer and access status entries for **Other**, **Uterine Activity** and **Fetal Heart rate**.

Some status panel entries might not be visible depending on the simulator's hardware configuration.

**Status** [Battery] [Signal Strength] [Expand]

**Breathing**  
Resp. Pattern: **normal**  
Rate: **13** bpm In: **33** %  
(O2-Sat : **98** %)  
(EtCO2 : **40** mmHg)

**Circulation**  
Cardiac Rhythm: **Sinus**  
Heart Rate: **80** bpm  
Blood Pressure: **120 / 80** mmHg

**Other**  
(Temperature **37.5** °C)  
Seizures: **none**

● Hemorrhage  
● Uterine Hemorrhage  
UP: **0** %

**Expanded Panel**    **Automatic Mode**

**Status** [ABC] [More] [Expand]

**Breathing**  
Resp. Pattern: **normal**  
Rate: **13** bpm In: **33** %  
(O2-Sat : **98** %)  
(EtCO2 : **40** mmHg)

**Circulation**  
Cardiac Rhythm: **Sinus**  
Heart Rate: **70** bpm  
Blood Pressure: **120 / 80** mmHg

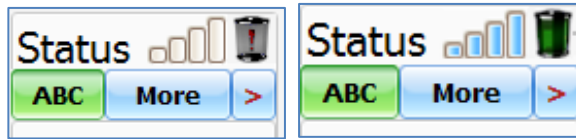
**Other**  
(Temperature **37.5** °C)  
Seizures: **none**  
● Hemorrhage  
● Uterine Hemorrhage  
UP: **0** %

**Parameter undergoing change** (points to yellow highlights in both panels)

**Active/inactive feature indicator** (points to Hemorrhage/Uterine Hemorrhage indicators in both panels)

## Communication indicator

The communication indicator shows the status of the link between the computer and the simulator. The indicator is clear when no attempts to communicate with the simulator are being made; for example when the wired communication module is not connected to the computer or the if system is in STAND-BY mode.




## Backup battery indicator

The battery status indicator updates as the backup battery in the simulator is used. The exclamation mark indicator is shown when there is no communication with the simulator and the program cannot retrieve battery information from the simulator.




## Volumes

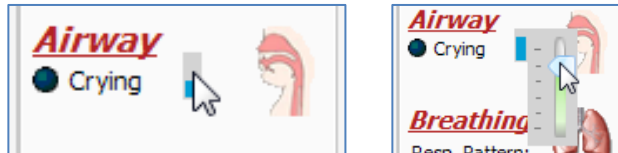The volume levels for sound features on the simulator are adjusted on the Status Panel. Click on the volume level indicator next to each sound level to adjust.



## Soft Power and Standby

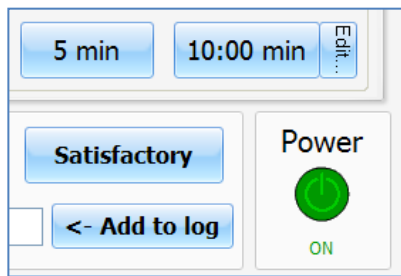The standby button is located on the bottom right corner of the GaumardUI software. Use the standby mode to save backup battery power between exercises.




The simulator will establish a connection less than one minute after resuming from standby. The connection bars will confirm the link between the controller computer and the simulator. **The wired communication module must be connected when resuming from stand by.**

## Details

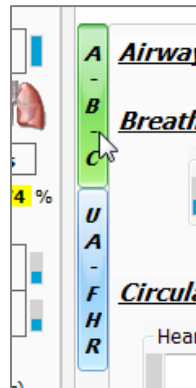From the **Details** tab, facilitators can manually control the simulator's vitals, enable/disable features and create **palette items** later used to build **scenarios**. Changing controls directly from the details page is the simplest form of control available to the facilitator.

Vital controls are divided in several categories. For each vital control on the details tab, there is a corresponding entry on the **Status Panel** displaying real time information.


Additional controls can be accessed by clicking on the **vertical tab** located on the left of the screen. Controls on the details tab vary from simulator to simulator. GaumardUI will automatically detect features and upgrades installed on the simulator and display the corresponding controls.




### Changing Vitals

Set the changes on the available fields, and use the **Apply panel** located at the bottom of the Details tab area to submit. GaumardUI only sends updated vital information to the simulator, **settings not specified will remain unchanged**. Outlined below are common control behaviors.


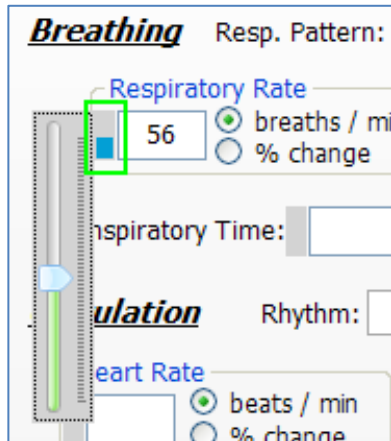
Click or tap **slider control** to quickly adjust numerical values using the track pad. Alternatively, use the keyboard to manually type a numerical value.



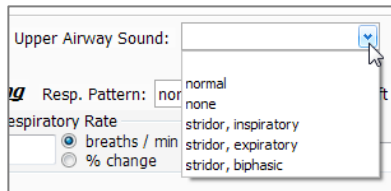Click on the **down arrow** to access sound types and rhythms.



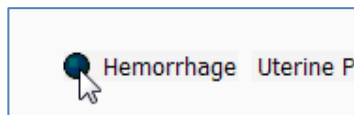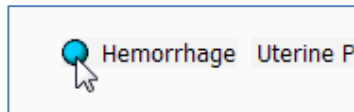
Click on the **feature control button** to enable (blue) or disable (black) features such as hemorrhage, chest rise, and crying. Then, click Apply to commit the state.



For more information about the simulator's features, navigate to the **Working with NOELLE/Newborn** section of this guide.

## Apply Panel

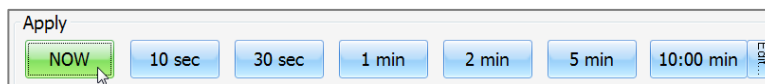
After the settings are selected in the details area, changes are submitted using any of the apply buttons located near the bottom of the page. Please note that **settings not specified will remain unchanged**.



Click the **NOW** button to change the simulator's condition instantly. Alternatively, click a **trending time** to gradually increase or decrease to the numerical value specified (e.g. heart rate, blood pressure) in the time allotted. The right-most button can be customized to any transition time by clicking the part of the button labeled "Edit".


Some settings are applied immediately, such as cardiac rhythm and breathing pattern, while numerical settings, such as heart rate and respiratory rate can trend.

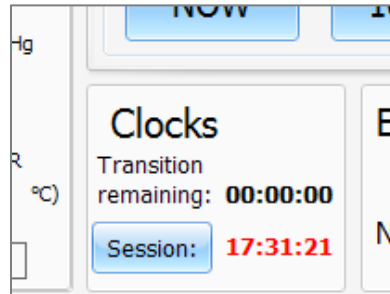As transitions are applied, the trend countdown is displayed in the **Clocks** panel at the bottom of the GaumardUI window. If there is already an ongoing transition at the moment you click an Apply button, it will stop, and a new transition will begin from the current physiological state.



## Ventricular Tachycardia

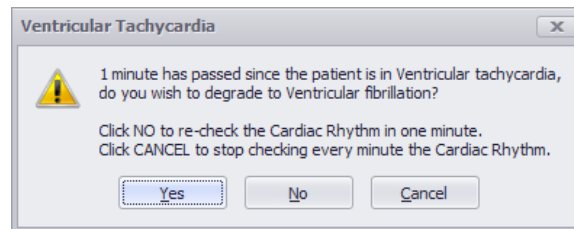
When cardiac rhythm is maintained at Ventricular Tachycardia for one minute, the following prompt dialog box will be displayed:



Click **Yes** to degrade the cardiac rhythm to **ventricular fibrillation**. Click **No** to maintain the cardiac rhythm and recheck after 1 more minute. Click **Cancel** to stop the software from checking the cardiac rhythm every minute.


## Uterine Activity/ Fetal Heart Rate

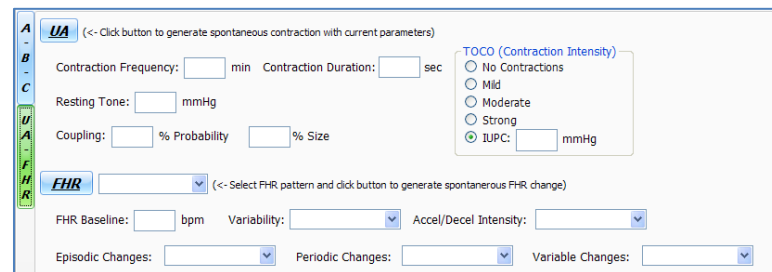Click on the **UA-FHR** vertical tab to access uterine activity and Fetal Heart Rate controls.



**UA** - The instructor can generate a uterine contraction at any time by clicking the **UA** button; the contraction will have the same duration and intensity settings applied on the status viewer.

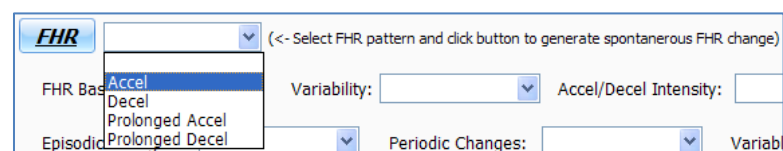**FHR** - Generate a fetal heart rate acceleration or deceleration at any time by selecting one of the options in the drop box below the button, and then clicking on the button.




## Palette Items


A **Palette** item stores several vital parameters in a single container. Use palette items to quickly change multiple parameters at once. For example, a facilitator may choose to create a palette item that stores all the vital parameters of a healthy patient. Then at any moment in the simulation quickly return to a healthy state by loading the palette item.

Palette items are also used to build scenarios as covered on page 30. As part of the scenario building process, a facilitator will first create a series of individual palette items which represent physiological states. The **scenario** feature then automates the loading of palettes in succession to simulate patient recovery or complications.

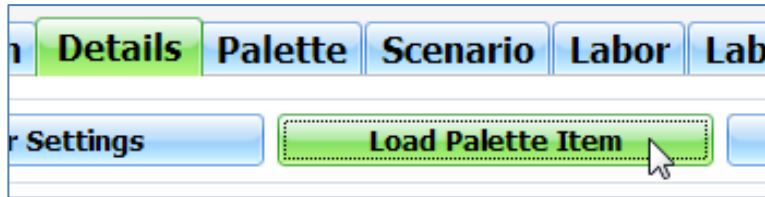A collection of palette items are included in the quick start profile. Click the **Load Palette** button to load a prebuilt palette item.




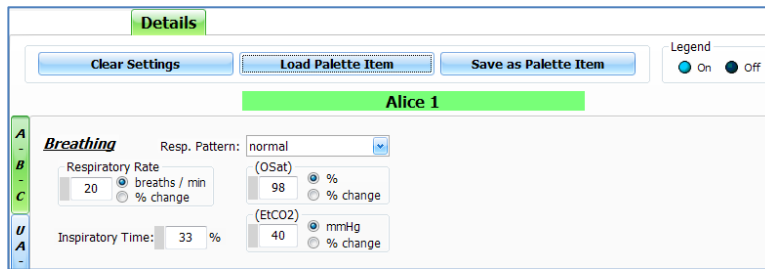Select a palette item from the list and click **Load**. Finally, click apply to submit the changes.




## Creating Palette Items

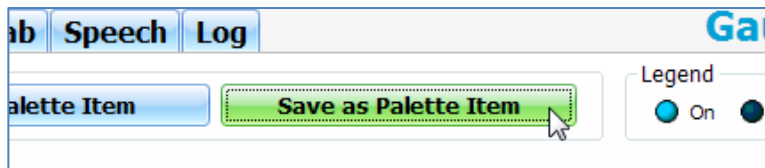To create a Palette Item, adjust the vital parameters and features on the Details tab and click the **Save as Palette Item** button near the top of the page.



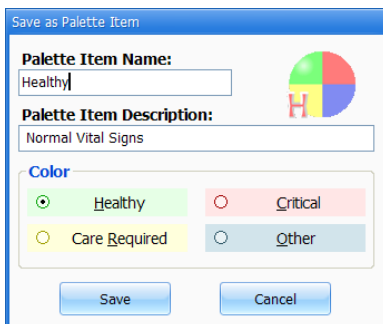Enter the palette's name, description, and color code and click **Save**.



## Palette

The second level of control is the Palette tab-page. Each item on the Palette represents a complete or partial physiological state. The Palette page displays all of the Palette Items in the active profile. Each profile has its own separately customizable Palette. Create Palette Items with the Details page, as described previously.

Name	Description
Alice 1	Admission
Alice 2	Transitional Labor
Alice 3	Second stage
Alice 4	Delivery
alice01	Stage 1 - Active
alice02	Stage 1 - Transition
alice03	Stage 2
alice04	Stage 3
Alice 1	Admission
Alice 2	Transitional Labor
Alice 3	Second stage
Alice 4	Delivery
Amy 1	Admission
Amy 2	Early Labor
Amy 3	Active Labor
Amy 4	Transition
Amy 5	Pushing

**View**

Sort by Name  
Healthy  
Care Required  
Critical  
Other

**Selected Item**

Edit...  
Delete  
Properties

**Apply**

NOW 10 sec 30 sec 1 min 2 min 5 min 10:00 min

Apply Palette Items using the buttons at the bottom of the page, exactly as changes to NOELLE's condition are applied on the Details page. Change NOELLE's vital signs and symptoms instantly by clicking to select a Palette Item and clicking the "NOW" button. Alternatively, create a gradual transition in physiological state with one of the other Apply buttons.

Palette Items can be sorted with the "View" buttons found on the right side of the page.

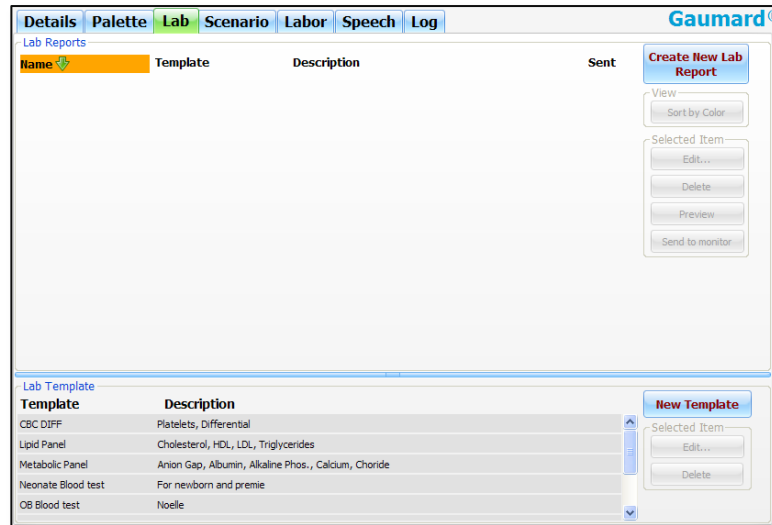
Editing existing Palette Items is simply a matter of selecting the item you wish to modify and clicking the Edit button. You will be automatically be taken to the Details page, and the settings that comprise the selected Palette Item will be displayed. Change them as desired, and click the "Save as Palette Item" button.

Many of the functions related to Palette Items are also available by clicking the second mouse button (usually the right button) while the pointer is positioned over an Item. Note that when using the tablet computer, this is best done by holding the stylus button while tapping the screen.

For more information on customizing the Palette, see the Tips on Palette Item and Scenario Creation section of the [Appendix](#) at the end of this guide.

## Lab


Laboratory tests are helpful tools for evaluating the health of a patient. To simulate this process, the Lab Tab allows the facilitator to create laboratory reports to aid providers during simulation. Once a laboratory test is created, the facilitator can display the results on the Virtual Monitor window for the provider to utilize.



The Gaumard Virtual monitors must be connected before using the file sharing feature. To verify the connection navigate to PAGE 91.

## Creating a Lab Template

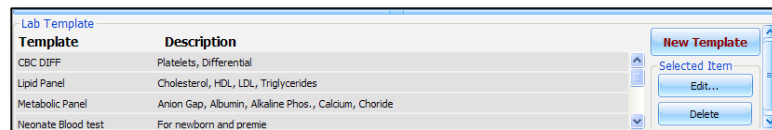To begin, navigate to the bottom of the tab and click on the **New Template** button on the right panel.





The **Edit Lab Results Template** window is used to create templates that will be later used to create lab reports. Enter a **name** for the lab template followed by a **description**.

Categories are used to group a series of tests in a lab report. From the **Category** box, click **Add** to name and create a category. Then, click **OK** to save the new category.

Use the **Single Entry** menu to add a test under the category previously created.

Variable	Normal Range	Value	Comment

Use the **Add/Edit Entry** menu to customize the different parameters on a specific test. Begin by providing the name of the test, unit and decimal precision. The **Normal Range** will be displayed on the lab report for the provider to use as guide while reading the results. Meanwhile, the **Allowed Range** restricts the minimum and maximum value a facilitator can input as a test result. Once the test entry is configured, click **OK** to add the new test.

**Add/Insert Entry...**

**Name:** HDL Cholesterol **Unit:** mg/dL

**Precision:** 2

**Normal Range**  
**Minimum:** 35.00 **Maximum:** 85.00

**Allowed Range**  
**Minimum:**  **Maximum:**

OK Cancel

Repeat the process to add more tests and categories using the Category and Single Entry menu. Please note that individual items can be moved, deleted or modified after they are created. After all the tests are entered, navigate to the bottom of the page to **Save** the new template.

**Edit Lab Results Template**

**Name:** Lipid Panel

**Description:** Cholesterol, HDL, LDL, Triglycerides

**Category** **Single Entry**

Add Remove Rename Insert Add Delete Edit Insert

**Lipid Panel - Female**

Variable	Normal Range	Value	Comment
HDL Cholesterol(mg/dL)	35.00--85.00		
Cholesterol/HDL(Ratio)	0.00--4.40		

**Lipid Panel - Male**

Variable	Normal Range	Value	Comment
HDL Cholesterol(mg/dL)	29.00--71.00		
Cholesterol/HDL Ratio(Ratio)	0.00--4.90		

**Lipid Panel - Both**

Variable	Normal Range	Value	Comment
Cholesterol(mg/dL)	0.00--199.00		
LDL Cholesterol(mg/dL)	0.00--129.00		
Triglycerides(mg/dL)	30.00--200.00		

\*Allowed value range is displayed in Value column. It won't be visible when entering data

Save Cancel

Once a new format is created, it will be listed on the Lab Format section at the bottom of the Lab tab. Use the buttons on the left panel to edit or delete lab formats.

**Lab Template**

Template	Description
CBC DIFF	Platelets, Differential
Lipid Panel	Cholesterol, HDL, LDL, Triglycerides
Metabolic Panel	Anion Gap, Albumin, Alkaline Phos., Calcium, Chloride
Neonate Blood test	For newborn and prenie

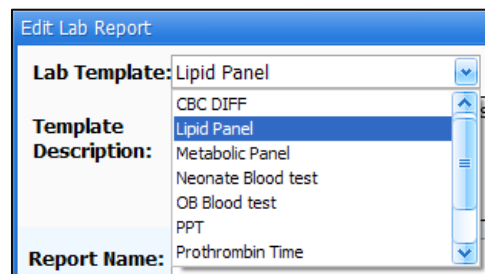
New Template Selected Item Edit... Delete

## Creating a Lab Report


To begin, click on **Create a New Lab Report** from the right panel.



The **Edit Lab Report** window is used to prepare new lab report. First, select a **Lab Template** from the drop down menu.



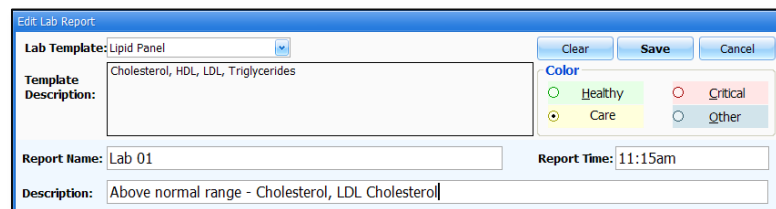Provide a **Report Name**, **Report Time** and **Description**. In addition, select a condition color tag for the lab report on the right panel. Color tags aid the sorting of lab reports on the report list window.



Input the results of the tests on the **Value** column. Values above the normal range specified will be displayed in bold. Include any comments associated with the test performed. Finally, click **Save** from the top right menu to create the lab report.

<b>Description:</b> Above normal range - Cholesterol, LDL Cholesterol			
<b>Lipid Panel - Female</b>			
Variable	Normal Range	Value	Comment
HDL Cholesterol(mg/dL)	35.00--85.00	85	
Cholesterol/HDL(Ratio)	0.00--4.40	5.1	Care required
<b>Lipid Panel - Male</b>			
Variable	Normal Range	Value	Comment
HDL Cholesterol(mg/dL)	29.00--71.00		
Cholesterol/HDL Ratio(Ratio)	0.00--4.90		
<b>Lipid Panel - Both</b>			
Variable	Normal Range	Value	Comment
Cholesterol(mg/dL)	0.00--199.00	190.00	
LDL Cholesterol(mg/dL)	0.00--129.00	135	
Triglycerides(mg/dL)	30.00--200.00	100	

The newly created lab report will be listed on the Lab Reports list. Lab reports can be sorted by name, template, and description, sent status or colors.

Name	Template	Description	Sent
Healthy 2	OB Blood test	Healthy Levels	
Healthy	Neonate Blood test	Healthy Levels	
James	Metabolic Panel	Session 1	
Susan	Lipid Panel	Session 4	
Lab 01	Lipid Panel	Above normal range - Cholesterol, LDL Cholesterol	
Lab 04	CBC DIFF	Simulation 1	
Lab 02	CBC DIFF	High levels	
Lab 03	Neonate Blood test	High levels	
Lab 05	OB Blood test	Session 2	
Lab 06	PPT	Session 2	

Select a lab report and click the **Preview** button on the right panel to view the report on the facilitator's screen. From the Lab Preview window, the facilitator can print as well as send the report to the virtual monitor computer. To make changes, click **Close** and then **edit**.

**Lab 01**

Report Time: 11:36am  
Description: Above normal range - Cholesterol, LDL Cholesterol

---

**Lipid Panel - Female**

Variable	Normal Range	Value	Comment
HDL Cholesterol(mg/dL)	35.00--85.00	85.00	
Cholesterol/HDL(Ratio)	0.00--4.40	5.10	Care required

---

**Lipid Panel - Both**

Variable	Normal Range	Value	Comment
Cholesterol(mg/dL)	0.00--199.00	190.00	
LDL Cholesterol(mg/dL)	0.00--129.00	135.00	
Triglycerides(mg/dL)	30.00--200.00	100.00	

The preview window also allows the printing of results for distribution and archiving. To make changes, click **Close** and then **edit**.

## Send to monitor

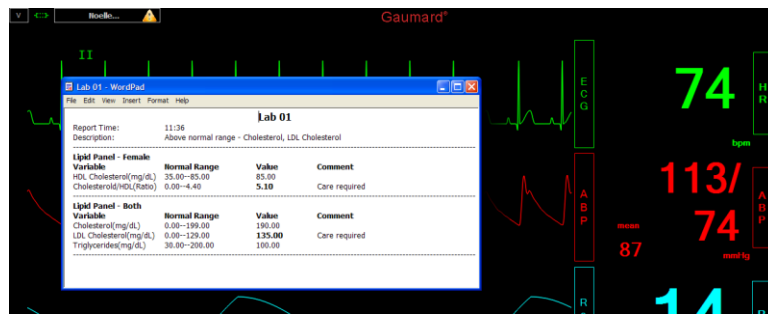
Begin by selecting the report from the lab reports list.

Click the **Send to Monitor** button to transfer the lab report to the Virtual Monitors.

On the Gaumard Monitor's window, an exclamation icon notifies the provider a file is ready for access. Click the NOELLE drop down menu to select the available lab report.



The lab report will open using the system's default application.



Once the report is sent, the letter Y will be present on the Sent column. Use the Stop Sharing button on the right panel to remove any items on the Gaumard Monitor file list.

Lab Reports				Create New Lab Report	
Name	Template	Description	Sent	View	Sort by Color
Healthy 2	OB Blood test	Healthy Levels	Y	Selected Item Edit... Delete Preview Stop sharing	
Healthy	Neonate Blood test	Healthy Levels			
James	Metabolic Panel	Session 1	Y		
Lab 01	Lipid Panel	Above normal range - Cholesterol, LDL Cholesterol	Y		
Lab 02	CBC DIFF	High levels			
Lab 03	Neonate Blood test	High levels	Y		
Lab 04	CBC DIFF	Simulation 1			
Lab 05	OB Blood test	Session 2			
Lab 06	PPT	Session 2			
Susan	Lipid Panel	Session 4			

An editable copy of the lab report is also copied onto the Gaumard\_UI folder on the tablet's home screen. For information on how to access other files from the Gaumard Monitor screen, navigate to page 89.

## Scenarios

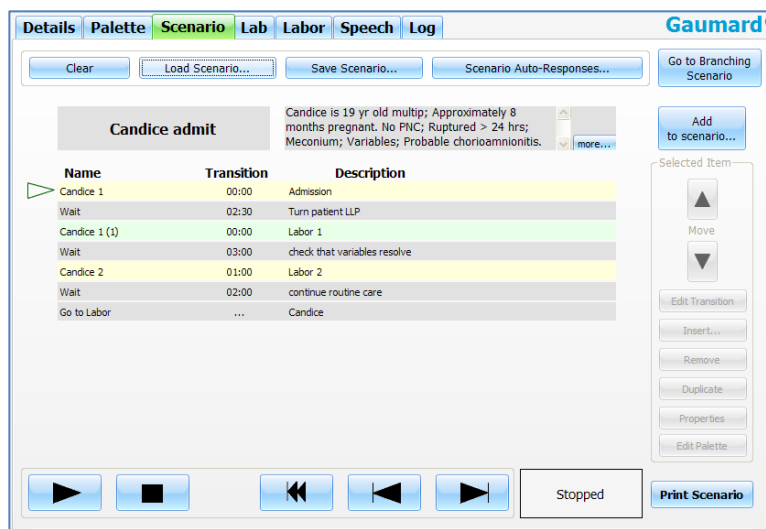
The most advanced method of controlling the system is to build a Scenario, a sequence of Palette Items and delay periods. Scenarios can be linear or branching as covered in the following section. Think of a scenario as a "playlist" of palette items. Consistent with this analogy, scenario controls at the bottom of the page look and behave just like traditional and software-based media players.

Scenarios let the facilitator automate most of the changes to the simulators 's condition, so their attention can remain on the providers' actions. The scenario system can also provide standardization of the patient's presentation of symptoms. For fair assessment of providers and any research application, such standardization is key. For tips on planning and creating scenarios, go to page 179.

NOELLE linear scenarios simulate pre partum and postpartum events. Scenarios that incorporate the labor mechanism are loaded using the Labor tab (page 46). In the example below, the simulation exercise begins during the admission process and continues onto the delivery when the scenario reaches the last item on the list.


## Linear Scenarios

Linear scenarios consist of palette items added in sequence with specific transition times as shown in the figure below.




## Using Factory Preset Scenarios

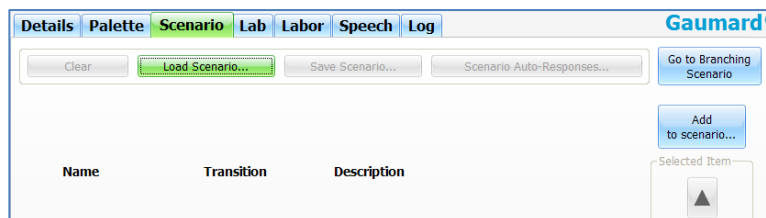
While operating in the NOELLE Advanced profile, go to the Scenario tab and Click on **Load Scenario**. On the Newborn software, load the Quick Start Profile.



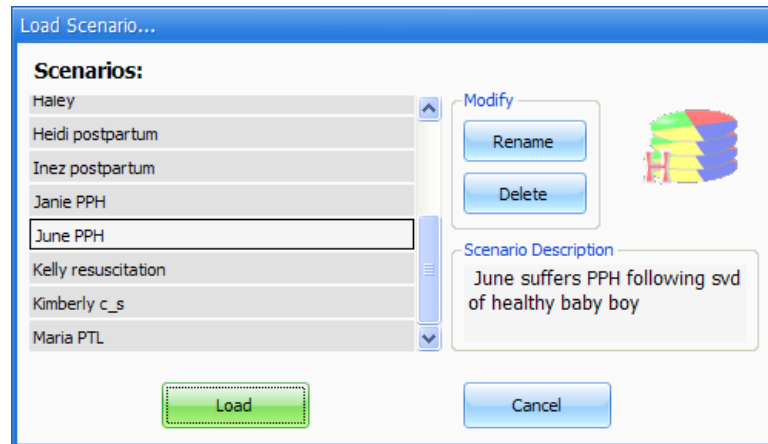Select the desired scenario and click on “Load”.



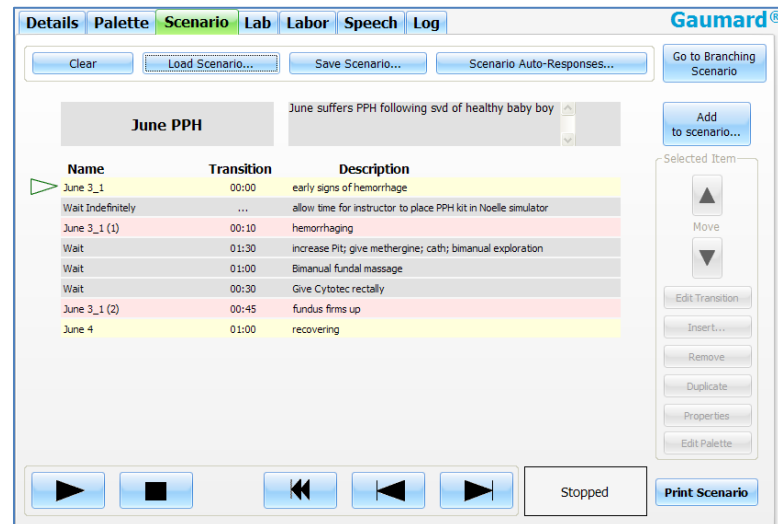After loading the desired scenario, click the **Play** button and watch the vital signs adjust according to the specifications of each palette.



## Scenario Controls

Scenarios are controls are located at bottom of the Scenario page. The same way a music player plays songs, the Scenario plays palette items. Intuitively, the facilitator can play, stop, pause, skip, or repeat items as appropriate. The Scenario Position Indicator points to the current item and shows the current status of the scenario. The following paragraphs describe in detail the behavior of each button and indicator.

### The Scenario Position Indicator



An unfilled triangle means that the scenario is stopped. When the Play button is clicked, the item pointed by the indicator is played.



A rapidly blinking triangle means that the scenario is playing the item to which the indicator is pointing.



A slowly blinking triangle means that the scenario is paused at the item to which the indicator is pointing.

### The Scenario Buttons



Plays the item to which the scenario position indicator is pointing. This button has two states: play or pause.



Pauses the scenario. This state of the play button is only active when the scenario is playing. It is disabled when a 'Wait indefinitely' item is playing because in such case the scenario is already paused.



The Stop button has 2 behaviors depending on when it is clicked. When clicked once, the Stop button halts the scenario at the end of the currently playing item. When clicked a second time, the scenario is stopped immediately. For example, if the currently playing item has a transition of 1:00 minute and the Stop button is pressed when it has 0:10 seconds left, the scenario will be halted at the end of the transition (i.e., in 10 seconds). If the Stop button is clicked again within those remaining 10 seconds, the scenario stops immediately.



The Next button advances the indicator to the next item on the scenario regardless if the scenario is playing, paused, or stopped. It can also be used to move the indicator to select an item before playing it.



Similar to the Next button, the Previous button returns the indicator to the previous item in the scenario.



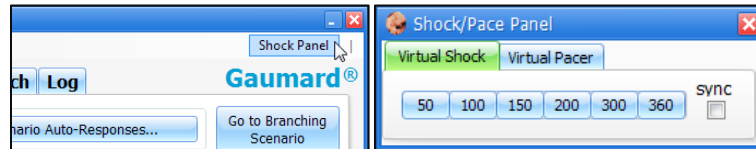
The Reset button stops the scenario immediately and returns the indicator to the first item in the scenario.




## Scenario Auto-Responses

The scenario auto-responses feature assists facilitators in automating electrical therapy responses during a scenario. First, click on the shock panel icon located near the top right of the screen to access the virtual **Shock /Pace** panel.

If the shock panel button is not displayed, go to **Setup>Options>General** and checkmark **Show defibrillation/cardioversion panel**.




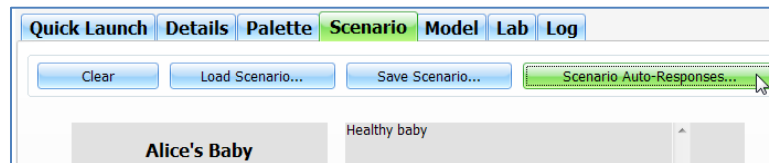Click on the **Scenario Auto-Responses...** button to access and configure the auto response feature.



During a wait indefinitely palette, the auto response feature will advance the scenario to the next palette item in response to the following interventions.

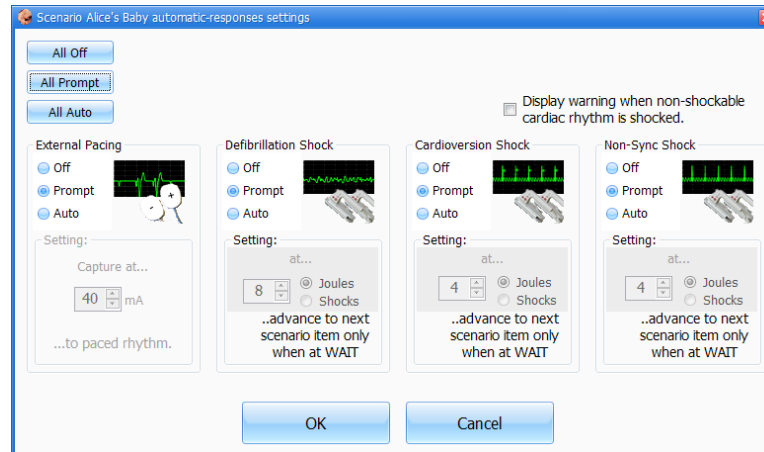- External Pacing -
- Defibrillation Shock – During shockable rhythms only
- Cardioversion Shock –
- Non-sync Shock -



The three auto response detection modes are defined below.


- **Off** - The software **does not** respond to the electric therapy.
- **Prompt** - The software detects the electrical therapy and prompts the facilitator if the scenario should proceed to the next palette.
- **Auto** - The software detects the electrical therapy and compares it to the joules or shock threshold configured by the facilitator. If the electrical therapy is equal to or greater than the threshold specified, the scenario will automatically continue without prompting the facilitator.

## Scenario Quick Launch

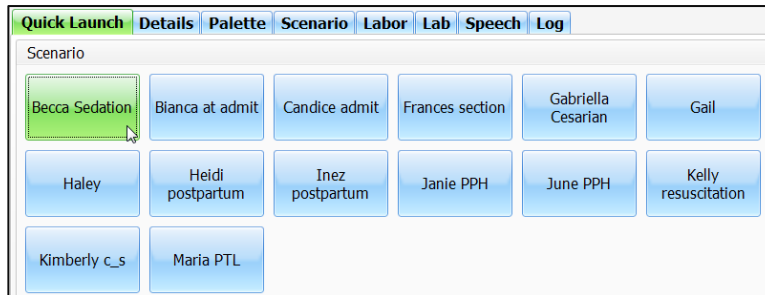
Use the scenario quick launch window to start scenarios with a single click. The Quick Launch window displays all the scenarios saved in the current profile.

### Linear Scenario Quick Launch

Click the scenario button to start the scenario.



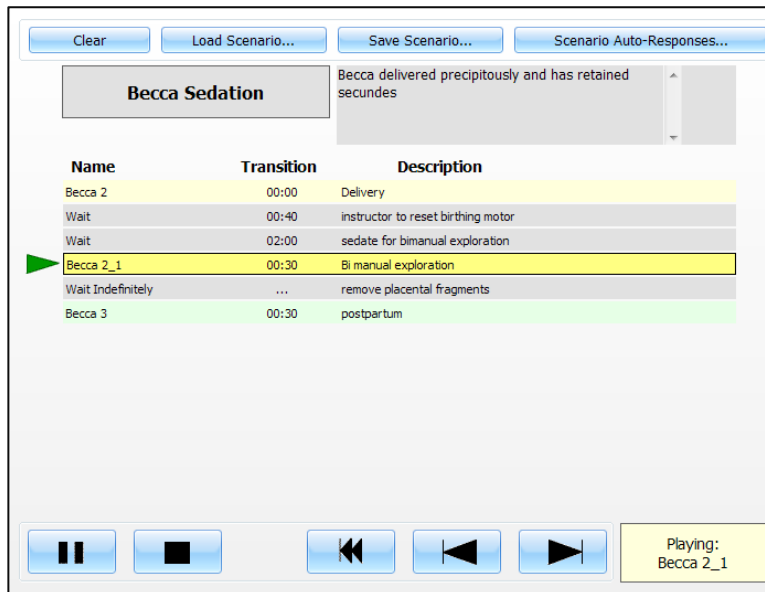
After the scenario button is clicked, GaumardUI will automatically start the scenario.



### Labor scenario Quick Launch

To start a labor scenario from the quick start tab, first specify the position of the birthing baby on the birthing mechanism by clicking on the corresponding labor position radio button.

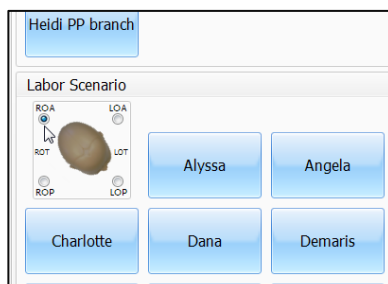
**Before using the labor quick launch feature for the first time, go to page 99 for important information on preparing NOELLE for delivery.**



After the position is set, click on the labor scenario to start.



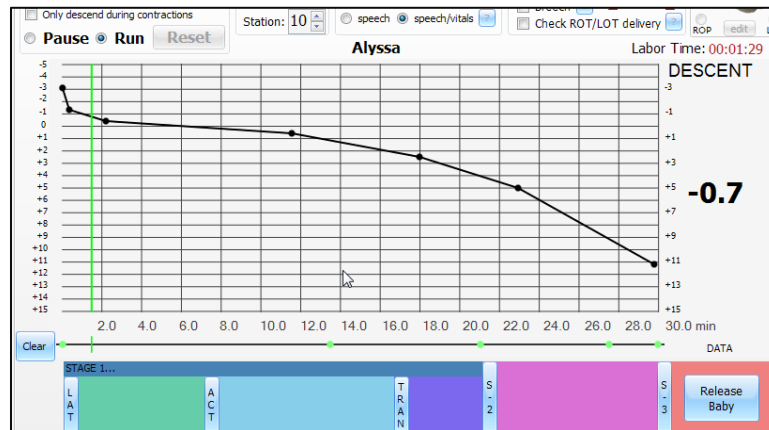
GaumardUI will automatically start the labor scenario and switch to the labor tab.




## Factory Preset Scenarios

A powerful way to supplement the labor feature is to load factory preset scenarios. To locate and load these scenarios follow the steps below.

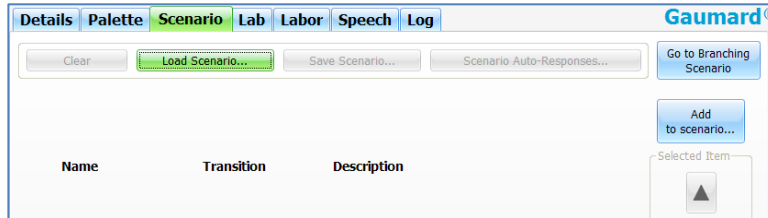
While operating under a Quick Start Profile, go to the Scenario tab and Click on **Load Scenario**.



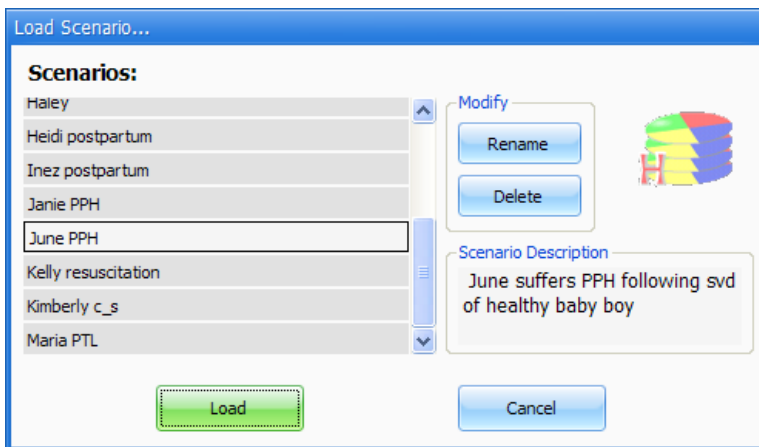Select the desired scenario and click on “Load”.



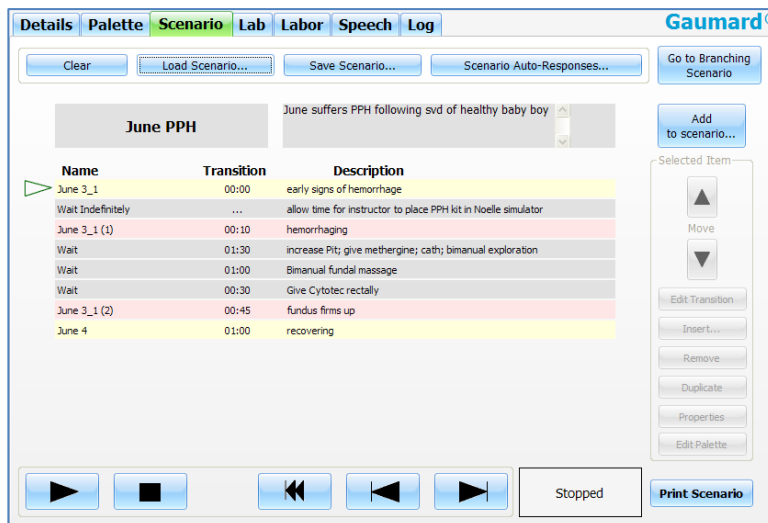After loading the desired scenario, click the **Play** button and watch the vital signs adjust according to the specifications of each palette.



## Creating your own Scenarios

First, create the palette items to be used in the scenario. To do so, go to the Details tab and change the controls that best describe the condition to be simulated. Not every field has to be populated in order to save a palette item

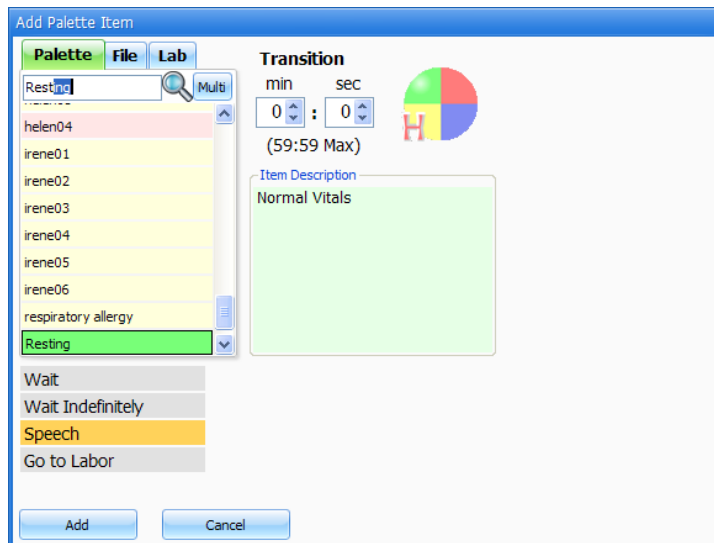
When operating in automatic mode, if inconsistent combinations or vitals are selected the model does not adjust to the specified values. The model is based on accurate physiologic principles, and therefore, choosing a combination of vitals that is inconstant with these principles, will not deliver the appointed results.

1. Save your palette(s) by clicking on the **Save as Palette Item** button on the upper right side of the page. Assign a **name** to the Palette Item and specify a brief description. Also, select a color that represents the palette's condition: Green for healthy, red for critical, yellow for care required, and blue for other. Then click **Save**.


2. After creating all your palettes, go to the **Scenario** tab.

3. Click on the **Add to Scenario** button. Select a single palette or enable **Multi** to select multiple palettes at the same time. Specify a transition time for the palette and click **Add**.

4. From the add menu, you may also choose the **Wait** item, which causes a delay of a specified duration, or a **Wait Indefinitely** item, which causes the scenario to pause until the facilitator **manually** advances to the next item.



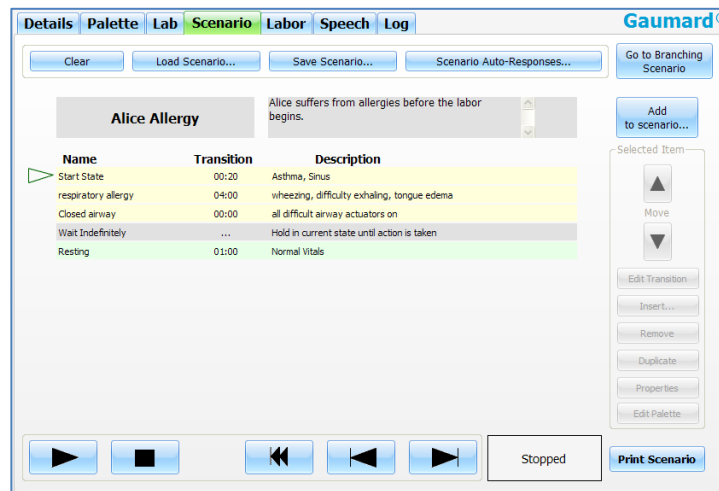5. Repeat the previous step to add more palettes, wait times, speech or labs.




6. Begin the scenario by pressing play. Watch the vital signs adjust according to the specifications of each palette item.

One can manipulate the scenario items with the **Selected Item** group of buttons on the right side of the page. Most of these functions are also available by right-clicking on scenario items. (When using the stylus, hold the stylus button and tap the screen to do a right-click).

### Auto-transition to Labor Scenario

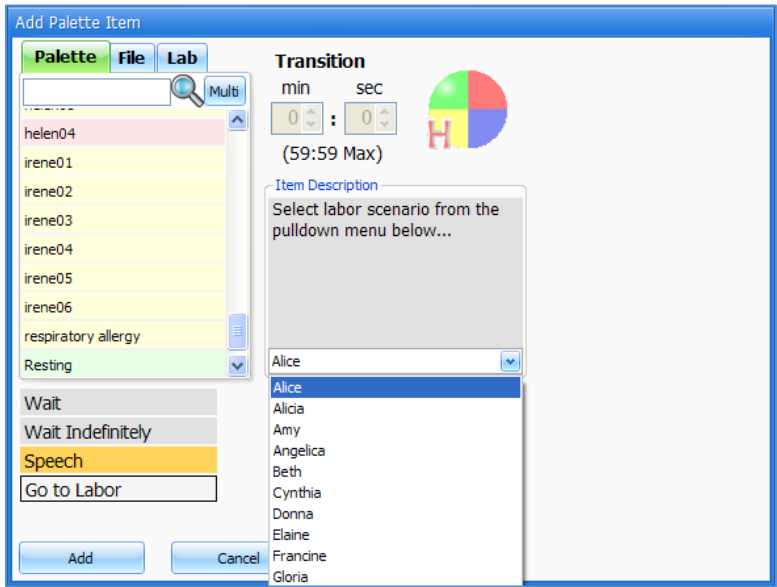To automatically transition from a linear scenario to a labor scenario, add a “Go to labor” palette at the end of the list. From the drop down, select the labor that will begin at the end of this scenario and click add.



Upon reaching the “Go to Labor” palette, the software will automatically begin the labor scenario named Alice.

Name	Transition	Description
Start State	00:20	Asthma, Sinus
respiratory allergy	04:00	wheezing, difficulty exhaling, tongue edema
Closed airway	00:00	all difficult airway actuators on
Wait Indefinitely	...	Hold in current state until action is taken
Resting	01:00	Normal Vitals
Go to Labor	...	Alice

**WARNING:** Prepare the simulator for delivery before using a linear scenario that automatically transitions to a labor scenario.

## Newborn Factory Preset Scenarios

These scenarios in the Quick Start Newborn profile when the Newborn software is activated. There are eleven linear scenarios.

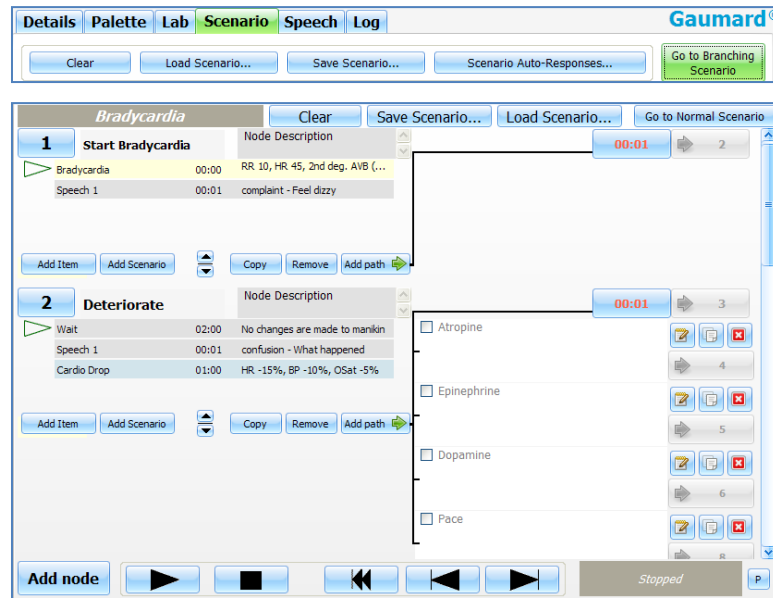
Scenario Name	Scenario Description
Linear	
Alice's Baby	Healthy baby
Asphyxia Scn	Baby has an asphyxia attack and the providers need to give ventilations to help bring back the vitals to a healthy state.
Beth's (Donna's) Baby	Baby is born with a mild asphyxia that needs attention. Once ventilations start the baby's vitals go to a healthy state.
Cynthia's Baby	Male infant with central cyanosis, limp, flaccid and requires immediate resuscitation. No spontaneous movement of right arm is noted. Stat CXR reveals a fractured right clavicle and right pneumothorax.
Elaine's Baby	This baby is born with moderate asphyxia, and will require CPR and oxygen to bring the vitals to a healthy state.
Francine's Baby	This baby was born through a C-Section and is responsive but need some attention, after a while all vitals go to a healthy state.
Gloria's Baby	This baby is born with mild asphyxia, but no matter how good the interventions are, this disastrous intrapartum complication results in neonatal death.
Helen's (Irene's) Baby	This baby is born with a severe asphyxia that has to be treated immediately, after ventilations and EPI have been given, the baby's vitals go towards a good outcome.
MAS	Meconium aspiration syndrome
RDS	Newborn with mild Respiratory Distress
TTN	Transient Tachypnea of the Newborn



## Branching Scenarios

The branching scenario is an advanced linear scenario editor. It allows the user to branch towards different scenarios/palettes depending on specific “Key Events” activated by the instructor.

To access the Branching window, click on the **Go to Branching Scenario** button on the Scenario tab. Similar to the linear scenario page, the buttons on the top panel can be used to clear, load and save a scenario, or to switch from branching to linear scenarios.




Go to the linear scenario page by clicking on the “Go to Normal Scenario” button.

## Adding Nodes

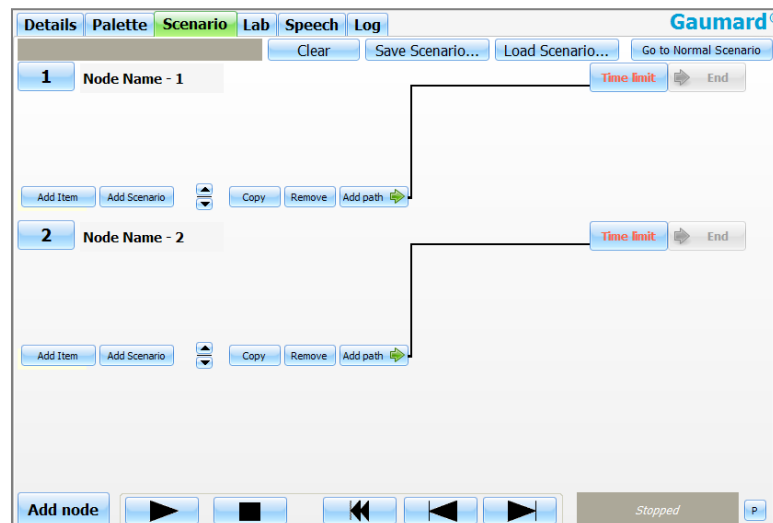A branching scenario will consist of several “Nodes” added by the facilitator. Each node is preconfigured to run a normal scenario or a series of palettes. The facilitator will then activate key events that will alter the trajectory of the nodes.

To add a node, click **Add node** near the bottom of the page.



To edit the node name and description, click the node's **Unique Identifier** number. Click **Save** to apply changes.

The 'Node Properties' dialog box shows the following details for Node 1:

- Unique Identifier:** 1
- Name:** Start Bradycardia
- Description:** Bradycardia

Buttons at the bottom: Save, Cancel.

In this example, the following nodes will be created: [1] Start Bradycardia, [2] Deteriorate, [3] Interventions, [4] Atropine, [5] Epinephrine, [6] Dopamine and [8] Pace. Each node has been programmed with specific palettes.

## Adding Palettes or Scenarios

Each node is configured with a set of palettes or scenarios.

The interface for Node 1 shows a list of palettes/scenarios. The first item is 'Bradycardia'. To the right of the list is a 'Time limit' icon and an 'End' button. At the bottom, there are buttons for 'Add Item', 'Add Scenario', 'Copy', 'Remove', and 'Add path'.

Click **Add Item** to add specific palette items or **Add Scenario** to add full scenarios to this node. Repeat the process and add palettes to the rest of the nodes.

The interface for Node 2 shows a list of palettes/scenarios. The first item is 'Wait' (02:00, No changes are made to manikin). The second item is 'Speech 1' (00:01, confusion - What happened). The third item is 'Cardio Drop' (01:00, HR -15%, BP -10%, OSat -5%). To the right of the list is a 'Time limit' icon and an 'End' button. At the bottom, there are buttons for 'Add Item', 'Add Scenario', 'Copy', 'Remove', and 'Add path'.

## Adding Paths

A path refers to the trajectory from one node to another after the last palette in a node expires. Click on the **Time Limit** icon to modify the **Go To** point for the default path. After the last palette expires, the scenario will move on to the node as indicated by the arrow.

The image shows a close-up of the 'Time limit' icon and the 'End' button. A blue arrow points from the text 'Path Go to point' to the 'End' button.

Configure the countdown timer and the “go to” point for the default path. Click **OK** to save.

The 'Set Time Limit' dialog box has a title bar 'Set Time Limit'. It contains two spinners for 'min' (set to 0) and 'sec' (set to 1), with '(59:59 Max)' below them. A label 'After...' is to the left. Below is a 'Go to' dropdown menu showing '2 - Deteriorate'. At the bottom are 'OK' and 'Cancel' buttons.

Node 1 is now configured to continue to Node 2 as indicated by the path's time limit.

The scenario editor shows Node 1 'Start Bradycardia'. The path list includes 'Bradycardia' (00:00) with description 'RR 10, HR 45, 2nd deg. AVB (...)' and 'Speech 1' (00:01) with description 'complaint - Feel dizzy'. A path line connects Node 1 to Node 2 with a timer of '00:01'. Buttons at the bottom include 'Add Item', 'Add Scenario', 'Copy', 'Remove', and 'Add path'.

## Key Events

Key events allow the facilitator to alter the trajectory of a branched scenario. This is done by assigning multiple paths to a single node, then selecting one of the paths when the provider completes a desired task. To add a key event to a node, first click the **Add path** button, then the **edit** button located on the right.

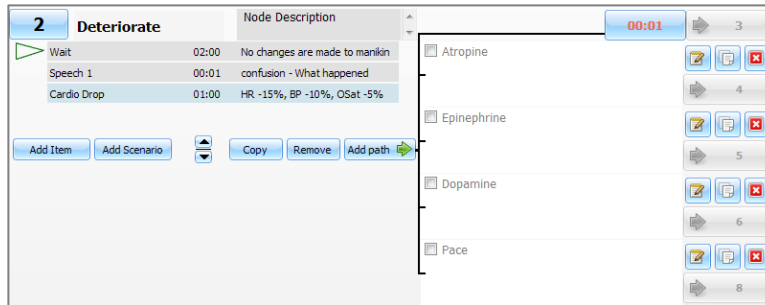
The scenario editor shows Node 2 'Deteriorate'. The path list includes 'Wait' (02:00) 'No changes are made to monitor', 'Speech 1' (00:01) 'confusion - What happened', and 'Cardio Drop' (01:00) 'HR -15%, BP -10%, OSat -5%'. A path line connects Node 2 to Node 3 with a timer of '00:01'. On the right, there are buttons for 'Edit', 'Delete', and 'Copy'.

Use the **Edit Path** window to name, sort and create **key events**. Begin by naming the *key event*, and then assigning a node destination from the **Go to** drop-down menu. Click **OK** to save changes.

*In this example, the provider will have the option to administer Atropine.*

The 'Edit Path' dialog box has a title bar 'Edit Path'. It contains a 'Key events:' list box. To the right are 'Selected', 'Move' (up/down arrows), 'Rename', and 'Delete' buttons. Below is a radio button for 'Action' with a text field containing 'Atropine' and an 'Add' button. There is also a radio button for 'Linked Action'. Under 'Linked Action', there are tabs for 'Electrical Therapy' and 'Drug Recognition'. Under 'Electrical Therapy', there are options for 'Defibrillation Shock' (200 Joules), 'Cardioversion Shock' (100 Joules), 'Non-Sync Shock' (100 Joules), and 'External Pacing' (40 mA). Under 'Drug Recognition', there are options for 'Shocks'. At the bottom is a 'Go to' dropdown menu showing '4 - Atropine' and 'OK' and 'Cancel' buttons. A list of nodes is shown at the very bottom: 1 - Start Bradycardia, 2 - Deteriorate, 3 - Interventions, 4 - Atropine (highlighted), 5 - Epinephrine, 6 - Dopamine, 8 - Pace, Previous, End.

The facilitator added a total of four paths to **Node 2 Deteriorate**: Atropine, Epinephrine, Dopamine and Pace. Once the provider performs an action (administration of Atropine or epinephrine) the facilitator clicks the key event to alter the path of the scenario.



## Linked Action Key Events

A linked action key event could either be a facilitator's checkmark, student interaction or automatically triggered actions.

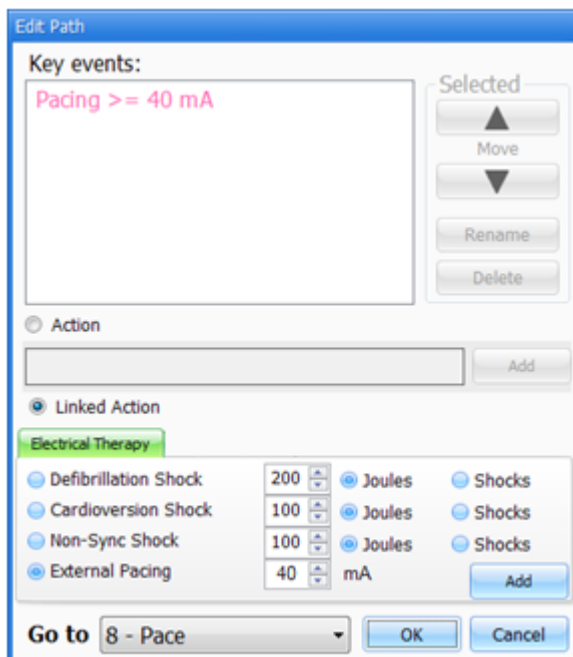
**Electrical Therapy** – The key event can be triggered by the electrical therapy applied by the students.

Configure the electrical therapy linked action menu to automatically trigger a key event or palette progression when electrical therapy is detected. The software responds to physical interaction and GaumardUI's software shock panel.


The settings for each Auto Response configured as a linked action event will override the general settings made in Setup, Auto Responses (auto, prompt, off).

Select from Defibrillation Shock, Cardioversion Shock, Non-Sync Shock or External Pacing. Set the parameters for response to either joules or number of shocks using the number fields and radio buttons.

Click **Add** to insert the electrical therapy into the Key Events field, and then **OK** to apply to the path. Electrical therapy can be identified by pink text in the node display.



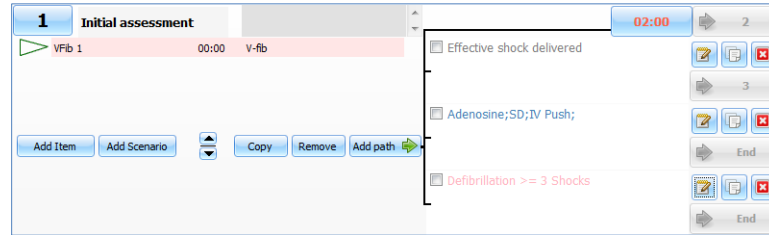*Node 1 of a Cardiac Arrest branching scenario with Model Medication and Electrical Therapy paths added:*



## Path Controls



Once a path is opened for any node, you can edit, copy or delete the path using the buttons below:

### Edit



This button allows you to add, remove, or change the order of actions in a specific path. Clicking on this button also enables you to select where the scenario should go after the path is finished playing.

### Copy



This button allows the user to copy the entire path. You can paste it to a different node, or paste it in the same node to duplicate it.

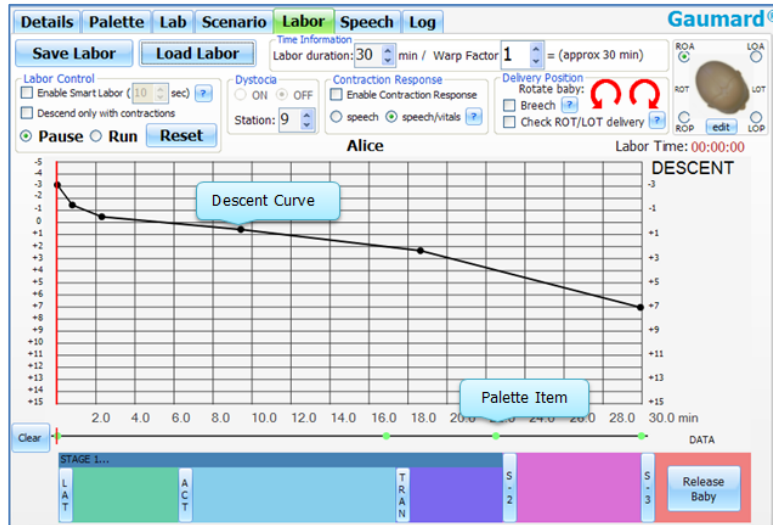
### Delete



This button allows the user to delete any undesired path.

## Labor


The most advanced method of controlling the NOELLE birthing system is to build a Labor Scenario, which is a sequence of Palette Items with delay periods corresponding to a labor curve.



The sections that follow describe in detail the function of the various controls found in the Labor page.

## Save Labor

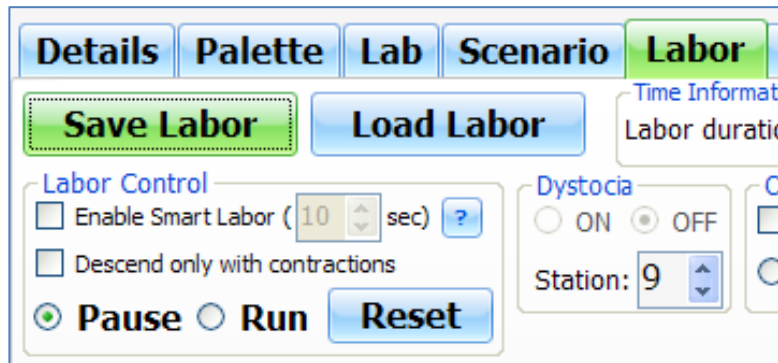Click on the **Save Labor** button to save the current configuration of the descent curve.



The window shown below opens. Use this window to input relevant information about the scenario for future reference. You **must** enter the Mother's name which will become the name of the scenario. All the remaining information can be left blank, as it is optional.

The Save Labor Scenario dialog box is divided into the various sections listed below:

- **Mother:** you can enter specific information about the mother including age, height, weight, among others. The only field that **must** be completed is the name of the Mother.
- **OB History:** Use this section to select specific details about the medical history of the patient.
- **Treatment History:** Use this section to check off any treatment histories that apply.
- **Labor type / Other:** Use this field to enter specific text regarding the labor. Entries may include notes such as "Postpartum Hemorrhage" or "Preterm Labor."
- **Stage Information:** this section contains a total of five tabs. Each tab corresponds to a labor stage. Use each tab to describe specific details about the patient during each stage of labor.
- **Post-Partum Scenario:** You can link a labor scenario to a postpartum scenario. This allows the software to automatically start the post-partum scenario upon delivery of the fetus and completion of the labor scenario.

**To view the Stage information while a labor is in progress, click on the label for the corresponding labor stage. Refer to Page 115 for more information.**

## Load Labor

To load a labor scenario, click the “Load Labor” button on the Labor tab and the Load screen will open.

The screenshot shows the 'Labor' tab selected in a top navigation bar. Below the navigation bar are three buttons: 'Save Labor', 'Load Labor' (highlighted with a green border), and 'Time Information'. Under 'Labor Control', there are checkboxes for 'Enable Smart Labor' (set to 10 sec) and 'Descend only with contractions'. There are also 'Pause', 'Run', and 'Reset' buttons. To the right, under 'Dystocia', there are 'ON' and 'OFF' radio buttons, and a 'Station' dropdown menu set to '9'.

All the previously saved files will appear on the left-hand side text box. A single click over a name displays the patient information for that labor scenario. To see the details for each stage, select each tab on the right hand side individually.

The screenshot shows the 'Load Labor Scenario...' dialog box. On the left is a 'Select' table with columns 'Name' and 'Last Modified'. The table lists several names: Alice, Alida, Amy, Angelica, Beth, Cynthia, Donna, Elaine, Francine, Gloria, Helen, and Irene. Below the table are 'Load', 'Cancel', and 'Delete' buttons. On the right is the 'Patient Information' section, which includes a text area with patient details and a 'Maternal Info' section with fields for 'Respiratory status', 'Cardiac status', 'Pain', 'Bleeding', 'Membranes', 'Dilation', 'Station', and 'Other'. There is also a 'Fetal Info' section with a checkbox for 'Decreased/No fetal movement' and an 'Other' field.

The preset labor scenarios are only found in the “Quick Start NOELLE or “NOELLE Advanced” profiles. There are no factory preset scenarios in the “Default Profile.”

## Time Information


**Labor duration:** This is the time of the delivery being simulated from stage one to stage three.

**Warp Factor:** By increasing this number the labor duration will be shortened so that the labor can be simulated faster. An approximation of the “Labor duration” will be given in parenthesis. In this way, a two-hour labor can be simulated in five minutes as an example.

The screenshot shows the 'Time Information' section in the 'Labor' tab. It features a 'Labor duration' field set to '30 min' and a 'Warp Factor' field set to '1'. To the right of the Warp Factor field is the text '= (approx 30 min)'.



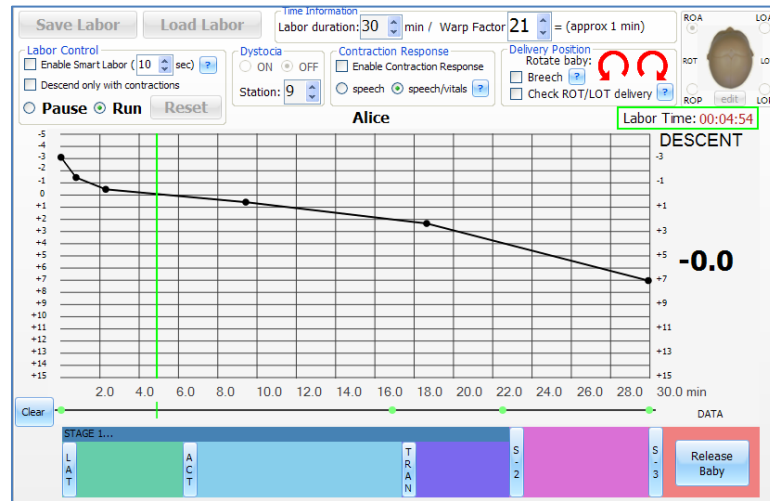
**Labor Time:** This timer represents the time on the labor curve.



## Delivery Position


NOELLE's birthing mechanism can be preset to four different initial configurations, if a regular vertex delivery is desired. These conditions are ROA, LOA, LOP and ROP. It is very important that once the motor has been reset, the instructor selects the desired delivery configuration in the software. Then the baby has to be connected representing this state.

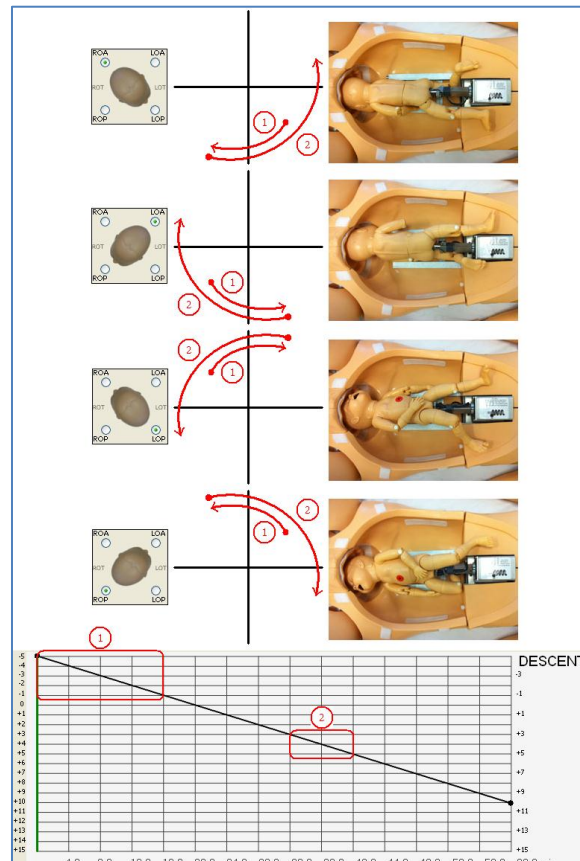**Select the initial position for the baby.** Make sure it represents the position at which the baby was inserted inside NOELLE.



Each of these states has a different rotation program that will represent the internal rotation and the external rotation. The rotation is also dependent on the labor curve. Notice that the internal rotation (marked as “1” on the diagram above) will be calculated between the stations -5 and -1, and the external rotation (marked as “2” on the diagram) will be calculated between the stations +3 and +5. If by any chance a labor curve is designed in such a way that the first point of the curve is below -1, the internal rotation will not be calculated.

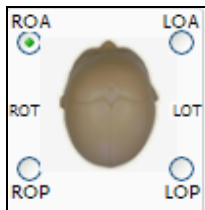The rotations may be edited by the instructor. Click on the “Edit” button at the bottom of the head rotation selection box to expand the fields. The first two rows allow you to set where in the delivery curve the two rotations occur. The software will not allow an overlap of the rotations – the stations specified for Rotation 2 will always be values greater than those specified for Rotation 1 (further down the curve)



Each head position has two pairs of fields, one pair each for Rotation 1 and Rotation 2. The provider can set the arc of the rotation to be as shallow or complete as desired. Please remember that the zero point in the arc will always be at the top of the head position rotation box, between ROA and LOA.

Below is an expanded view of the “Edit” control.

The 'Edit' control interface consists of two side-by-side panels. Each panel shows a fetal head diagram with rotation points labeled ROA, LOA, ROT, LOT, ROP, and LOP. Below the diagram, there are two rows of rotation settings. The first row is for Rotation 1, with 'Station' set to -5 to -1, and 'ROA (135)' set to 135 to 195. The second row is for Rotation 2, with 'Station' set to 3 to 5, and 'ROP (45)' set to 45 to -15. Below these are two more rows for LOA (225) and LOP (315). At the bottom of each panel are 'Apply', 'Cancel', and 'Reset' buttons. In the right panel, the 'Apply' button is highlighted with a green border.

These settings will be retained until either altered again, or the “Reset” button is pressed. “Reset” will return all stations and arcs to factory defaults.

Make sure the labor curve covers all the points from -5 to +5 if both rotations are desired. Also notice that the rotations are somewhat exaggerated — this is done because the head swivels a couple of degrees separately from the alignment of the torso. The torso is the portion of the baby that’s being rotated by the motor, so the head lags behind by a small amount.

**Rotate baby:** Both of these arrows allow the instructor to rotate the baby to any desired position while the delivery is in progress.




**Breech:** Use this control when planning a breech delivery. When selected, it disables the motor rotation.


**Check ROT/LOT delivery:** this control limits the release mechanism to trigger only when the neonate's shoulders are positioned vertically. If the shoulders are not aligned vertically, the mechanism will not release. The only way to detach the baby then is by manually clicking on the "Release Baby" button.

## Labor Control

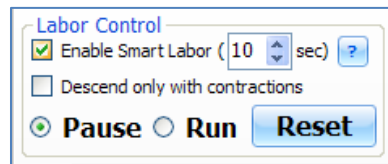**Enable Smart Labor:** this option allows you to run all non-speech palettes in real time. You can specify for how long the palettes will run in real time. The software defaults to ten seconds.



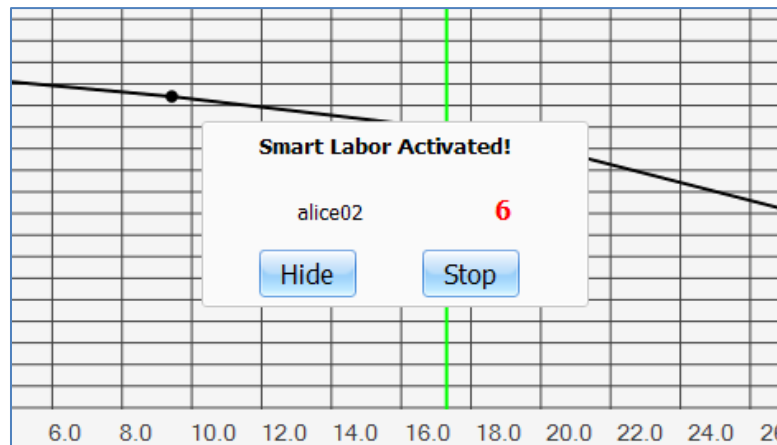When the Smart Labor is activated, the window circled in red below is displayed:




Hide the message by clicking the "Hide" button or "Stop" that palette from running in real time; otherwise, the message will display for ten seconds. The number in red indicates how many seconds are left before switching back to warped speed. The name of the palette being applied in real time is also shown in this dialog box.

**Descend only with contractions:** If checked, this option configures the motor to descend ONLY when there is a contraction.

**Run:** Begins a labor and **initializes the birthing mechanism**. The red bar on the left of the labor screen turns green when labor scenario is running.

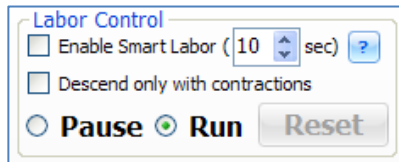**WARNING:** Do not initialize the birthing mechanism until reading **Working with NOELLE** and **Care and Cautions** sections of this guide for important information on preparing NOELLE for delivery.




**Pause:** This button pauses the labor.

**Reset:** This button will reset the motor to the initial position and also reset the labor scenario to its initial conditions. The reset button is enabled only when a labor is paused. If a labor is running, this button will be grayed out. To troubleshoot the motor mechanism, navigate to page 181.


## Dystocia

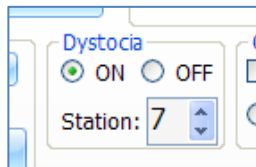The dystocia controls will be activated only if the last point on the descent curve is lower than the dystocia threshold set on the “Options” or on the “Dystocia Control Box”. Dystocia can be used during any scenario by switching the control on in the Dystocia control box. When the dystocia is set to “On” and the labor starts, the warp factors for labor and perinatal monitor will change to real time (warp factor 1) as the descent curve passes the preset dystocia threshold. While in dystocia mode, a “Turtle Sign” will occur with each uterine contraction.




The dystocia controls are disabled when the last point on the descent curve is less than the indicated dystocia station.

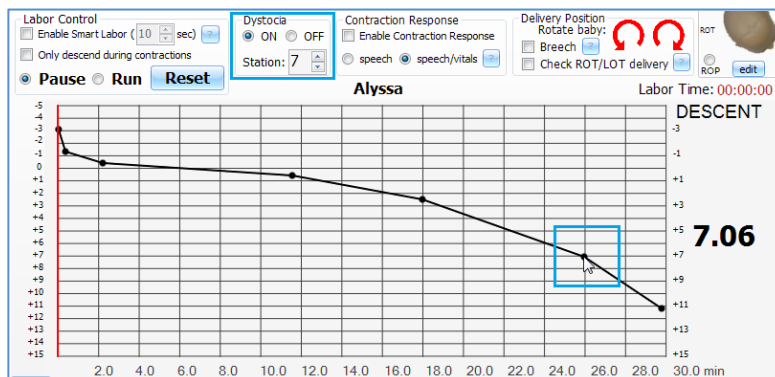In the example below, dystocia controls are disabled because the last point is less than **seven** (the current dystocia station).




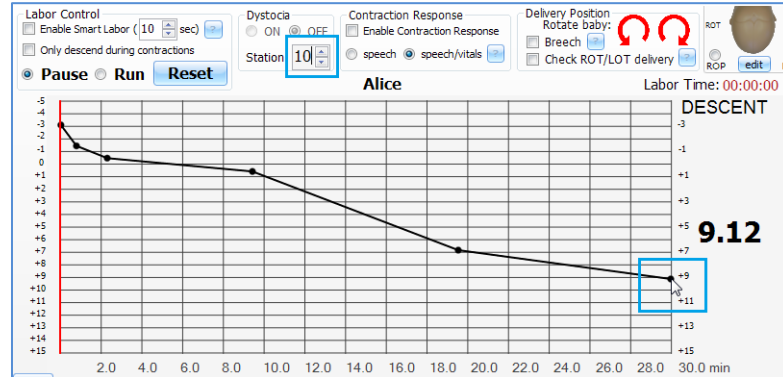In the example below, shoulder dystocia and turtle signs begin station 10.



## Contraction Response

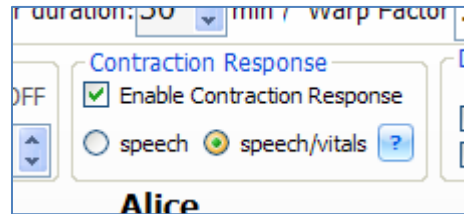The contraction response feature makes your labor simulation more realistic. When enabled, this feature can automatically make changes to the maternal vital signs or add speech palettes during the contraction.



There are two kinds of patients that you can simulate: Pain-controlled patients and patients with uncontrolled pain. The changes made are different depending on the selection speech or speech/vitals.

Selection	Maternal Vital Signs Changes	Audio Palettes	Notes
Speech	No Changes	"I think I'm having contractions"	Simulates a patient that is pain-controlled
Speech/vitals	Increase HR, BP and RR	"Ouch," "Ahhh," "Please, I need something for the pain," "Please give me an epidural."	Simulates a patient that is NOT pain-controlled

The changes are applied only during the length of the contraction. At the end of a contraction, the maternal vital signs-if changed- adjust back to the previous setting.

When a contraction response is applied, the changes are recorded in the Log page.

Pain controlled patient:

```
00:00:44 Contraction Response: Speech - Contractions ;
00:00:56 Applied (00:00): Palette: alice02: {Details: ContractionFreq. 4 min; ContractionDurat. 40 sec; }
00:00:59 Contraction Response: Speech - Epidural ;
00:01:11 Contraction Response: Speech - Aahh ;
00:01:12 Applied (00:00): Palette: alice03: {Details: RR 23; OSat 97%; ContractionFreq. 3 min; ContractionDurat. 70 sec; ContractionI
00:01:21 Contraction Response: Speech - Ouch ;
00:01:30 Contraction Response: Speech - Epidural ;
00:01:33 Applied (00:00): Palette: alice04: {Details: HR 70; RR 19; }
```


No pain controlled patient:

```
00:00:41 Contraction Response: Decrease HR, BP, RR 9% during contraction.
00:00:42 Contraction Response: Decrease HR, BP, RR 9% during contraction.
00:00:55 Applied (00:00): Palette: alice02: {Details: ContractionFreq. 4 min; ContractionDurat. 40 sec; }
00:00:57 Contraction Response: Speech - Epidural ; Increase HR, BP, RR 12% during contraction.
00:00:58 Contraction Response: Decrease HR, BP, RR 9% during contraction.
00:01:07 Contraction Response: Speech - Aahh ; Increase HR, BP, RR 12% during contraction.
00:01:08 Contraction Response: Decrease HR, BP, RR 9% during contraction.
00:01:11 Applied (00:00): Palette: alice03: {Details: RR 23; OSat 97%; ContractionFreq. 3 min; ContractionDurat. 70 sec; ContractionI
```

## Descent Curve

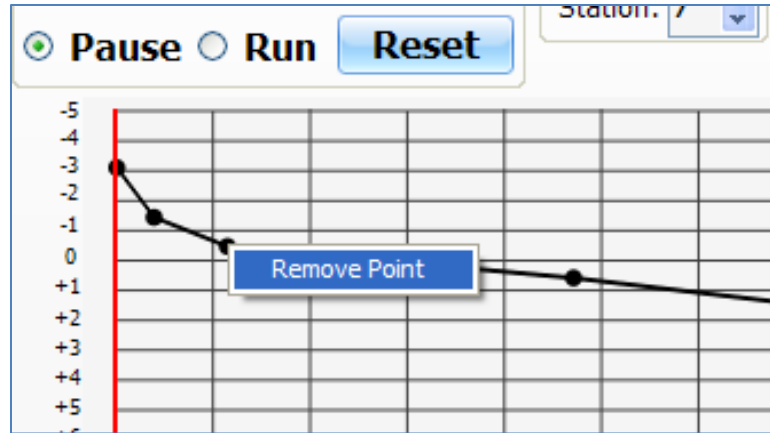
The descent curve represents the position of the fetal head relative to the ischial spines in centimeters, from -5 cm to +10 cm (station). The user can manipulate this curve by clicking over a point on the line and dragging it to the desired location on the grid. In order to add more points, right-click on the line between existing points and select "Add Point". To remove points, right-click over the point and select "Remove Point".





## Time Line

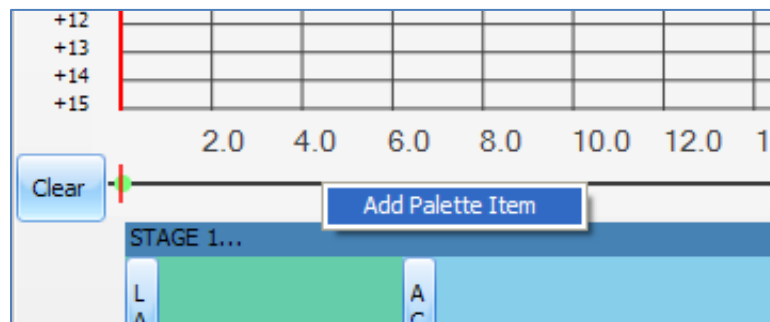
The time line is where Palette Items are added to build the labor scenario. As the time indicator line passes over each point, the corresponding Palette Item is implemented. To add a palette item, right-click on the time line and select "Add Palette". This will open a window with all the palette items available in the current profile.



Select a palette item or speech and specify a transition time, then click **Add**.

Once the point has been added, you can touch it with the stylus and drag it to the desired location on the time line. Double-click over the point to bring up a window displaying the properties of the Palette Item represented by the point.

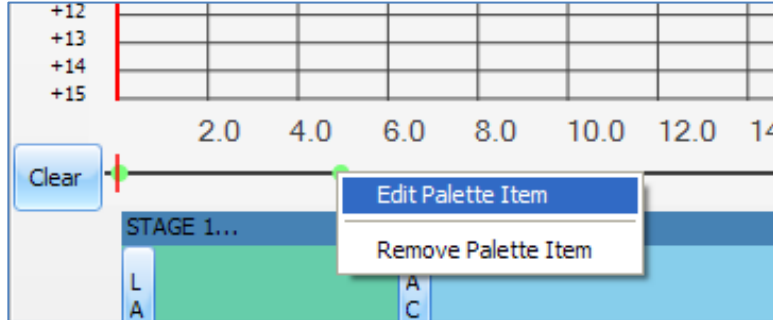
To remove or edit a palette item from the time line, right click over the point and select the desired option.




The user can remove all palette items at once by pressing the "Clear" button located to the left of the timeline. When all palette items are removed the user is left with a clean slate to either rebuild or to completely change the scenario.

## Stages

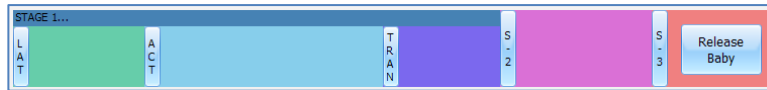The stages (shown in colored blocks under the timeline) are also dependent on the descent curve.



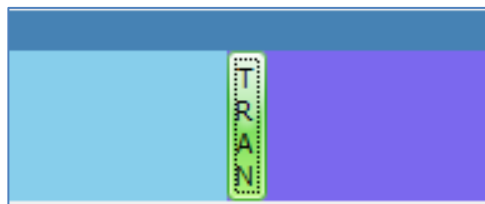The blocks, representing different stages of labor, change size as the points on the descent curve are moved. Click on each stage label to view the details.





This information was assigned when the scenario was created.

**Stage Information**

S1 - Latent | S1 - Active | **S1 - Transition** | Stage 2 | Stage 3

**Maternal Info:**

Respiratory status: Normal

Cardiac status: Adequate

Pain: 3

Bleeding:

Membranes:

Dilation:

Station:

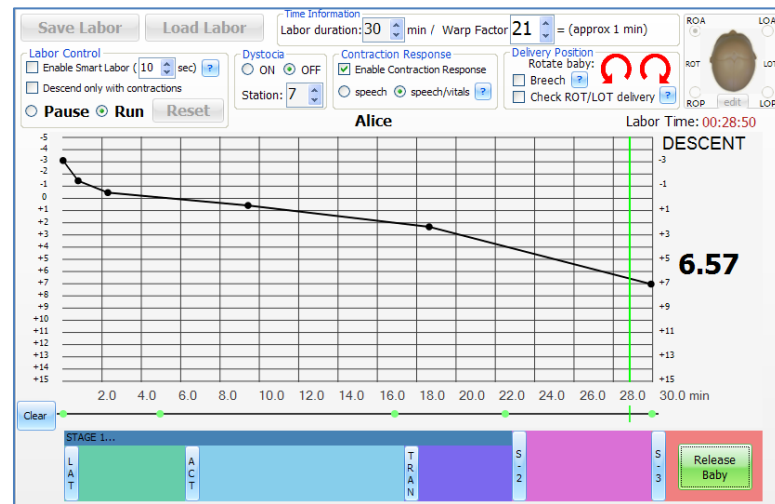
Other:

**Fetal Info:**


☐ Decreased/No fetal movement

Other:

**Release Baby:** This button activates the release mechanism on the baby, which allows either the providers to disengage the baby or the instructor to connect the baby.



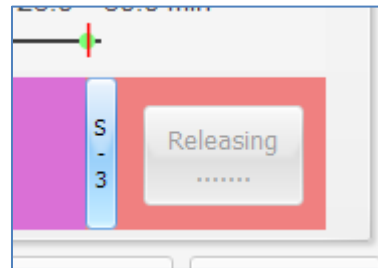After clicking on the Release Baby button, as shown below, the system will take a few seconds to send the release commands to the simulator.



## Locking Mechanism

Once the baby is plugged into the motor arm, the locking mechanism acts twice to secure the birthing baby (a clicking noise can be heard coming from inside the baby). The user then must hold the baby until the noise stops, and make sure that it is locked by pulling the baby to check that it is held fast. For more information on preparing the simulator for delivery, navigate to page 99.

Follow the guidelines below when releasing the baby from the birthing mechanism. To troubleshoot any additional issues, please navigate to page 181.

### If shoulder dystocia turned OFF:

- The release mechanism is automatically actuated ten times in one-second intervals once vertical progress bar reaches **STAGE 3**.
- The baby has descended more than 80% of its entire translation (baby's head is out) and user pulls the baby with more than three pounds of force.
- By pressing the **Release Baby** button inside the "Labor" tab.

### If shoulder dystocia is turned on:

- When user pulls the baby with more than 35 lbs., the baby is released as a safety feature in order to avoid damage
- By pressing the **Release Baby** button inside the "Labor" tab.
- Allow at least 30 seconds after the baby has been initially locked before trying to release it.
- The labor progress line reaches Stage 3 and dystocia is turned off.

**WARNING:** Guide the baby out of the birth canal by gently pulling in line with the birthing mechanism. Pulling the baby upward or downward in contrast to the birthing mechanism's linear trajectory may bend the motor arm and cause damage to the birthing mechanism.

## Labor Activity

As labor progresses, the manikin will detect manipulations performed on the fetus by the care providers. Activity on the fetus prompts the "Activity on Fetus" window to appear on the tablet screen, providing the instructor with feedback on pulling force, uterine contractions, torque on the baby (twisting force) and shoulder position. The instructor can evaluate the providers in terms of pulling in, or out of, sync with uterine contractions. The labor activity is a beneficial tool for the instructor as excessive forces applied to the fetus by care providers during delivery can lead to brachial plexus injuries. All of these values will be given in real time and they will be recorded by the graphical interface for further analysis.

**Peak Force** - measurement of how hard the provider is pulling the fetus. (Can be displayed in lbs. or Newton)

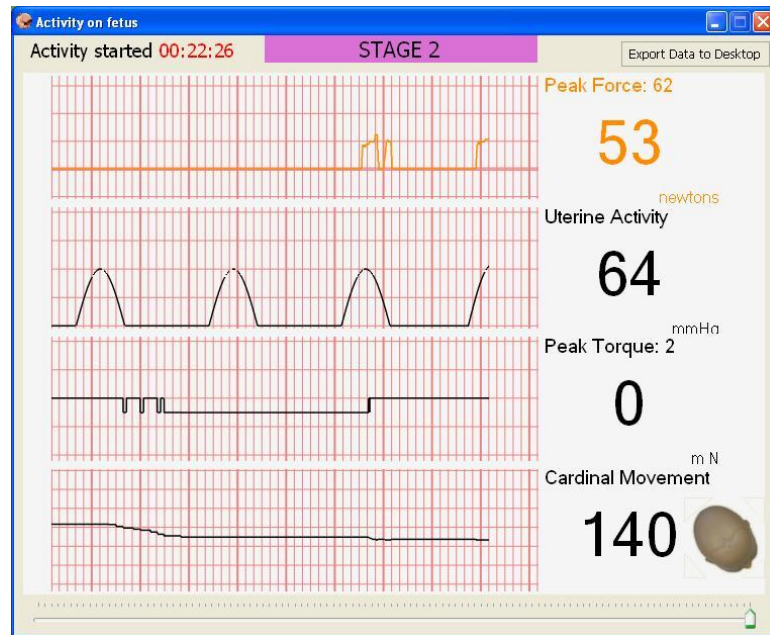
**Uterine Activity** - contractions generated by NOELLE.

**Peak Torque** - measurement of the torque induced as the provider rotates the fetus; negative readings represent counter clockwise rotation, and positive values represent clockwise rotation.

**Cardinal Movement** - the angular position of the baby is measured in degrees, taking into account that "0 degrees" corresponds to the fetus facing upwards with the shoulders horizontal.

**Export Data to Desktop** - exports labor force data to an Excel file. Use the data to make graphs to study trends and performance patterns.

**Slider** - The slider at the bottom of this screen lets you rewind the strip to see the previous graph data.



## Factory Preset Labor Scenarios

NOELLE has a total of thirty-six factory preset labor scenarios, which were designed in conjunction with healthcare professionals. The scenarios are divided into two separate profiles. Below is information on the scenarios available under the quick start and NOELLE Advanced scenarios.

### Quick Start Scenarios

Scenario Name	Labor Type	History	Overview
<b>Alice</b>	Normal	Alice is a 24 year old gravida 2/1 at 39 weeks. She weighs 170 pounds. She has had prenatal care. She has not been using medications of any kind.	Runs for 30 minutes. Labor progresses normally and fetal heart tones remain within normal limits. The normal male infant earns good APGAR scores.
<b>Alicia</b>	Variations on Normal	Alicia is a 24 year old gravida 2/1 at 39 weeks. She weighs 160 pounds. She has had prenatal care. She has not been using medications of any kind.	Runs for 20 minutes. Labor progresses normally and fetal heart tones remain within normal limits. The normal infant earns good APGAR scores.
<b>Amy</b>	Variations on Normal	Amy is 19 years old 1/0 at 40 weeks. She weighs 160 lbs.	Runs for 30 minutes. Labor progresses normally and fetal heart tones remain within normal limits. The baby earns good APGAR scores.
<b>Angelica</b>	Variations on Normal	Angelica is a 31 year old gravida 5/3 at 41 weeks. She weighs 160 lbs. She has experienced no prenatal complications and has a history of fast labors.	Runs for 20 minutes. Labor progresses normally and fetal heart tones remain within normal limits. The normal infant earns good APGAR scores.
<b>Beth</b>	Variations on Normal	Beth is a 16 year old gravida 2/0 at 37 weeks. She has had one elective abortion. She has had prenatal care.	Runs for 10 minutes. Fetal descent is rapid. Fetal baseline is maintained at 150. Nuchal chord is evident and moderate bleeding is noted immediately following delivery from a second degree perineal laceration. The normal female infant is limp, dusky and does not cry spontaneously. Baby is hypothermic and tachypneic.
<b>Cynthia</b>	Shoulder Dystocia	Cynthia is a 31 year old gravida 3/1 at 41 weeks. She weighs 170 lbs.	Runs for 30 minutes. Labor progresses normally and fetal heart tone baseline remains within normal limits. Patient is unable to fully "crown." Vacuum extractor is required to deliver the head. Shoulder dystocia is encountered and McRobert's and suprapubic pressure maneuvers are required. Male infant is centrally cyanotic, limp, and flaccid and requires immediate resuscitation. Stat CXR reveals a fractured right clavicle and right pneumothorax.
<b>Donna</b>	Breech	Donna is a 20 year old gravida 4/2 at 31 weeks. She weighs 180 lbs. She has had one elective abortion. She has had prenatal care.	Runs for 20 minutes. Labor progresses quickly and breech is delivered by the nurse. Meconium is noted. The preterm female baby cries weakly with stimulation but color and tone are poor. She is transferred to the nursery for stabilization and continuing care.
<b>Elaine</b>	Preeclampsia	Elaine is a 23 year old gravida 1/0 at 37 weeks. She weighs 140 lbs. She has had prenatal care. She complains of mind frontal headache. 3+ tibial edema and 4+ DTRs with 2 beats clonus are	Runs for 40 minutes. Progressive cervical change and fetal descent are noted during 9 hour induction. Fetal descent continues with little active pushing. FHTs show decreased variability and mild to moderate decelerations. Delivery is accompanied with outlet forceps. Female infant

Scenario Name	Labor Type	History	Overview
		noted.	is dusky, limp and does not breathe spontaneously at delivery. Meconium is noted and a small amount is observed below vocal cords. Baby is suctioned and is eventually transferred to NICU for continued observation.
<b>Francine</b>	Cesarean Delivery	Francine is a 19 year old female gravida 2/1 at 37 weeks. She weighs 145 lbs. She has had prenatal care. She has STD, Herpes.	Runs for 10 minutes. Delivery of male infant is accomplished through a low transverse uterine incision. The infant exhibits good tone and cries spontaneously at delivery, peripheral.
<b>Gloria</b>	Cord Prolapse	Gloria is a 34 years old gravida 1/0 at 25 weeks. She weighs 190 lbs. She has had prenatal care.	Runs for 10 minutes. Gloria arrives at the hospital with ruptured membranes and an obviously prolapsed cord. Profound fetal bradycardia is noted. Delivery occurs almost immediately after the patient is moved to the delivery room. C&S is obtained and sent to pathology. The male infant is placed on infant warmer. His earns APGAR score
<b>Helen</b>	Hemorrhage	Helen is a 25 year old gravida 1/0 at 35 weeks. She weighs 180 lbs. She has had prenatal care.	Runs for 30 minutes. FHT remain WNL. Midline episiotomy is performed and the delivery is spontaneous. The placenta delivers but is not intact. Bimanual uterine exploration removes small amount of placental tissue. The male infant is pale and flaccid. He is covered with vernix and blood. Neonatal resuscitation is begun immediately.
<b>Irene</b>	Cesarean Delivery	Irene is a 19 year old gravida 2/0 at 29 weeks. She has had one spontaneous abortion.	Runs for 45 minutes. Sterile speculum exam indicates a shortened cervix with a cerclage in place. Fluid is noted in the vaginal vault and ferning is positive. Breakthrough contractions occur and sterile speculum reveals cervical change. Cerclage is removed. Fetal heart tones remain stable throughout the short labor. Spontaneous delivery occurs. The female infant born earns good APGARs.

## NOELLE Advanced

Patient Name	Labor Type	History	Overview
<b>Alyssa</b>	Normal Labor	Alyssa is a 23 year old primip at term. She has received prenatal care in a multidisciplinary practice. Her general health is good and she has experienced no prenatal complications. She wants to have an unmedicated natural birth.	Runs for approximately 30 minutes. Simulates 8 hour labor without complications resulting in SVD over intact perineum. Baby is vigorous and earns good APGARS.
<b>Angela</b>	Normal Labor	Angela is a 31 year old grand multip. Even though she has had only a few contractions she came in anyway because she has a history of rapid labors.	Runs for approximately 30 minutes. Simulates normal labor and delivery of grand multip resulting in shorter labor duration.
<b>Becca</b>	Variations on Normal	Becca is a young pregnant teen who is living on the streets. She is a heavy smoker and drug user. She was seen twice in the Adolescent Clinic and referred to Social Services, but she only saw the social worker once and did not go to the follow-up appointment. The extremely precipitous delivery and patient's agitated state should point participants to order tox screen to check for illicit drug use.	Runs for approximately 18-22 minutes. Following ROM, nurse must support fetal head to prevent explosive delivery. Patient has retained secundes therefore bimanual exam is required. <b>Safety Note:</b> Instructor is required to reset motor before allowing learner to perform bimanual exam.
<b>Bianca</b>	Variations on Normal	Bianca is a 16 year old gravida 2/1. She had an elective abortion at age 13. She lives with her 17 year old boyfriend and has no contact with her family. She has been inconsistent with appointments at the teen clinic due to transportation issues.	Runs for approximately 30 minutes. La3-25 <b>Note:</b> To run full scenario, load Bianca at Admit from scenario tab. Instructor has to set up fetus with Nuchal chord prior to running scenario.
<b>Candice</b>	Shoulder Dystocia	Candice is a 19 year old multip. She and her boyfriend and 3-year-old daughter are homeless and currently living in a car. She has not seen a doctor, but believes that she is about 8 months pregnant.	Full scenario runs for approximately 23-25 minutes. Dystocia drill runs for 15 min. Labor progresses normally for about 6 hours, but after 45 minutes of pushing, patient is unable to bring vertex to perineum. A vacuum extractor is necessary to bring the head to the perineum. Patient is still unable to deliver. McRobert's, suprapubic pressure, Woods and Rubin maneuvers fail. Baby is finally delivered using Gaskin maneuver. <b>Note:</b> To run full scenario, load Candice Admit from scenario tab. To run dystocia drill only, load Candice from labor tab.
<b>Charlotte</b>	Shoulder Dystocia	Charlotte is a 31 year old gravida 3/1 at 41+5/7 weeks. Her physician stripped her membranes yesterday and she began contracting during the night. She is admitted in active labor.	Scenario runs for approximately 40 minutes. Baby is delivered after basic shoulder dystocia maneuvers such as McRoberts and stiff arm suprapubic pressure.

Patient Name	Labor Type	History	Overview
<b>Dana</b>	Breech Presentation	Dana is a 24 yr. old multip @ 29 weeks who was admitted because she began contracting. Upon V/E physician discovers that she is 4-5cm with bulging membranes. She was given Terbutaline subQ and then transported to the regional medical center.	Runs for approximately 45 minutes. Instructor must set up baby for double footling breech delivery. To skip resolution of contractions portion of scenario, it is suggested that the instructor warp to around 20 minute mark on the timeline.
<b>Demaris</b>	Breech Presentation	Demaris is a young Hispanic teen who has received prenatal care in the Adolescent OB clinic. She kept the pregnancy a secret as long as was possible and did not attend any childbirth classes. Her plan is to return to high school while her mother cares for the baby. The baby's father will not accept any responsibility and does not wish to be involved.	Runs for approximately 30 minutes. Frank breech position has been confirmed by real time U/S. Patient refuses Cesarean delivery. Pinard maneuver must be utilized to bring the legs down. MLE is performed and baby is delivered.
<b>Eleanor</b>	Preeclampsia	The ER is notified that EMS is about 4 minutes away with a 19 year old pregnant, post-ictal patient named Eleanor. Her aunt found her convulsing in the bathroom and called 911. The aunt told the EMS providers that Eleanor was 8½ months pregnant with her first baby and that it was a difficult family situation. She added that Eleanor had just moved in with them and had not yet seen a doctor. The paramedic reports to the ER physician by radio that the patient is responsive only to pain. Her initial blood pressure is 180/120. The EMS crew applies a C-collar and move her onto a backboard for transport. The paramedic initiates ECG monitoring and does a genital exam before they move her. She notes a small amount of vaginal bleeding. The crew moves her to the ambulance and the paramedic starts an IV of LR and initiates oxygen @ 10L by non-rebreather mask. She monitors the ECG and VS during transport, and also tilts the backboard about 15° to the left with a blanket roll to decrease vena caval compression. Eleanor has no further seizure activity during transport, but her BP remains consistently 180/110. The ER notifies L&D about the patient and asks that an OB nurse come to the ER to assist.	Runs for approximately 20 minutes. Instructor must place C-collar on NOELLE simulator before beginning scenario. During delivery, patient has tonic-clonic seizure followed by tetanic contractions. Fetal baseline drops to approximately 60 bpm. SVD occurs very quickly. Infant is cyanotic and limp, and no respiratory effort is evident. Baby is pronounced dead after 20 minutes of resuscitation.
<b>Erin</b>	Preeclampsia	Erin is a 28 yr old @ 38 weeks admitted by her physician for preeclampsia. She is started on Magnesium sulfate per protocol, induced with pitocin and her membranes are ruptured. Light meconium staining is noted following AROM.	Runs for approximately 60 minutes. Upon induction of pitocin and rupture of membranes, instructor should tell participants to switch audio output on fetal monitors to FSE tones.

Patient Name	Labor Type	History	Overview
<b>Faye</b>	Cord Prolapse	Faye is a 34 year old gravida 1/0 @ 25 weeks' gestation. She had been involved with a married man and this unexpected and unwanted pregnancy caused a great deal of stress in her life. After much emotional upheaval, she decided to have the baby. The affair ended and she is no longer involved with the father of the baby. None of her family or friends are aware of the pregnancy. She began cramping about 3 hours ago and decided to drive herself to the hospital when she began leaking clear fluid. An admitting clerk helps her into a wheelchair and takes her to L&D.	Runs for approximately 20 minutes. Pro-lapsed cord is protruding into the vaginal opening. Faye is put into deep Trendelenburg position. There is no palpable cord pulsation and FHTs are not audible with Doppler. Upon VE, fetus presents as double footling breech and is already partially into birth canal. Therefore, baby has to be delivered vaginally. Fetus is non viable.
<b>Frances</b>	Cord Prolapse	Frances is admitted into a small town hospital due to regular contractions @ 4 minutes apart and bloody show. She labors without problems for about 4 hours and then the fetus starts to brady down after SROM. A V/E reveals a prolapsed cord in the vagina.	Runs for approximately 22-27 minutes. Instructor must set up prolapse cord prior to beginning scenario. Instructor may disconnect NOELLE simulator from power outlet and continue running scenario while transferring the NOELLE birthing simulator to the OR.
<b>Gabriella</b>	Uterine Rupture	Gabriella is a young Hispanic woman who presents to a small hospital just across the Mexican border. She is alone and speaks little English. Her nurse is fluent in Spanish, but Gabriella is very quiet and gives little information.	Runs for approximately 25 minutes. Shortly after admission, patient clutches her abdomen and fetal baseline bradys down to the 80s. Patient becomes very pale and diaphoretic. BP drops to 80/60 and pulse is 120. She suffers heavy vaginal bleeding. Patient is rushed for emergency C-Section. Baby is limp and severely depressed, needs to be intubated and ventilated.
<b>Gail</b>	Uterine Rupture	Gail is a 29 year old primip @ 35 weeks. She was admitted to L&D from the ER after being involved in car accident. Both she and her husband, Alan, were seriously injured and she is on a backboard wearing a c-collar to stabilize the spine. Her right humerus is fractured and seat belt marks are visible across the abdomen.	Runs for approximately 181-20 minutes. Use file sharing feature to display pictures of C-spine x-ray as patient was involved in car crash.
<b>Haley</b>	Peripartum Hemorrhage Previa	Haley is a 33yr old G2 @ 35 weeks. - Previous U/S revealed a low lying placenta and this is the 5 <sup>th</sup> time in 11 weeks she been admitted for bleeding. This time the bleeding is heavier and is not resolving. Her OB is on the way to the hospital; bimanual palpation shows the uterus to be soft and non-tender.	Runs for approximately 15 minutes. Instructor must set up for partial placenta previa. Instructor must fill hemorrhage kit with fluids before running the scenario
<b>Heidi</b>	Peripartum Hemorrhage Previa	Heidi is a 25 yr old primip @ 35 weeks. - She has experienced several mild bleeding episodes during pregnancy and is known to have a low lying placenta. She arrives in L&D complaining of abdominal cramps and has bright red vaginal bleeding.	Runs for 35-45 minutes. This scenario utilizes a vaginal delivery even though there is a low lying placenta. Instructor has option to run a linear or a branching postpartum section of the scenario. Instructor must fill hemorrhage kit with fluids before running the scenario

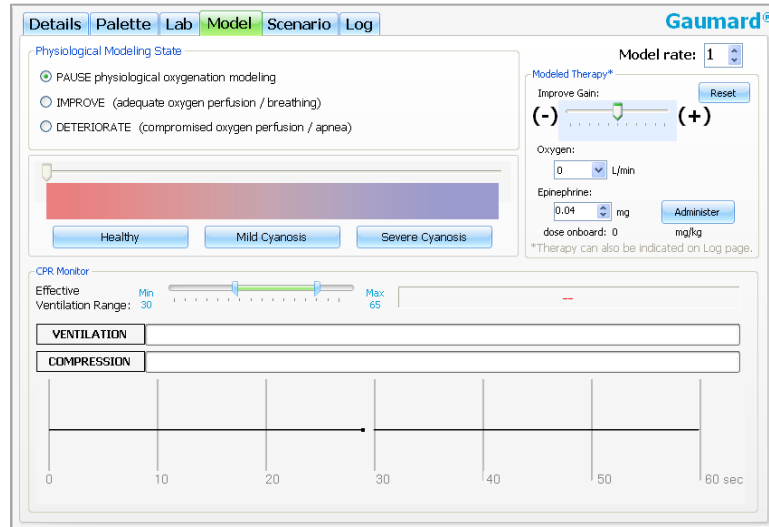


Patient Name	Labor Type	History	Overview
India	Peripartum Hemorrhage Abruption	India is a 19 yr old gravida 2 @ 37 weeks. She arrives at hospital with her husband who says she fell down the stairs and she has been cramping and bleeding for about an hour. During admitting interview husband answers all the questions and India doesn't make eye contact. The nurse palpates uterus, initiates fetal monitoring and starts a pad count.	Runs for approximately 25 minutes. This scenario is a suspected spousal abuse case. Patient has a central abruption of the placenta. Instructor must fill hemorrhage kit with fluids before running the scenario.
Inez	Peripartum Hemorrhage Abruption	Inez is a 27 yr old primip @ 35 weeks. She arrives at hospital one evening crying and doubled over in pain. She is admitted to a birthing room and the nurse notices bright red blood on Inez's panties. She is having very intense and close contractions.	Runs for approximately 18 minutes. This scenario involves a precipitous delivery with heavy bleeding due to placental abruption. Abruption was brought on by use of alcohol and cocaine. Instructor must fill hemorrhage kit with fluids before running the scenario.
Janie	Peripartum Hemorrhage/PPH	Janie is a 23 yr old G 2 @ 38 weeks. She has experienced several bleeding episodes due to a low lying placenta. She has been counseled about the potential for postpartum hemorrhage. Her religious beliefs prohibit the administration of any blood products.	Runs for approximately 25 minutes. Even though patient suffers heavy blood loss, blood products are not to be used so alternative methods must be found. Instructor must fill hemorrhage kit with fluids before running the scenario.
June	Peripartum Hemorrhage/PPH	June is a 31 year old who is about to deliver her 5th baby. She has had a normal pregnancy and plans natural childbirth and breastfeeding. Her husband and oldest daughter attended the CBE refresher course and she has arranged for the 11 year old to attend the delivery.	Full scenario runs approximately for 22-25 minutes. PPH runs for 8-10 minutes. Labor and delivery progress uneventfully. Heavy bleeding begins immediately upon delivery of placenta. Uterus remains atonic despite Fundal massage and rapid pitocin infusion. Patient's BP drops to 80/50 and uterus begins to clamp down following Cytotec rectal insertion. <b>Note:</b> Running full scenario requires retraction of motor mechanism and insertion of PPH kit following delivery. To run PPH only, load June PPH from the scenario tab.
Kelly	Amniotic Fluid Embolism	Kelly is a 34 year old gravida 5/2 @ 38 weeks. She is scheduled for an induction due to problems with her last pregnancy. That baby weighed almost 10 pounds, and she had experienced a severe shoulder dystocia with the delivery. She has gained 43 pounds with this pregnancy and her glucose tolerance test (GTT) is borderline. An ultrasound a few days ago estimated current fetal weight at 3800 to 4000 grams. Kelly's physician feels that her cervix is favorable, so he has recommended an elective induction.	Runs approximately for 35-40 minutes. Shortly following SROM, Kelly begins having mild variable decels that are associated with contractions. Over the next few minutes, variables become more severe. Kelly vomits and begins gasping then suddenly becomes unresponsive. Decels now have late characteristics. She is given SQ Terbutaline and fetal baseline drops to 90. She is intubated and very difficult to ventilate. Patient develops VFib and then arrests. Baby is removed by emergency C-Section. Kelly remains on a ventilator. Baby is limp and unresponsive and requires resuscitation, and suffers repetitive seizure activity.

Patient Name	Labor Type	History	Overview
<b>Kimberly</b>	Amniotic Fluid Embolism	Kimberly is a 27 yr old multip @ 42 weeks. She began having contractions at home and now they are becoming stronger. She is excited to be finally going into labor. By the time the nurse completes admission Kimberly is requesting pain meds as her labor is progressing quickly.	Runs for approximately 30 minutes. In this scenario, patient has an AFE and resuscitation efforts are in vain. Fetus is delivered via perimortem C-section.
<b>Madonna</b>	Preterm Labor	Madonna is a 41 year old multip @ 31 weeks. She has experienced difficult pregnancies in the past and has one Down's Syndrome baby, so she is very apprehensive. An early U/S and genetic studies showed this baby to be a normal female. She has had several episodes of preterm contractions that resolved with LLP bed rest and oral hydration. She began contracting again 3 hours ago and again tried LLP bed rest and oral hydration, but the contractions have continued & become stronger.	Runs for approximately 35 minutes. Patient starts having contractions and is given SQ Terbutaline to stop contractions. She begins to have unfavorable reactions to the drug, and contractions break through about 30 minutes later. Another dose of Terbutaline is given, and her reaction is worse, and contractions break through yet again. Eventually, she's given Procardia which resolves the contractions and does not give adverse reactions.
<b>Maria</b>	Preterm Labor	Maria is a 30 yr old multip @ 27 weeks. She has an 11 yr old and has been trying for more children. She has had 2 miscarriages in the last 4 years and she lost both due to an incompetent cervix. This time a McDonalds suture was placed @ 14 weeks.	A run for approximately 15 minutes and it is a tocolysis emergency. Patient has reaction to overmedication of magnesium sulfate.

## Model (Newborn)

The physiologic cyanosis model adjusts vital signs and skin color in response to ventilation effectiveness. Set the physiological modeling state to deteriorate and monitor the provider's intervention effectiveness using the CPR monitor window. The model speed can be controlled through the Model warp-factor, which goes from 1 to 5 (1 representing real-time).



### Physiological Modeling State

**Pause:** Model will pause at the current state.

**Improve:** Model trend to a healthy state. Once the model reaches the complete Healthy state, the model will go to *Pause* mode.

**Deteriorate:** Model will trend to a severe cyanotic state. If ventilations given to the neonate are of correct depth and between 40 and 60 per minute, the vital signs will improve. Otherwise, they will continue to deteriorate.

### Cyanosis Levels

The facilitator can quickly jump to any of the three points in the state of the model.

**Healthy:** Pedi is pink with adequate oxygenation.

**Mild Cyanosis:** Pedi is bluish and vital signs are starting to deteriorate.

**Severe Cyanosis:** Pedi is blue, apneic and vital signs are rapidly worsening.

### Modeled Therapy

**Improve Gain:** Moving this slider will help increase or decrease the cyanotic response to ventilations.

**Oxygen:** By selecting an oxygen rate, the baby will improve faster with proper ventilation. Flow On must be selected to activate the oxygen response.

**Epinephrine:** Select the desired epinephrine dose and then select "Administer". Immediately the heart rate of the Newborn should rise and the dose on board should start diminishing over time. The dose should be enough to increase the heart rate a small percentage in order to help the oxygen delivery in the system, therefore helping improve the neonate with proper ventilation.

**Reset:** By selecting "Reset" the oxygen flow and the epinephrine dose onboard will be eliminated.

# Speech

## Prerecorded Sounds

NOELLE has over 90 pre-recorded expressions, which can be initiated with a single click on the Speech page. The collection of speech and other sounds was chosen to cover a wide range of simulated emergencies.

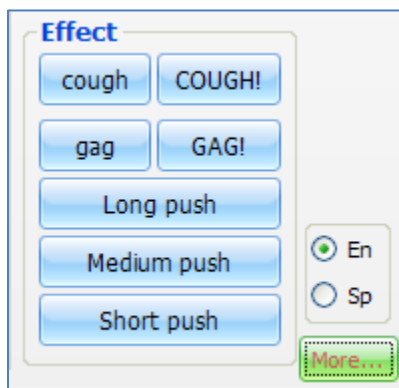


For ease of use, NOELLE's expressions are divided into categories: complaints, history, description, confusion, reply, location, occurrence, urgency, effect and numbers. These categories are labeled in blue on the speech page. Any of the speech items listed on the Speech page can be incorporated into labor and post-partum scenarios.

## Streaming Audio

Streaming audio makes simulation even more realistic. It allows the facilitator to hear everything the providers are discussing around the simulator. At the same time, the facilitator can interact as the patients voice for the provider t. The instructor will also be able to record his own speech phrases that can be used at any given time or within a scenario.


The streaming audio controls can be accessed by clicking on the "More" button in the lower right corner of the "Speech" tab.





Please verify that the simulator is connected by serial number on the Setup>Options>Environment FIXED field to access the more button.

## Streaming Voice Control window

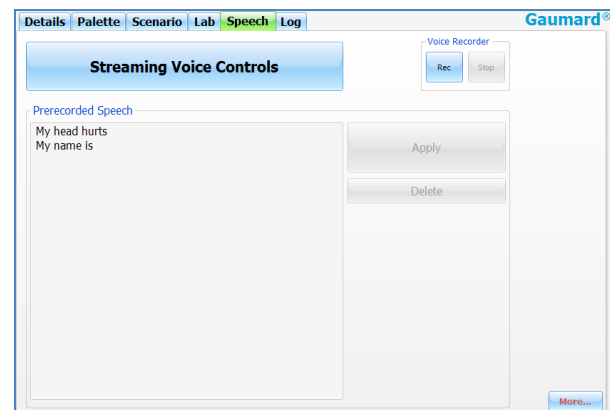**Streaming Voice Controls:** This button opens a new dialog box that is available to the user at all times. Selections on that dialog box include "Mute," "Listen," and "Talk". Select "Mute" to stop the communication; "Listen" to hear what providers are saying, or "Talk" to speak to the providers as the manikin's voice.



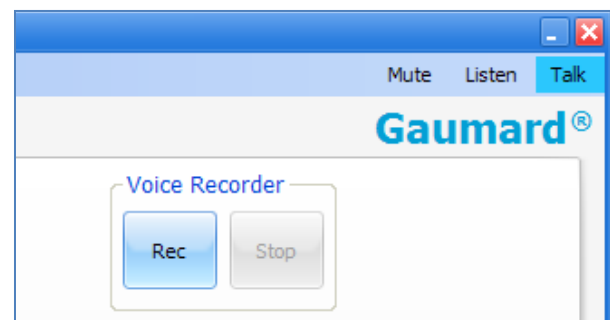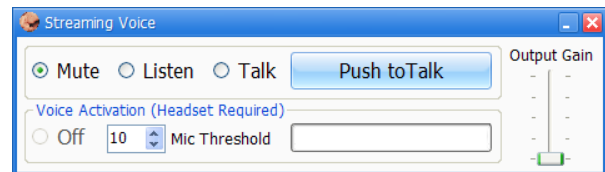
Mute, Listen and Talk controls are always available on the top right corner of the user interface.




It is recommended that the instructor uses a headset to improve sound quality. The headset also allows the user to use the "Voice Activation" so that at any time the instructor wishes to speak, it is sent directly to the manikin without user intervention.

The voice activation threshold can be adjusted. The "Mic Threshold" is used to adjust how sensitive the microphone is to the user's voice. The higher the threshold, the less sensitive the microphone is; and vice versa. For instance, if the threshold is set to high, users must speak loudly for the microphone to detect the audio.

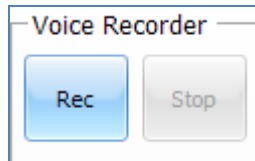**Voice Recorder:** The instructor is able to record his/her own speech phrases at any time. Once the instructor clicks on the “Rec” button the software automatically starts capturing everything that is said into the microphone.




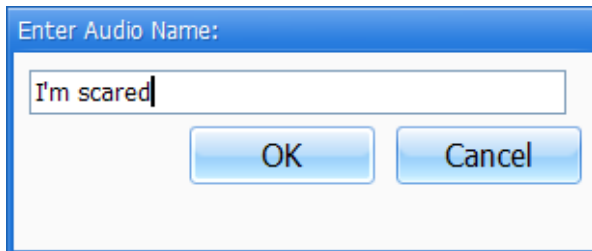Press “Stop” to finish recording. Use the Enter Audio Name window to name the speech phrase.



Enter a name and click “OK.” The phrase will now be accessible under the “Prerecorded Speech” menu.

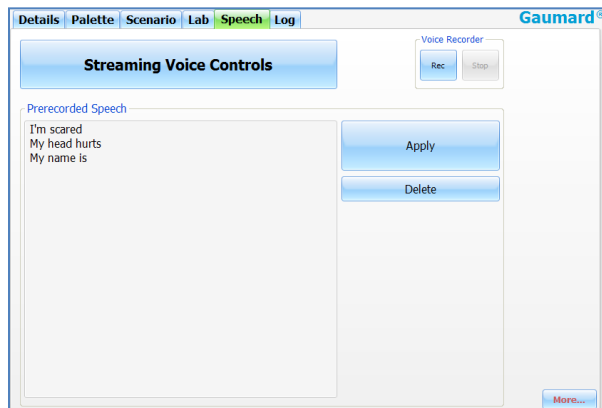**Prerecorded Speech:** On this section the instructor is able to play any of the custom phrases by either typing on the text box or clicking on the menu option. Once a phrase is selected, it can be played as the manikin’s voice or deleted.



**More...:** If the instructor wishes to play any of the Simulator’s prerecorded speech phrases, he will need to click on the “More...” button located on the bottom right corner of the Speech page.

- At any time the user can change the sound volume for sending and receiving. To change volume at simulator’s end, change the tablet’s microphone gain. To change volume on the tablet, use the speaker volume control.
- See troubleshooting section in the manual for more info.

## Log

The Log page allows the facilitator (Instructor or tablet operator) to keep track of every event during a session. It automatically creates an entry whenever a detected event occurs as well as every change in the condition of the patient. In addition, the facilitator can enter observed provider actions to the log with a simple click.

The Log page consists of four different areas (from bottom to top): the text log, provider action buttons, team logging buttons, and session info.

## Text Log

This is the large panel at the bottom of the Log Page, containing all of the time-stamped text entries. Every event that occurs in a session is recorded as an entry in the Text Log. The types of entries recorded by the log are categorized as follows: Actions, Applied Changes, Detected Events, Evaluations, Speech, and Notes.

### Actions

The term actions refers to tasks performed, by one or more of the providers, on the manikin during the session. The facilitator can quickly log actions from the Provider Actions section of the log page. To assign the entry to a particular provider or to the team in general the instructor may utilize the Team Logging feature. The following is an example of an unassigned Action entry:

"00:07:24 Action (Assess responsiveness)"

### Applied Changes

An "Applied" log entry occurs automatically every time a change is applied to the physiological condition of the manikin. In other words, each time changes are applied to the manikin from the Details page, the Palette page, or from a Scenario a log entry similar to the following is created:

"00:04:01 Applied (00:30): Details: Rhythm Sinus; Cardiac event 0; HR 80;"

## Detected Events

Each time one of the various sensors within the manikin detects a provider action, it is automatically logged as a “Detected” entry. These actions include intubation, BP cuff placement, artificial ventilations, chest compressions, and electrical therapy (pacing, defibrillation, cardioversion, inappropriate shock). The following example shows a detected log entry after a provider attempts to defibrillate:

"00:03:26 Detected (defibrillation): Shock # 2 - 300 Joules."

## Evaluations

Evaluations are added by the facilitator clicking on the “Satisfactory” or “Unsatisfactory” buttons on the Evaluation panel. The Evaluation panel is present at the bottom of the screen next to the clocks panel and is accessible at all times. Team Logging allows the facilitator to evaluate individual providers with a single click. For example, if provider Chris performed a procedure satisfactorily, the Evaluation entry would be:

"00:07:43 [Chris] Evaluation (Care Provided): SATISFACTORY"

## Speech

When the facilitator makes NOELLE speak by pressing buttons on the Speech page, an entry into the text log is automatically generated:

"00:18:10 Speech (Urgency): "Don't touch me"

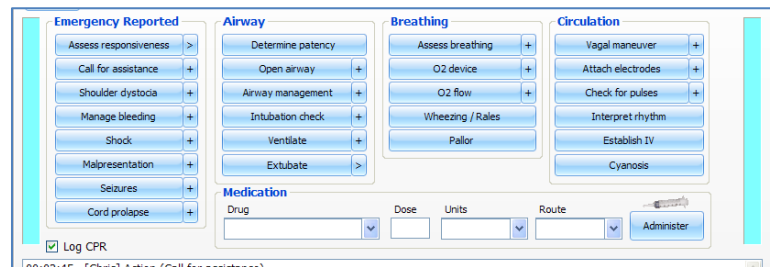
## Notes

Notes can be entered directly from the Evaluation panel or by right-clicking on the text log at any time by the facilitator. Note entries display each and every character the facilitator types into the text box. The following is an example of a Note entry:

"00:10:10 Note: provider took too long to assess patient."

## Provider Actions

The Provider Actions section refers to the collection of buttons in the middle of the log page. These buttons permit the facilitator to accurately and quickly track common provider actions.



The buttons are categorized into five groups: Emergency Reported, Airway, Breathing, Circulation, and Medication. Any time the facilitator clicks one of the buttons, a time-stamped log entry is generated documenting the action. For example, if the “Assess responsiveness” button is clicked when the session clock reads 00:07:24, the following entry is automatically generated:

"00:07:24 Action (Assess Responsiveness)"


## Special Buttons

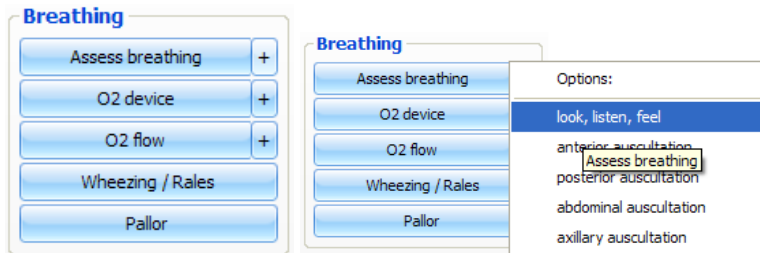Some provider-action buttons are accompanied by a special option button.

The first special button, “+”, lets the facilitator be a log actions in more detail. For example, if the button “Assess breathing” is clicked, the following entry is created:

"00:01:28 Action (Assess breathing)"


On the other hand, if the "+" button next to "Ventilate" is clicked, a list of additional options appears. The facilitator can be more specific and choose, for example, "look, listen, feel"...



...and the following entry is added:

"00:01:28 Action (Ventilate): look, listen, feel"

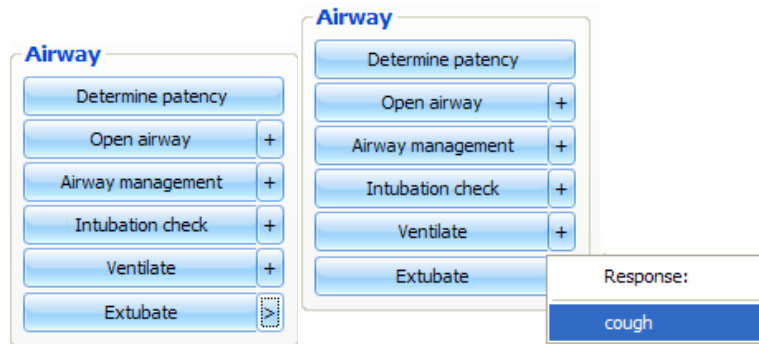The second special button, ">", allows pre-programming common responses to specific actions. For example, the facilitator can pre-program normal respiratory sounds and re-enable the lungs when the provider performs a needle decompression.




When the "extubate" button is, the following entry is created:


"00:01:28 Action (Extubate): cough"

## Medications

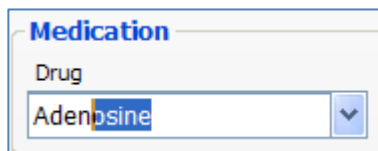The Medications section allows for quick and easy logging of drug administration, including dosage and route. The software comes preloaded with a set of commonly used drugs. Each of these drugs has a default dosage unit and a default route for administration (which can be overwritten by just typing over it). For example, Adenosine has the default dosage unit of "mg" and the default route is "IVP" (intravenous push). In order to enter, an administered dose 6 mg of adenosine via IVP, the facilitator need only enter the text "ad", which prompts the software to automatically search the drug list and display the best match (if any).



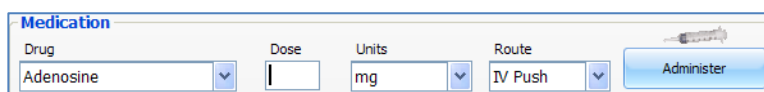After a drug has been selected, the action of clicking on the "dose" text field, the units and route fields with the default values for that particular drug are filled automatically.





The facilitator then enters the dose and clicks on the “Administered” button prompting the event to be recorded in the log. Following the example, suppose the dose entered was “6”:

"00:05:43 Action (Medication Administered): Adenosine, 6 mg, IVP"

## Team Logging

The Team Logging feature allows the facilitator to designate which member of the team performed a particular action. The Team Logging section is right above the Provider Actions section on the Log page.

Before beginning a session, the facilitator can add the names of all providers in the team to the team log. This is done by clicking on the **Add** button and filling in the **Add Provider** name field.

As shown in the Log Page image, a colored button is inserted onto the Team Logging region for the provider just added. The software allows up to six different providers, each with a corresponding button to be entered. Each time one of the provider buttons is clicked, the indicated person becomes the active provider in the evaluation window.

The colored vertical bars on either side of the log text box will match the color that is coded to the chosen provider. On the Log Page image, above, **Chris** is the active provider, so the vertical bars are turquoise. Notice that Patricia, who has a green colored button, is not the active provider. While there is an active provider, every time a Provider Action or Evaluation log entry is created it will have the name of the provider prefixed to it as follows:

"00:07:41 [Chris] Action (Check for pulses): radial"

To deactivate this feature, deselect the active provider, and return to general logging, click the "Team" button and the vertical bars will return to neutral color. All provider buttons can be edited or deleted by right-clicking them and selecting an option from the menu that appears.

## Session Information

The session info area contains the "Session Title" and "Facilitator" fields at the top of the page. At the beginning of each training session the session title and facilitator fields can be filled in and the information contained in them is stored with the text file when the log is saved or printed.

## Evaluation

The Evaluation panel, always visible at the bottom of the GaumardUI window, allows the facilitator to insert standard evaluations or other relevant notes into the log. The stylus device and hand-writing recognition technology makes annotation in real-time rapid and convenient.

Standard evaluations (satisfactory or not) are given context by their position in the log relative to detected and observed provider actions. The following example illustrates this idea.

The evaluation panel is part of the team-logging system, described previously. When a particular provider is selected, log entries generated via the Evaluation panel will be prefaced with the provider's name. For more information on Team Logging, see the previous section of this guide on Logging.

## Evaluation Form

The evaluation tool assists facilitators in reporting and assessing provider interaction using a questionnaire form. A completed evaluation form can then be stored as a digital document or printed for distribution.

### Using built in evaluation templates

Several evaluation templates are built-in to GaumardUI's evaluation tool. Each template includes a set of multiple choice questions related to the type of assessment being performed.

To begin, select an evaluation template from the **Load Template** dropdown. New templates created in the edit mode will also be listed.

Enter the name of the facilitator performing the review in the **On-Site Reviewer** field.

Complete the form by first entering the written response in the “Problems Identified” field and selecting the applicable multiple choice responses throughout.

	Yes	No	N/A
Calls for additional help (Nursing or Physician)	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Calls for anesthesia.	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Providers turn or assist in turning patient to side during seizure.	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Inserts tongue blade into mouth.	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Places supplemental oxygen on patient.	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Places pulse-oximeter on patient.	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>

Once the evaluation is completed, click the **Print** button located at the top right corner of the screen.

The 'Form' window displays the following information:

- plate:** Eclamptic Seizure (with a dropdown arrow)
- Clear** button
- Print** button
- Go To Edit Mode** button
- on Name:** Eclamptic Seizure
- Reviewer:** James C
- Date:** 08/08/2011

To save the finished evaluation as digital document, select **Microsoft XPS Document writer** and click **Print**. It is recommended that documents are first saved as XPS files before being printed into hard copies.

The 'Print' dialog box shows the following settings:

- General** tab selected.
- Select Printer:** Microsoft XPS Document Writer
- Status:** Ready
- Location:**
- Comment:**
- Print to file:** ☐
- Preferences** button
- Find Printer...** button
- Page Range:**
  - All** (selected)
  - Selection**
  - Current Page**
  - Pages:** [ ]
- Number of copies:** 1
- Collate:** ☒ (with page sequence icons 1, 2, 3)
- Print**, **Cancel**, and **Apply** buttons at the bottom.

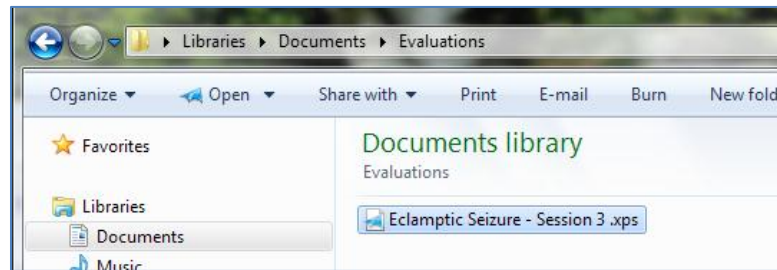
Enter the name of the evaluation, Select a location and click **Save**.

The 'Save the file as' dialog box shows the following details:


- Save the file as** title bar.
- Location: Documents > Evaluations
- Documents library** header.
- Message: No items match your search.
- File name:** Eclamptic Seizure - Session 3 .xps
- Save as type:** XPS Document (\*.xps)
- Save** and **Cancel** buttons at the bottom right.

## Printing an evaluation


If a print device is connected to the tablet, first select and open the evaluation document saved in the previous step.



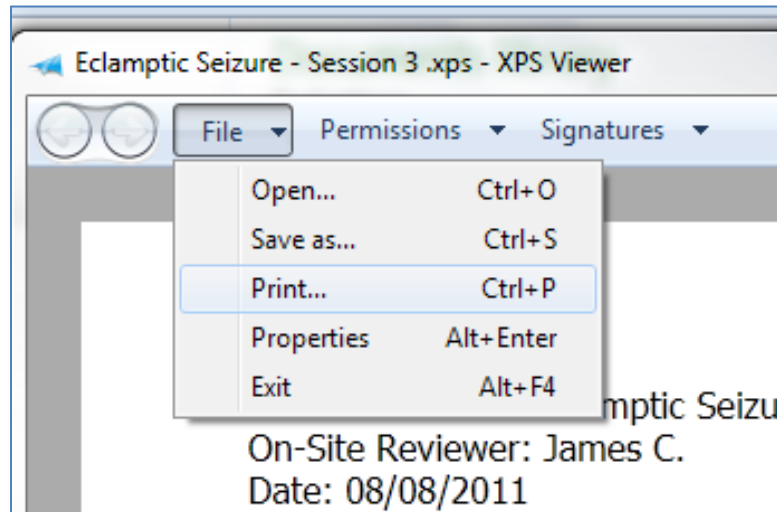In the XPS Viewer, navigate to the file menu and select print.



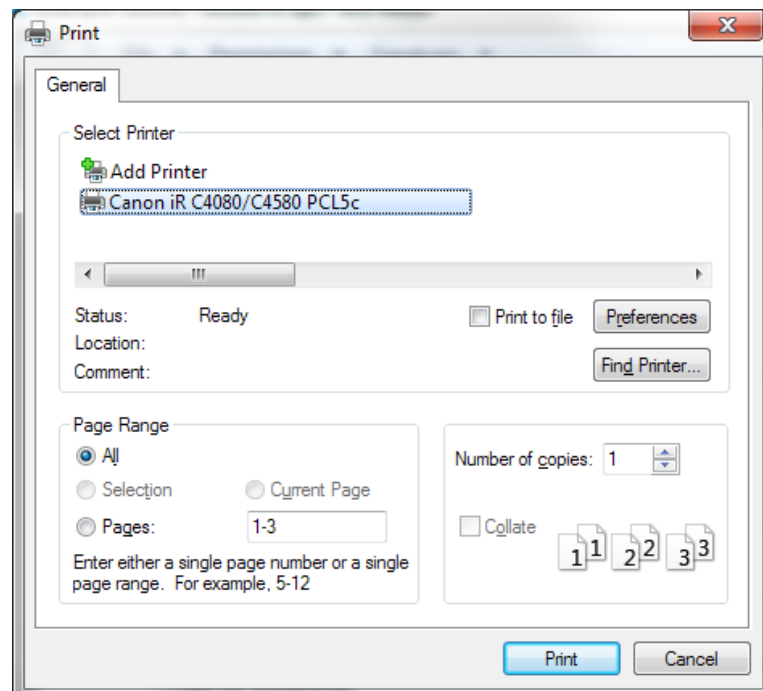Select the printer device from the list box and click **print**.



## Creating new evaluation templates

Each evaluation is based on an evaluation template. Facilitators can create new evaluation templates to tailor specific assessments. Template design and creation is done in the evaluation form **edit mode**. To enter the edit mode, toggle the Go To button located on the top right of the evaluation form window.

In edit mode, enter the name of the evaluation template in the **Evaluation Name** field.

From the **Category** menu, click **Add**.

Select the **Question with answer region** category type and enter the category information title. To save the changes, click **OK**.

The newly created category is shown.

Click, **Add** to include another category. Select the **Multiple choice** category type, enter the category information title and click **OK**.

Category Options...

1. Select Category Type:

- ☐ Question with answer region.
- ☒ Multiple choice: Yes - No
- ☐ Multiple choice: Disagree - Neither Agree or Disagree - Agree
- ☐ Multiple choice: Unacceptable - Poor - Average - Good - Perfect

2. Category Information:

Provider Tasks and performance

OK Cancel

Highlight the **Provider tasks and performance** category and then click **Add** from the **Multiple Choice Entry** menu.

Evaluation Name:

On-Site Reviewer:

Category:  Add Remove Insert ▲ ▼

Multiple Choice Entry (2):  Add Delete Insert ▲ ▼

1. Problem Description

2. Provider Tasks and performance

Yes No N/A

Enter the multiple choice information in the entry fields.

Category:  Add Remove Insert ▲ ▼

Multiple Choice Entry (2):  Add Delete Insert ▲ ▼

1. Problem Description

2. Provider Tasks and performance

	Yes	No	N/A
Physician/Provider introduces themselves to the patient.	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Positions baby with neck slightly extended.	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>
	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>

Repeat the previous steps to add more categories, questions and multiple choice options.

**On-Site Reviewer:**

Category: [Add] [Remove] [Insert] [▲] [▼]

Multiple Choice Entry (2): [Add] [Delete] [Insert] [▲] [▼]

How would you rate teamwork during this delivery/emergency?

**5. Communication**

	Unacceptable	Poor	Average	Good	Perfect	N/A
Overall Communication Rating						
Orient new members (SBAR)						
Transparent thinking						
Directed communication						
Closed loop communication						

**6. Situational Awareness**

	Unacceptable	Poor	Average	Good	Perfect	N/A
Overall Situational Awareness Rating						

After the evaluation template design is complete, click **Save** at the top of the window.

**Evaluation Form**

Load template: [▼] [Save] [Clear] [Delete] [Go To Evaluation Mode]

Evaluation Name: Vacuum Delivery

On-Site Reviewer: [Text Field]

Category: [Add] [Remove] [Insert] [▲] [▼]

Multiple Choice Entry (6): [Add] [Delete] [Insert] [▲] [▼]

Finally, select the new template from the Load template drop down to use the completed evaluation form.

**Evaluation Form**

Load template: [▼] [Save] [Clear] [Delete] [Go To Evaluation Mode]

Evaluation Name: [Text Field]

On-Site Reviewer: [Text Field]

Category: [Add] [Remove] [Insert] [▲] [▼]

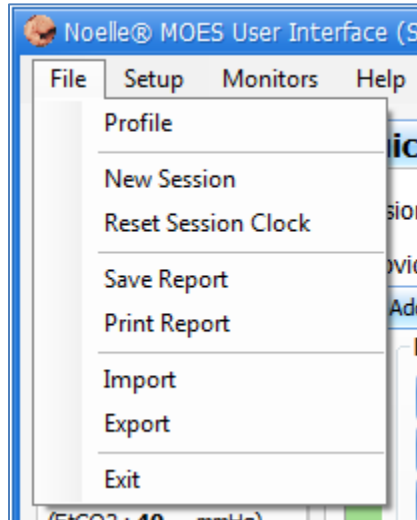
Multiple Choice Entry (6): [Add] [Delete] [Insert] [▲] [▼]

- Breech Vaginal Delivery
- Eclamptic Seizure
- Forceps Delivery
- Postpartum Hemorrhage
- Shoulder Dystocia
- Umbilical Cord Prolapse
- Vacuum Delivery**



# Menus

## File



## Profile

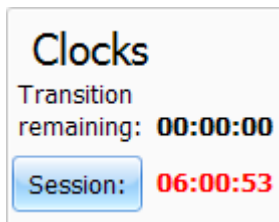
This option allows you to change your current profile. The profiles dialog box displays the available profiles. A software restart is not necessary to switch between profiles.

## New Session

Clicking New Session in the file menu will:

- Clear any loaded/playing scenario
- Clear any loaded/playing palette
- Reset vital signs to normal values
- Clear out log page
- Restart the session clock.


The session clock is located at the bottom of the dialog box.



The shortcut key for starting a new session is: **Ctrl + N**


## Reset Session Clock

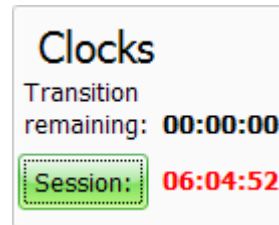Clicking on Reset Session Clock resets the clock back to zero. It does not have any effect on the transition time remaining on a scenario; it does not reset the vital signs, or clear out loaded scenarios. The facilitator can also reset the session clock by clicking on the Session button next to the session time.




## Save Report

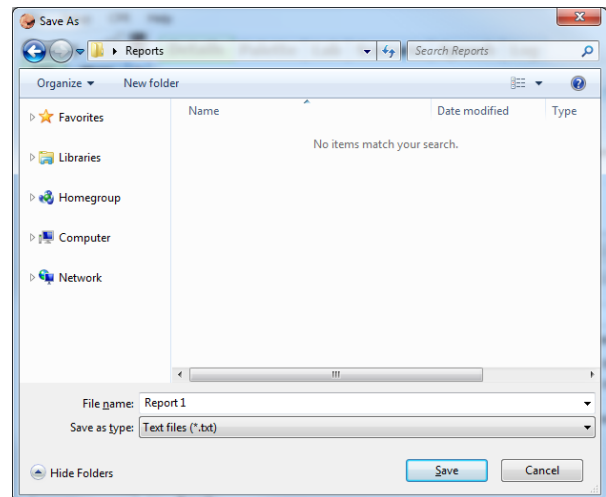This option allows you to save all the information recorded in the log page as a text file. Clicking on it brings up the “Save As” dialog box:



Select the desired name and path, and click “Save”.

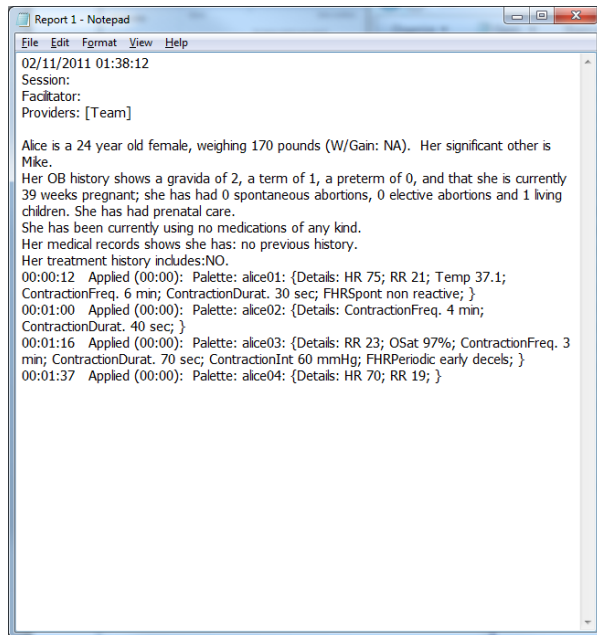The shortcut key for saving a report is Ctrl + S. For a sample report, look at the figure below:




## Print Report

This option allows you to print a text file containing all the information in the log for the latest session. Clicking on “Print Report” brings up the Print dialog box. The shortcut key for this option is **Ctrl + P**.

## Import


This tool allows the import of palettes, scenarios or modeling patients that may have been created on another tablet PC.

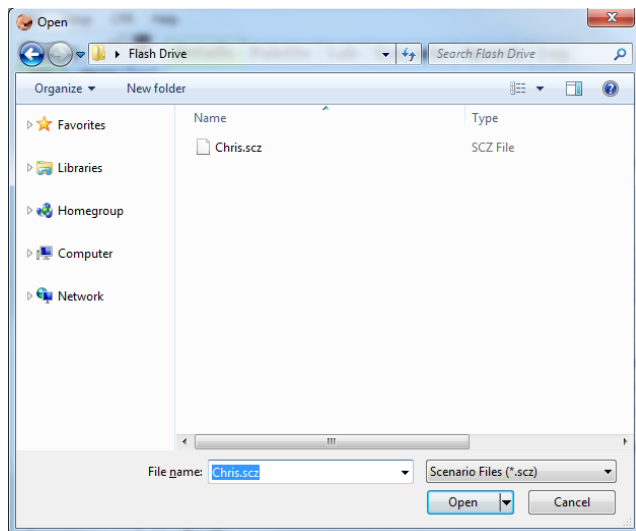When Import is clicked, the “Open” dialog box is displayed:




Browse to the location where the palette, scenario, or patient files have been saved and open it. They are automatically brought into the GaumardUI.

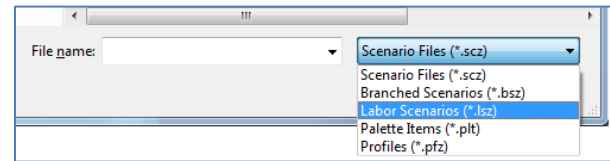Make sure that you have the correct file type selected:




## Export

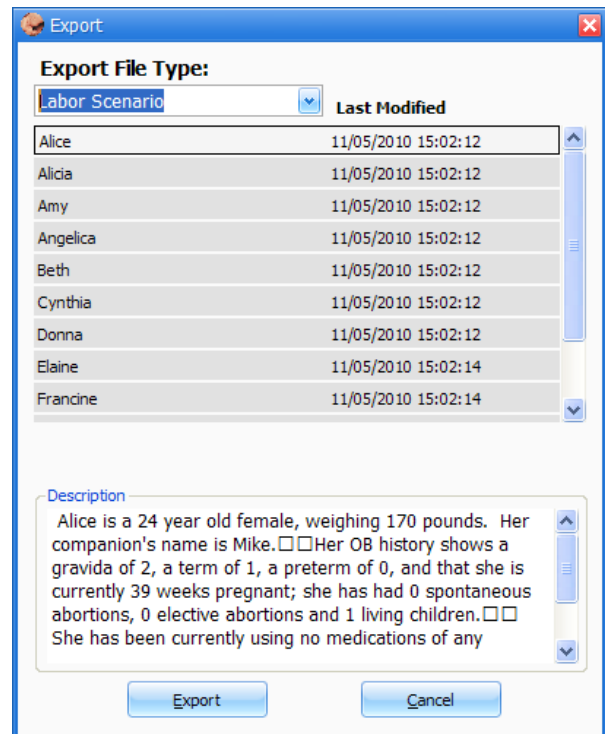You can export palettes, scenarios (branched or linear), and model patients. After selecting the kind of file to be exported, the following dialog box is displayed.



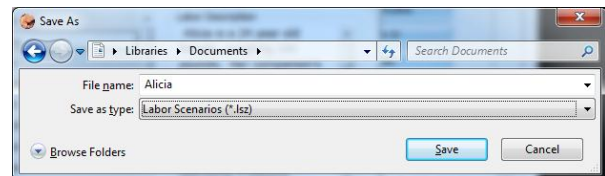Make a selection and click “Export”. The “Save As” window is then displayed.




Once the files are saved on to the flash drive, plug the drive to the computer where the files will be imported. From the GaumardUI, select import from the file menu.

Navigate to the location where the file was saved on the flash drive and click open. GaumardUI will copy the scenario to the computer during the import process.

## Exit

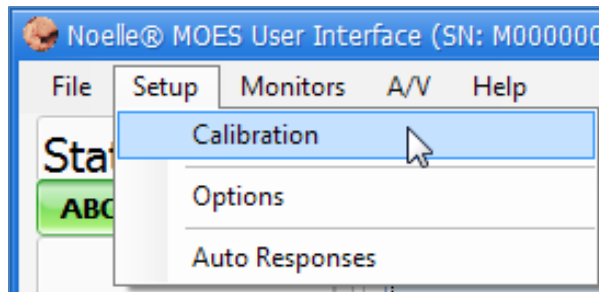
You can exit the software at any time by going to File, Exit or by clicking on the “x” button at the top right corner of the user interface.

## Setup



## Calibration

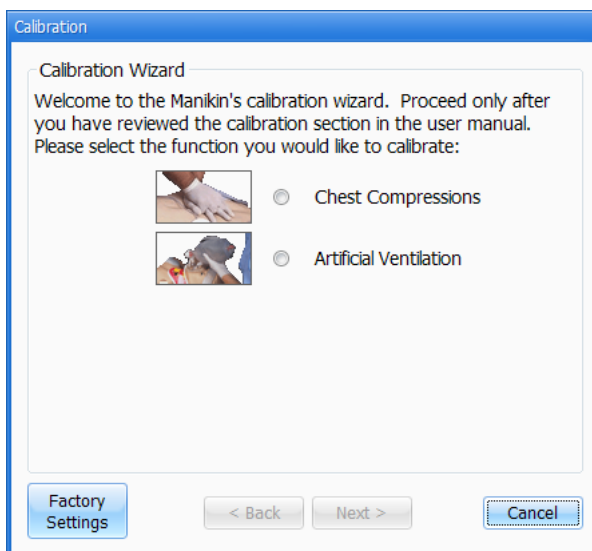
Use the calibration window to calibrate features and to reset sensors the simulator.



## Chest Compressions/Artificial Ventilations


This tool helps you calibrate the chest compressions and the artificial ventilations to your specific criteria. That is, you will be telling the system what a correct chest compression is and/or what a correct artificial ventilation is. Providers will be evaluated by the system based on this criteria.

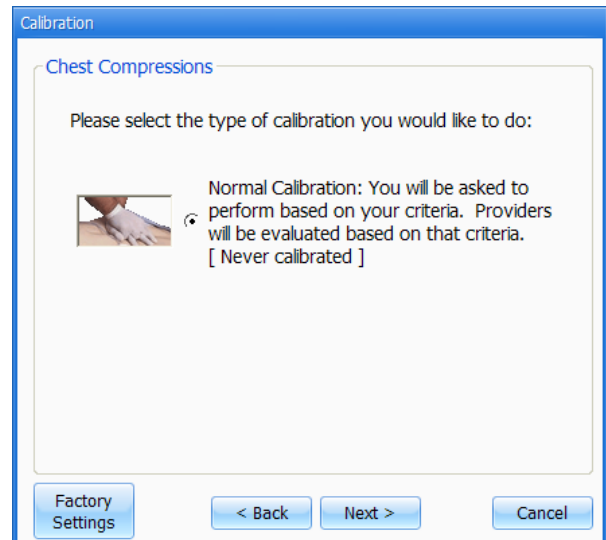The chest compressions and ventilations are calibrated the same way. After making a selection, this dialog box is displayed:



Click next to proceed with the calibration.

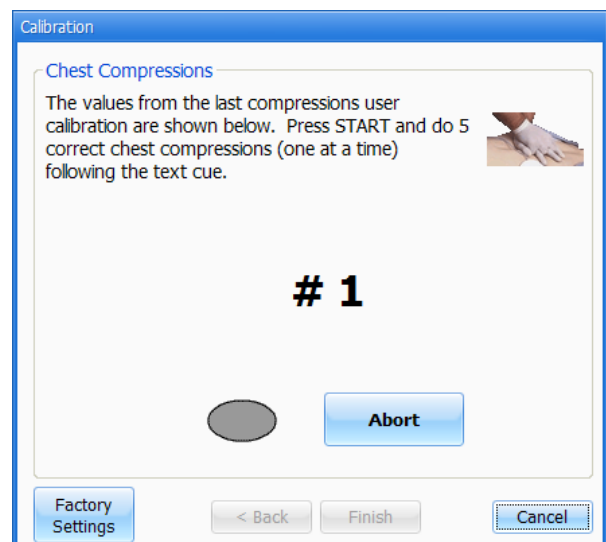The software will now ask you to perform a number of “correct” chest compressions or artificial ventilations, depending on what you are calibrating.



The facilitator should follow the text cue on the screen to perform just ONE compression or ventilation at a time, until prompted for the next one.

For example, if calibrating chest compressions:



1. The wizard prompts you with a “#1”.
2. Perform one correct chest compression.
3. A green filled oval indicates that the chest compression was successfully recorded.
4. The wizard prompts you with a “#2”.
5. Perform a second correct chest compression.


6. A green filled oval indicates that the chest compression was successfully recorded....
7. Continue the same process to finish the calibration.

At the end of the calibrating session, the wizard shows the average peak, depth, and duration values for the procedure. If you feel you performed the procedures correctly, click the "Save" button. Otherwise, press the "Back" button to repeat the calibration.

Notice that you can go back, abort or cancel at any time during the procedure.

## Factory Settings

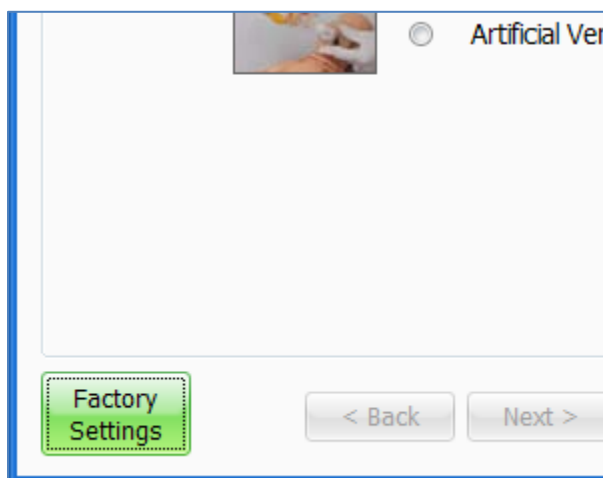Use the factory settings menu to restore sensors to the factory calibration. A factory reset will over-ride any calibrations performed by the facilitator.



Make sure that when you are restoring the sensors to the factory settings that no one is practicing chest compressions or ventilations. Any of these actions may interfere with the reset. Each time that one of these sensors is clicked a message will appear at the bottom left of the screen notifying the user of the status of the reset (OK, or TRY AGAIN). Should the sensor not respond, please refer to the [troubleshooting guide](#) or contact Customer Support.

At the end of calibrating a function, the Calibration Wizard resets the simulator for the changes to take effect and displays the message "Done". If the wizard displays the message "Can't reset", it simply means that the new calibration values will take effect next time you start the software.

If the changes need to take immediate effect, simply close the GaumardUI software, wait one minute (for the simulator to turn off), and then start the GaumardUI software again.

## NOELLE specific sensors

The options below are specific to the NOELLE simulator. Remove the birthing baby from the birthing mechanism before continuing.

**Reset Force sensor** – Reset the force sensor if the labor activity is displaying incorrect force readings.

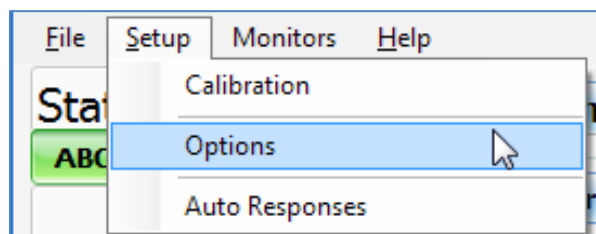
First remove the baby from the birthing arm and allow the motor arm to remain free for 30 seconds. While the birthing arm is free, click reset the force sensor button to record the current state as zero force.

**Reset Uterine Pressure** – Disconnect the uterus from the pressure port and click reset to restore the calibration to factory default.

**Reset Labor Motor** – Recalibrate the position of the birthing mechanism on the track. Always remove the birthing baby prior to recalibrating the birthing mechanism position.

## Options

The GaumardUI has several options that can be preset by the instructor.




NOELLE specific options include: NOELLE Features.

Newborn specific options include: Newborn Add-Ons and Neonate Features.

## Environment

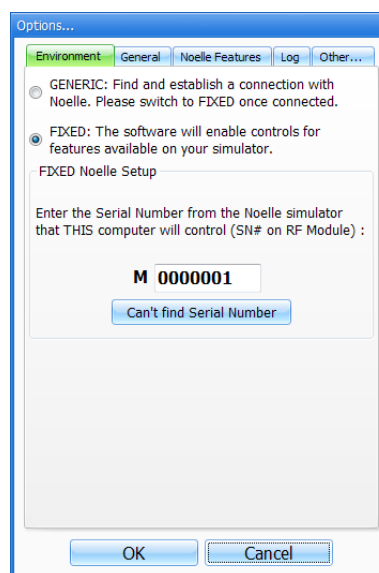If the **FIXED** button is selected, the simulator's serial number must be entered in the text box.

If you cannot find the serial number, first connect to the simulator using **GENERIC** and then press on the button that says "Can't find Serial Number" and follow the instructions.




## General


This tab allows the facilitator to:

- Enable auto saving of the log.
- Save your current log report.
- Enable stand-by warning.
- Select units (SI or English).
- Enable electrical therapy shock panel

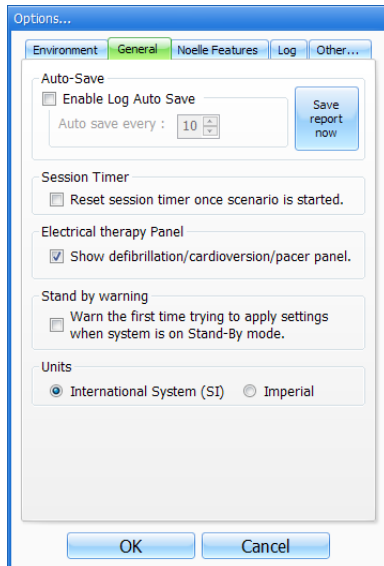
The shock panel is a floating window used for simulating electrical therapy. It can also be used in conjunction with “auto-responses”. For more information, navigate to page 87.



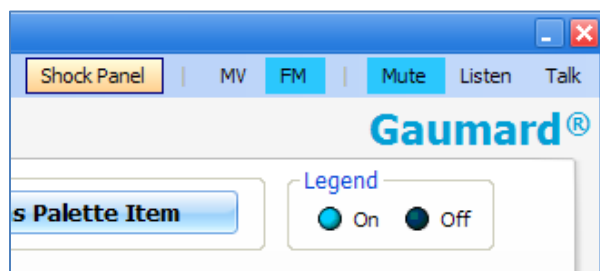If the Electrical therapy panel is enabled, a new control will be accessible from the top right of the screen.



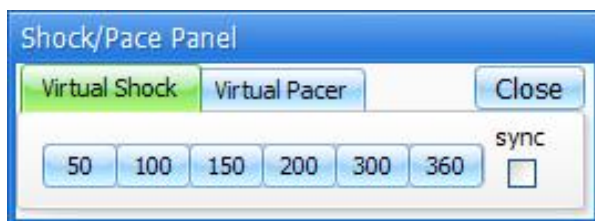
Use this control to bring focus to the floating shock/pace control window.



At any moment the shock/pace window can be moved or closed, while remaining accessible by clicking the button above.

## Log

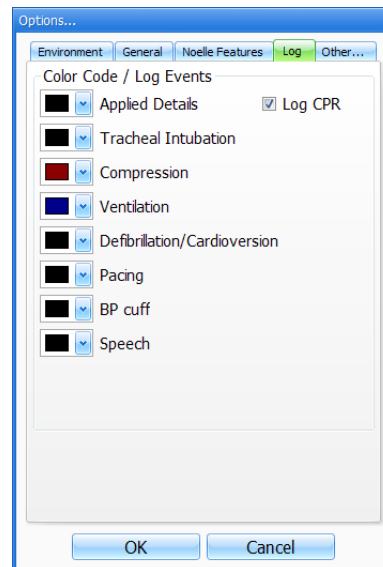Assign a color code for each of the log entry categories and click OK to save.




## CPR Options

In this tab you can:


- Select the number of desired compressions per minute.
- Specify the compression/ventilation ratio

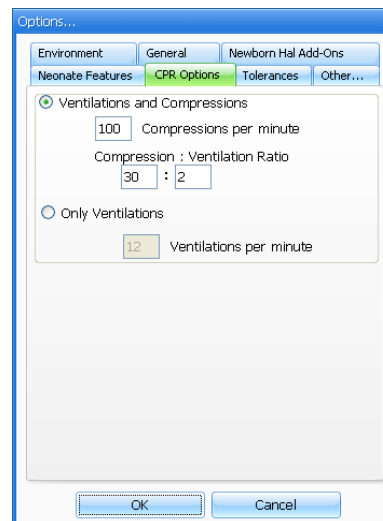Select number of ventilations per minute (if the “Only Ventilations” button is selected).



## Other

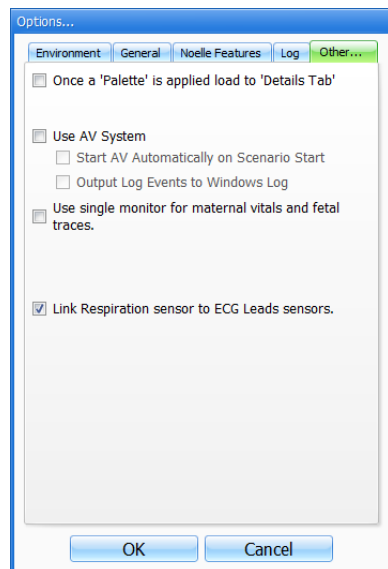The first option "Once a 'Palette' is applied load to 'Details Tab'" is useful for facilitators who want to keep track of the latest parameters that were updated using the Details Tab.




## AV System

GaumardUI is capable of interfacing with a number of third-party A/V recording systems. Checkmark the Use AV System to display the AV drop down as shown on page 94.

## NOELLE specific options

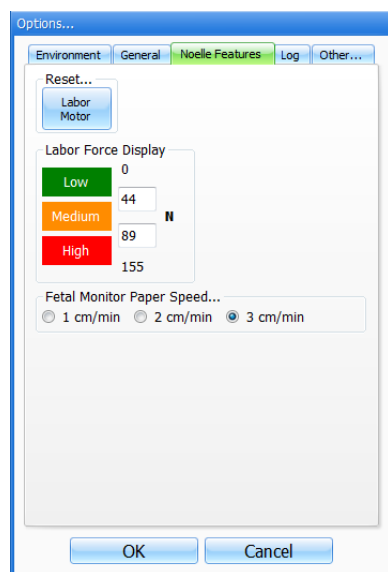Configure NOELLE specific features.



## Reset the labor motor

If the motor is not moving properly through the track during a labor, remove the baby and click Reset Labor Motor. The motor will move down the track and calibrate the proper starting position.

## Labor force display

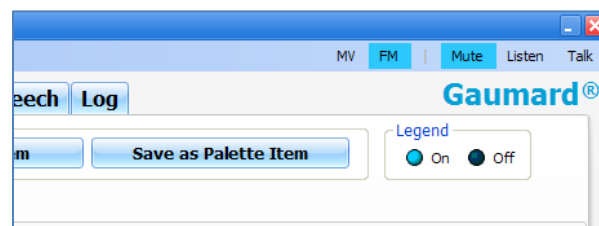
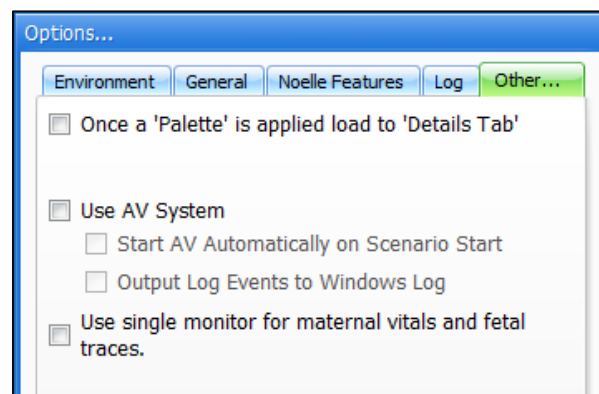
Set the threshold for the labor force display. To read more about the labor force display, go to Page 58.

## Fetal Monitor paper Speed

Use this feature to adjust the speed of the fetal monitor paper on the FHR monitor screen.

## Single Vital Monitors

You can display Gaumard Monitors using one monitor. To do so, select the checkbox "use single monitor for maternal vitals and fetal traces." When you enable this checkbox, the Maternal Vitals/Fetal Monitor control is shown in the GUI. These controls enable you to change the display on the vital signs monitor.




Select FM to display the fetal monitor and MV for maternal vitals. You can also have both screens auto change from one to the other. To do so, select the box labeled auto-flip and specify how often you will like the switch to take place.

## Options (Newborn specific)

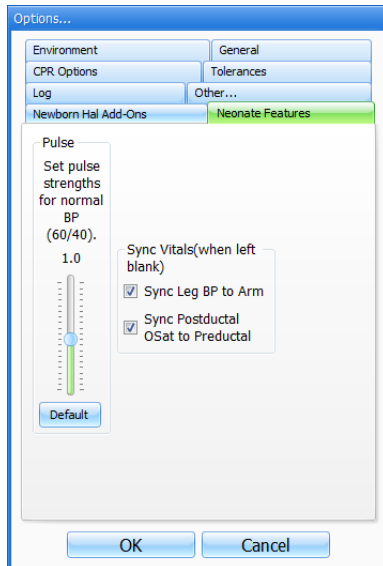The following options are specific to the Newborn control software.


### Neonate Features

Configure Newborn features available in the Newborn software.




## Tolerances

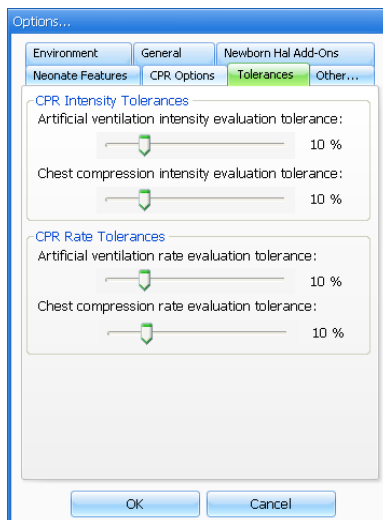This tab is used to select the tolerance and intensity of both chest compressions and ventilations.




## Auto Responses

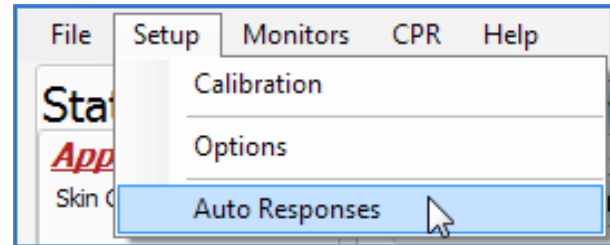The GaumardUI Module software can be configured to respond to electrical therapy when changes are applied using the detail list.



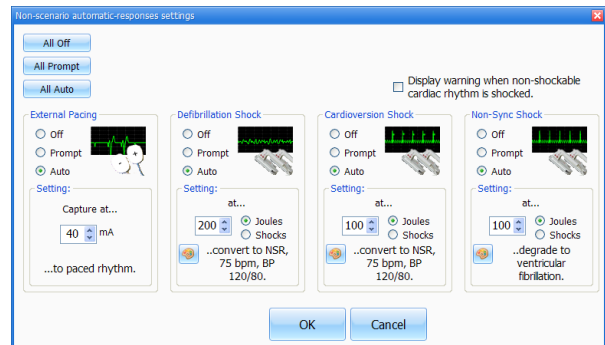Auto responses configured from the set-up menu respond to electric therapy administered when a scenario is not in progress. This type of auto response is referred to as a Non-Scenario Automatic Responses.




These three states are defined below:


- **Off** - The software **does not** respond to the electric therapy.
- **Prompt** - The software detects the electrical therapy and prompts the user if they would want to change the simulator's vitals to some preset healthy vitals.
- **Auto** - The software automatically detects the electrical therapy and compares it to a threshold selected by the provider, and once this threshold is accomplished the vitals automatically change to a healthy vital state.

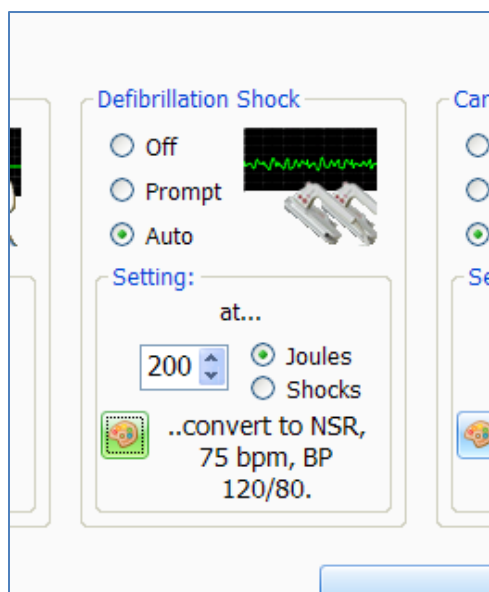Click on the palette button to program a specific palette to be applied after the electrical therapy.






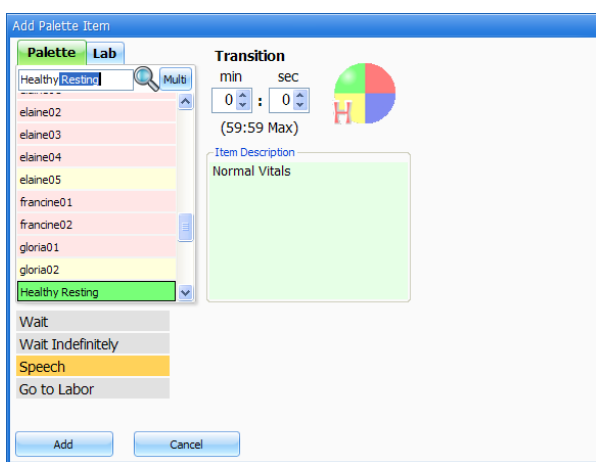The “Load Palette Item” window is displayed. Highlight the desired palette and click “Load”.



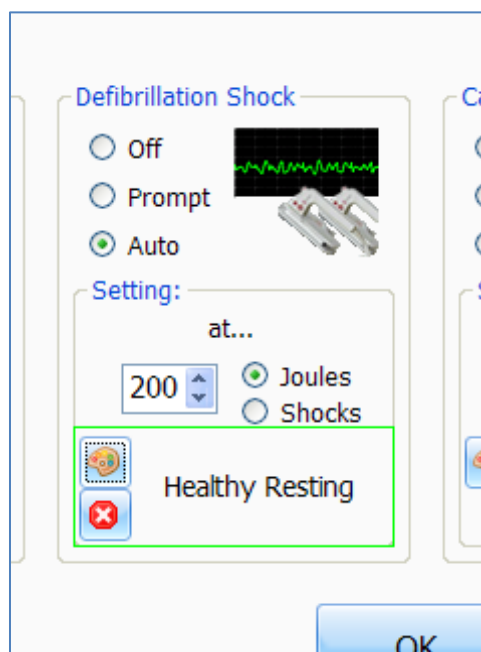
The desired palette is now displayed in the “Setting” section. You can delete the palette by clicking the “X” button. Deleting the palette defaults the electrical therapy to NSR, 75 bpm, BP 120/80.





## Monitors

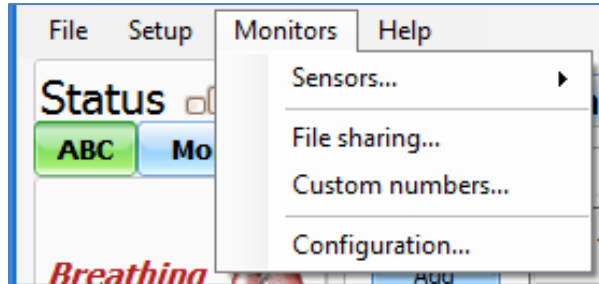
The optional Virtual Signs Monitor computer simulates a vital signs monitor attached to the patient. The vital signs information is sent through the Wi-Fi wireless network from the facilitator's tablet to the computer running the Gaumard Virtual Monitor software.



Use the Monitors drop down menu to enable/disable sensors on the virtual monitor screen, share files, program custom scalars and configure the connection between the GaumardUI and Gaumard Virtual Monitor software. For detailed instructions on how to configure and connect the Gaumard Vital Signs Monitors, refer to page **Error! Bookmark not defined.**

## Sensors

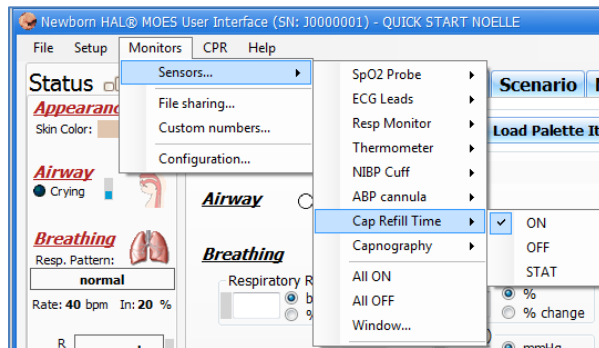This tool allows you disable any of the waveforms present in the vital signs monitor.




Select any of the waves that you will like to display and select "ON". If you want to turn any of them off, click "OFF". The vital signs monitor defaults to "All On."

Some sensors, such as NIBP and Thermometer are equipped with a STAT control that will allow the facilitator to activate readings on the virtual monitors from the controller software.

Another way to control the sensors is to go to Monitors, Sensors,

## Window

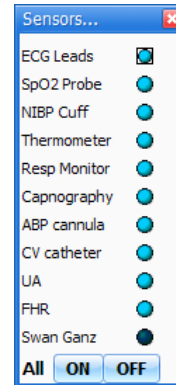This option brings up a floating window that can be viewed from any of the tabs in GUI. It can also remain opened as the users work in different scenarios. The sensors dialog box will look like the example below:




Light blue buttons indicate that a particular sensor is turned on and dark blue buttons indicate that a sensor is turned off. In the example above, all the sensors are turned on, except the thermometer and the ABP cannula.


## File Sharing

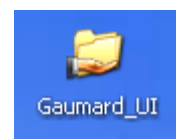File sharing is only available when using the Gaumard vital signs monitor. To use this tool you must first locate the "GaumardUI" folder on the desktop (of the tablet). Make sure you enable it for file sharing. This can be done by doing a right click over the folder, selecting 'Properties' and then enable sharing.


Add to this folder any kind of files that you wish to share with your students or providers.



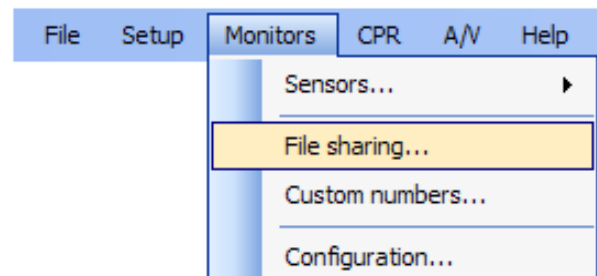
The File Sharing Setup menu is used to manage shared files.



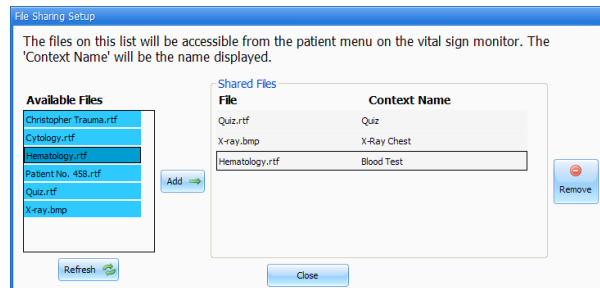Files in the Gaumard\_UI folder will be listed on the **Available Files** panel located on the left. To share a file, click on the **Add** button in the middle of the screen. Enter a context name on the pop-up menu and click **OK**. The shared file will appear on the right list box. Remove individual files by using the remove button on the right.




A yellow icon will be shown on the top left of the screen once a file is shared. This will inform the provider that a file is available for viewing.

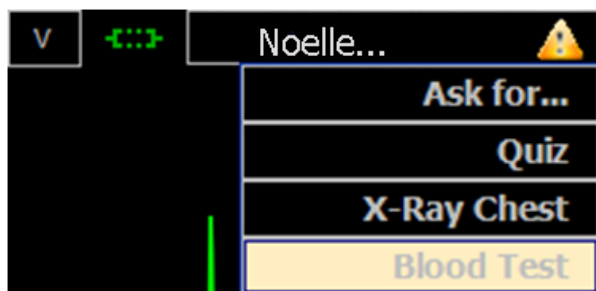


Click on the NOELLE button to bring down the selection of available files. Once a file is selected, it will automatically open on the Gaumard Monitor screen.

## Custom Numbers

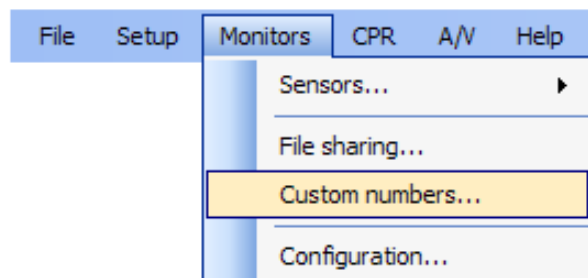
Use the custom numbers tool to add a new parameter, such as glucose level or platelet count, to NOELLE's virtual monitor.



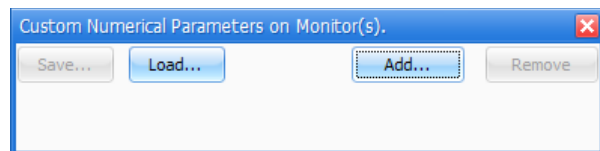
On clicking the "Custom numbers" option, the following window is displayed:



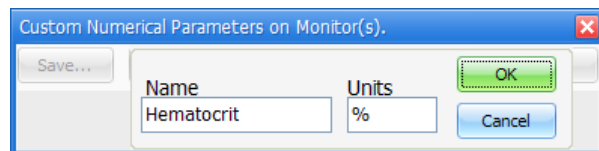Click Add to enter a new parameter for display in the virtual monitor. This dialog box is displayed:




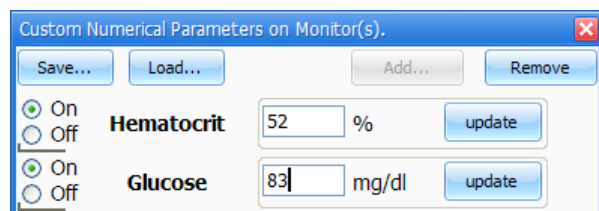Enter the name and units of the new parameter and click "OK". The new parameter is displayed. Several custom numbers may be entered at a time.




Enter the amount to be displayed in the virtual monitors and click "update".

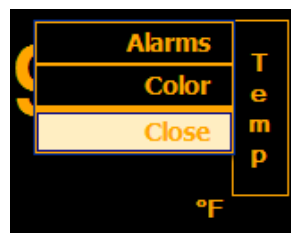
Have the student or provider close one of the parameters currently displayed by the virtual monitor by clicking on the button of the value to be removed, and selecting "Close".



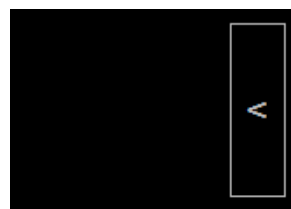Now the new custom number has a display slot.



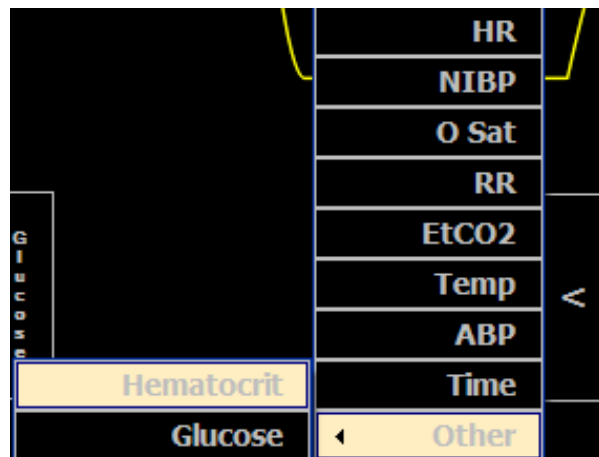
Click on the button and select "Other". Choose the value to display.




Do the same for as many new values as preferred for display. The figure below shows three new values: Glucose level, hematocrit, and TSH levels along the bottom of the display.



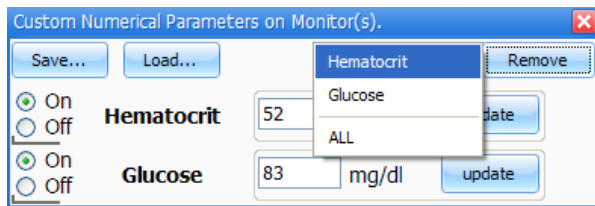Delete any parameter by clicking “Remove”.



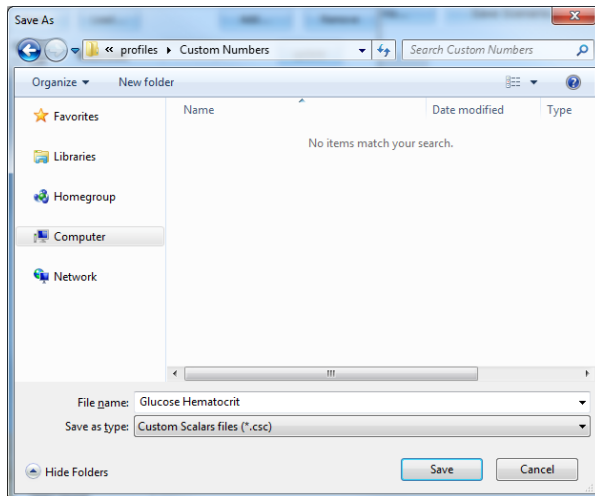
Save any list of added custom numbers by clicking on the “Save...” button. After clicking this button, the “Save As” window is displayed:



Type a file name, and click “Save”. You can load any of the pre-saved combination of custom numbers by clicking on the “Load” button.

## Configuration

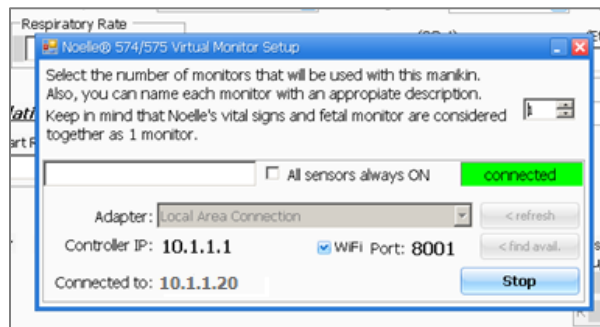
The configuration button is used to properly connect the virtual monitor to the tablet. Clicking on it, displays the following dialog box:



To properly configure this window and the virtual monitors software, refer to the Appendix on page 184.


## CPR

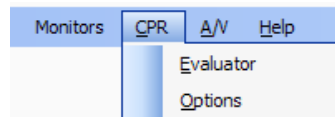GaumardUI features a CPR performance evaluator and trainer. From the menu bar, click on CPR and select Evaluator to access the CPR Evaluator window.



**Chest compressions will not provide feedback to the instructor if the cardiac rhythm is set to a healthy state, and ventilations will only be reported if the respiration rate is set to zero.**

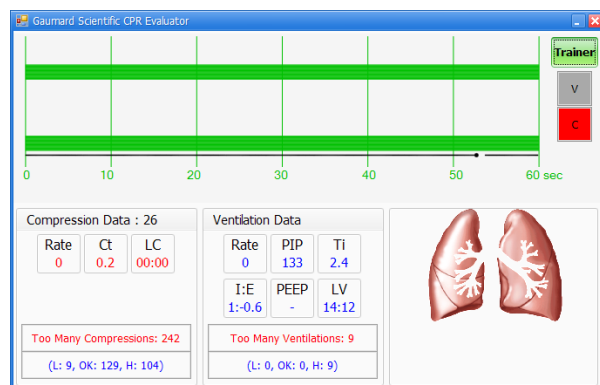The CPR evaluator feature provides real time feedback on the provider's compression and ventilation performance.




The provider performance indicator boxes are located on the right. The V (ventilation) and C (compression) box fill color changes between the following states:



- **Grey** – No intervention was detected.
- **Yellow** – Compression was too shallow. Ventilation was too weak.
- **Green** – Compression/ventilation was performed correctly.
- **Red** – Compression was too deep. Ventilation was too strong.

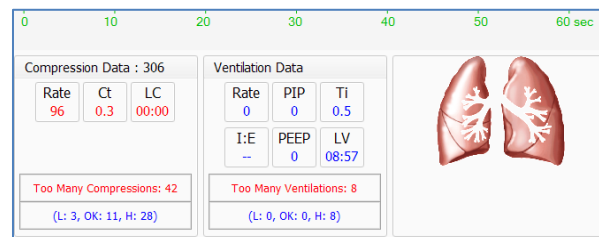Compression and ventilation data is displayed at the bottom of the window as CPR is performed by the provider.



### Compression Data

**Rate** – Rate of compressions in real time.

**Ct** (Compression time) – Average length of each compression in seconds.

**LC** (Last Compression) – Time elapsed since the last compression performed.

### Ventilation Data

**Rate** – Ventilation rate in real time.

**PIP** – (approx.) Peak Inspiratory Pressure

**Ti** – Time Inspiration

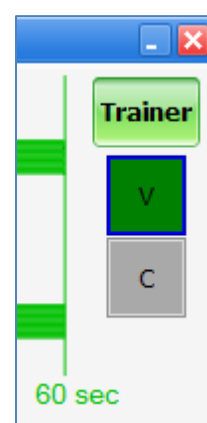
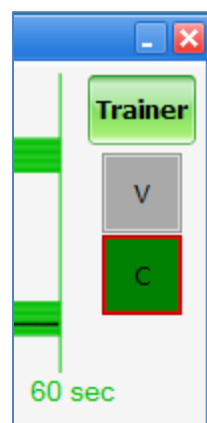
**I:E** – Inspiratory: Expiratory Ratio

**PEEP** - (approx.) Positive end-expiratory pressure.

**LV** (Last Ventilation) – Time elapsed since the last ventilation performed.


### Trainer

The CPR trainer features generates a visual queue of the compression to ventilation ratio programmed in the CPR Options menu. When the **Trainer** button is clicked, the V (ventilations) and C (compressions) box borders blink to indicate the correct reference CPR rate.




## CPR Options

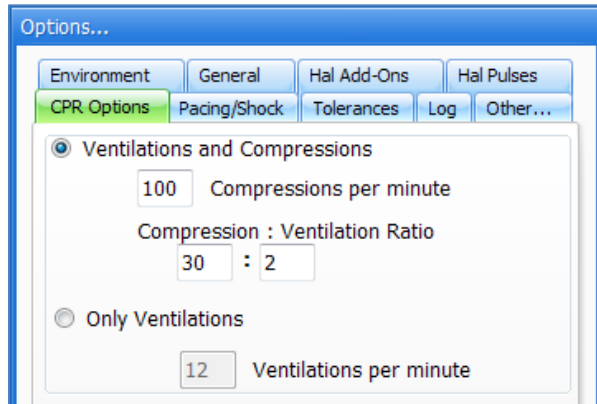By default, the trainer is configured to blink the reference borders at a 30:2 compression to ventilation ratio. To change the ratio, navigate to the menu bar and click the CPR dropdown and CPR Options.




### Performance Examples

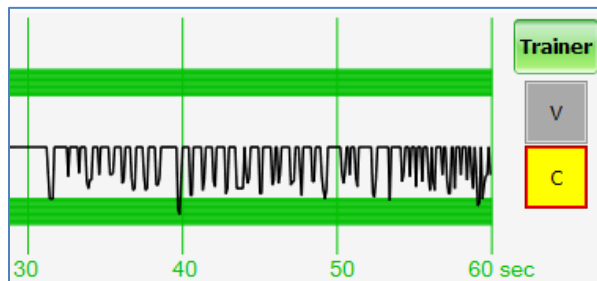
Compressions are too shallow. Waveforms mostly do not reach the green zone. Compression indicator is yellow.



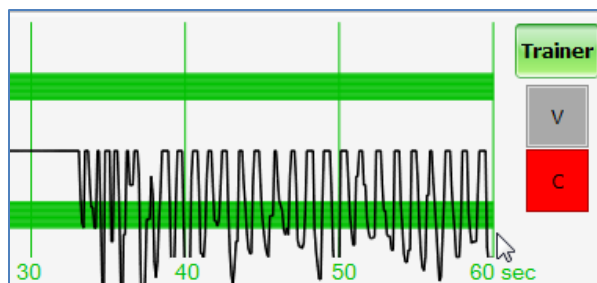
Compressions are too deep. Waveforms mostly exceed the green zone. Compression indicator is red.



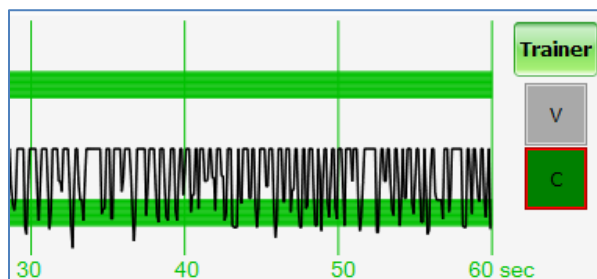
Compressions are performed correctly. Waveform peaks are mostly inside the green zone.



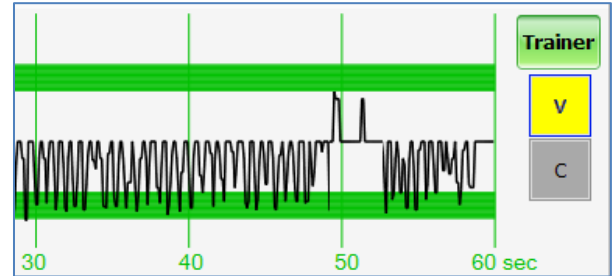
Ventilations are too shallow. Waveform peaks do not reach the green zone.



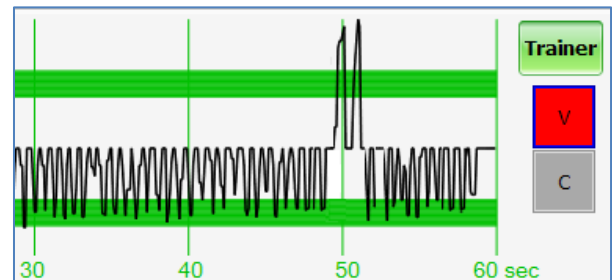Ventilations are too strong. Waveform peaks exceed the green zone.




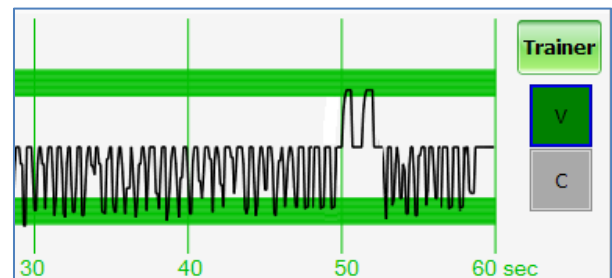Ventilation was performed correctly. Waveform peak is inside the green zone.



## Options

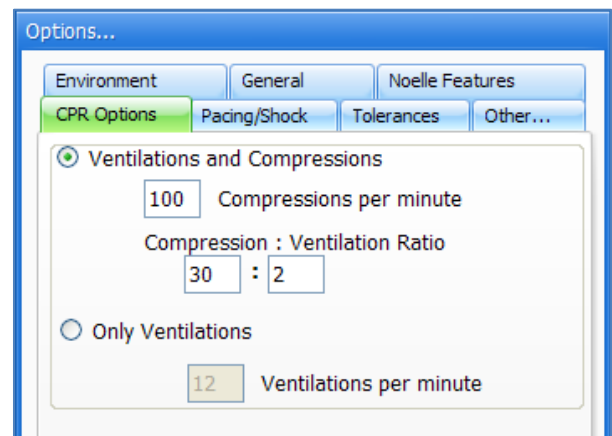
Navigate to CPR>Options to modify the ventilation and compression ratio.



## A/V

GaumardUI is capable of interfacing with a number of third-party A/V recording systems that enable the capture of Audio and Video interlaced with the events recorded in the software log.

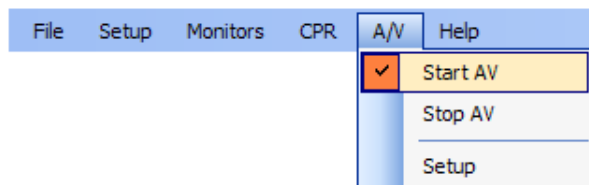Enabling the AV Link displays the following A/V menu:




Click “Start AV” to enable the Audio Video messages in the log system. GaumardUI also permits automatic sending of a “Start Record” message to the A/V Unit.

**Because it is possible to extend a simulation session beyond the last step in a scenario, the “Stop Recording” message does not have an “automatically stop” option.**
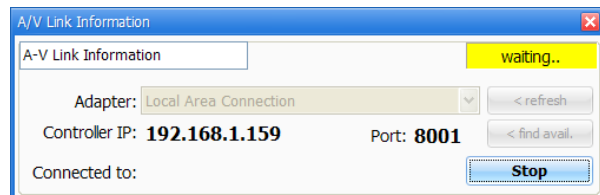
## Setup A/V

Clicking on A/V, Setup displays the following dialog box: This menu permits sending Start and Stop messages to the recorder, as well as displaying the connection status. In order to set up the connection on the A/V System side, please consult your A/V System's documentation.



## Help

The help menu has four options: Instruction Manual, About GaumardUI, Check for Updates and Diagnostics.

## Instruction Manual

Instruction Manual allows you to view a soft copy of the entire simulator help manual.

## About GUI

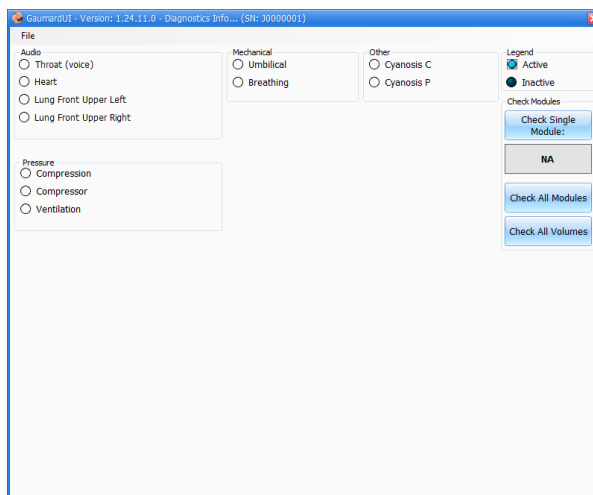
This window gives you the version of the software you are currently using and also the firmware version.



If updates are not available, the “Install” button is disabled and the following dialog box is displayed.

## Diagnostics

Use this feature as a troubleshooting tool.



For more information, go to the Appendix.

## Working with NOELLE

---



# Airway

## Nasal and Oral Intubation

Airway management techniques can be practiced on NOELLE including BVM, nasal/oral intubation, and suctioning. Endotracheal tubes, NG tubes and LMAs can be used.

Procedure	Recommended Device Size
Intubation(Blade size)	Miller 4 or MAC 3.5
LMA	Size 4
Nasal Intubation	8 mm outer diameter max
Oral Intubation	ETT Fr 7 or 7.5

### WARNING

Do not introduce liquids when performing nasal and oral intubation. Doing so can permanently damage the system.

Always lubricate tubing, airway and nasal opening prior to performing any nasal or oral intubation.

# Breathing

Software controlled breathing patterns: Kussmaul's, Cheyne-Stokes, Biot's, Apneustic, apnea, and normal.

## Pulmonary Ventilation

The airway contains nominal landmarks permitting either BVM or intubation exercises including the use of a LMA. The trachea extends to the bronchi and lungs.

## CPR

Use a normal size adult BVM which will seal around the mouth and nose. The ribs have normal anatomic landmarks and the lungs permit an adequate chest rise. Normal CPR procedures can be followed with aid of GaumardUI's CPR trainer.

# Circulation

## Bilateral IV arms

NOELLE has a bilateral IV training arms that can be used for bolus or intravenous infusions as well as for drawing fluids.



### WARNING

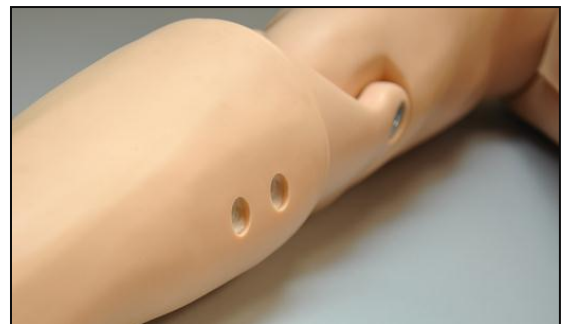
Do not attempt to fill IV system without the drain connector in place.

Always leave the drain port connected when injecting fluids into the system.

Use only Gaumard's provided simulated blood. Any other simulated blood brand containing sugar or any additive may cause blockage and/or interruption of the vasculature system.

Always purge the system with 70% isopropyl alcohol solution at the end of every simulation.

1. First, locate the fill syringe with tubing and the drain tube with pinch-clamp. Fill the syringe with the desired fluid -- water or simulated blood.



2. Connect the syringe with tubing to one port and the drain tube with clamp to the other port as shown.

**WARNING:** A drug recognition arm is equipped with a black drainage port. Reversing the fill and drain connections on a drug recognition arm will damage the system and void the warranty.



3. Leave the drain tube clamp opened and depress the syringe until all air has been pushed from the IV system and fluid runs from the drain.



4. To simulate a patient with no accessible peripheral IV sites, connect only the syringe. Pull the plunger to create suction, which will collapse the veins. Disconnect the syringe tube from the arm port while maintaining suction. The port will seal, and the veins will remain collapsed.

# Obstetrics

## Normal Labor and Delivery

### Birth Canal Maintenance:

- Ball point pens, ink and markers permanently stain the birth canal insert.
- Do not wrap this or any other Gaumard product in newsprint.
- The birth canal insert can be cleaned by wiping with a mild solution of soap and water. After cleaning, dust with talcum powder.
- Store the unit in a cool, dry place.
- After exercise is completed, **DO NOT** leave birthing baby in contact with the birth canal.
- Always lubricate the birth canal prior to delivery.

### Preparing for a delivery

1. Lubricate the fetal head and shoulders, plus the inside of the birth canal insert, using the silicone oil provided.


**WARNING:** Always lubricate the fetus and the birth canal before every delivery. Failure to do so will result in damage to the birthing mechanism and the birth canal.




2. Lubricate the head and shoulders of the fetus.



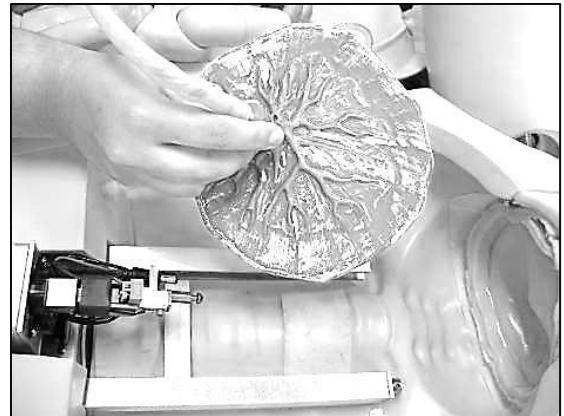3. Attach the umbilicus to the placenta.



4. Attach placenta to either side of the abdominal wall. Orienting the Velcro patches in parallel causes segments of the placenta to be retained, if secundes are reversed; orienting them at right angles causes the placenta to release with modest traction.

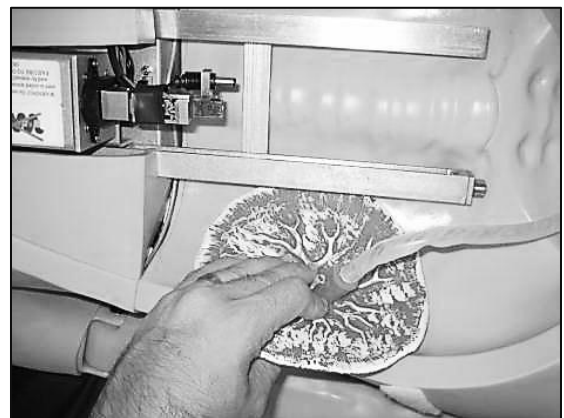

5. Attach the umbilical cord to the baby, route the cord so it does not bind in the mechanism and attach the placenta to the pelvic cavity using the Velcro® fastener. Note that the fetus has 2 receptacles at the perineum into which the matching pins located on the birthing mechanism are inserted.



6. Position the fetal arms and legs as shown.



7. Umbilical cord can be wrapped around the neck, demonstrating a nuchal cord.



## Locking Mechanism

The birthing mechanism arm attaches and locks into the fetus attachment port. The motor arm low voltage cable feeds power to the locking mechanism and fetal heart tones speaker inside the fetus.

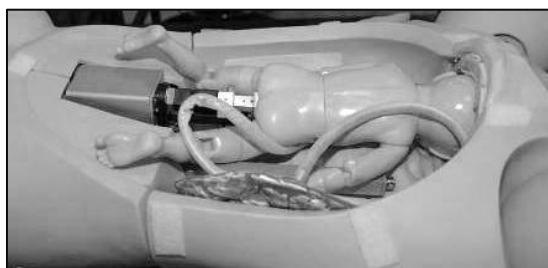
There will be one actuation of the locking mechanism when the articulating baby is inserted, which is normal.

**WARNING:** Never operate the birthing mechanism without the tummy cover in place.


1. Position the baby so that its face is upward (anterior). Connect the baby to the birthing mechanism while the manikin is "ON" so that the electromechanical mechanism allows the baby to be locked into place.



2. Once the baby is connected and locked onto the mechanism, position the baby to its initial birthing position. The user has the choice of ROP, LOP, ROA, LOA.



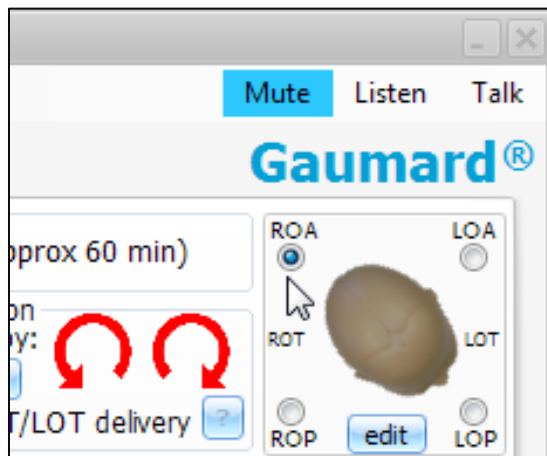
3. On the labor tab, select the fetal position that matches the birthing baby's initial position on the simulator.




## Loading a Labor Scenario

Several prebuilt labor scenarios are included in the Quick Start NOELLE Profile. Click the **Load Labor** button to load a labor scenario.

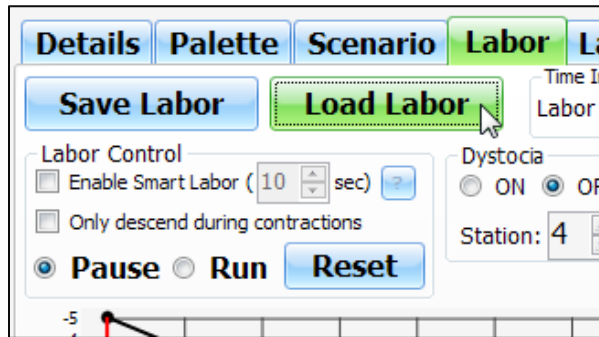**Tip:** Navigate to File>Profile from the menu bar to quickly change between profiles.



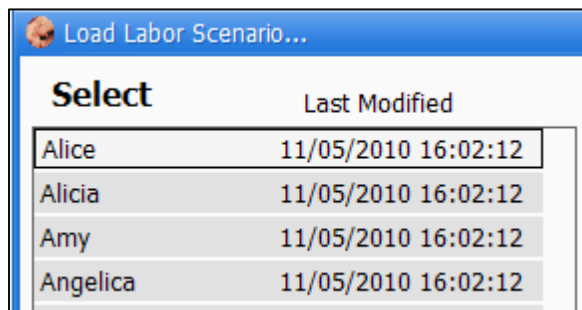On the **Load Labor Scenario...** window, select the **Alice** labor from and click **Load**.




## Starting the Labor

The labor process is started from the labor control panel. Increase the **Warp Factor** to simulate a 30 minute labor in a fraction of the time.

**WARNING:** Always place the tummy cover on the birthing cavity prior to starting the labor mechanism.

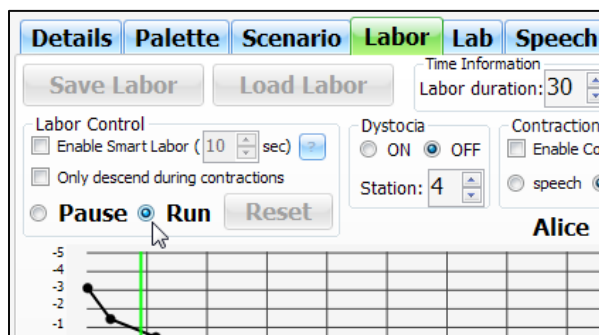Click the **Run** to begin the labor process.



The labor mechanism will descend as indicated by the green vertical progress bar and descent curve.




During delivery, fetal heart tones can be heard by placing the bell of a conventional stethoscope on the abdomen. Move it around until the tones are clearly heard. Tones are supplied via a small speaker located in the fetus.

## Completing a Delivery

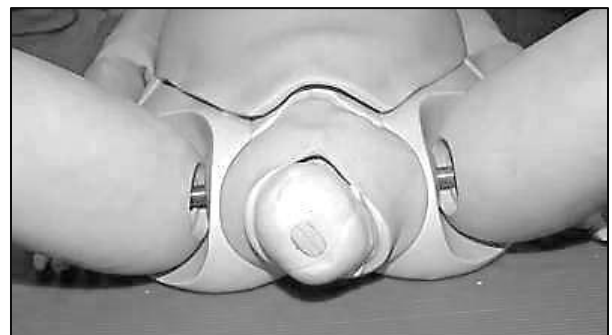
The first few centimeters of movement normally take about half the total delivery time. The baby rotates internally as it moves forward, after the head is delivered and before the shoulders are delivered.

The fetus is turning and the head is crowning.






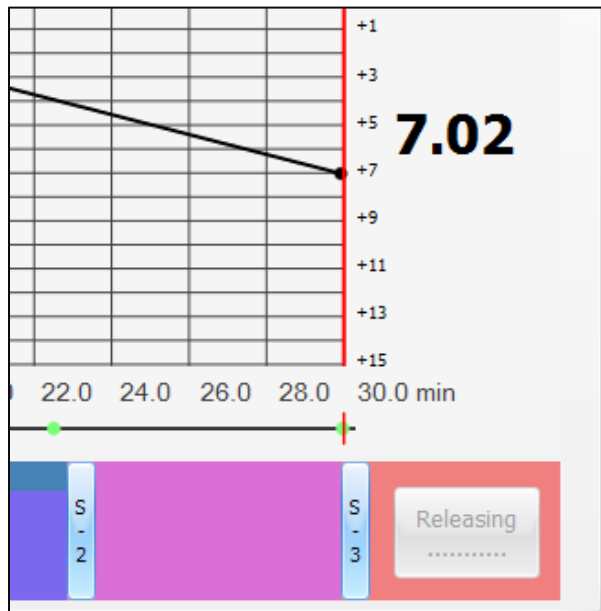The student or instructor may help the fetal head and shoulders through the vulva just as in real life. However, the fetus will **automatically unlock** only after the vertical progress bar reaches **Stage 3** of the delivery process. The automatic unlock process is indicated by the **Releasing...** status.



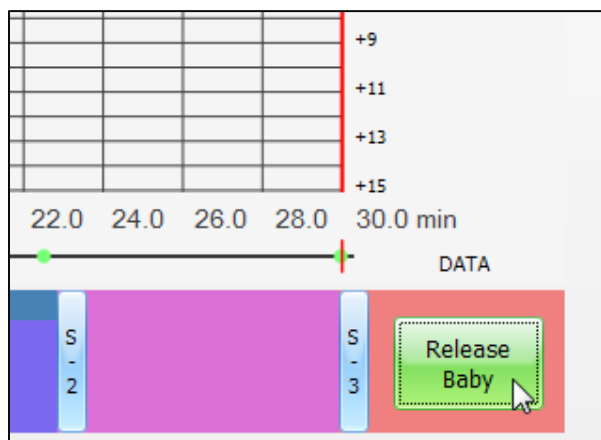Click the **Release Button** to manually disengage the baby.




**WARNING:** Do not pull the baby upward in contrast to the birthing mechanism's linear trajectory. Doing so can bend the motor arm and cause damage to the birthing mechanism.



## Resetting the Delivery Mechanism

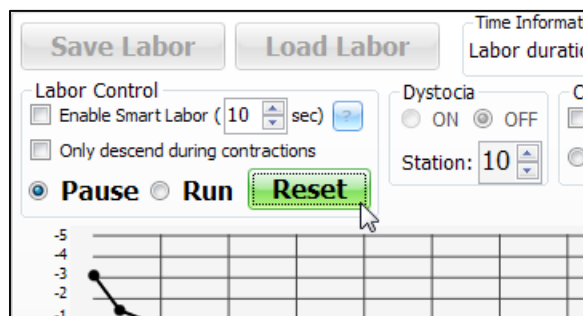After the delivery is complete, click the reset button to return the delivery mechanism to the initial position.

**WARNING:** Do not turn off the simulator until the birthing mechanism has fully retracted to its initial position.



## Vacuum-Assisted Delivery

Vacuum-assisted delivery is a technique for the management of arrest during the second stage of labor. Criteria for successful delivery include: (1) cervical dilation is complete; (2) cephalic presentation is confirmed; (3) the fetal head is no more than 1/5 palpable above the pubic bone; (4) effective uterine contractions continue; (5) maternal expulsive efforts continue. A soft skin scalp cover for the articulating baby is provided for vacuum-assisted deliveries.

**WARNING:** The soft skin scalp cover is meant to be used during the process of a vacuum delivery, and not resting against the cervix during non-simulation times. If pressure is left on the cervix from the scalp cover for lengthy times, both pieces will be damaged.

Always remove and store the soft skin after simulation is complete.

Vacuum-assisted delivery may be practiced with the NOELLE simulator using a vacuum cup available from a number of suppliers. Vacuum-assist device attaches to fetal scalp between fontanelles:



Await the next contraction that may be simulated by asking NOELLE to bear down, and have the student apply steady traction perpendicular to the plane of the cup. Some vacuum-assisted delivery devices are equipped with a means for measuring the amount of traction, which may be on the order of about fifteen pounds.

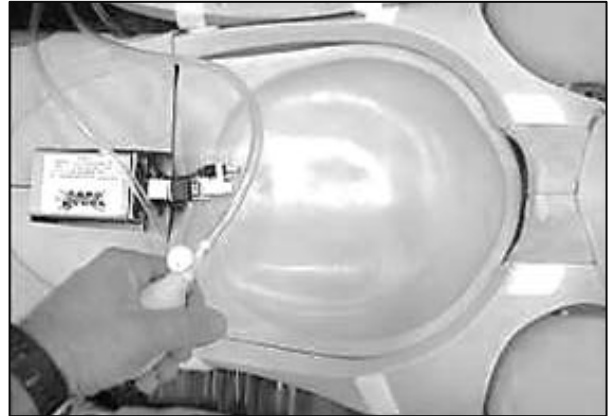
The software will also give the instructor a graph that indicates applied force, as the delivery mechanism is equipped with a strain gauge. The student must stop traction when the simulated contraction ceases. Repeat this procedure of waiting for the simulated contraction and providing traction during the contraction if and only if the fetus is moving down the birth canal with each contraction.

## Leopold Maneuver

To perform Leopold Maneuvers, retract the birthing mechanism fully and remove the articulating birthing baby. Place the elevating cushion within the birthing torso. Route the inflation bulb outside NOELLE through any space open on the left side. Place the birthing baby in the elevating cushion in the vertex, breech, or transverse positions. Install the “tummy cover”. Inflate the elevating cushion until the fetus can be felt under the abdomen cover.

**WARNING:** Do not enable the motorized birthing mechanism while performing this exercise.

Placing elevating pillow within simulator:



Place fetus onto elevating pillow and lift fetus anteriorly using the squeeze bulb:



Set abdominal cover into place:



Lift fetus anteriorly using squeeze bulb until it can be felt under the cover.



Conduct the four Leopold Maneuvers.



## Shoulder Dystocia

Shoulder dystocia is a dangerous condition defined in the NOELLE Guide as the “arrest of delivery of the fetal body after the successful delivery of the fetal head”. It may be characterized by the so-called “turtle-sign” wherein the fetal head moves forward and then retracts.

During dystocia, the fetal shoulders become wedged behind the symphysis pubis. NOELLE may be used to practice the resolution of dystocia using episiotomy techniques, the McRobert’s maneuver, suprapubic pressure, posterior arm sweep, or elbow-knee delivery.

The McRobert’s maneuver causes pelvic tilt that helps release the fetal shoulder from behind the pubic bone:



Suprapubic pressure may also release the fetal shoulder:

To demonstrate shoulder dystocia, place the fetal baby in the ROA position. Activate the delivery mechanism moving the fetus down the birth canal until the fetal head is delivered. Simulate dystocia by clicking the “Turn ON Dystocia” button on the Labor Tab. Once the dystocia mode is active, the fetal traces will automatically convert to real-time mode, and with each subsequent contraction there will be a “Turtle Sign”.

Students must use the various maneuvers including fetal manipulation to deliver the baby. Once the students perform all of the appropriate maneuvers required by the instructor, the dystocia mode can be deactivated by clicking on the “Turn OFF Dystocia” button. Once the dystocia mode is off, the fetal traces and labor warp factor will adjust to the previous settings, and the labor will continue resulting in delivery of the baby.

It is *very important* that the students are aware of the “Turtle Signs”. If, for any reason, the students fail to recognize the proper corrective procedures required, the labor can be stopped completely and set up again. The instructor can then add an unsatisfactory note to the log of the providers’ actions.

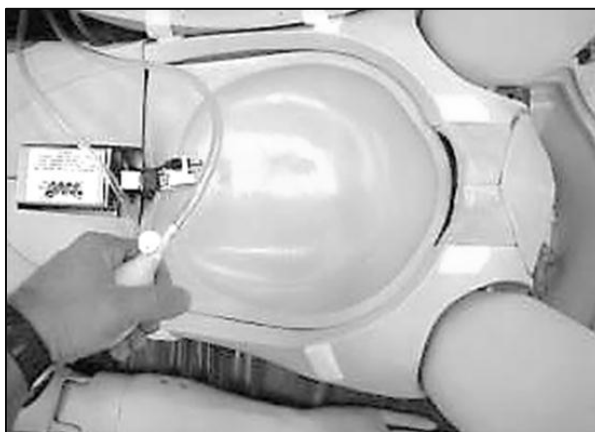
Normally, the fetus is retained by the delivery mechanism so that it can be rotated but not delivered. The facilitator can either turn off the dystocia mode to allow the baby to be delivered in the normal fashion, or press the “Release Baby” button on the Labor tab. This action unlocks the baby allowing students to pull the baby through the birth canal.

## Cesarean Delivery

Cesarean birth is the delivery of the fetus through an abdominal and uterine incision. A Cesarean delivery, also called a C-section, may be performed as a result of breech presentation, pre-term or dysfunctional labor, fetal distress, prolapsed umbilical cord, placenta previa, placental abruption, or a variety of other abnormalities.

Demonstrate a C-section using NOELLE by unfastening the snaps just above the pubic bone and birthing the baby between the tummy cover and the pubic bone. An optional abdominal cover is available if the Instructor wishes to demonstrate midline or “bikini” incisions.

Delivery mechanism fully retracted and inflatable cushion inserted:





NOELLE C-section delivery using abdominal cover with “bikini” incision. P/N S575.029



**WARNING:** Do not enable the motorized birthing mechanism while performing this exercise

## Prolapse of the Umbilical Cord

Prolapse of the umbilical cord is a dangerous complication which involves the presence of the umbilical cord in the birth canal in front of the presenting fetal part. This condition may occur as a result of breech presentation, transverse lies, a small fetus, an overly long cord, a placenta placed low in the uterus, or other abnormalities.



If the cord is observed in the birth canal ahead of the presenting part, gloved fingers should be inserted and the presenting part lifted off the cord to relieve pressure from the cord. This procedure must be maintained until the prolapse has been resolved, either by termination of the compression of the cord, or until delivery of the fetus by C-section.

## Placenta Previa

Placenta previa is a condition in which the placenta is in the lower half of the uterus, located near to or covering the cervical os. There are three types of placenta previa: Total, partial and marginal.

- Total placenta previa is when the placenta completely covers the cervical os.
- Partial placenta previa is when the cervical os is partially covered by the placenta.
- Marginal placenta previa is when the edge of the placenta extends to the internal os where the uterus opens into the vaginal canal.

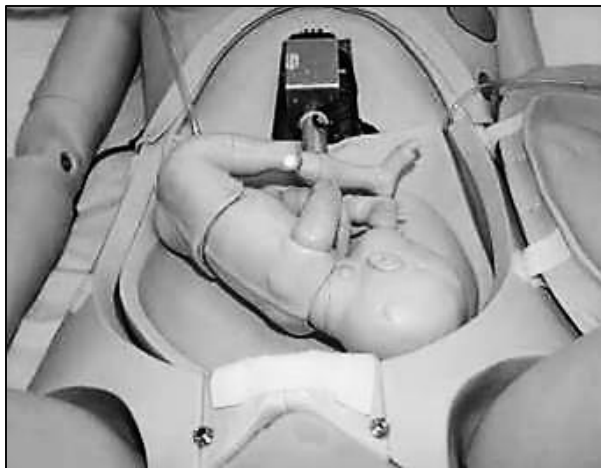
To simulate placenta previa with NOELLE, place the placenta in the desired position to simulate the condition with the maternal side against the uterine wall, or the cervical os. Then position the fetus within the uterine cavity with the presenting part closest to the placenta.



## External Version

Version may be attempted by the care provider to rotate the fetus from a breech position into one permitting normal vertex presentation. To practice “version” remove the abdominal cover and the fetus, retract the delivery mechanism fully and insert the inflatable cushion. Thoroughly lubricate the inside surface of the abdominal cover, the fetus, and the inflatable cushion.

Place the lubricated fetus onto the lubricated inflatable cushion and snap the lubricated abdominal cover into place. Inflate the cushion lifting the fetus anteriorly. Inflate the cushion at the base of the pelvic cavity to position fetus.



Confirm the breech position and attempt to manually turn the fetus within the uterus by trans-abdominal manipulation.



## Breech Birth

Breech birth occurs when either the buttocks or lower extremities of the fetus are the presenting part. There are three types of breech birth: frank, complete and incomplete, or footling.

- Frank breech occurs when the buttocks are the presenting part and the legs of the fetus are extended up toward the baby's head.
- Complete breech occurs when the buttocks are the presenting part and the baby's legs are flexed along the lower torso.
- Footling or incomplete breech occurs when one or both of the legs are the presenting part.

There are many differences in labor between the breech presentation and the vertex presentation. During the descent, the posterior hip encounters the pelvic floor and internal rotation takes place, allowing the anterior hip to move beneath the pubic arch. The anterior hip then delivers, followed by the posterior hip, the legs and the feet. External rotation allows the shoulders to move into the maternal pelvis and internal rotation allows the shoulders to deliver. Downward traction allows the delivery of the anterior shoulder, with a finger inserted into the birth canal to free the arm. Upward traction allows the posterior shoulder to deliver and the posterior arm is freed in the same manner. After the delivery of the shoulders, the fetal head delivers in a flexed or heads up position.

Although it is possible for a vaginal delivery of breech presentations, once a breech presentation has been confirmed, a Cesarean is often performed to lower the risk of infant mortality due to cord prolapse or birth asphyxia.

To simulate breech presentations with the NOELLE, retract the birthing mechanism fully, remove the cover in the fetal head, insert the birthing mechanism into the fetal head and place the fetal legs in either an extended position to simulate "footling" delivery or retract the legs for a "frank" delivery.

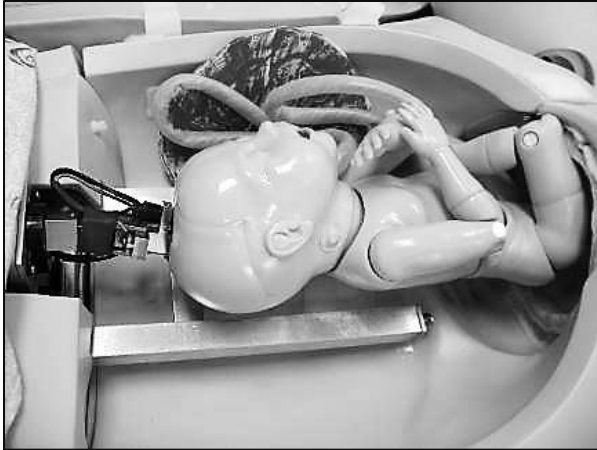
Remove plug in fetal head for breech delivery:



Removing plug reveals aperture for birthing mechanism:



Attach the fetal head to the birthing mechanism:



### Frank delivery



The Pinard or leg-flip maneuver frees one leg then another.



The fetal arms may also require a similar maneuver during delivery.



The fetal arms are delivered and the fetus rotated anteriorly to birth the head.

### Placenta Delivery

The placenta supplied with NOELLE may be positioned so that it births spontaneously, or requires either modest cord traction or manual removal. In addition, note that the placenta is designed with two removable placental fragments. These fragments are attached to the body of the placenta with Velcro. You may reverse one or both fragments causing one or both to birth with the placenta or remain affixed to the uterine wall.



Students must carefully inspect the birthed placenta to make sure it is complete and that no fragments remain internally. If retained fragments are noted, the student must retrieve them using a gloved hand under appropriate sterile conditions.

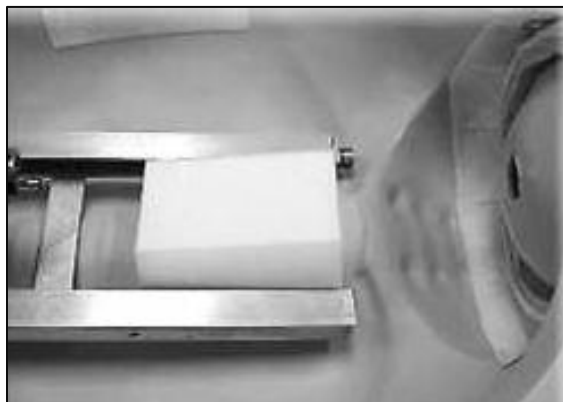
## Postpartum Activity

After delivery the uterus normally contracts reducing postpartum bleeding. Under certain conditions contraction does not occur and extensive bleeding may continue. If this condition is not recognized and treated in a timely manner the new mother may go into shock and die.

Inadequate uterine contraction may present as a “boggy” or soft uterus assessed through abdominal palpation. Uterine contraction may be augmented using certain drugs and/or uterine massage.


To simulate postpartum activity, follow these simple steps:

1. Connect the white Styrofoam block between the rails of the birth mechanism.



2. Place the uterus over the Styrofoam block lining up the Velcro attachment so that the uterus does not move around, and then insert the tip of the uterus inside the cervix.

**The cervix should fit on the ring on the tip of the uterus.**

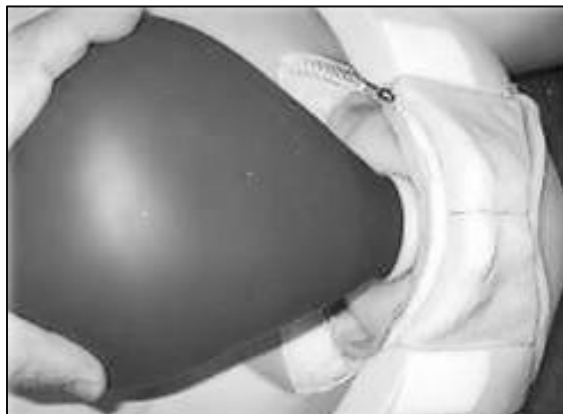

3. At the top of the uterus, connect the center hose to the air pressure port.

Should the instructor wish to use the tummy with the inflatable bladders the tummy can be connected to the extra port on the air hose.

4. Connect the lateral hose to the hemorrhage port inside the abdomen.

Once the uterus is in place, the instructor will be able to change the uterus pressure using the “Uterine Pressure” option inside the Details tab, or a palette item that has this setting saved.

Use bimanual massage to shrink the “boggy” uterus into a smaller and firmer condition.



## Post-Partum Hemorrhage

Noelle can hemorrhage real fluid from the cervix or birth canal. Follow the guide below enable the post-partum hemorrhaging feature.

**WARNING:** Always position the simulator so post-partum hemorrhage fluid flows away from the birth canal and the simulator itself.

**Do not allow PPH fluid to puddle beneath the simulator or reach the lower back.**

## Filling the PPH Reservoir

The post-partum hemorrhage fluid reservoir port is located behind the right knee. Fill the PPH reservoir with water or Gaumard simulated blood mix using the **PPH Fill Kit Syringe**.

The hemorrhage reservoir has a maximum capacity of 850 mL or approximately 15 fill kit syringes.

**WARNING:** Only use Gaumard provided simulated blood. Any other simulated blood brand containing sugar or any additive may cause blockage and/or interruption of the vasculature system.

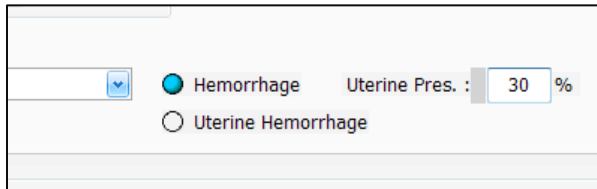
At the end of every simulation, always flush the system with distilled water to prevent clogging.

Do not overfill the reservoir.



## Enabling post-partum hemorrhage


To manually enable the PPH feature, go to the Details tab and click the hemorrhage feature icon **blue** (active).



Then, click the Apply **NOW** button to start the hemorrhage.

To stop the hemorrhage, click the control icon to black (disabled) and then click apply **NOW**.

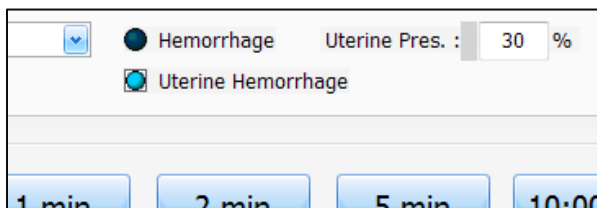Recall that when the uterus insert is connected to the hemorrhage port inside the abdomen, the instructor can simulate hemorrhaging from the cervix.



Several built-in labors are preprogrammed with post-partum hemorrhage activity. Go to page 60 for more information.

## Episiotomy Repair

Remove the fully dilated vulva used during delivery and select one of the three episiotomy repair modules. Snap a repair module into place. Use a "000" size suture and small curved needle to repair the surgical incision or repair.



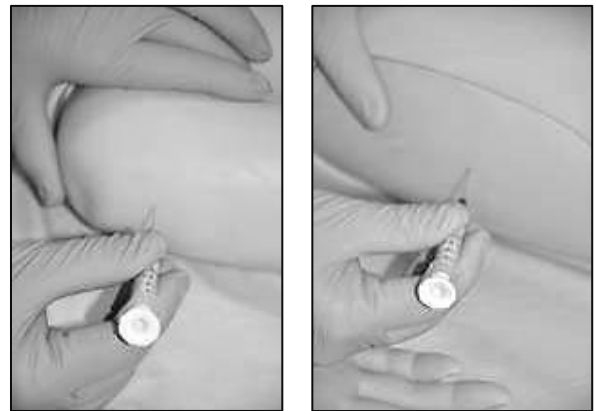
Episiotomy repair modules snap into birth canal:

"000" sized sutures are recommended to extend the life of the repair modules.

## Systemic

### Intramuscular injection sites

IM sites for placement exercises are located on both deltoids and quadriceps.




**WARNING:** Do not inject fluids into the IM sites.

## Other

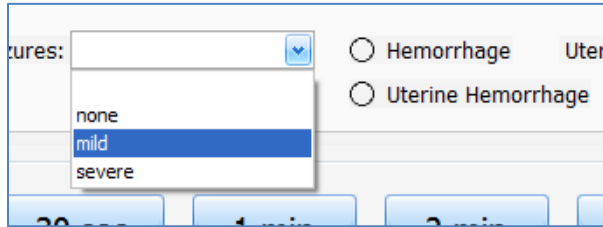
### Seizures

NOELLE is outfitted with a convulsion mechanism that can be used in conjunction with scenarios, preeclampsia for example.



## Vital Signs Monitor

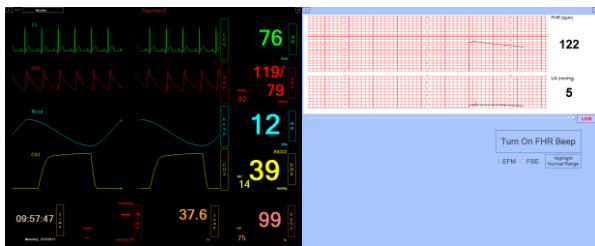The vital signs monitor simulates a vital signs monitor attached to the simulated patient. The vital signs are synchronized through a wireless network between the facilitator's tablet and the computer running the monitor. You can customize each trace independently of each other; users can set alarms, time scales, boundaries and grid options.



For information on how to setup Gaumard Monitors with GUI, please refer to the Appendix.

### Micro + (optional)

The Micro+ system is an all-in-one debriefing and simulator control solution for facilitators working in a lab or mobile environment. The system combines GaumardUI and powerful audio/video recording software in one control laptop.



## Working with Newborn

---



## Airway

Newborn's airway can be intubated orally using LMA or endotracheal tubes and nasally using a nasogastric tube.

Procedure	Recommended Device Size
Intubation (Blade size)	Miller 0
LMA	Size 1
Nasal Intubation	8 Fr catheter
Oral Intubation	ETT 3.0 no cuff, 6 Fr suction catheter

**WARNING:** Always lubricate tubing, airway and nasal opening prior to performing any nasal or oral intubation. Failure to do so will make intubation very difficult and is likely to result in damage.

Do not place silicone oil directly into the mouth and airway.

The nasogastric intubation feature is used for placement techniques only. Fluids cannot be inserted through nasogastric tubes as internal damage will result.

## Breathing

Control the breathing and lung sounds synchronized with selectable breathing patterns.

Use a BVM which will seal around the mouth and nose. The ribs have normal anatomic landmarks and the lungs permit an adequate chest rise. Normal CPR procedures can be followed with aid of GaumardUI's CPR trainer.

## Cardiac

Newborn is equipped with several realistic heart sounds which are tied to heart rates and selectable rhythms.

## Circulation

Newborn's software controlled umbilical pulse is blood pressure and heart rate dependent.

Newborn's umbilical cord can be catheterized or injected with up to 2ccs. The umbilical cord can be trimmed or clipped for simulations and is therefore considered a consumable item.

To fill the umbilical cord with fluid, inject any of the three blood vessels with 2 cc of water using the fill syringe. For catheterization, use a 6 Fr urethral round tip catheter lubricated with silicon oil.



To remove the umbilical cord, first turn the simulator OFF (by clicking on File, Exit on the GaumardUI software), or set the manikin on STAND-BY mode. The replacement of the umbilical cords should also be done while the manikin is turned off or set on STAND-BY mode.

**WARNING:** Do not remove umbilical cord while Newborn is in operation.

## Cephalic

Newborn displays central and peripheral cyanosis at various intensities – healthy, mild, and severe.

## Systemic

Intramuscular Injection sites are located on both deltoids and quadriceps for placement and technique exercises.





# Other

## Vital Signs Monitor

The Virtual Signs Monitor simulates a vital signs monitor attached to Newborn. The vital signs are synchronized through a wireless network between the facilitator's tablet and the computer running the monitor. Each trace can be customized independently of each other; users can set alarms, time scales, boundaries and grid options. In addition, it allows the facilitator to display lab reports, x-rays, and other files on the virtual monitor screen for use by the provider.




## Appendix

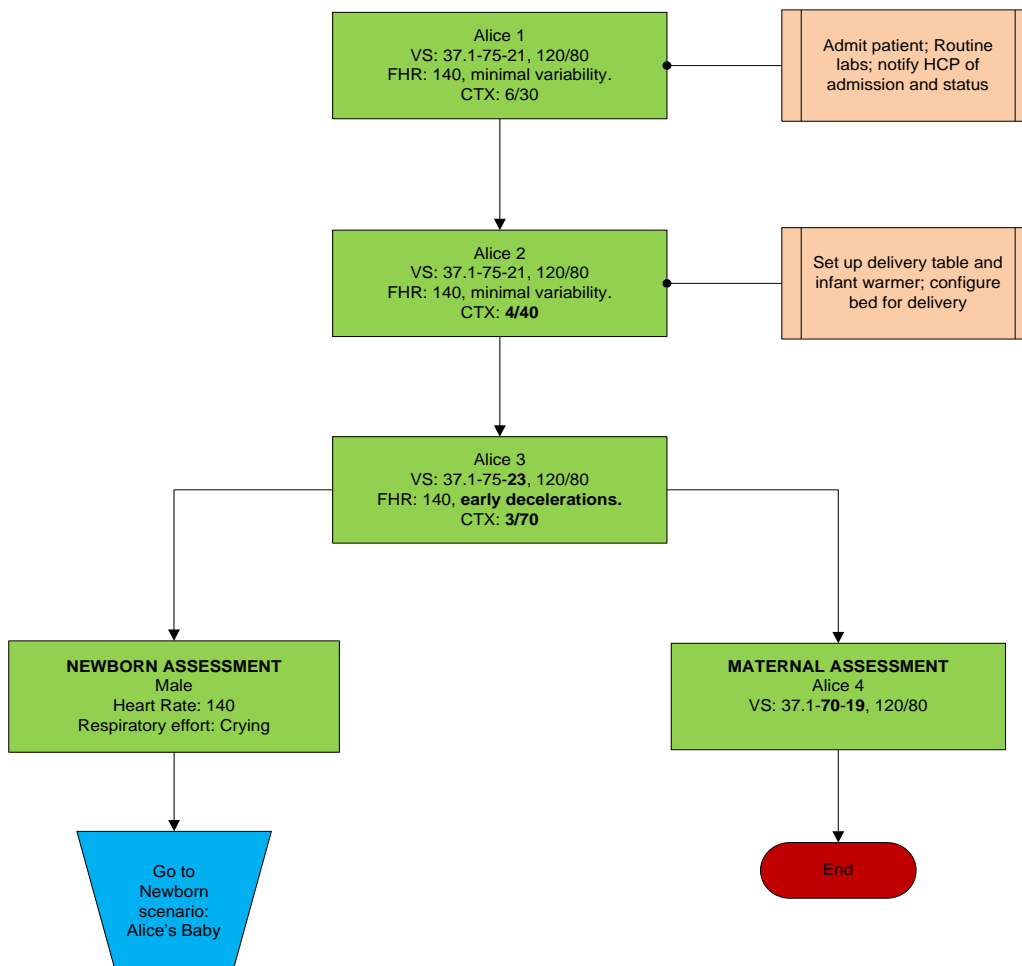
---


# Factory Preset Labor Scenarios


## Quick start Scenarios


	Scenario Name	Labor Type
1	Alice	Normal
2	Alicia	Variations on Normal
3	Amy	Variations on Normal
4	Angelica	Variations on Normal
5	Beth	Variations on Normal
6	Cynthia	Shoulder Dystocia
7	Donna	Breech
8	Elaine	Preeclampsia
9	Francine	Cesarean Delivery
10	Gloria	Cord Prolapse
11	Helen	Hemorrhage
12	Irene	Cesarean Delivery


 <div data-bbox="300 231 508 283"><b>Gaumard®</b> Simulators for Health Care Education</div>	<p>Noelle® - Labor Scenario</p> <p><b>Alice</b></p> <p>Normal</p>
<p>Alice is a 24 year old female, weighing 170 pounds. Her OB history shows a gravida of 2 and a term of 1. She is currently 39 weeks pregnant and has one living child. She has had prenatal care. She has not been using medications of any kind. Labor duration: 30 minutes.</p>	



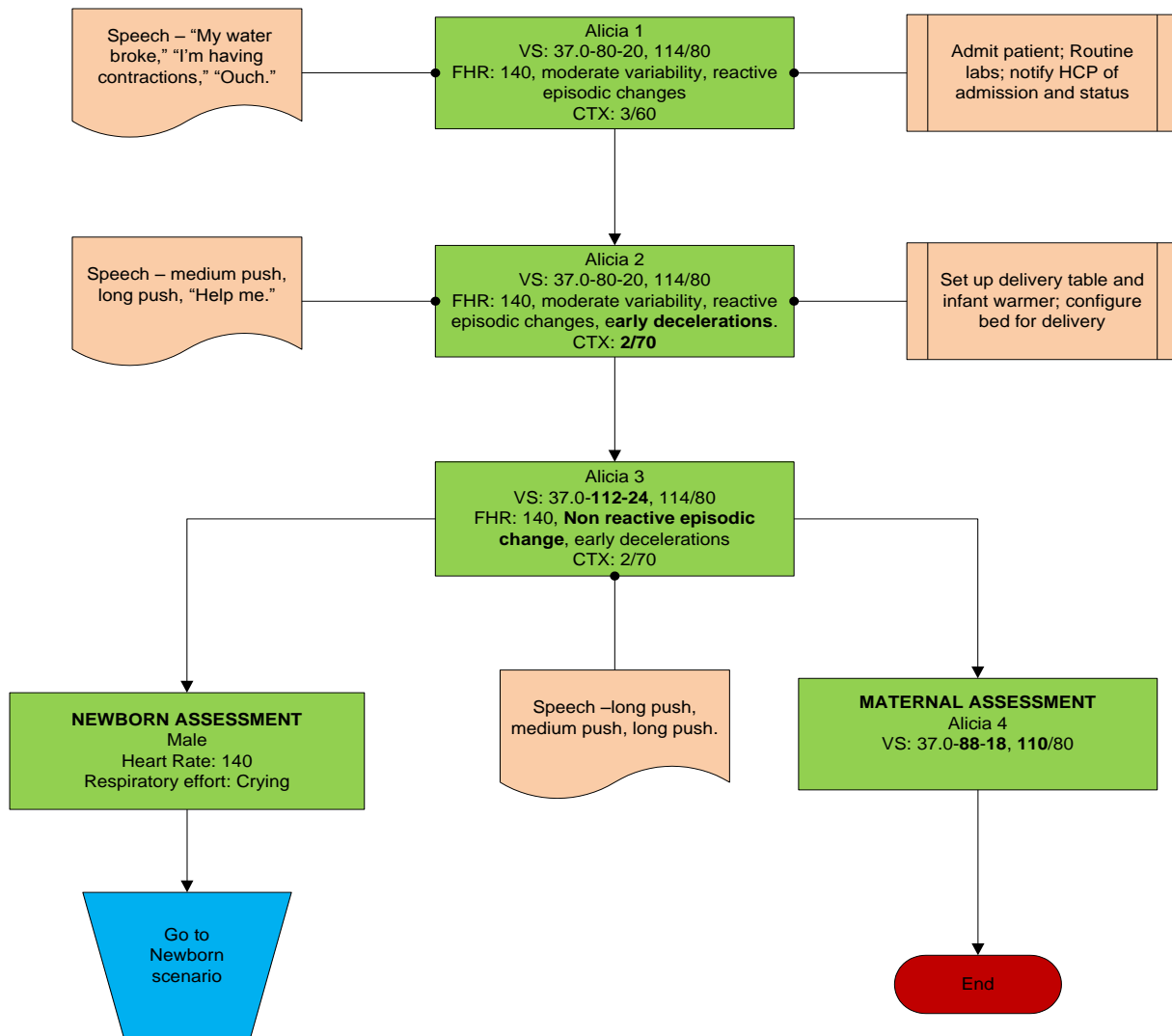	<p>Noelle® - Labor Scenario  <b>Alicia</b>          Variations on Normal</p>
<p>Alicia is a 24 year old gravida 2/1 at 39 weeks. She weighs 160 pounds. She has had prenatal care. She has not been using medications of any kind. Labor duration: 20 minutes.</p>	



 <b>Gaumard®</b> Simulators for Health Care Education	Noelle® - Labor Scenario <b>Amy</b> Variations on Normal
Amy is 19 years and she weighs 160 lbs. Her OB history shows a gravida of 1. She is currently 40 weeks pregnant. She enters LD accompanied by her mother. Labor duration: 30 minutes.	

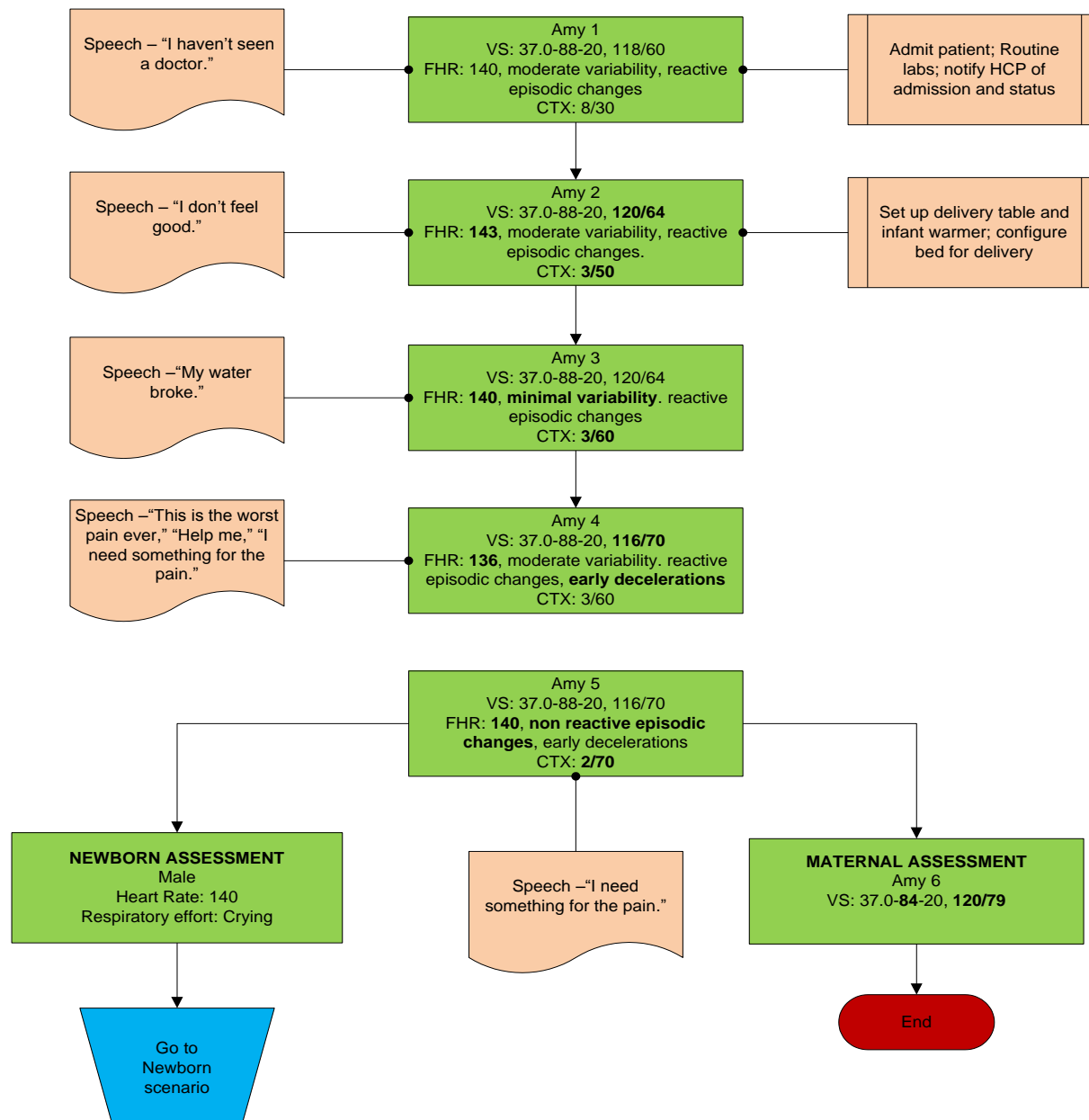




**Gaumard®**  
Simulators for Health Care Education


Noelle® - Labor Scenario


## Angelica

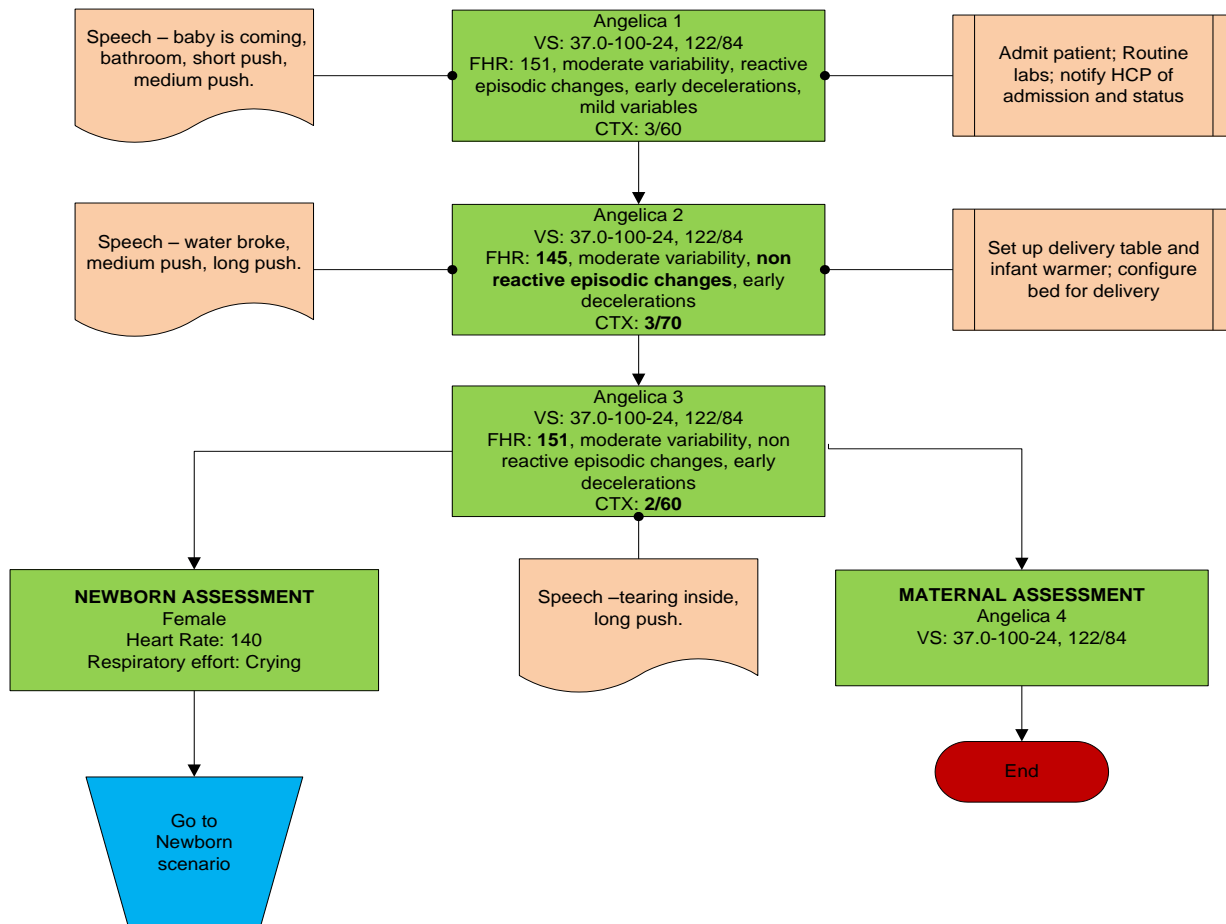Variations on Normal


Angelica is a 31 year old gravida 5/3 at 41 weeks. She weighs 160 lbs. She has experienced no prenatal complications and has a history of fast labors. Labor duration: 20 minutes.




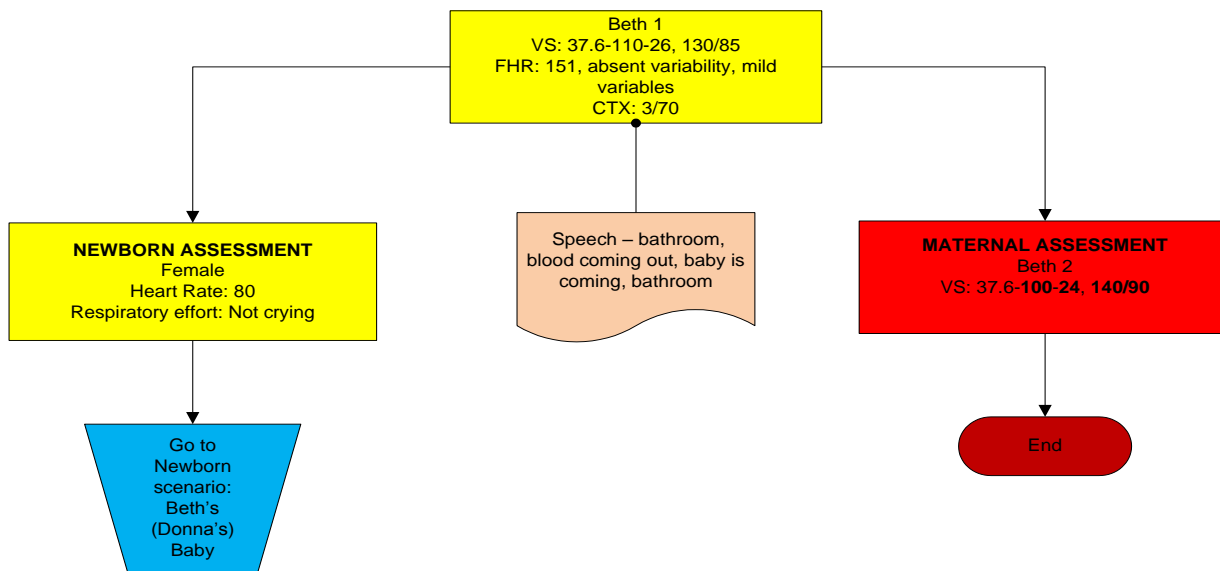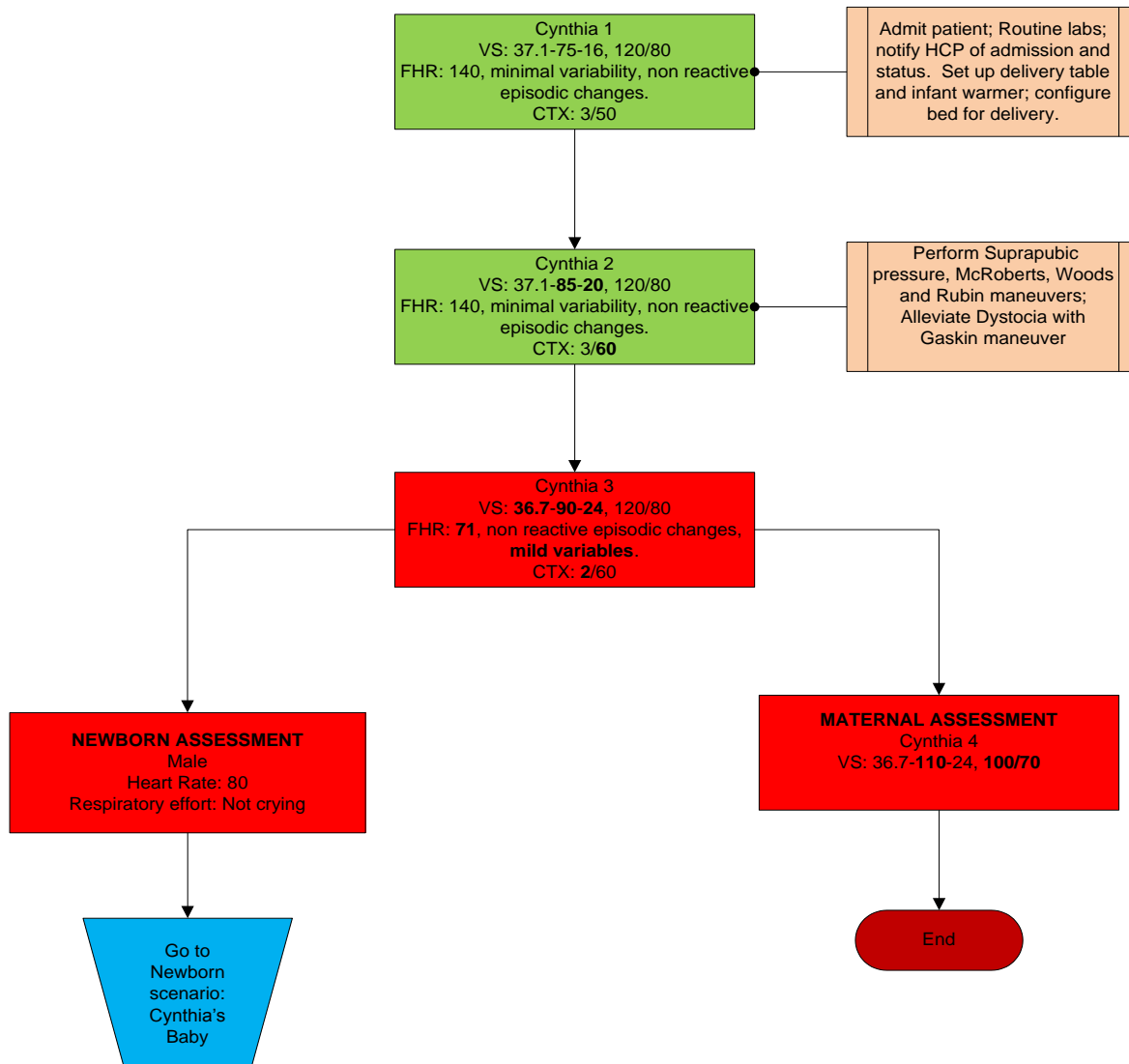
	<p>Noelle® - Labor Scenario <b>Beth</b> Variations on Normal</p>
<p>Beth is a 16 year old gravida 2/0 at 37 weeks. She has had one elective abortion. She has had prenatal care. Labor duration: 10 minutes.</p>	

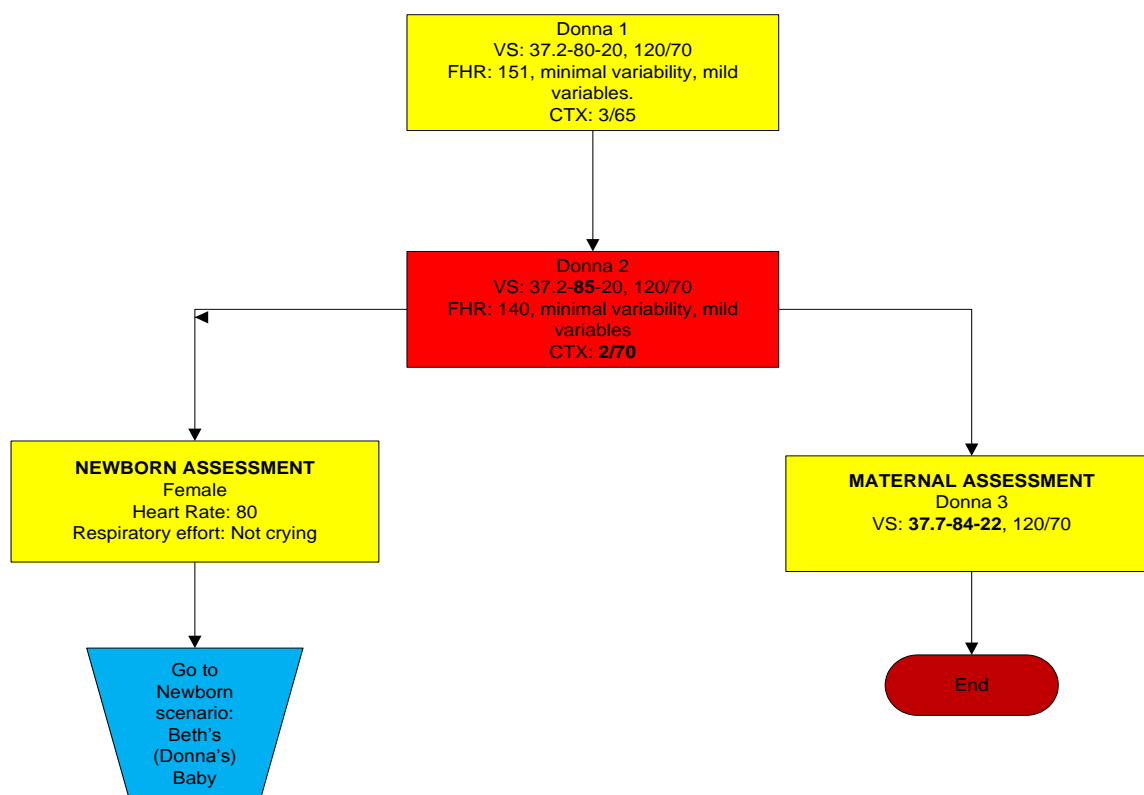





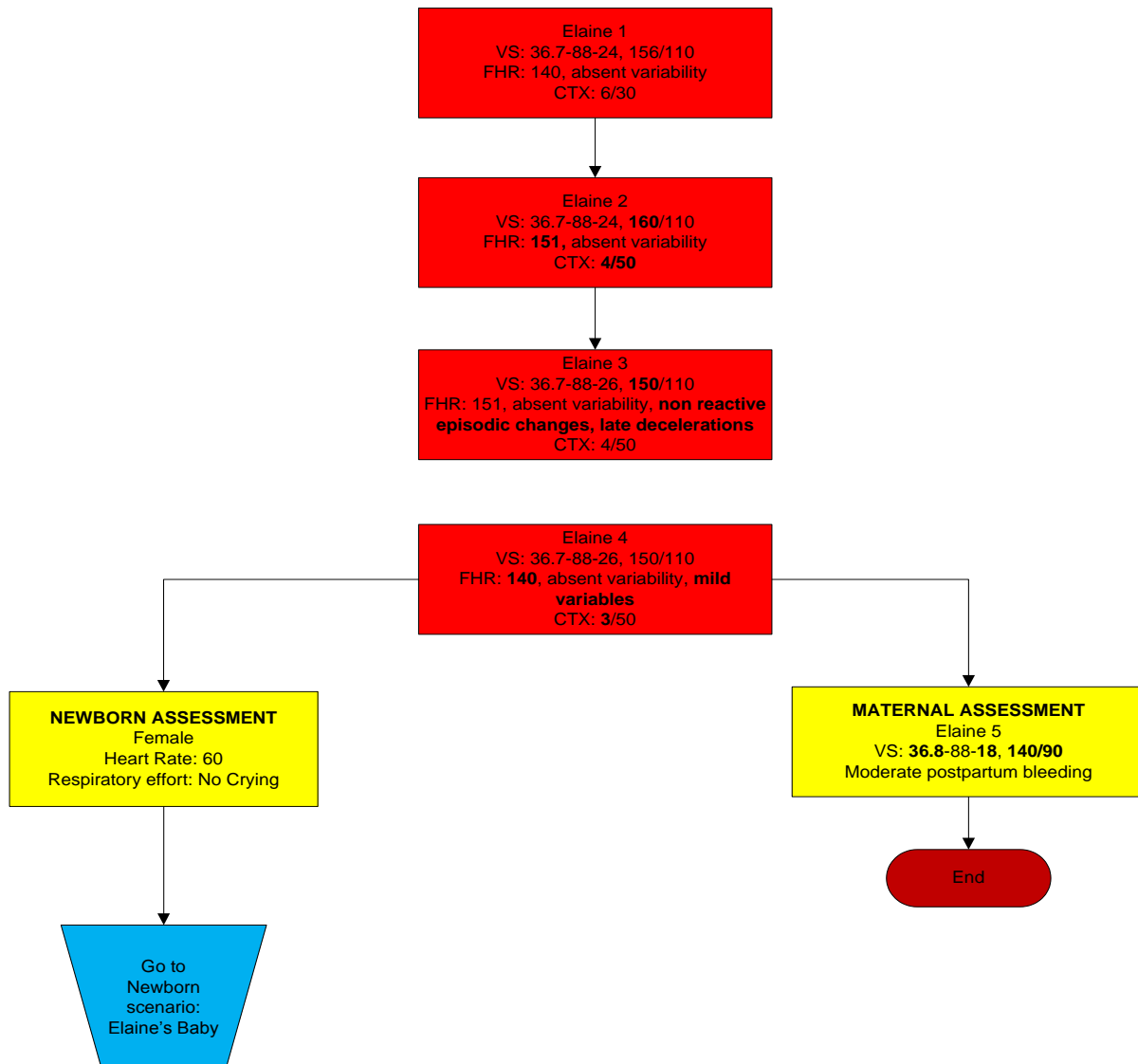
	<p>Noelle® - Labor Scenario</p> <p><b>Cynthia</b></p> <p>Shoulder Dystocia</p>
<p>Cynthia is a 31 year old gravida 3/1 at 41 weeks. She weighs 170 lbs. Labor duration: 30 minutes.</p>	




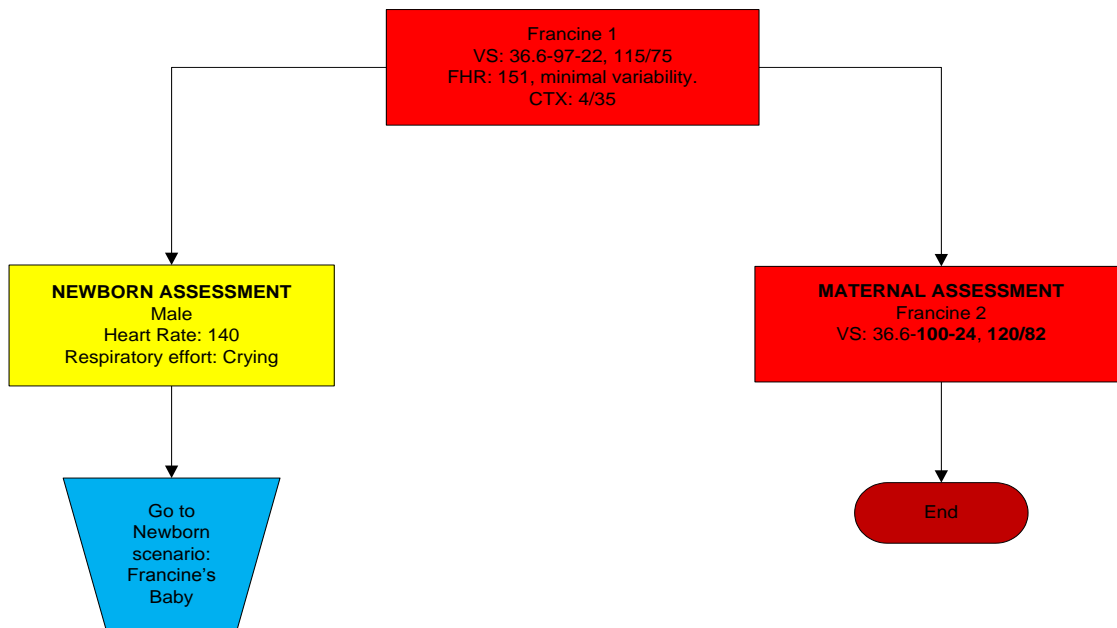
 <b>Gaumard®</b> Simulators for Health Care Education	Noelle® - Labor Scenario <b>Donna</b> Breech
Donna is a 20 year old gravida 4/2 at 31 weeks. She weighs 180 lbs. She has had one elective abortion. She has had prenatal care. Labor duration: 20 minutes.	




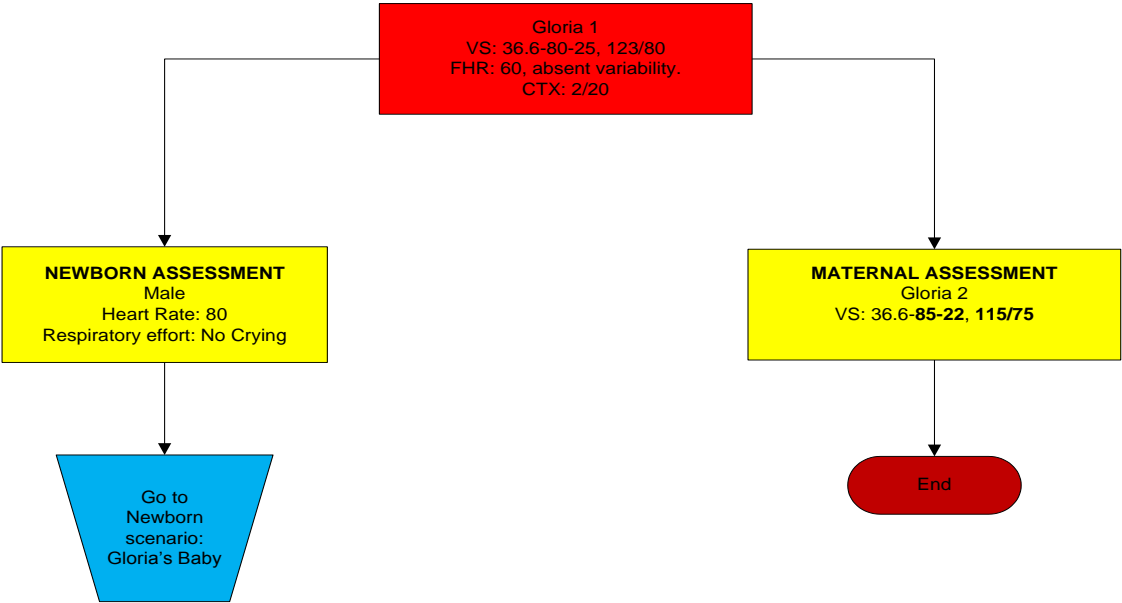
 <b>Gaumard®</b> Simulators for Health Care Education	Noelle® - Labor Scenario <b>Elaine</b> Preeclampsia
Elaine is a 23 year old gravida 1/0 at 37 weeks. She weighs 140 lbs. She has had prenatal care. She complains of mind frontal headache. 3+tibial edema and 4+ DTRs with 2 beats clonus are noted. Labor duration: 40 minutes.	




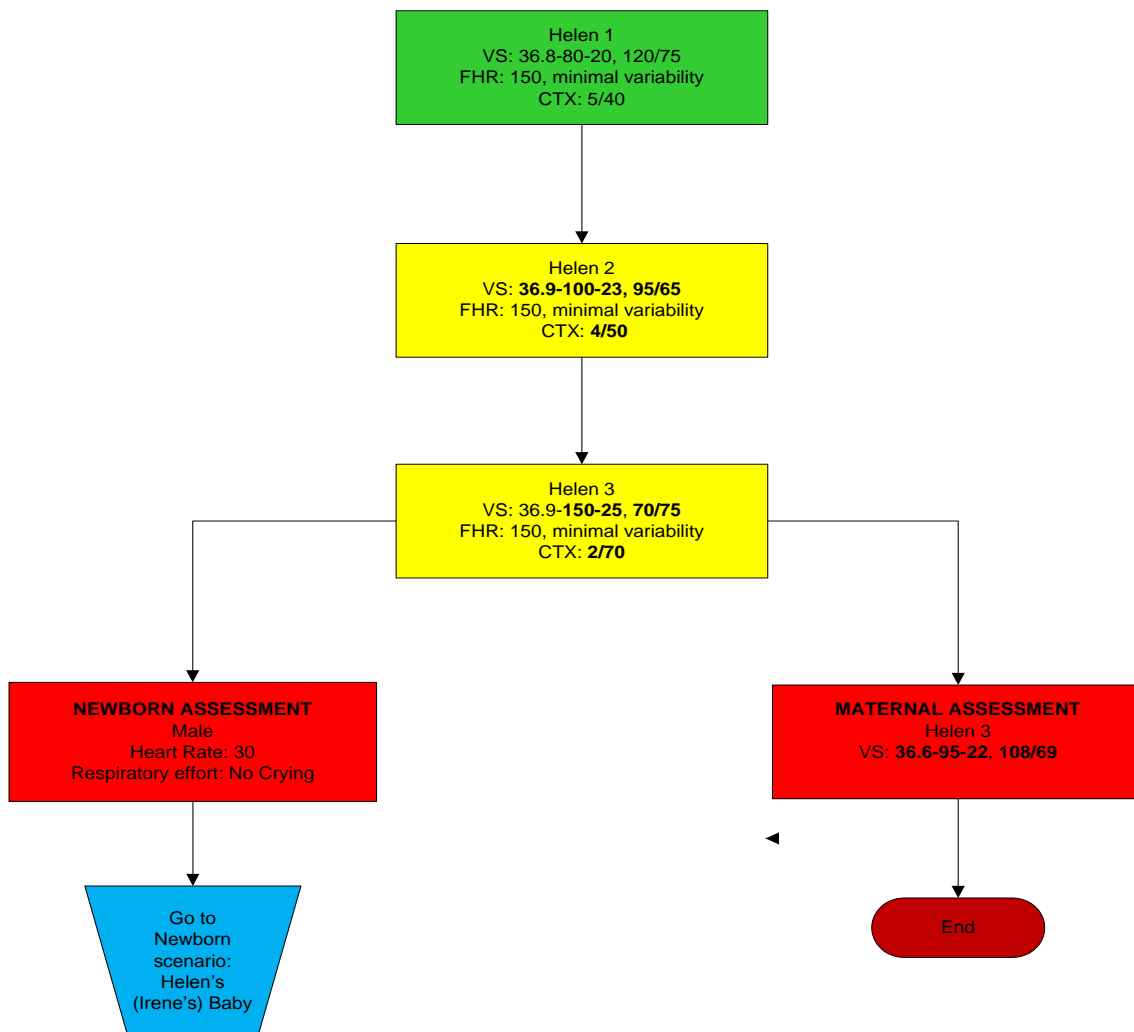
	<p>Noelle® - Labor Scenario  <b>Francine</b>  Cesarean Delivery</p>
<p>Francine is a 19 year old female gravida 2/1 at 37 weeks. She weighs 145 lbs. She has had prenatal care. She has STD, Herpes. Labor duration: 10 minutes.</p>	




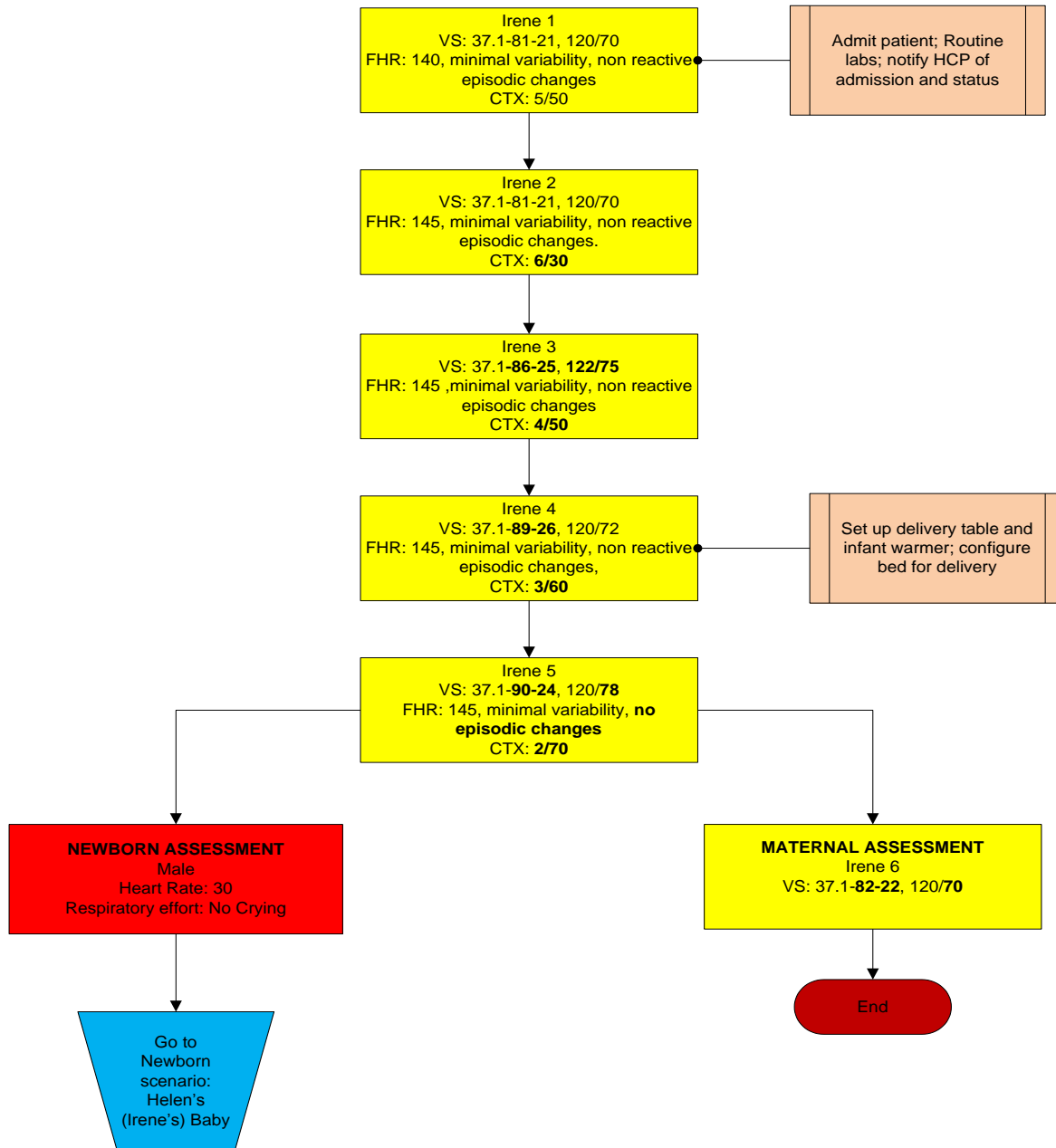
	<p>Noelle® - Labor Scenario <b>Gloria</b> Cord Prolapse</p>
<p>Gloria is a 34 years old gravida 1/0 at 25 weeks. She weighs 190 lbs. She has had prenatal care. Labor duration: 10 minutes.</p>	



 <div><b>Gaumard®</b> Simulators for Health Care Education</div>	Noelle® - Labor Scenario <b>Helen</b> Hemorrhage
Helen is a 25 year old gravida 1/0 at 35 weeks. She weighs 180 lbs. She has had prenatal care. Labor duration: 30 minutes.	




 <b>Gaumard®</b> Simulators for Health Care Education	Noelle® - Labor Scenario <b>Irene</b> Cesarean Delivery
Irene is a 19 year old gravida 2/0 at 29 weeks. She has had one spontaneous abortion. Labor duration: 45 minutes.	





## NOELLE Advanced

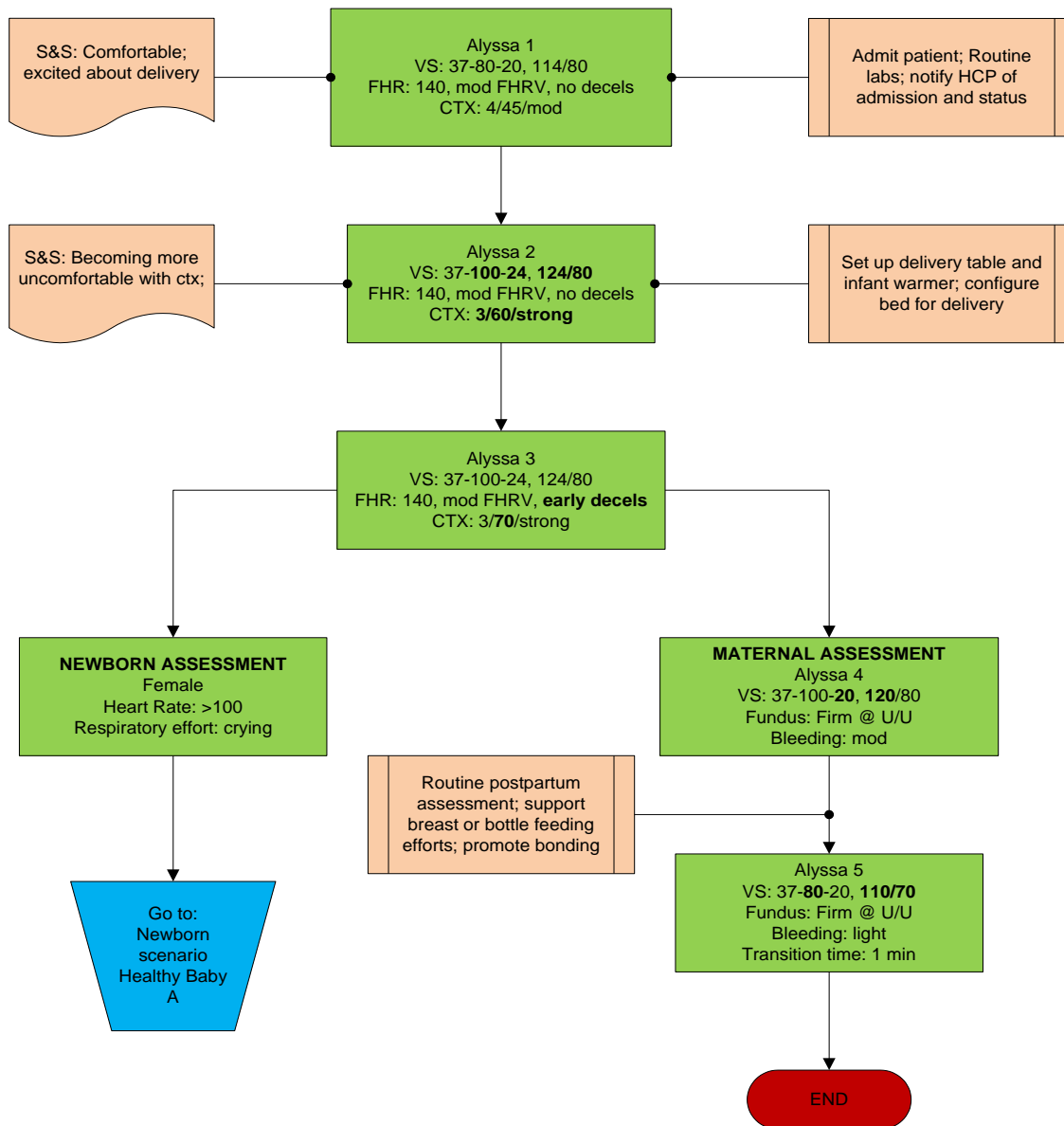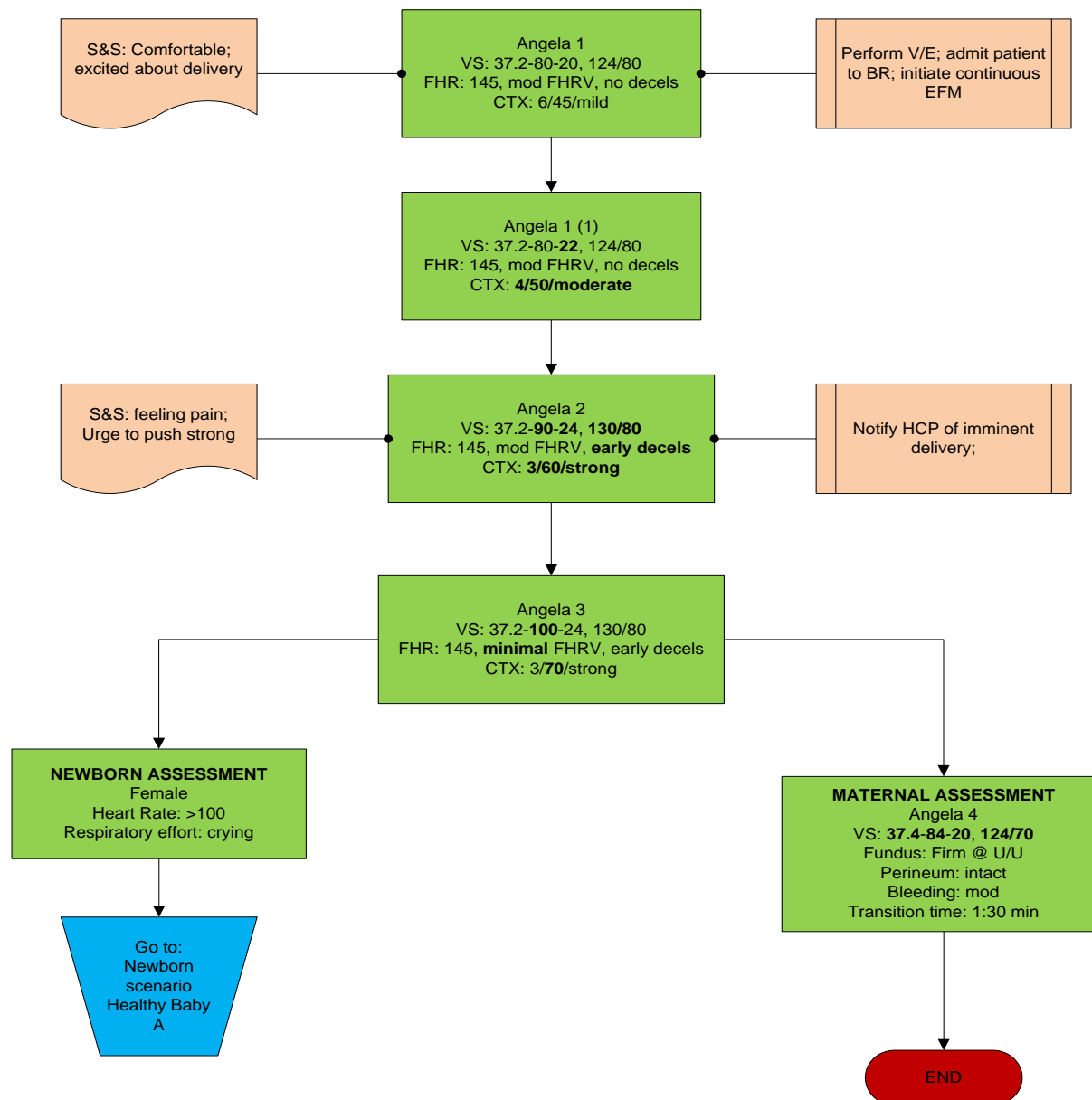
	Scenario Name	Labor Type
1	Alyssa	Normal Labor
2	Angela	Normal Labor
3	Becca	Variations on Normal
4	Bianca	Variations on Normal
5	Candice	Shoulder Dystocia
6	Charlotte	Shoulder Dystocia
7	Dana	Breech Presentation
8	Demaris	Breech Presentation
9	Eleanor	Preeclampsia
10	Erin	Preeclampsia
11	Faye	Cord Prolapse
12	Frances	Cord Prolapse
13	Gabriella	Uterine Rupture
14	Gail	Uterine Rupture
15	Heidi	Peripartum Hemorrhage - Previa
16	Haley	Peripartum Hemorrhage - Previa
17	India	Peripartum Hemorrhage - Abruptio
18	Inez	Peripartum Hemorrhage - Abruptio
19	Janie	Peripartum Hemorrhage/PPH
20	June	Peripartum Hemorrhage/PPH
21	Kelly	Amniotic Fluid Embolism
22	Kimberly	Amniotic Fluid Embolism
23	Madonna	Preterm Labor
24	Maria	Preterm Labor




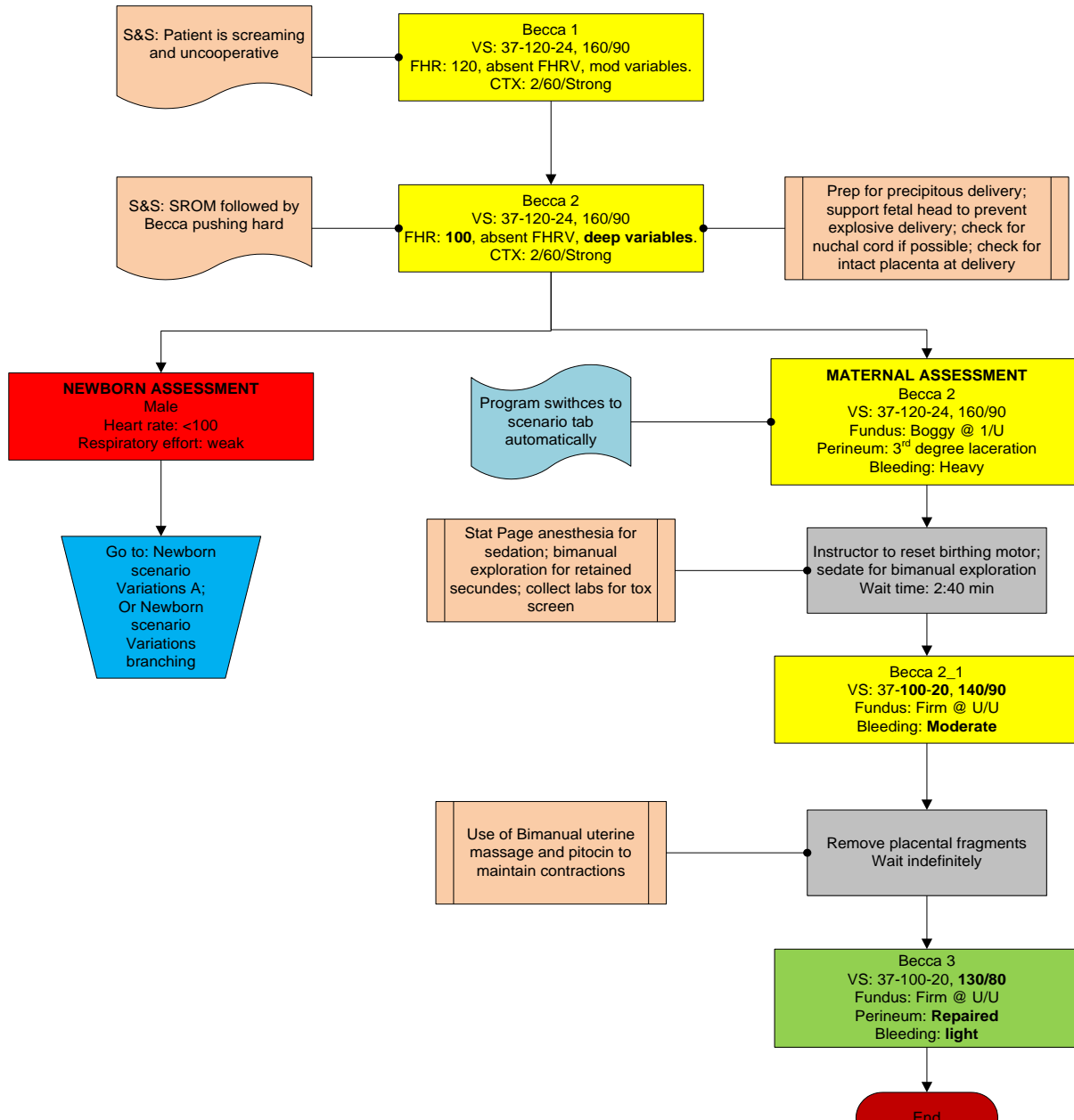
	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Alyssa</b></p> <p>Normal Labor</p>
<p>Alyssa is a 23 yr old primip at term. Her health is generally good and she has experienced no prenatal complications. She wishes to receive no medications and will have the CNM attending her delivery. Labor duration: 30 minutes.</p>	



 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® - Labor Scenario <h2 style="text-align: center;">Angela</h2> Normal Labor
Angela is a 31 yr old grand multip @ term. She tells the triage nurse that even though she has had few contractions she came in because she has a history of rapid labors. Her general health is good and she has had no problems during this pregnancy. V/E shows the cervix to be paper thin and Angela is admitted to birthing room. Labor duration: 30 minutes.	



	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Becca</b></p> <p>Variations on Normal</p>
<p>Becca is a young pregnant teen who is living on the streets. She is a heavy smoker and drug user. She was seen twice in the Adolescent Clinic and referred to Social Services, but she only saw the social worker once and did not go to the follow-up appointment. Labor duration: 18-22 minutes.</p>	






**Gaumard®**  
Simulators for Health Care Education

Noelle S574-575® - Labor Scenario

## Bianca

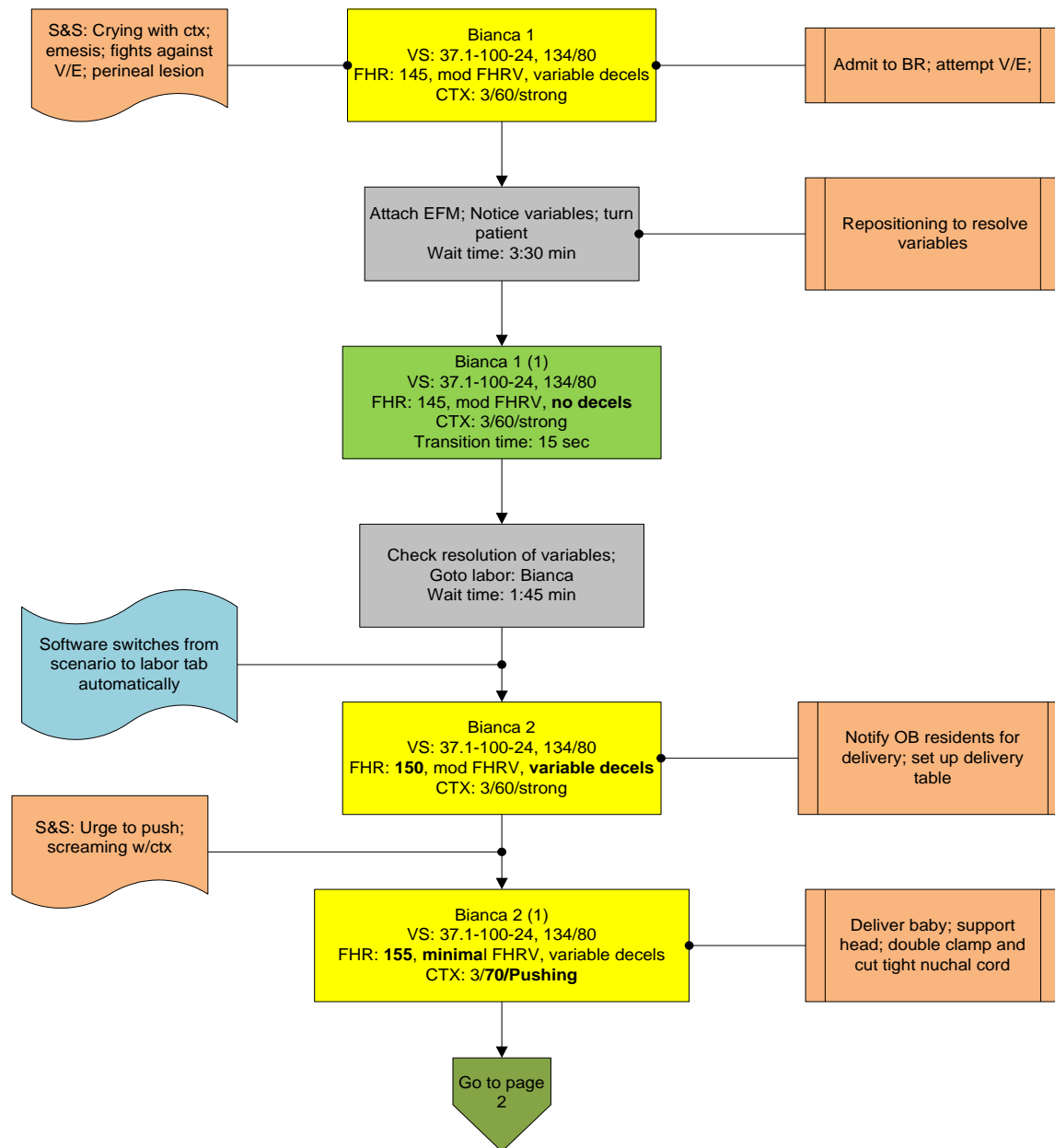Variations on Normal

Bianca is a 16 yr old @ 38 weeks who shows up in L&D in active labor. She had a previous elective AB as a result of incest at age 13. She now lives with her boyfriend and his mother who are both with her at the hospital. Her prenatal visits have been irregular due to transportation issues. She is leaking light meconium fluid and she vomits as she is undressing. Labor duration: 25 minutes. Scenario duration: 30 minutes.






**Gaumard®**  
Simulators for Health Care Education

Noelle S574-575® - Labor Scenario

**Bianca**

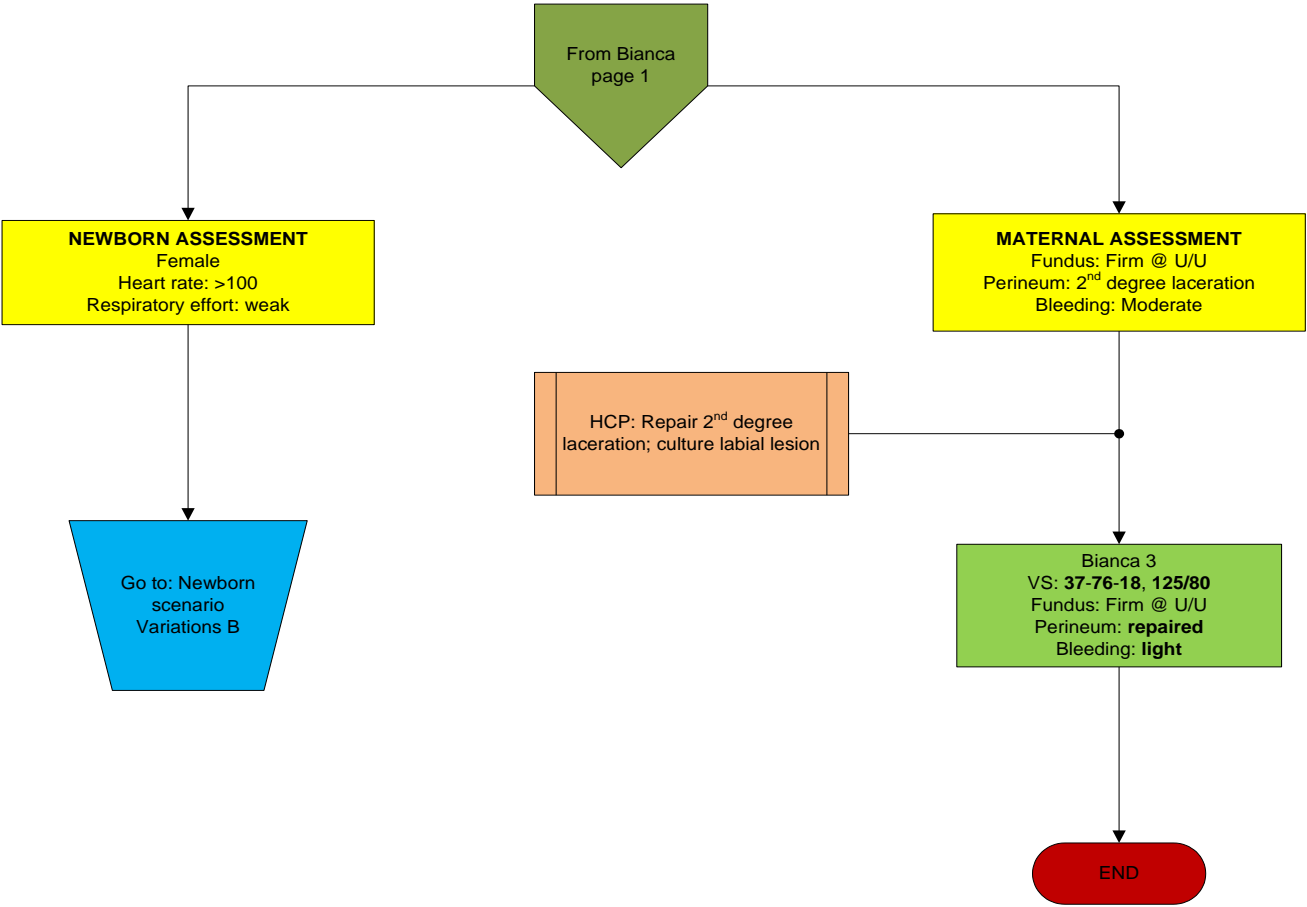Variations on Normal






Gaumard®  
Simulators for Health Care Education


Noelle S574-575® - Labor Scenario


## Candice

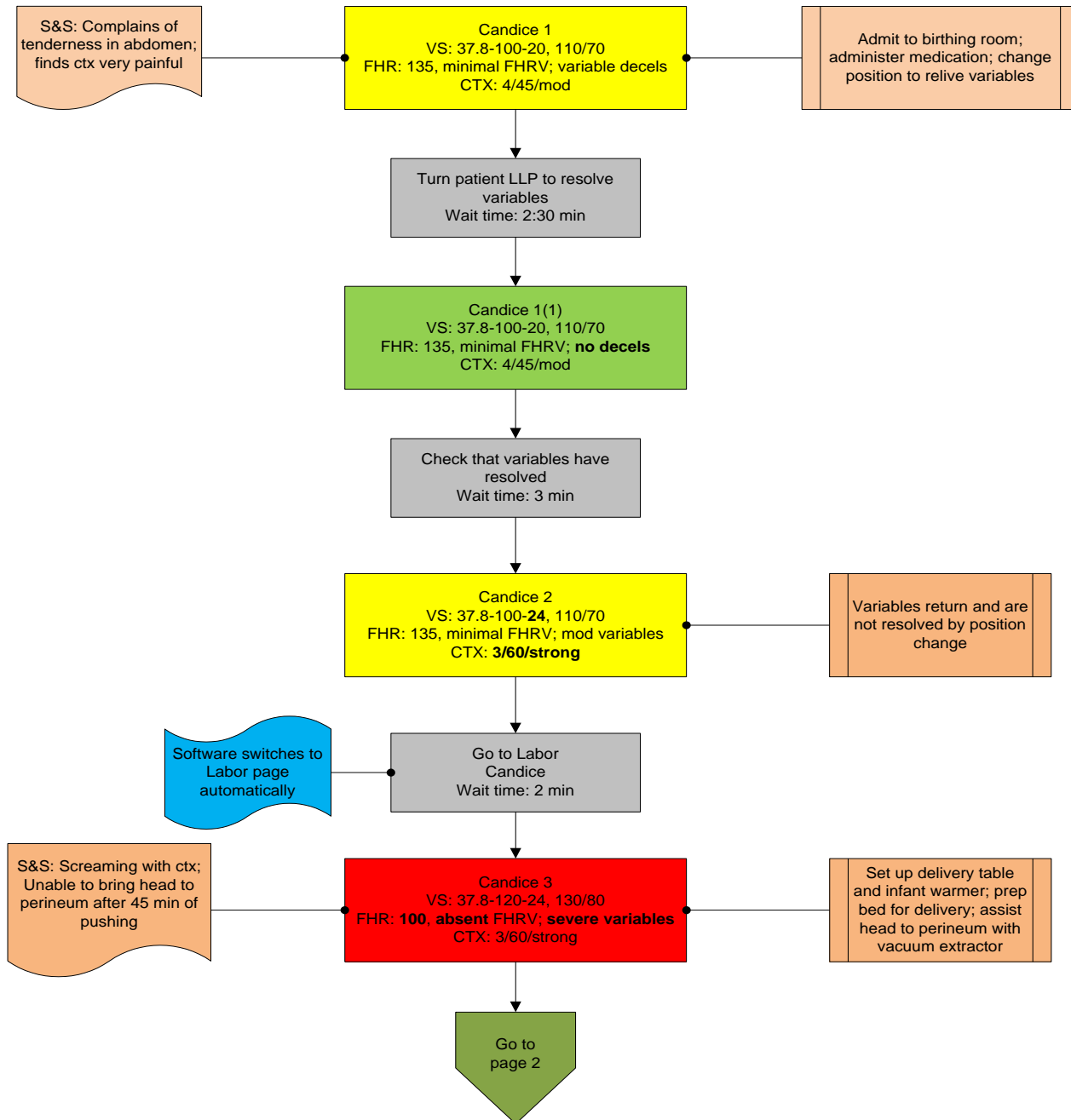
### Shoulder Dystocia


Candice is a 19 r old multip. She, her boyfriend and their 3 yr daughter are homeless living in a car. She has not seen a doctor and believes that she is about 8 months pregnant. Her water broke yesterday and she is leaking moderately thick meconium fluid. An ultrasound is performed to determine position and gestational age. Labor duration: 15 minutes. Full scenario duration: 23-25 minutes.




 <p>Gaumard® Simulators for Health Care Education</p>	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Candice</b></p> <p>Shoulder Dystocia</p>
--	---

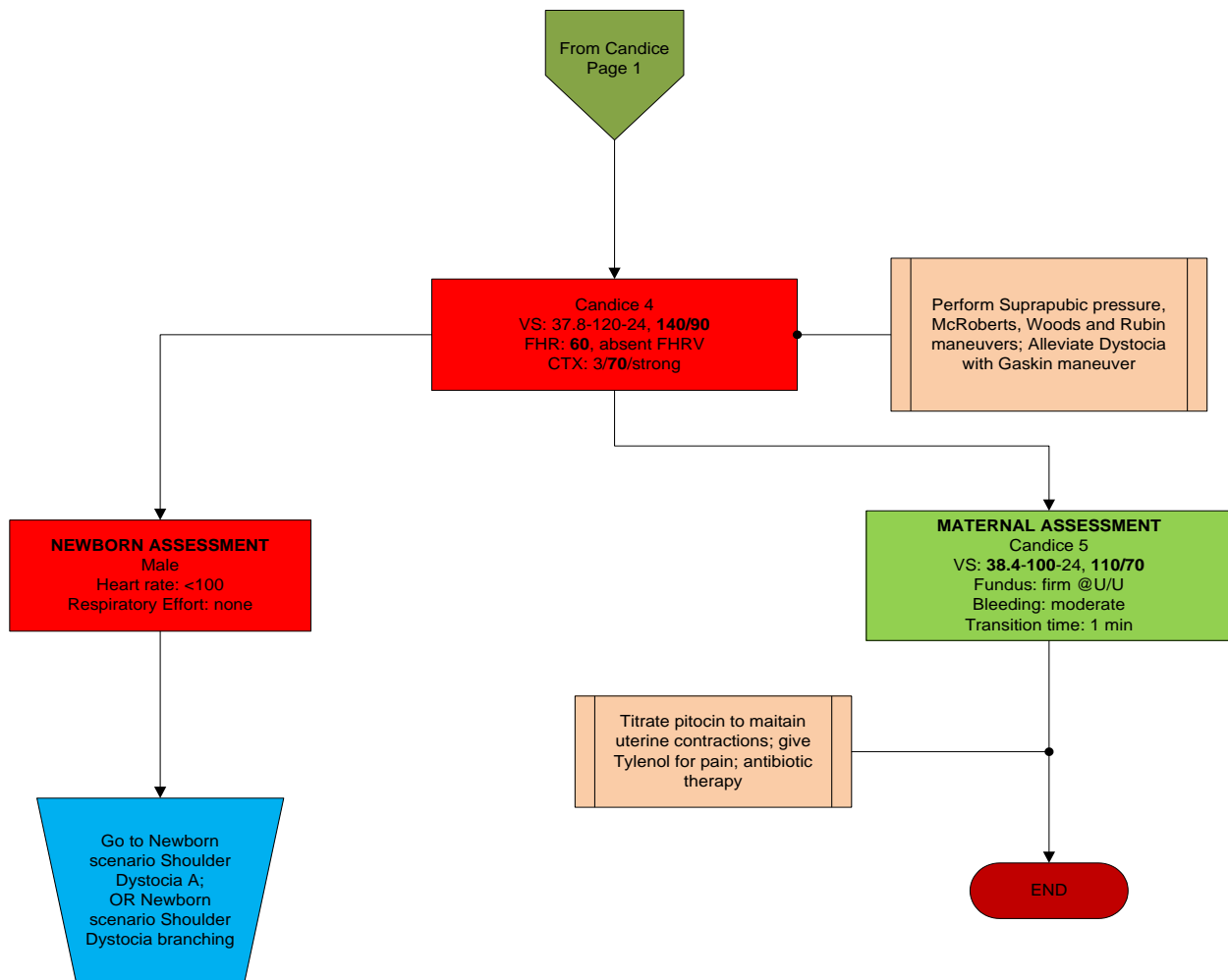



 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® - Labor Scenario <b>Charlotte</b> Shoulder Dystocia
Charlotte is a 31 yr old gravida 3/1 @ 41+ 5/7 weeks. Her physician stripped her membranes yesterday and she began contracting during the night. She is admitted in active labor. Labor duration: 40 minutes.	

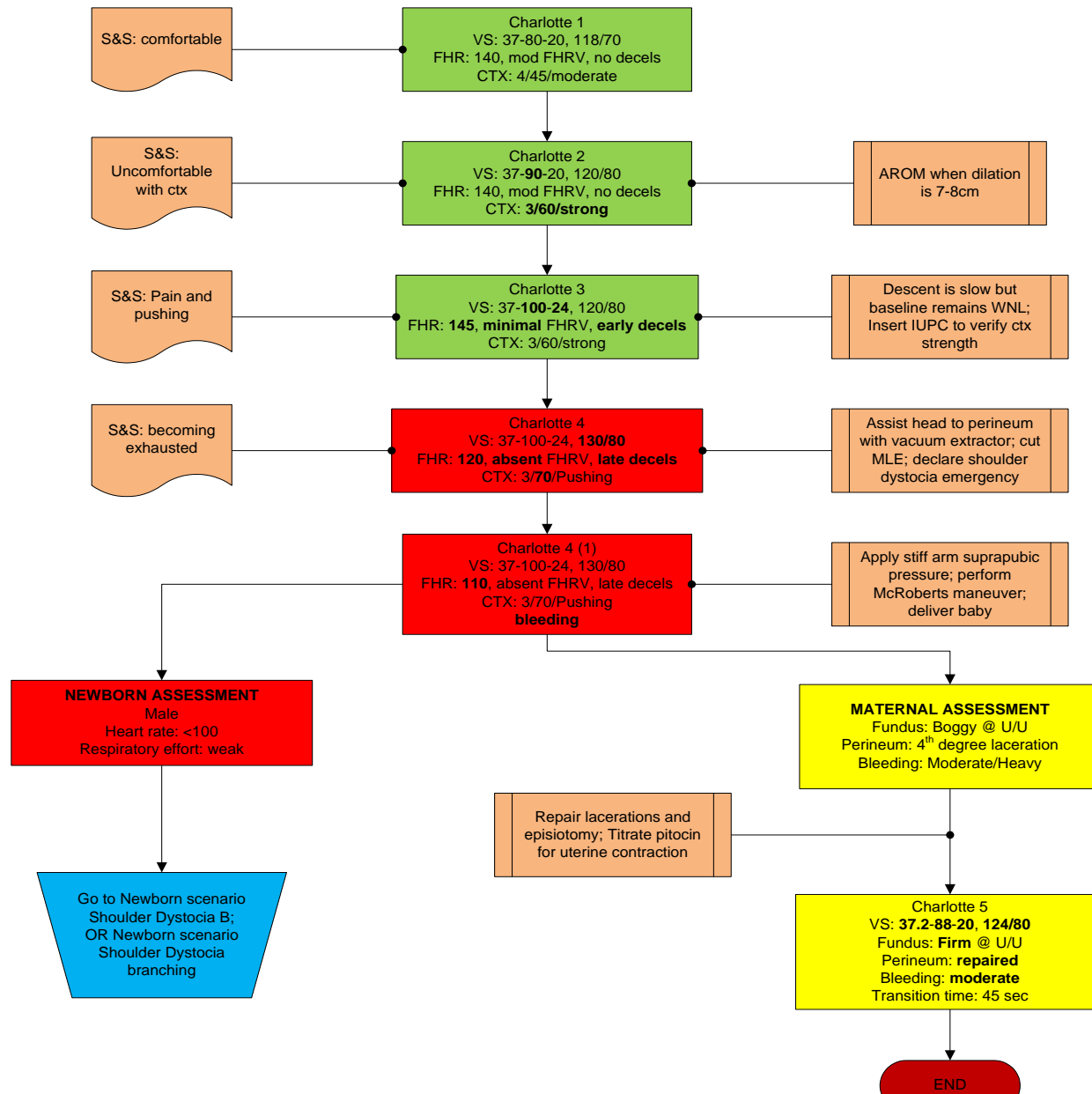




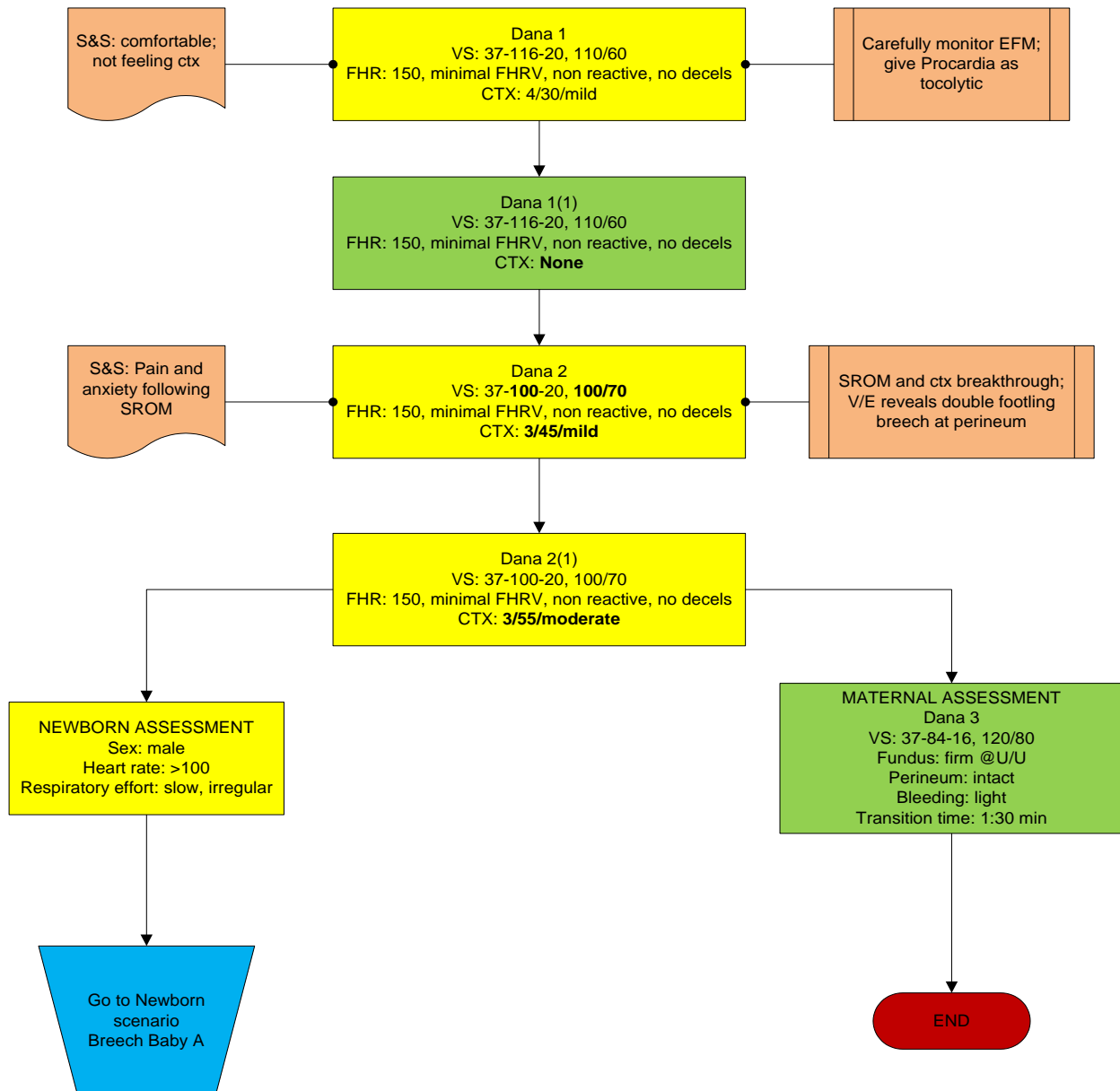	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Dana</b></p> <p>Breech Presentation</p>
<p>Dana is a 24 yr old multip @ 29 weeks who was admitted because she began contracting. Upon V/E physician discovers that she is 4-5cm with bulging membranes. She was given Terbutaline subQ and then transported to the regional medical center. Labor duration: 45 minutes.</p>	



	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Demaris</b></p> <p>Breech Presentation</p>
<p>Demaris is a young Hispanic teen who has received prenatal care in the Adolescent OB clinic. She kept the pregnancy a secret as long as was possible and did not attend any childbirth classes. Her plan is to return to high school while her mother cares for the baby. The baby's father will not accept any responsibility and does not wish to be involved. Labor duration: 30 minutes.</p>	

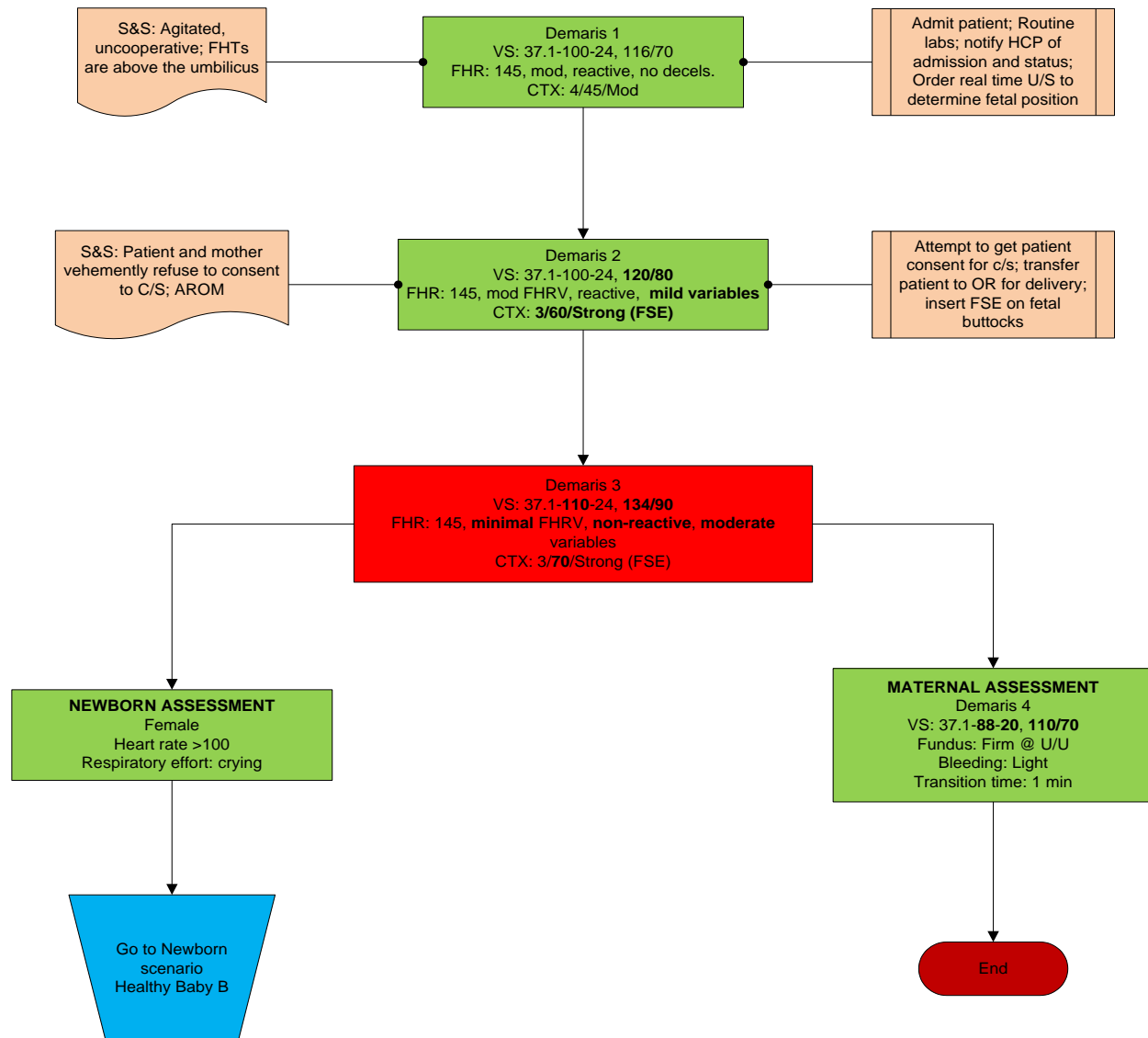




**Gaumard®**  
Simulators for Health Care Education


Noelle S574-575® - Labor Scenario


## Eleanor

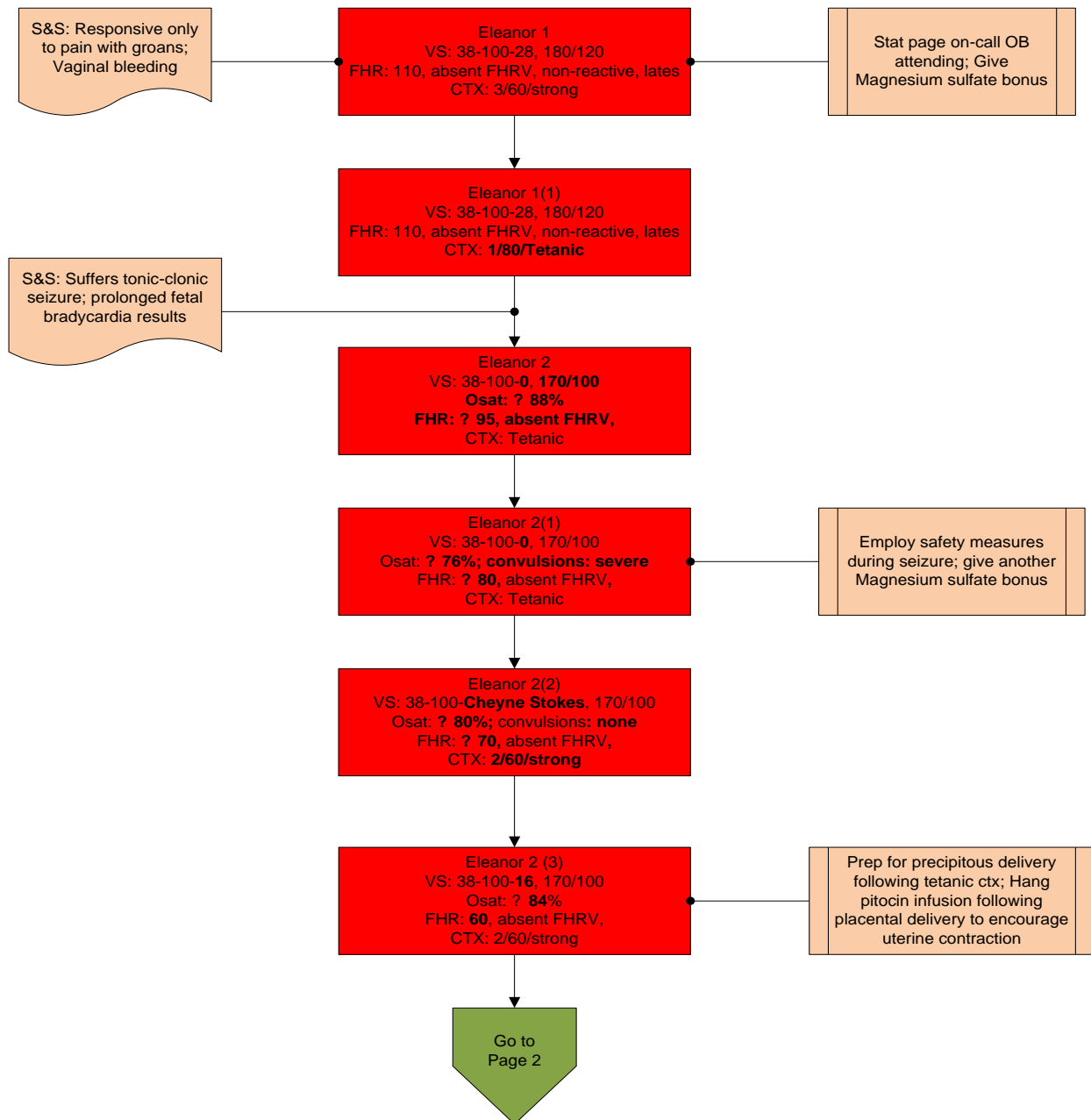
Variations on Normal


Eleanor is a 19 yr old post-ictal patient being transferred to the ER by EMS. She was found convulsing in the bathroom. According to relatives she is 8½ months pregnant with her first baby. She has been on an IV during transport and her BP is 180/120. Labor duration: 20 minutes.



 <p><b>Gaumard®</b> Simulators for Health Care Education</p>	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Eleanor</b></p> <p>Variations on Normal</p>
---	--

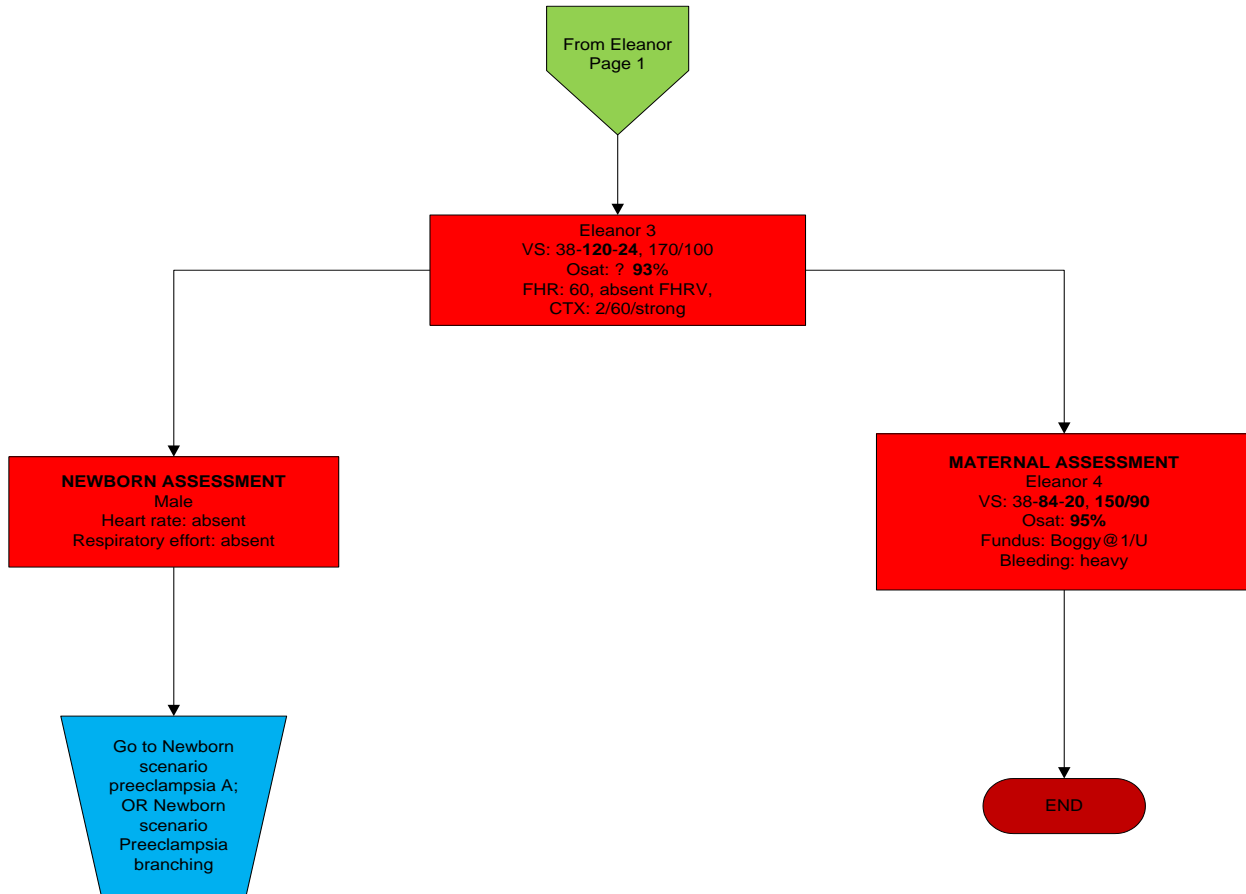

 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574- 575® - Labor Scenario <h2 style="text-align: center;">Erin</h2> <h3 style="text-align: center;">Preeclampsia</h3>
Erin is a 28 yr old @ 38 weeks admitted by her physician for preeclampsia. She is started on Magnesium sulfate per protocol, induced with pitocin and her membranes are ruptured. Light meconium staining is noted following AROM. Labor duration: 60 minutes.	

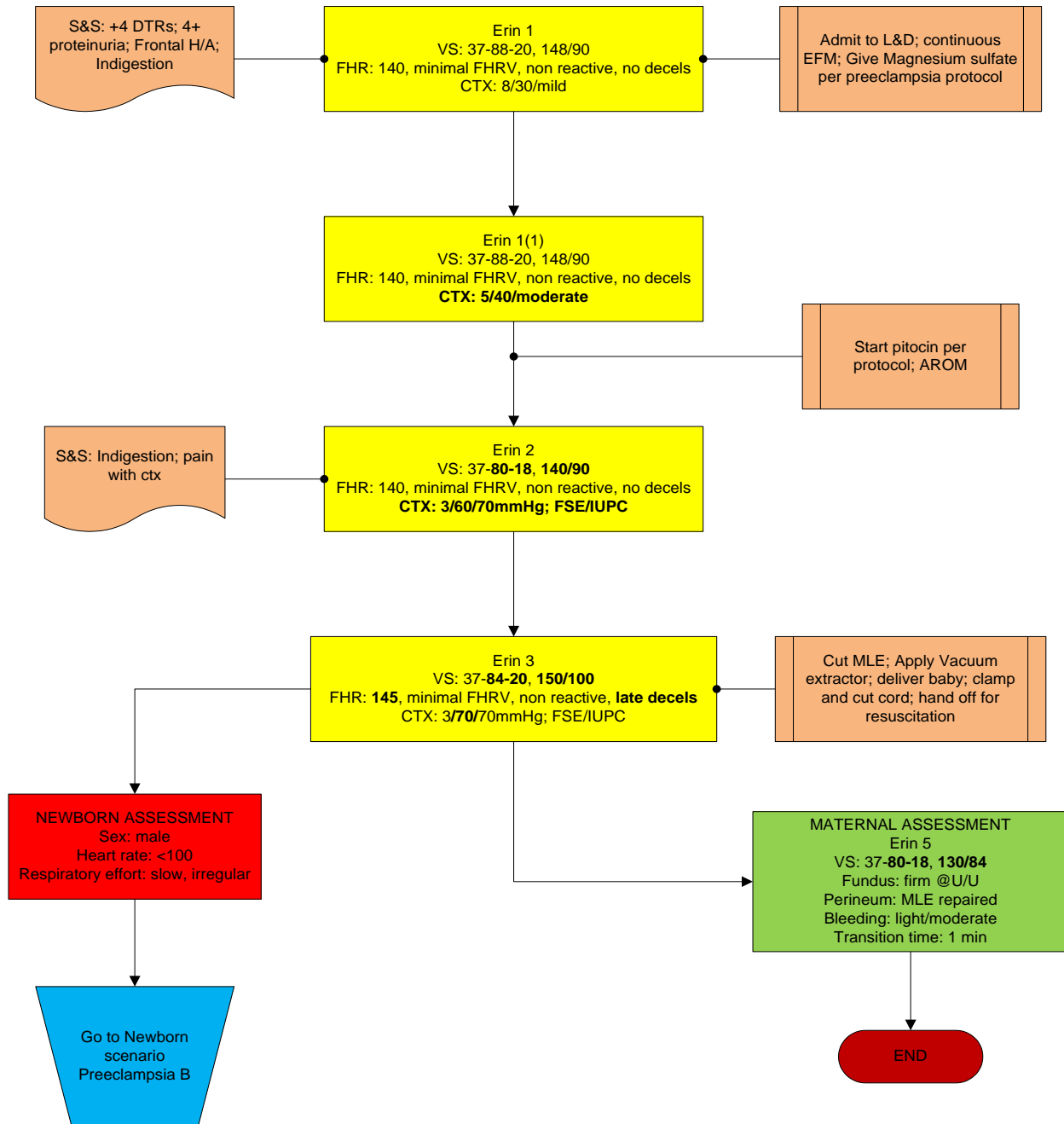





**Gaumard®**  
Simulators for Health Care Education


Noelle S574-575® - Labor Scenario


## **Faye** Cord Prolapse

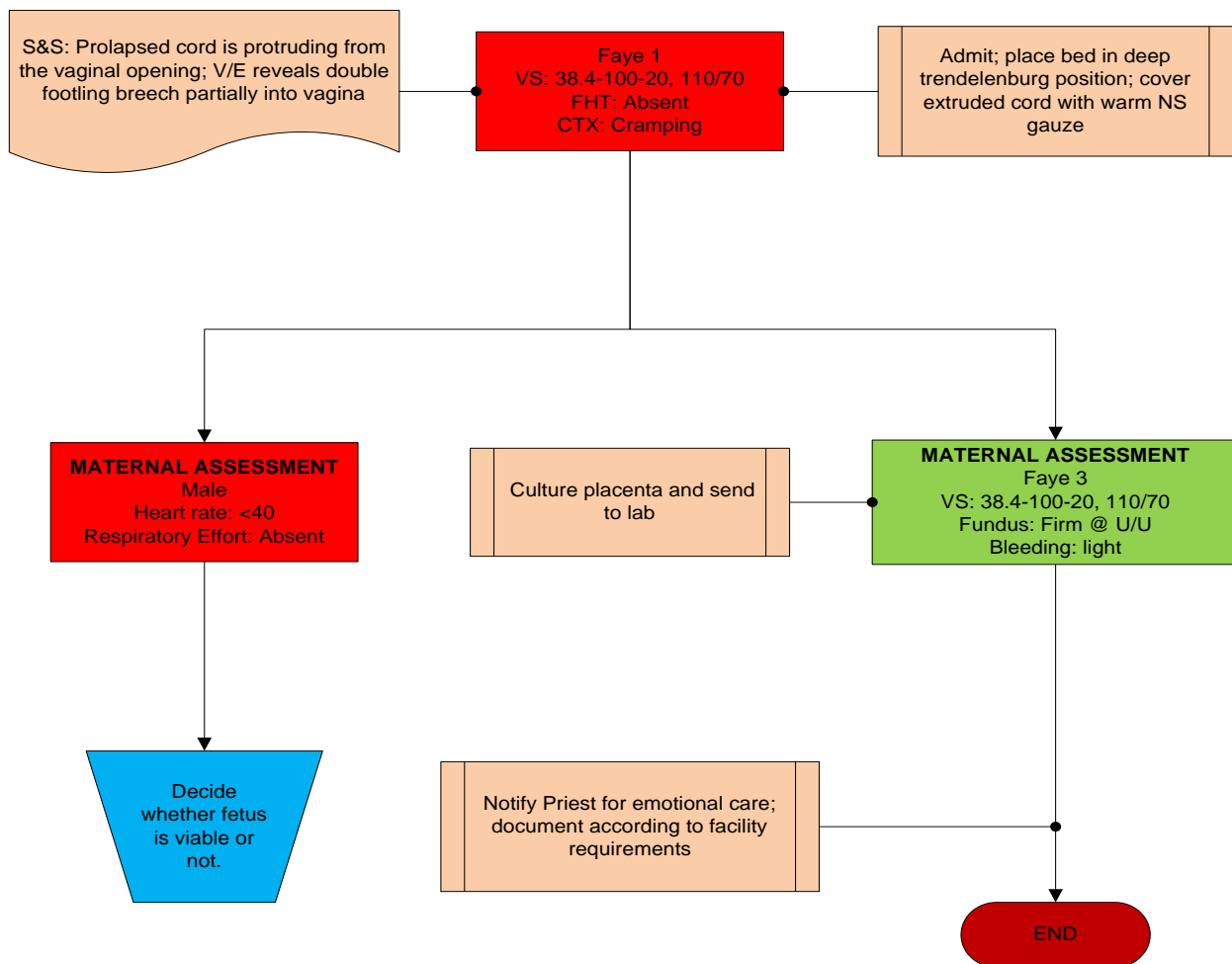
Faye is a 34 yr old gravida 1 @ 25 weeks' gestation. She began cramping about 3 hours ago and decided to drive herself to the hospital. She began leaking clear fluid on the way. An admitting clerk helps her into a wheel chair and takes her to L&D. Labor duration: 20 minutes.




 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® - Labor Scenario <h2 style="text-align: center;">Frances</h2> <h3 style="text-align: center;">Prolapsed Cord</h3>
<p>Frances is admitted into a small town hospital due to regular contractions @ 4 minutes apart and bloody show. She labors without problems for about 4 hours and then the fetus starts to brady down after SROM. A V/E reveals a prolapsed coed in the vagina. Labor duration: 20 minutes. Scenario duration: 22-27 minutes.</p>	

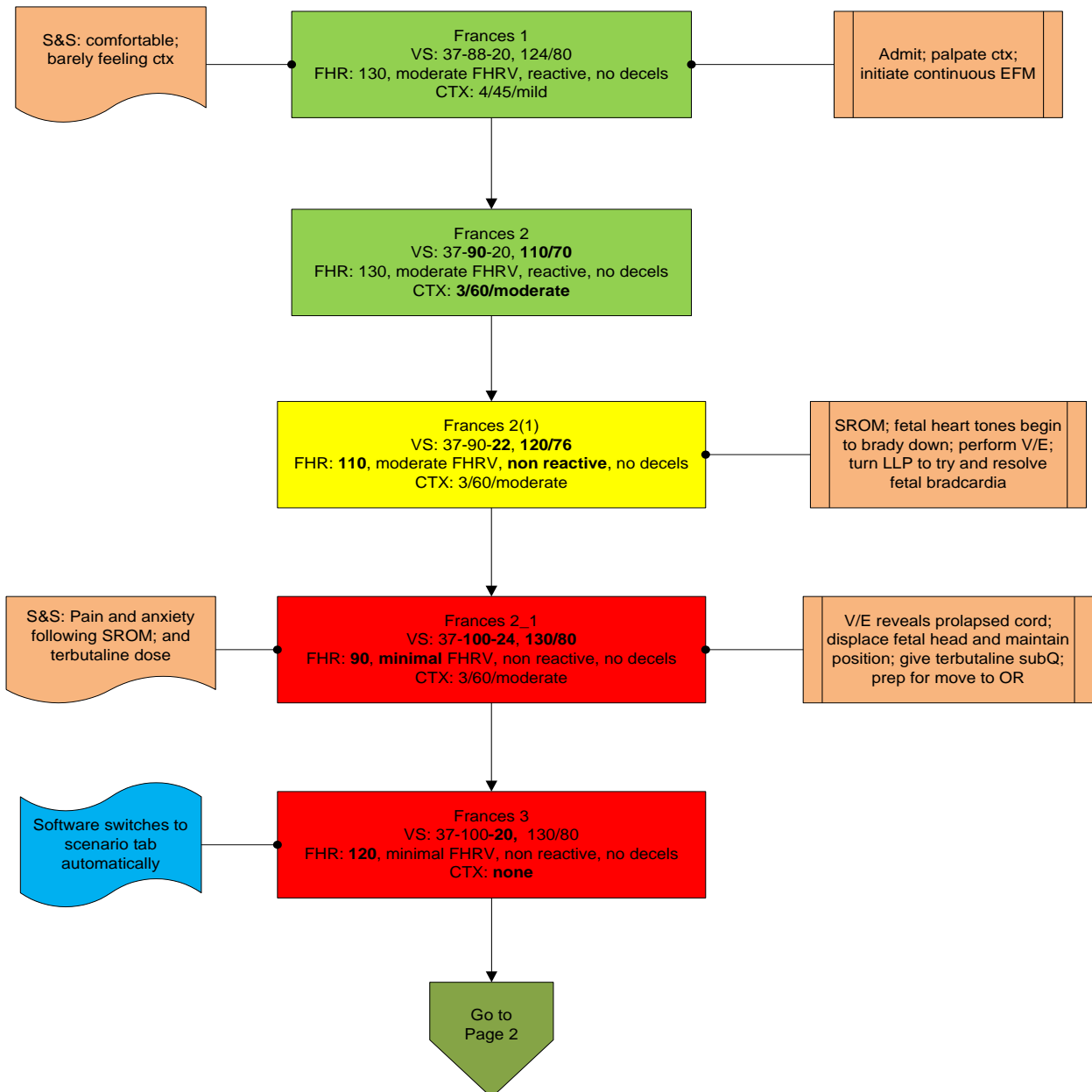


 <p><b>Gaumard®</b> Simulators for Health Care Education</p>	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Frances</b></p> <p>Prolapsed Cord</p>
---	--

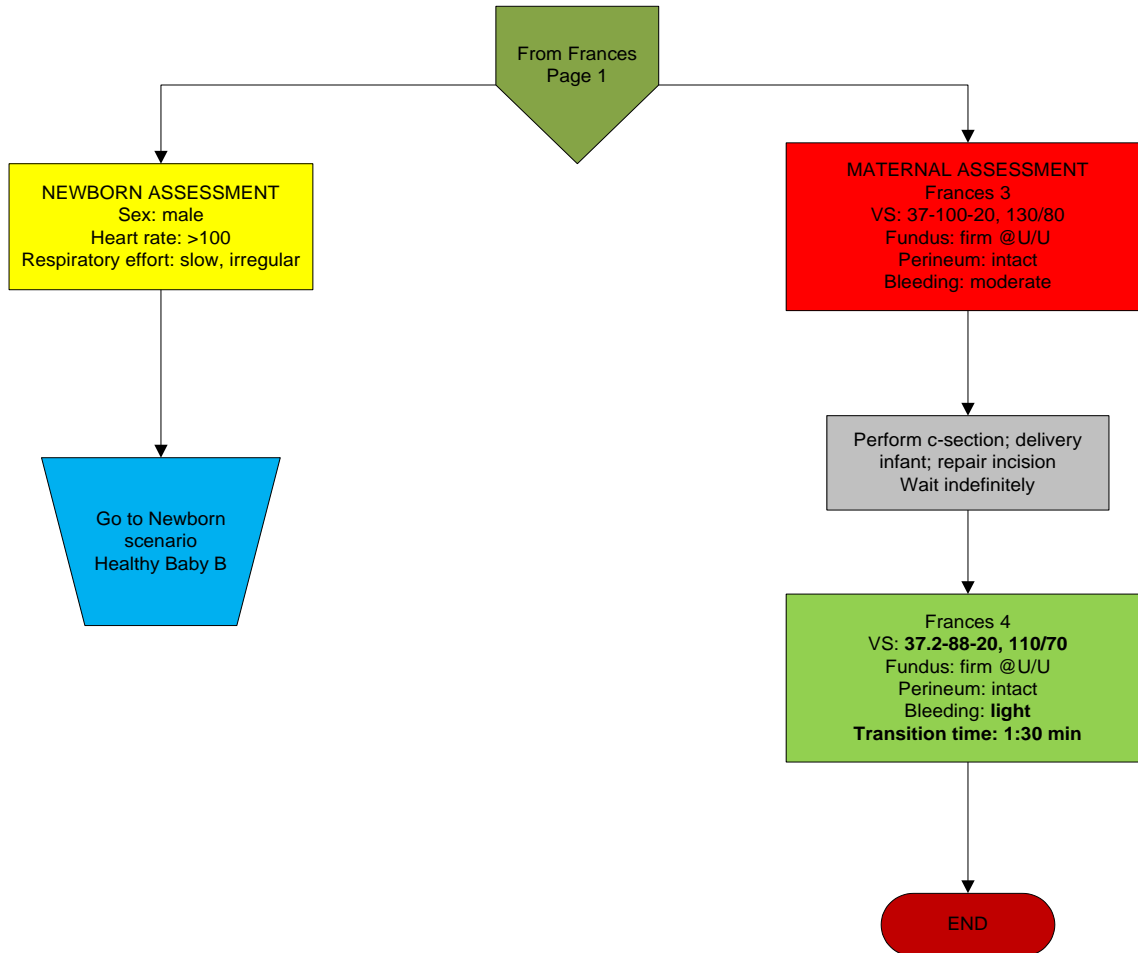




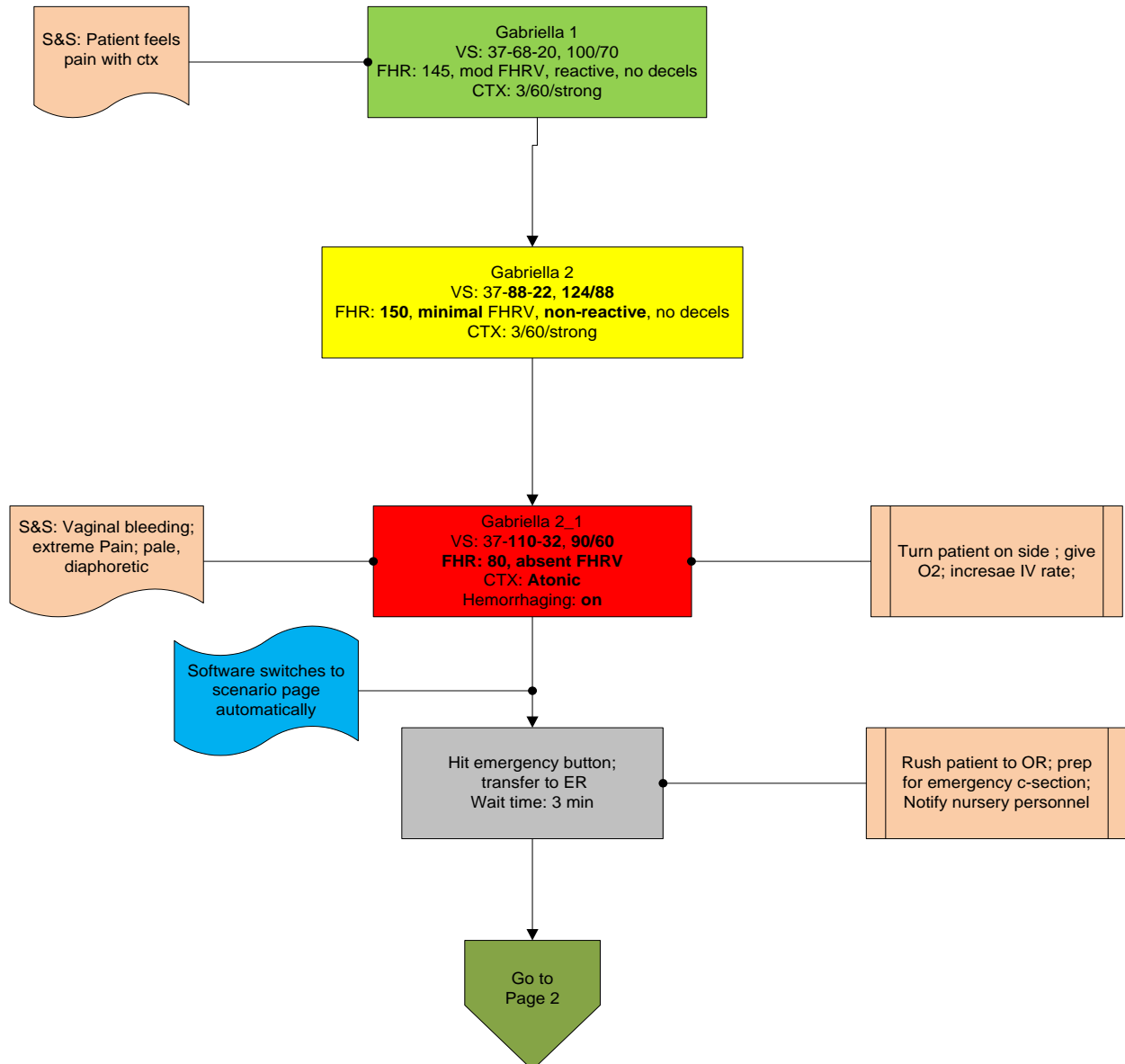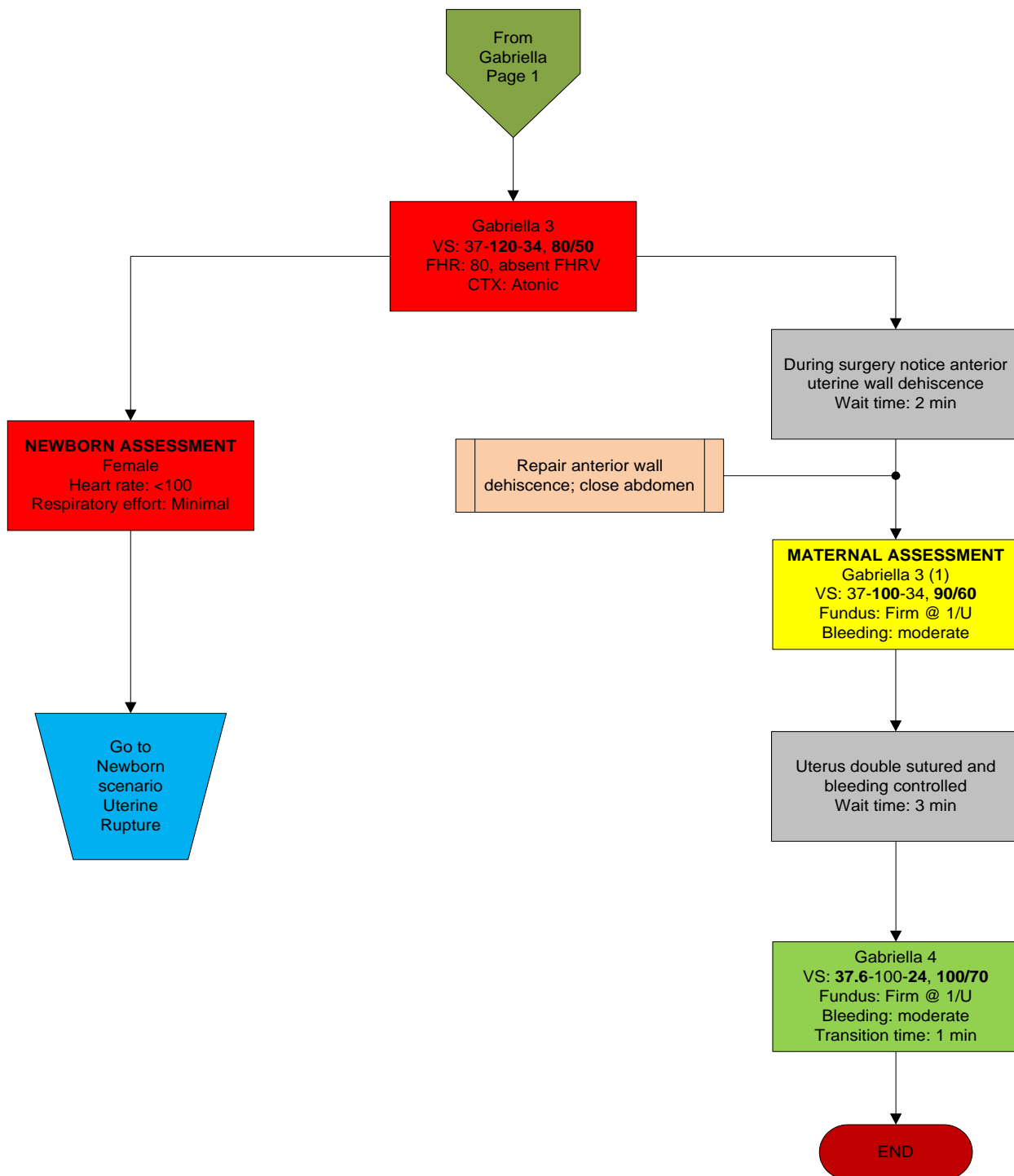
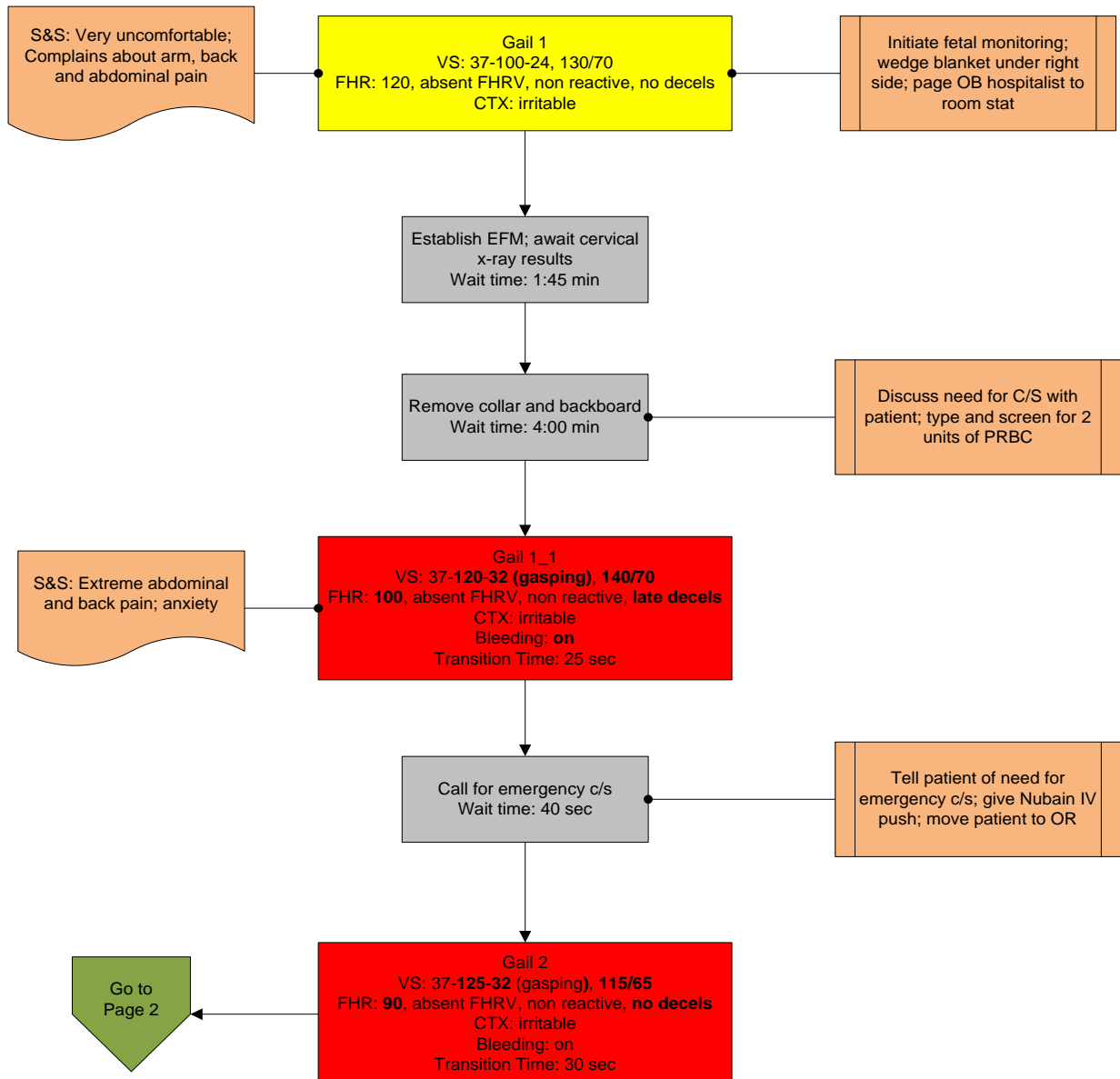
	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Gabriella</b></p> <p>Uterine Rupture</p>
<p>Gabriella is a young Hispanic woman who presents at a small hospital just across the Mexican border. She appears to be in late pregnancy and in active labor. As the nurse helps her to bed she notices a midline abdominal scar. Gabriella has had one previous child in Mexico, but shares no more information. Labor duration: 15 minutes. Scenario duration: 25 minutes.</p>	





	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Gabriella</b></p> <p>Uterine Rupture</p>
---	---




 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® - Labor Scenario <b>Gail</b> Uterine Rupture
<p>Gail is a 29 year old primip @ 35 weeks. She was admitted to L&amp;D from the ER after being involved in car accident. Both she and her husband, Alan, were seriously injured and she is on a backboard wearing a c-collar to stabilize the spine. Her right humerus is fractured and seat belt marks are visible across the abdomen. Labor duration: 18-20 minutes.</p>	



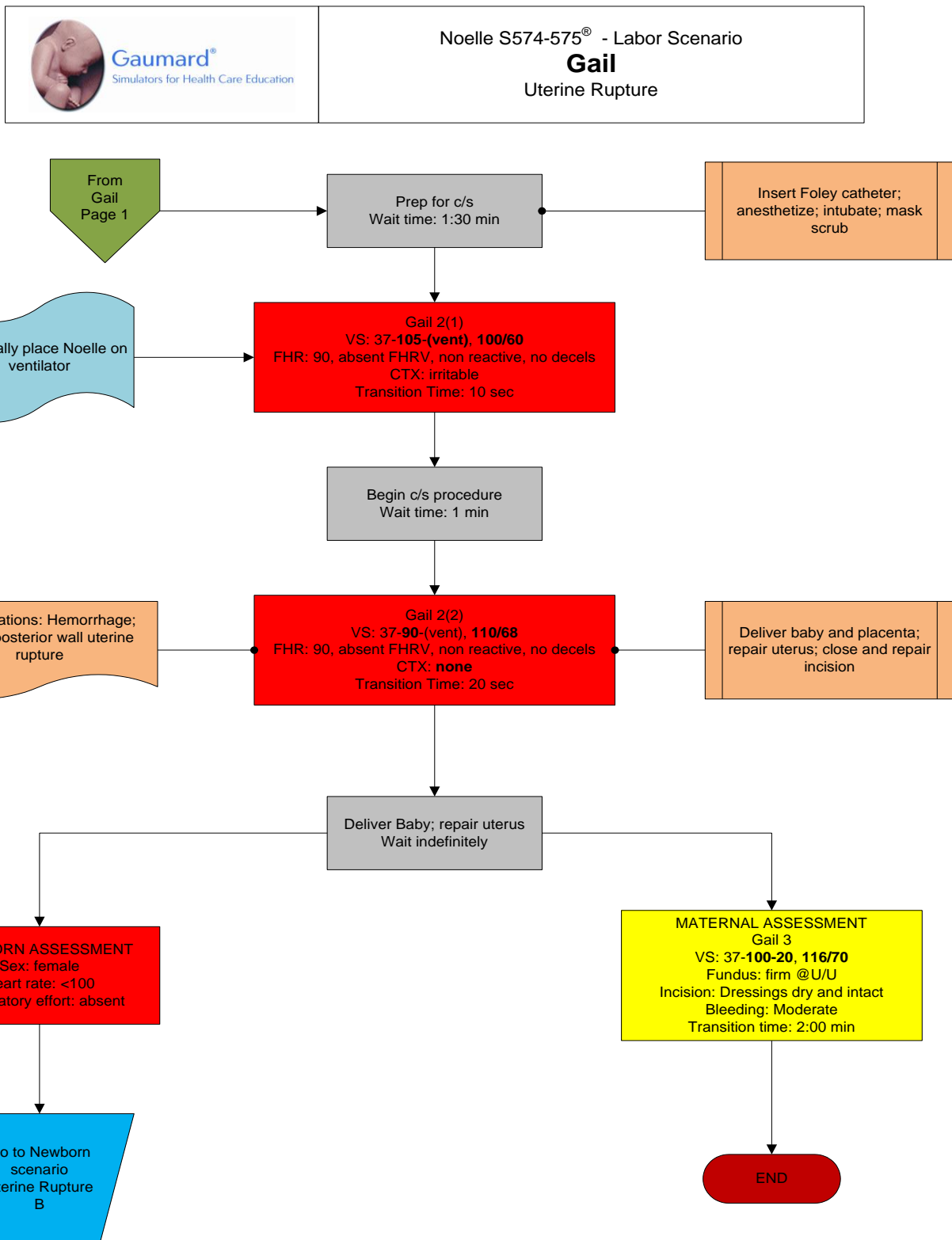



 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® Labor Scenario <b>Heidi</b> Peripartum Hemorrhage/Previa
Heidi is a 25 yr old primip @ 35 weeks. She has experienced several mild bleeding episodes during pregnancy and is known to have a low lying placenta. She arrives in L&D complaining of abdominal cramps and has bright red vaginal bleeding. Labor duration: 30 minutes. Labor duration: 35-45 minutes.	

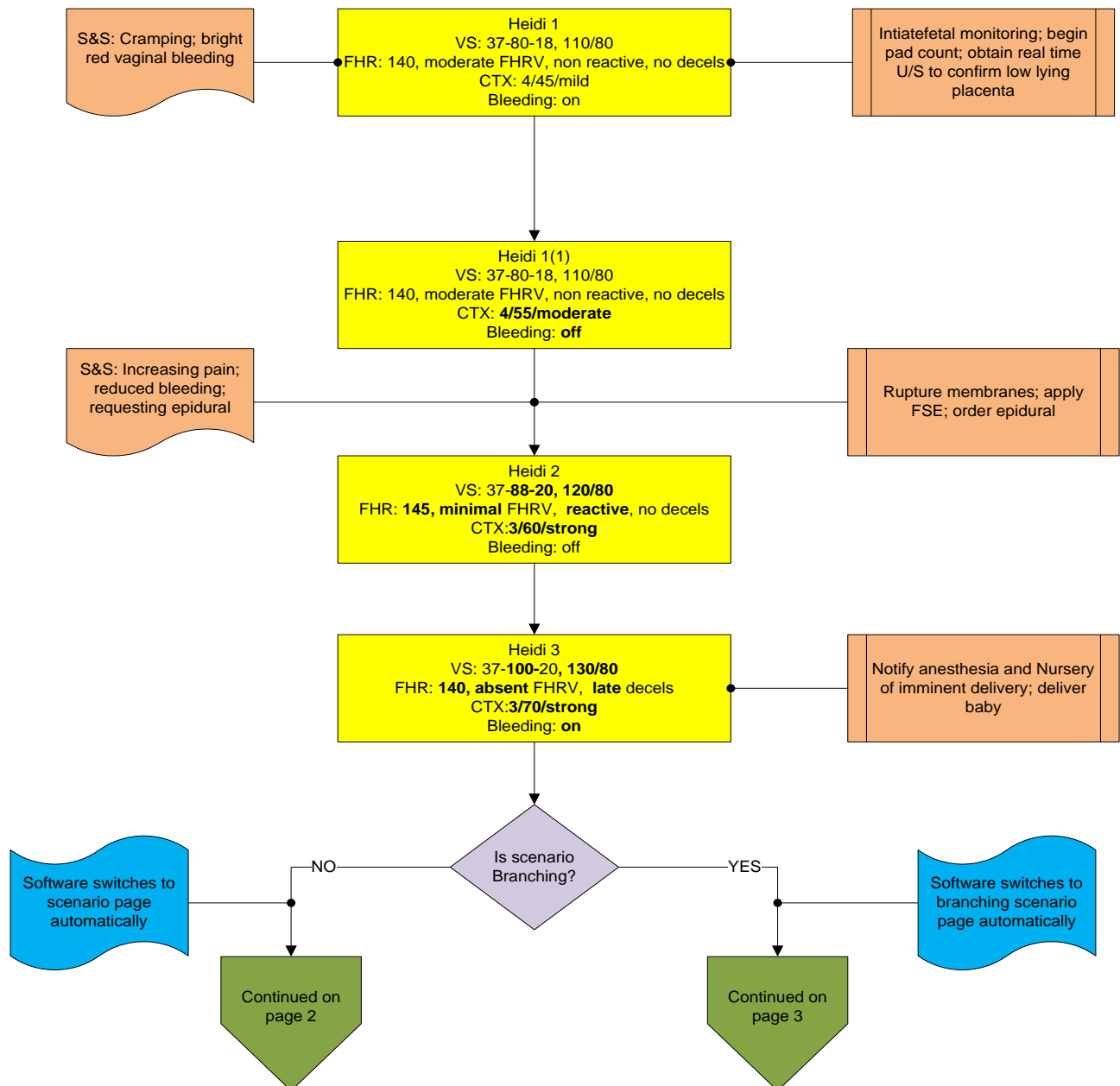





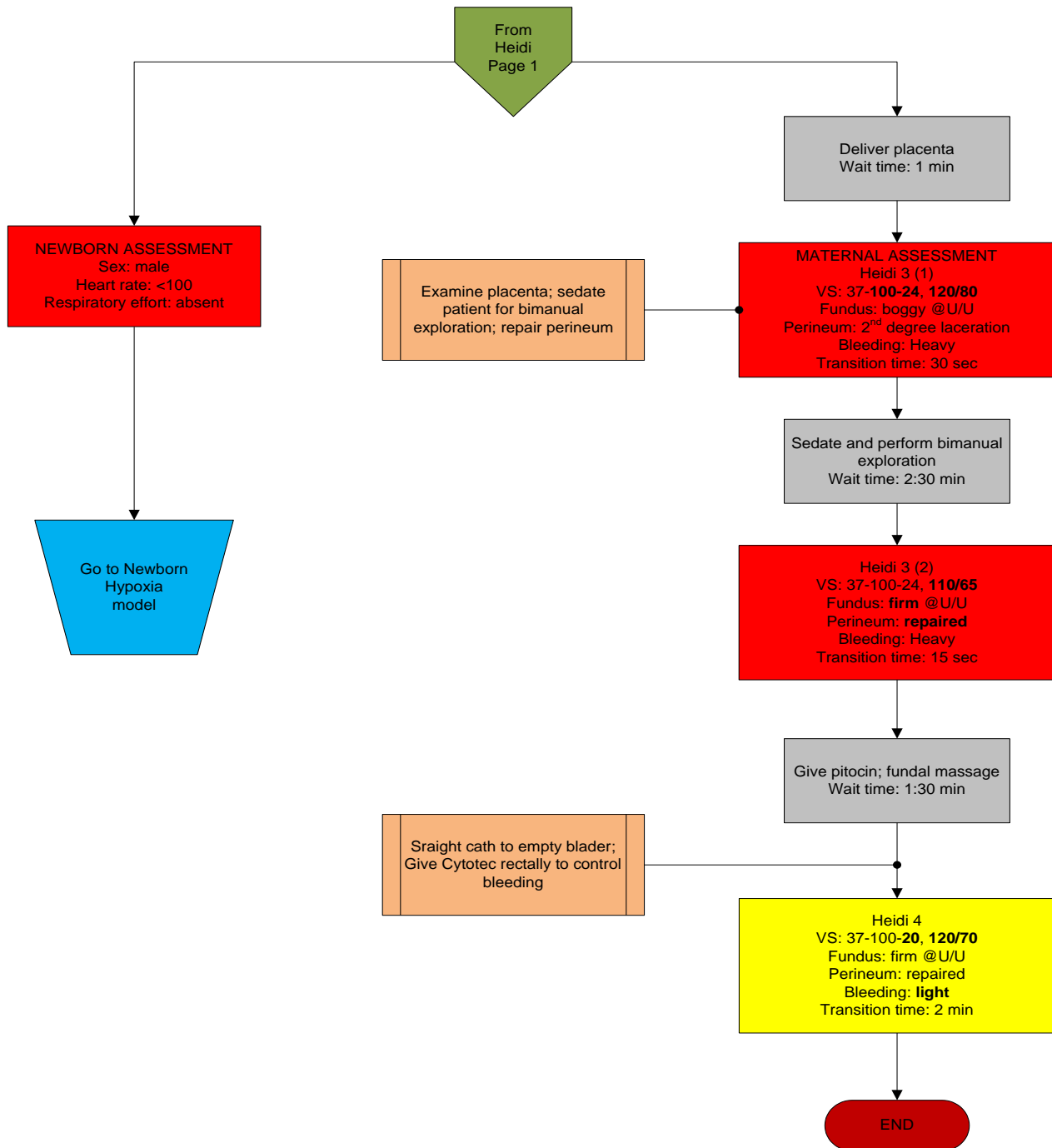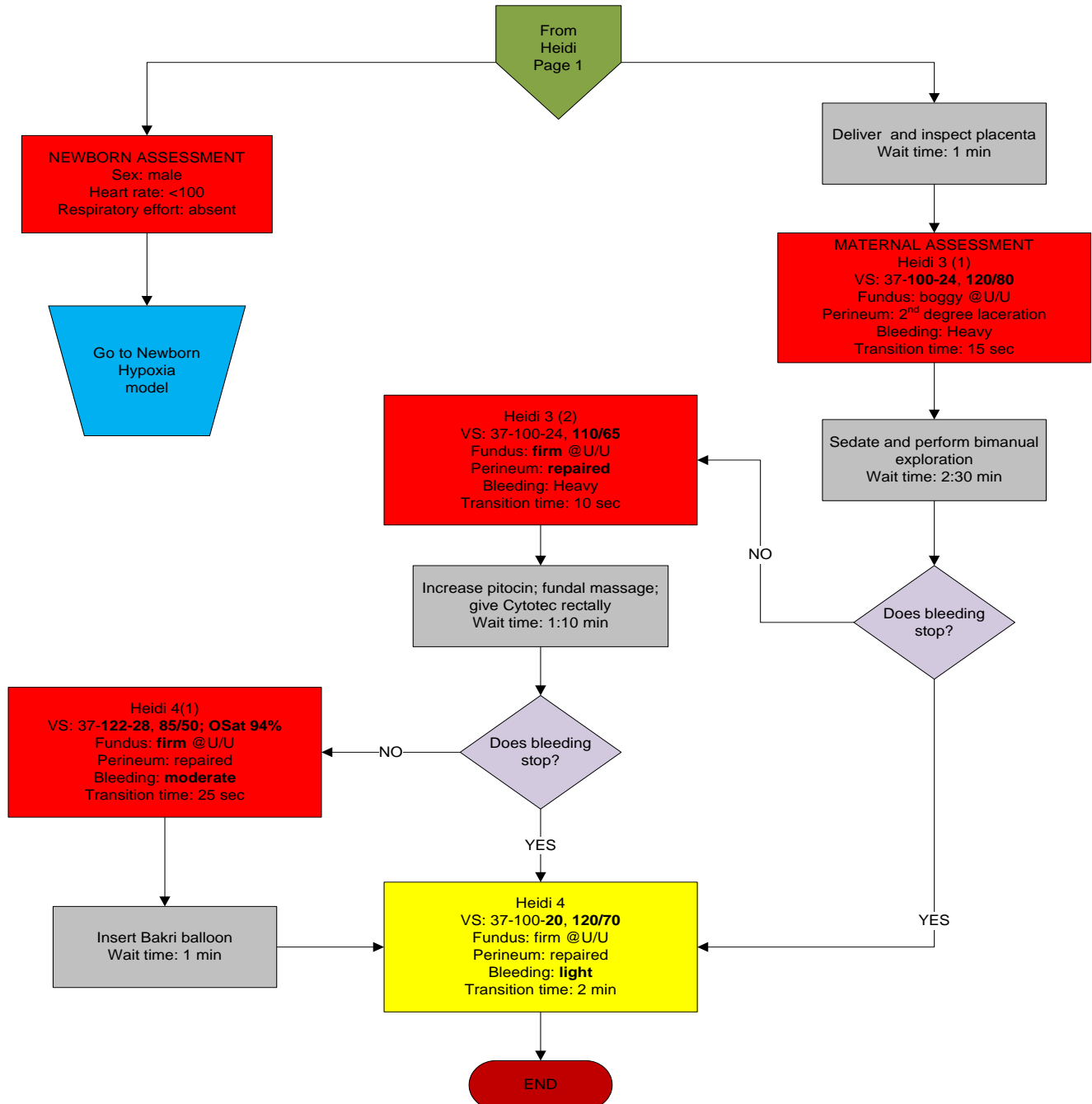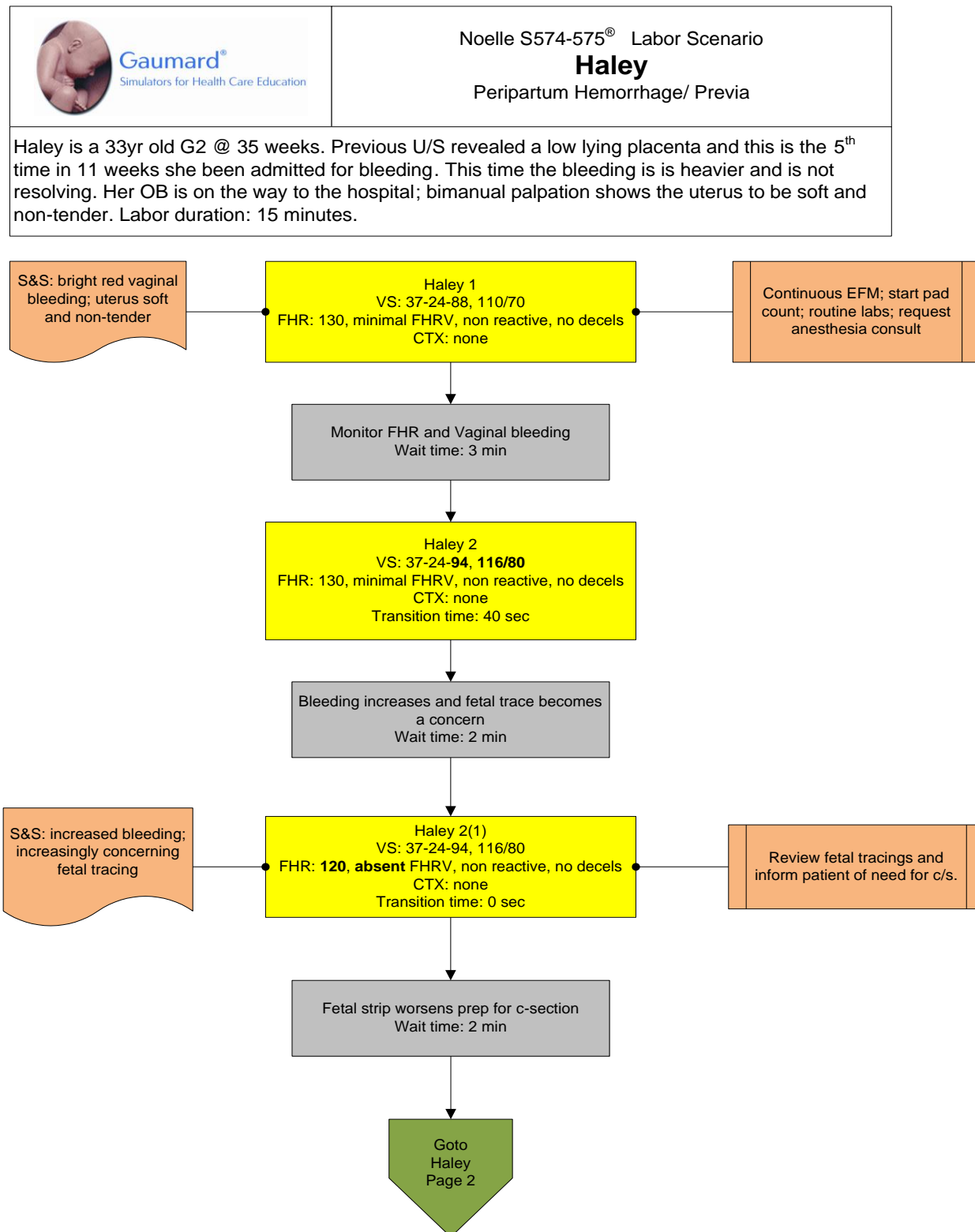Noelle S574-575® Labor Scenario  
**Heidi**  
 Peripartum Hemorrhage/Previa




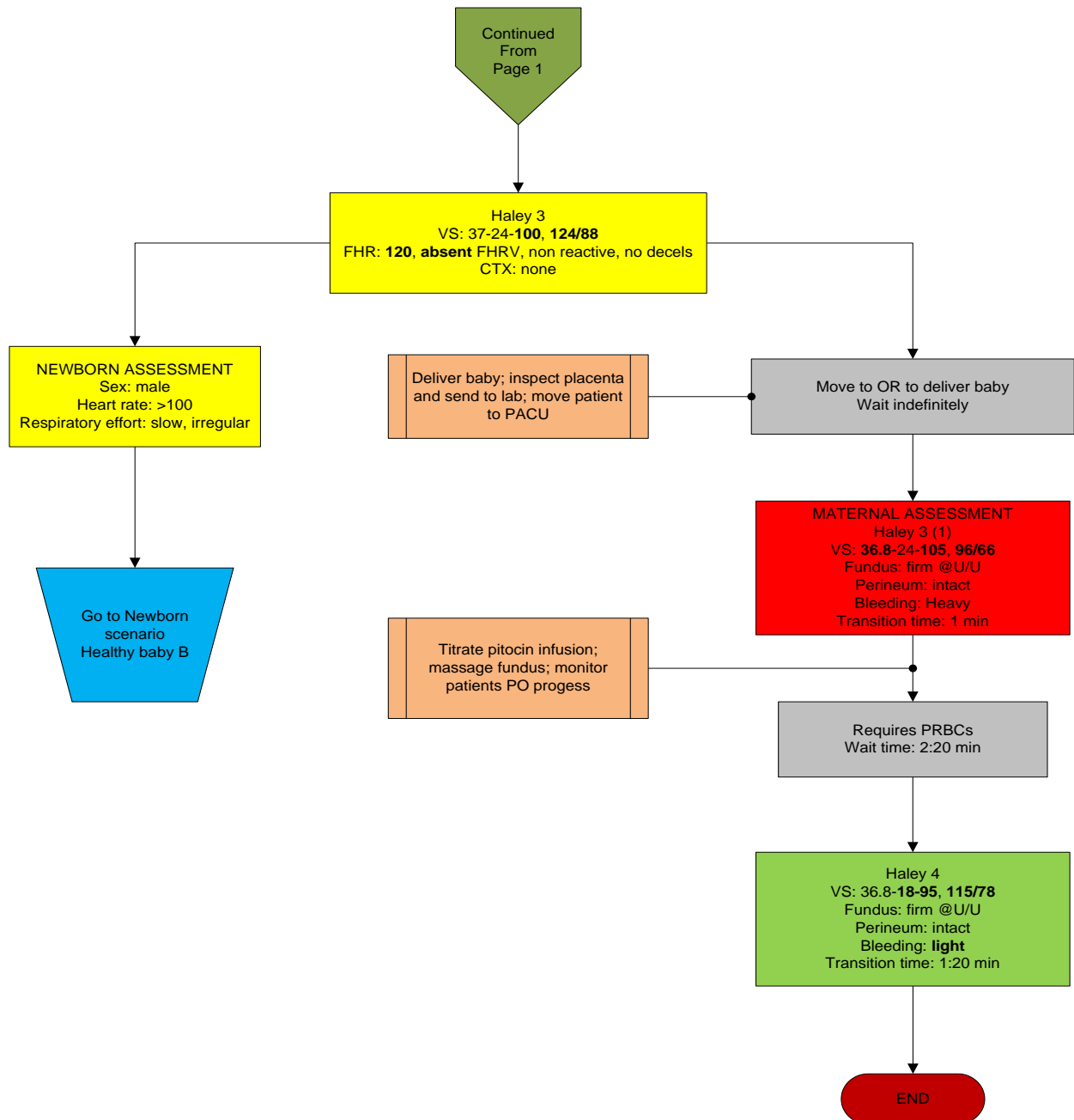
 <p>Gaumard® Simulators for Health Care Education</p>	<p>Noelle S574-575® Labor Scenario <b>Heidi</b> Peripartum Hemorrhage/Previa</p>
--	--

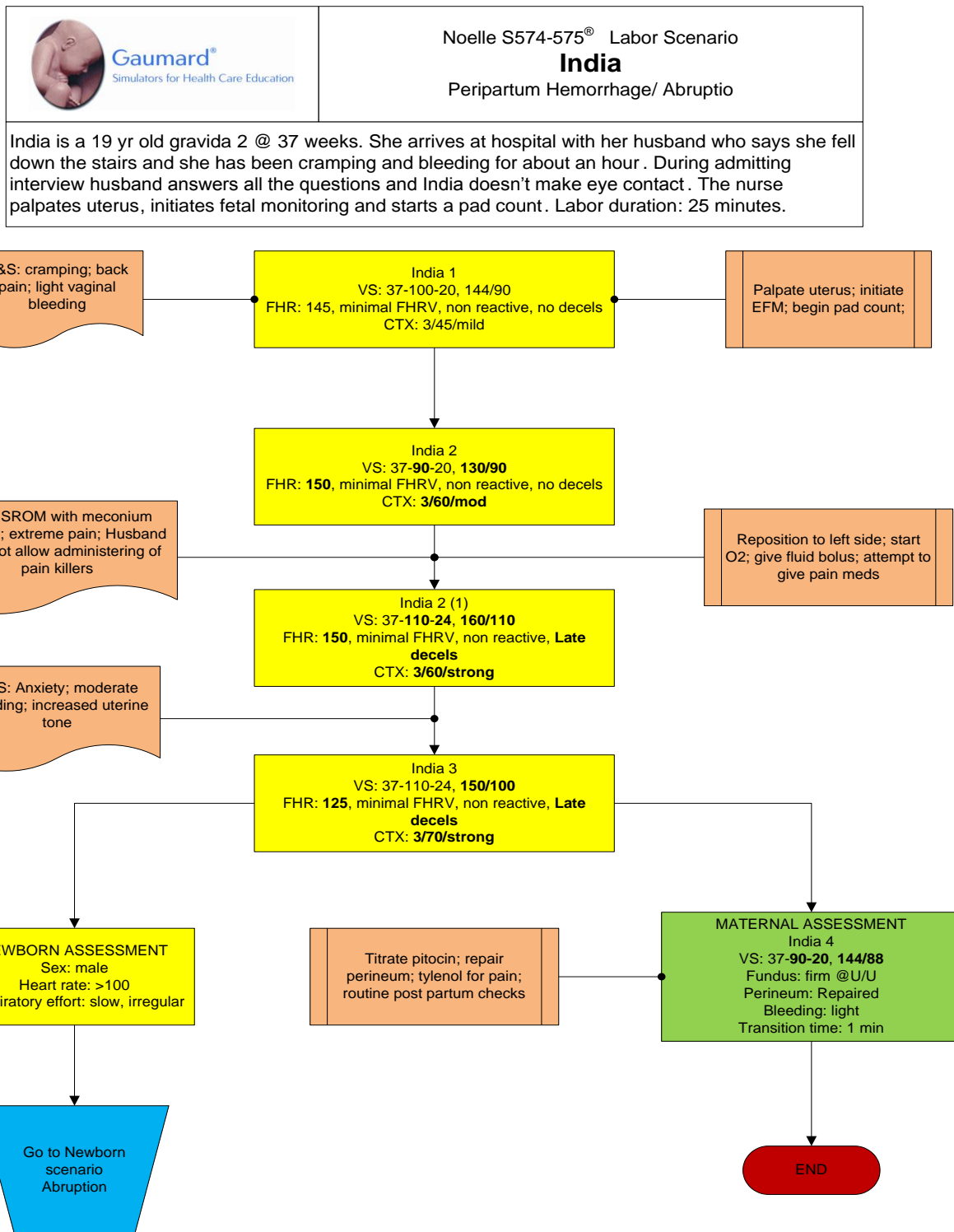





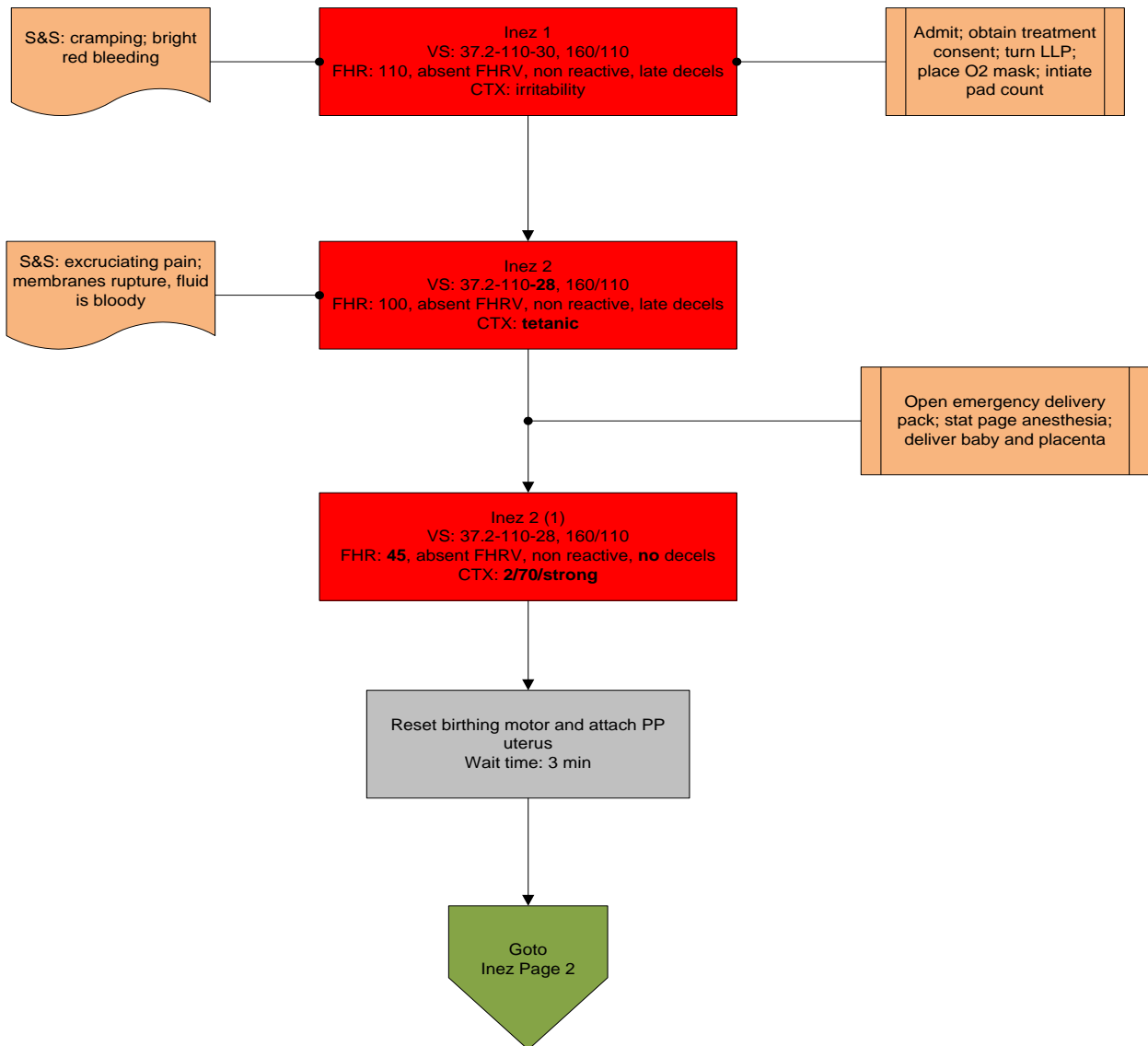



	<p>Noelle S574-575® Labor Scenario</p> <p><b>Haley</b></p> <p>Peripartum Hemorrhage/ Previa</p>
---	---

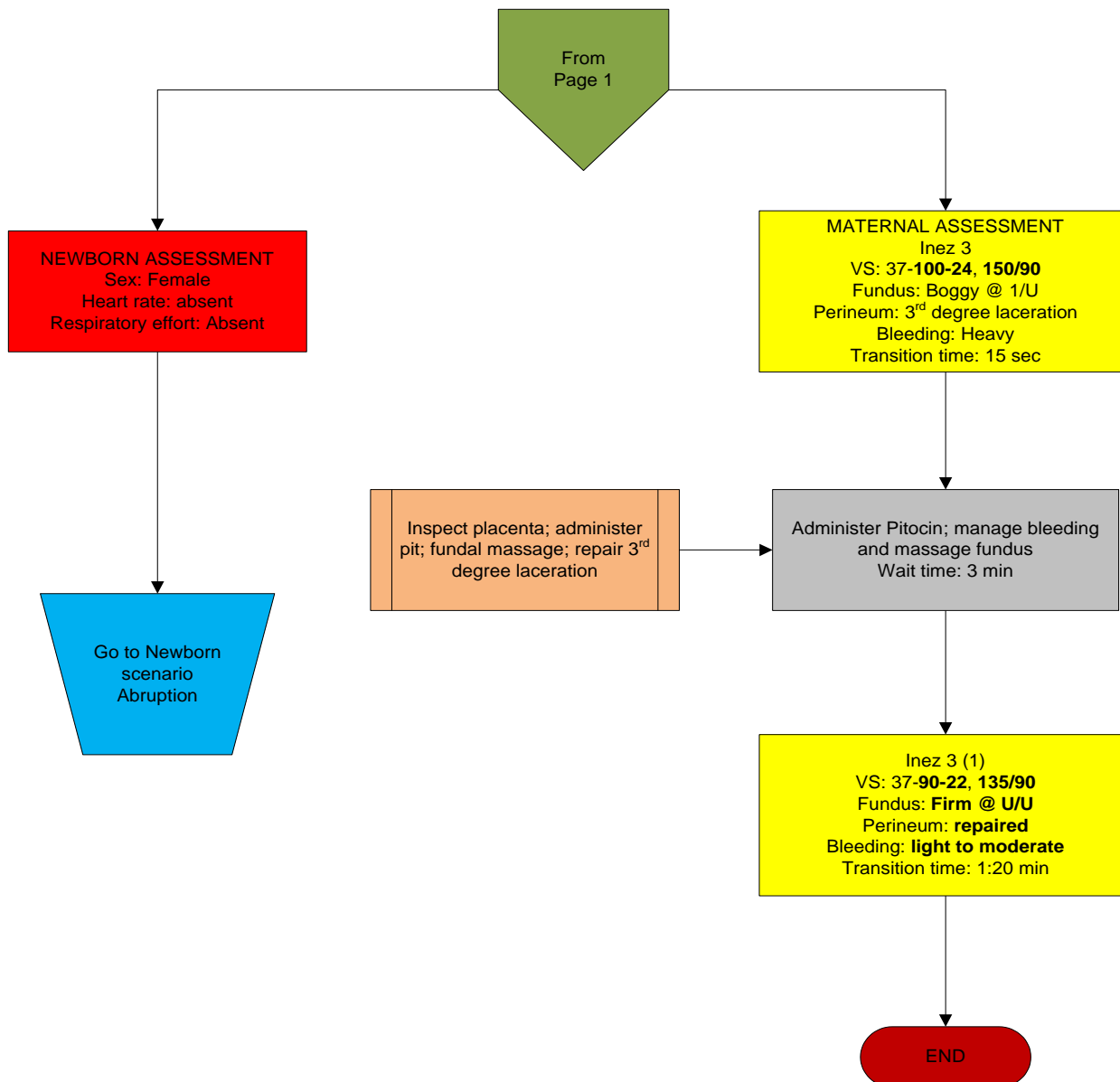





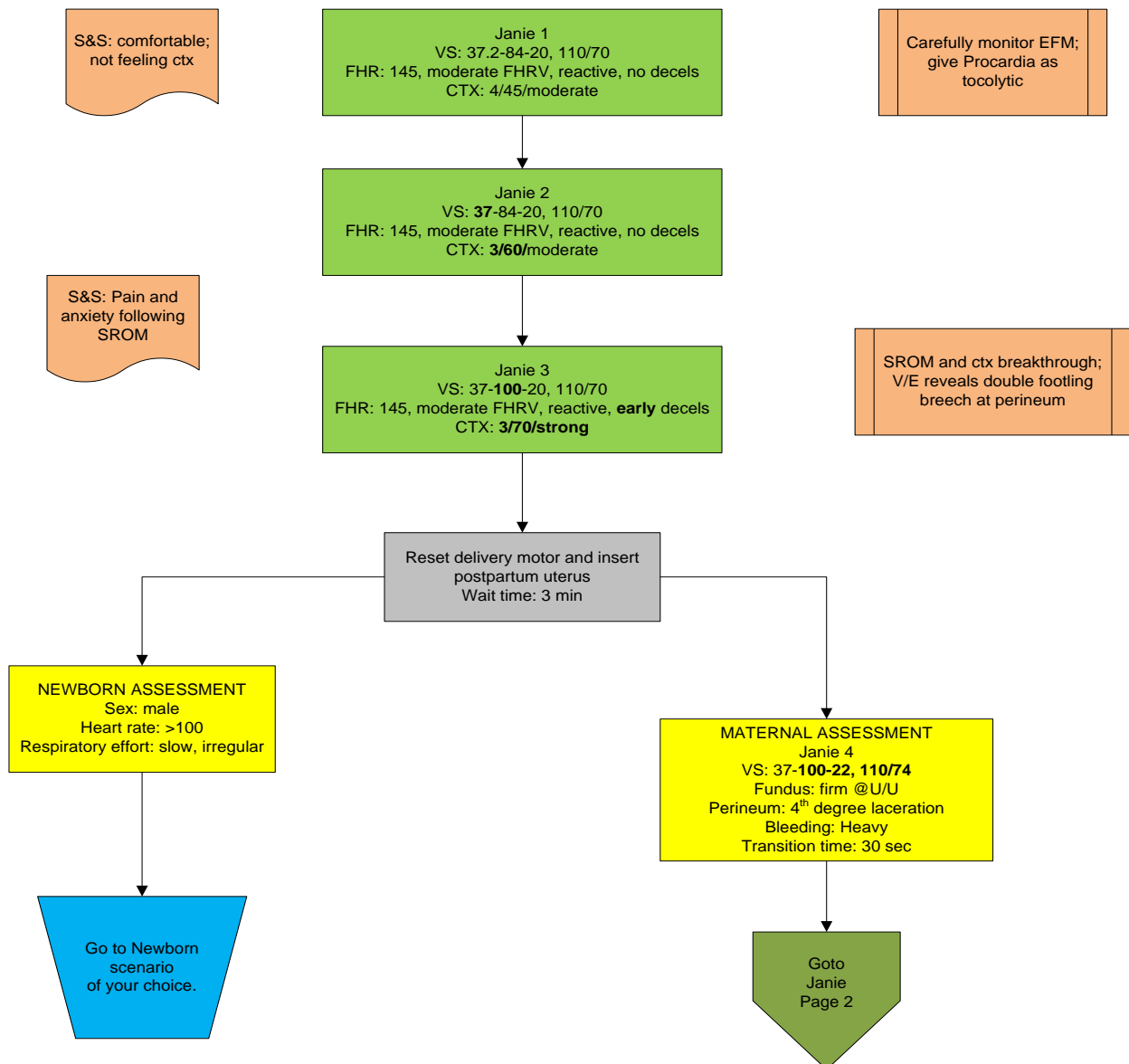
	<p>Noelle S574-575® Labor Scenario</p> <p><b>Inez</b></p> <p>Peripartum Hemorrhage/Abruption</p>
<p>Inez is a 27 yr old primip @ 35 weeks. She arrives at hospital one evening crying and doubled over in pain. She is admitted to a birthing room and the nurse notices bright red blood on Inez's panties. She is having very intense and close contractions. Labor duration: 10 minutes. Scenario duration: 18 minutes.</p>	




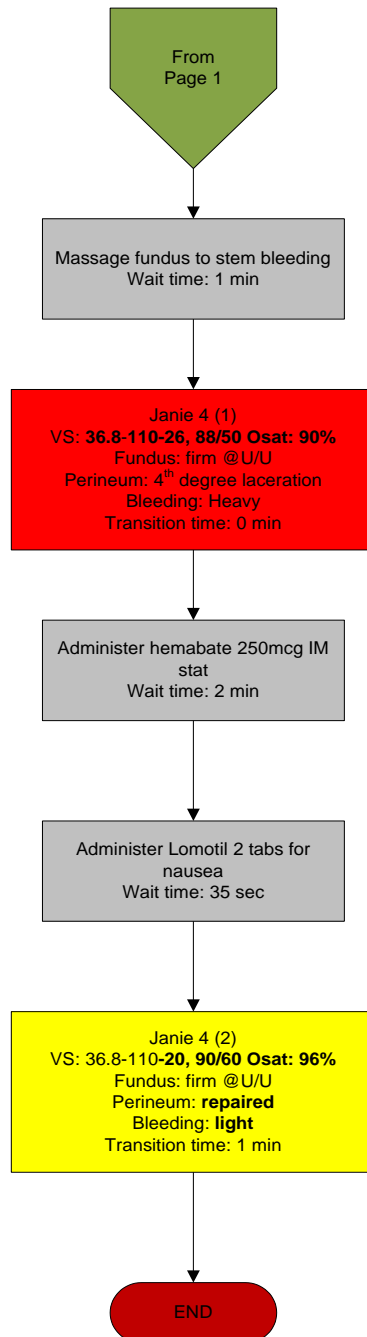
 <p><b>Gaumard®</b> Simulators for Health Care Education</p>	<p>Noelle S574-575® Labor Scenario  <b>Inez</b>                      Peripartum Hemorrhage/Abruption</p>
---	--




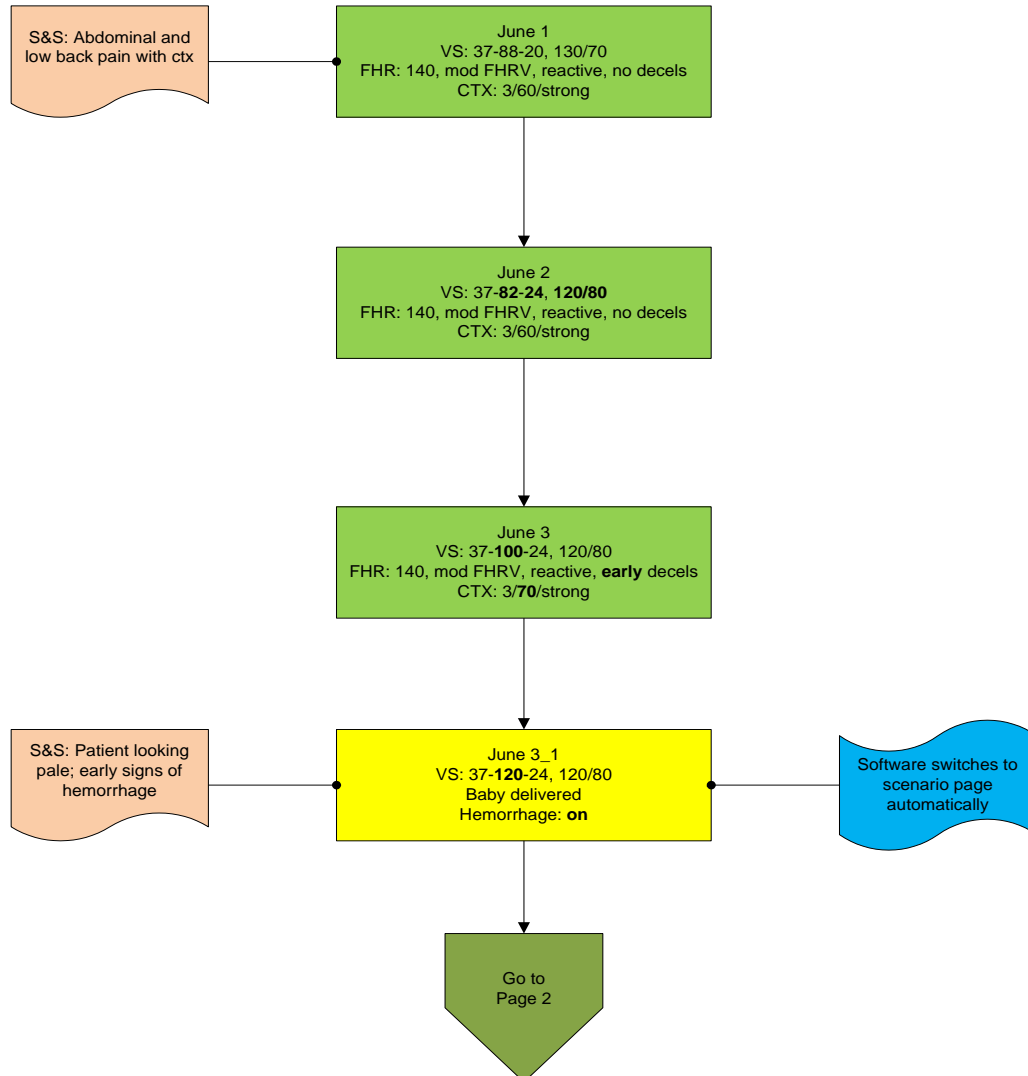
 <b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® Labor Scenario <b>Janie</b> Peripartum Hemorrhage/ PPH
Janie is a 23 yr old G 2 @ 38 weeks. She has experienced several bleeding episodes due to a low lying placenta. She has been counseled about the potential for postpartum hemorrhage . Her religious beliefs prohibit the administration of any blood products . Labor duration: 15 minutes. Scenario duration: 25 minutes.	




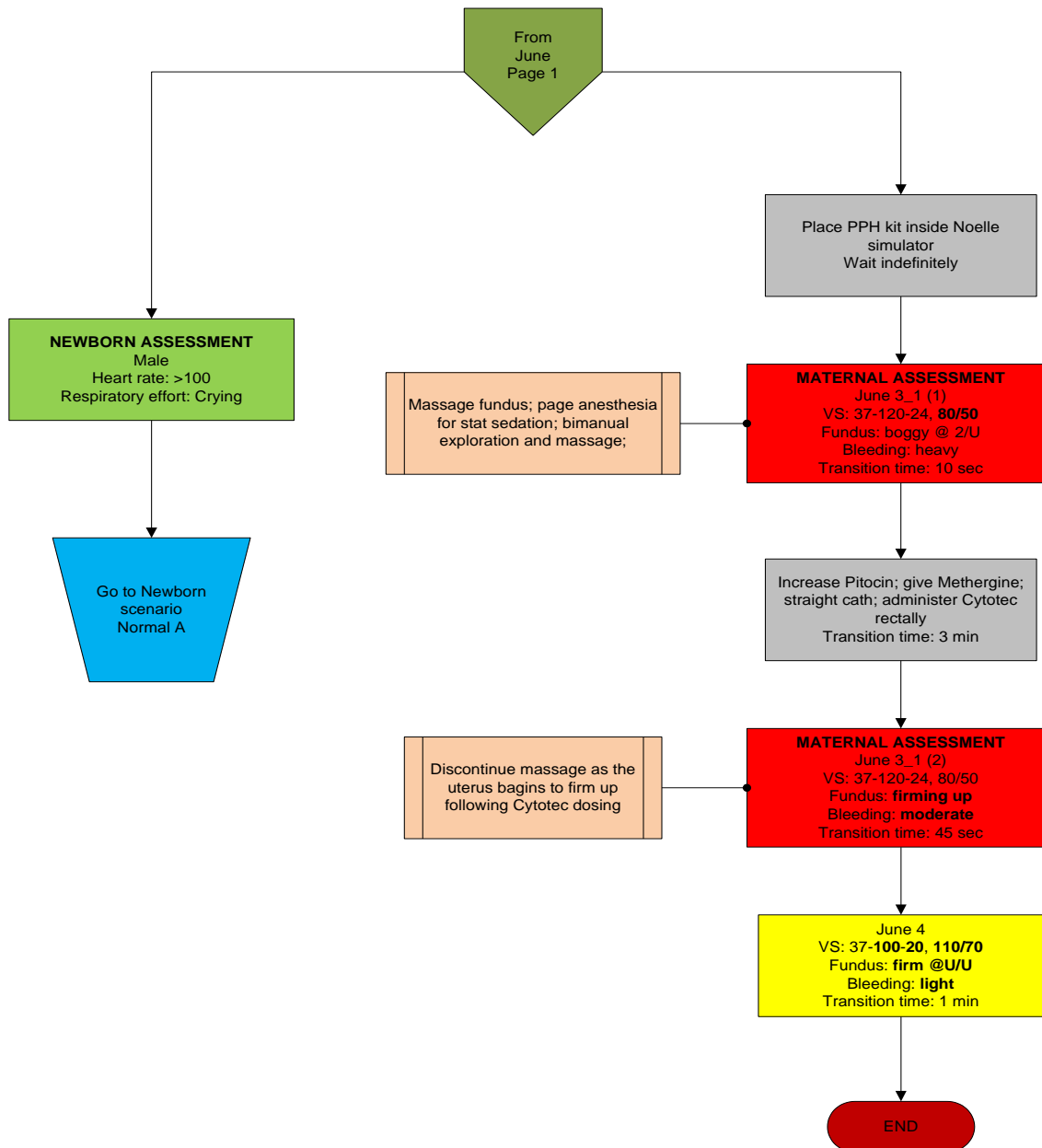
 <p>Gaumard® Simulators for Health Care Education</p>	<p>Noelle S574-575® Labor Scenario <b>Janie</b> Peripartum Hemorrhage/ PPH</p>
--	--




	<p>Noelle S574-575® - Labor Scenario</p> <p><b>June</b></p> <p>Peripartum Hemorrhage/PPH</p>
<p>June is a 31 yr old multip about to have her 5<sup>th</sup> baby. She has had a normal pregnancy and she is planning natural childbirth. She enters the hospital in active labor. The family is very excited as they know this baby is a boy. Labor duration: 15 minutes. Scenario duration: 22-25 minutes.</p>	




	<p>Noelle S574-575® - Labor Scenario</p> <p><b>June</b></p> <p>Peripartum Hemorrhage/PPH</p>
---	--

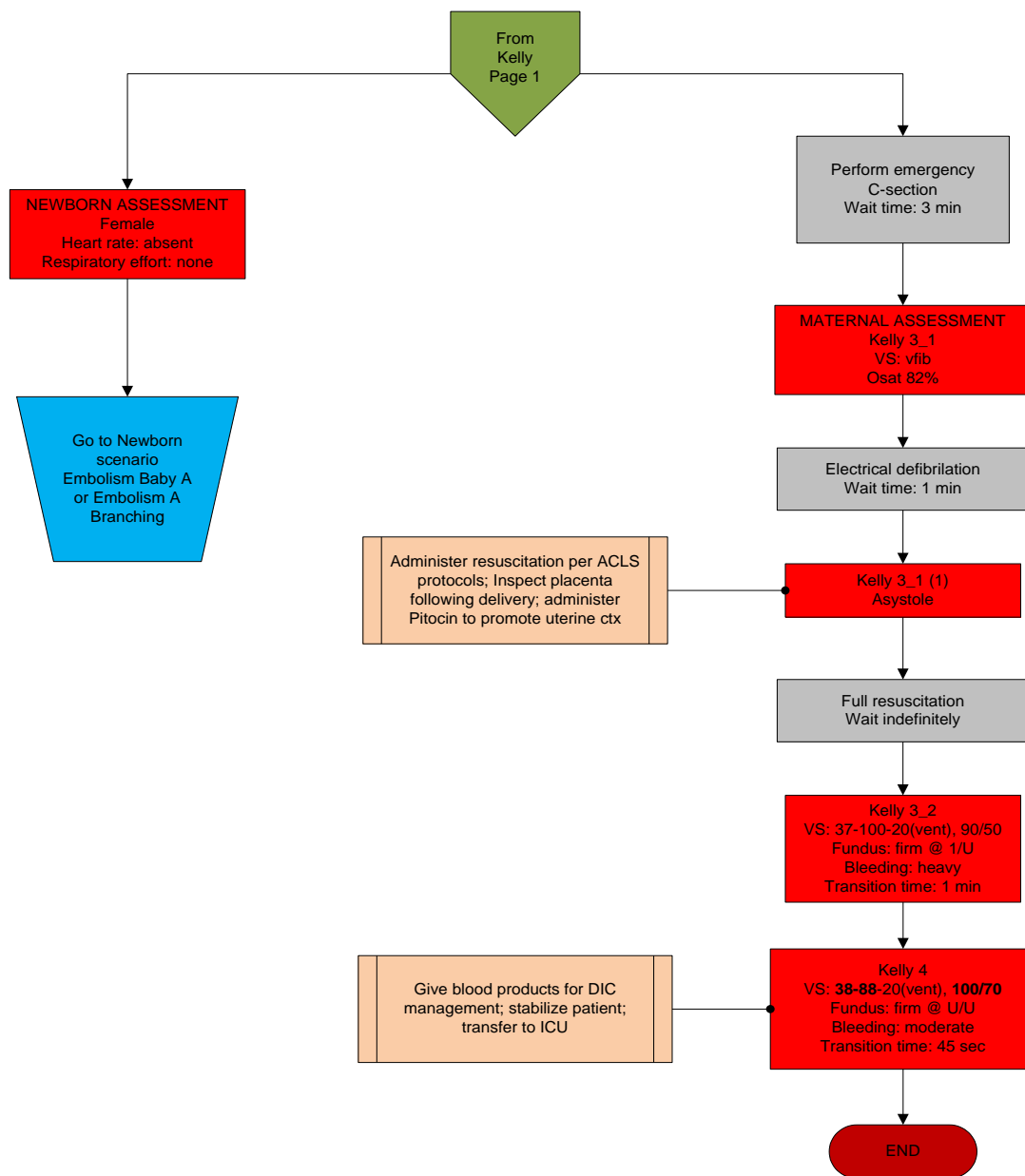





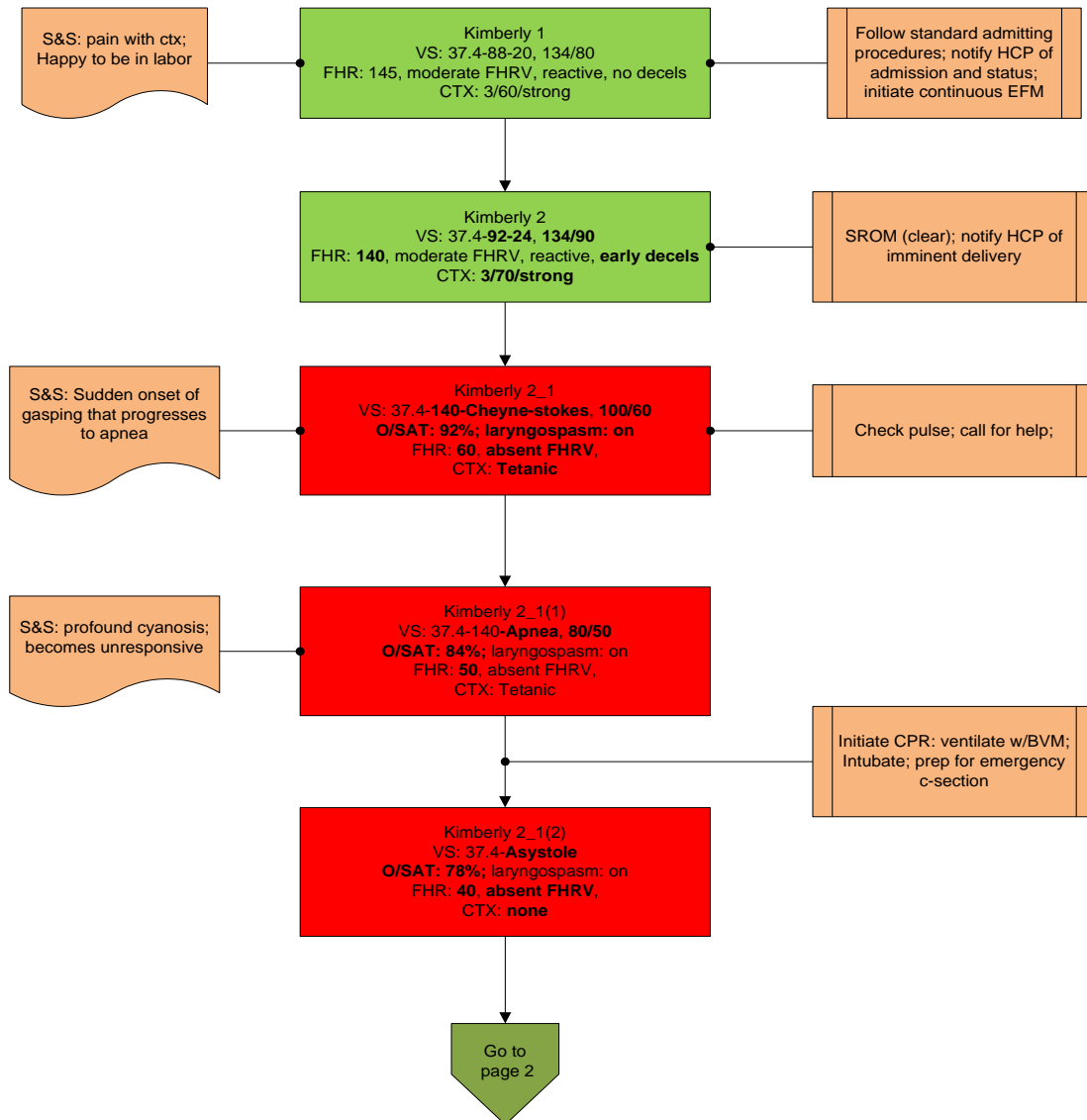
	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Kelly</b></p> <p>Amniotic Fluid Embolism</p>
<p>Kelly is a 34 yr old gravida 5/2 @ 38 weeks. She is scheduled for induction as her last baby weighed almost 10lbs and she experienced a severe shoulder dystocia with that delivery. She has gained 43lbs with this pregnancy and her GTT is borderline. Labor duration: 25 minutes. Scenario duration: 35-40 minutes.</p>	




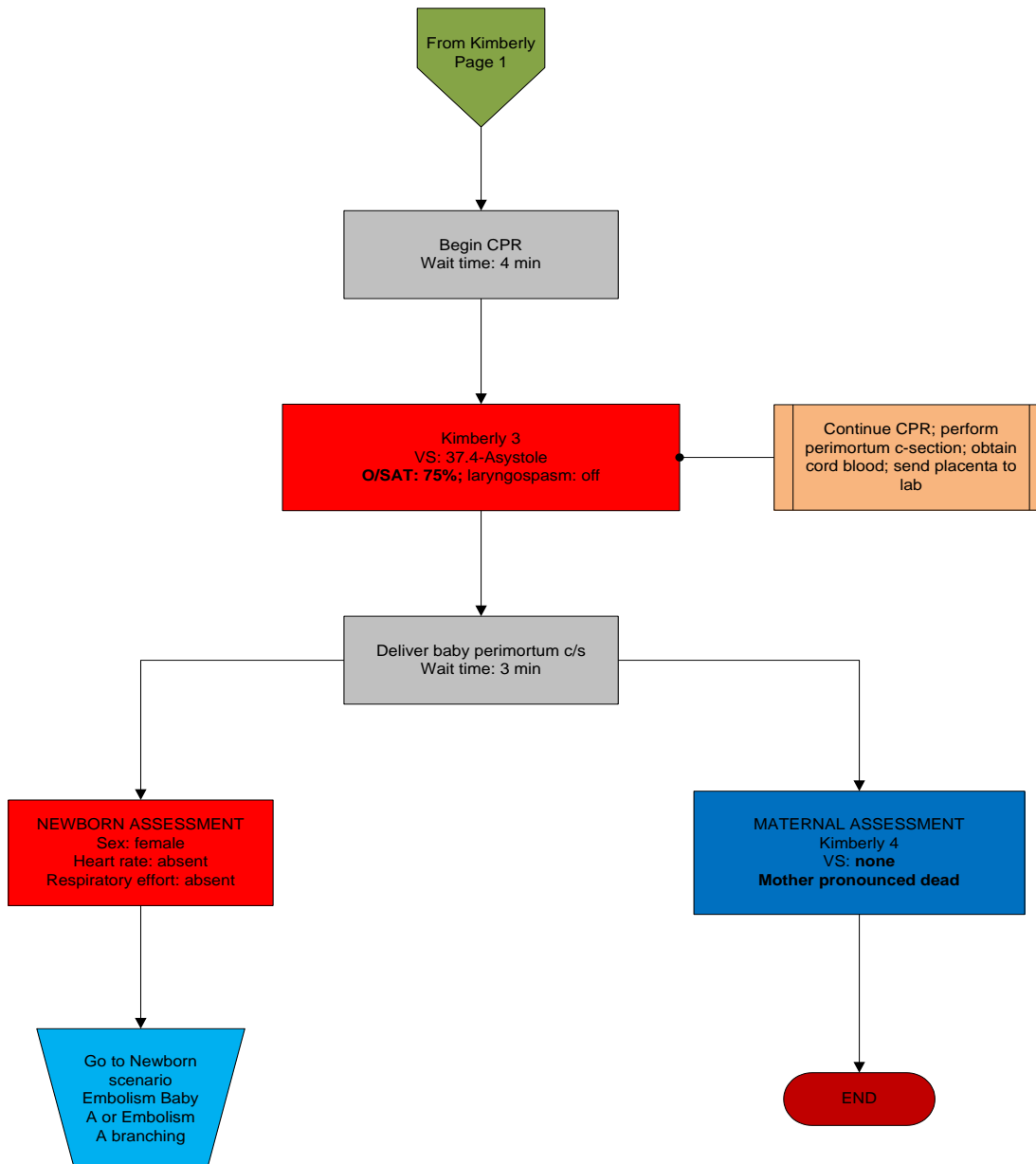
	<p>Noelle S574-575® - Labor Scenario</p> <p><b>Kelly</b></p> <p>Amniotic Fluid Embolism</p>
---	---



	<p>Noelle S574-575® Labor Scenario</p> <p><b>Kimberly</b></p> <p>AFE</p>
<p>Kimberly is a 27 yr old multip @ 42 weeks. She began having contractions at home and now they are becoming stronger. She is excited to be finally going into labor. By the time the nurse completes admission Kimberly is requesting pain meds as her labor is progressing quickly. Labor duration: 25 minutes. Scenario duration: 30 minutes.</p>	



	<b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® Labor Scenario <b>Kimberly</b> AFE
---	---	---






**Gaumard®**  
Simulators for Health Care Education


Noelle S574-575® - Labor Scenario


## **Madonna** Preterm Labor

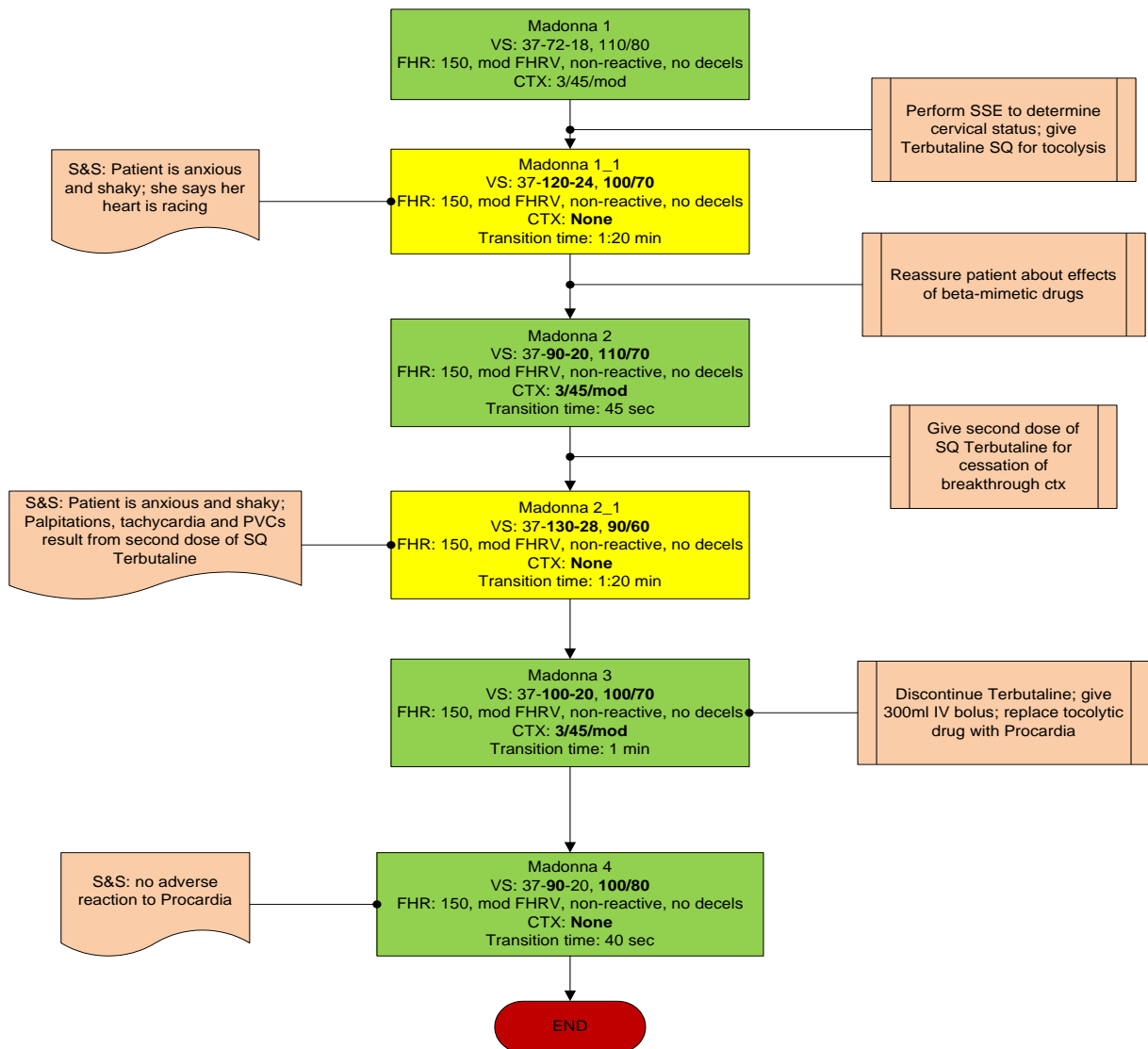Madonna is a 41 yr old multip @ 31 weeks. She has experienced difficult pregnancies and has one Downs Syndrome baby. She has had several episodes of preterm contractions that resolved with LLP bed rest. This time the bed rest and oral hydration are not resolving the contractions; in fact, they seem to be getting worse. Labor duration: 35 minutes.



© 2009 Gaumard Scientific  
All rights reserved


	<b>Gaumard®</b> Simulators for Health Care Education	Noelle S574-575® Labor Scenario <b>Maria</b> Preterm Labor
Maria is a 30 yr old multip @ 27 weeks. She has an 11 yr old and has been trying for more children. She has had 2 miscarriages in the last 4 years and she lost both due to an incompetent cervix. This time a McDonalds suture was placed @ 14 weeks. Labor duration: 15 minutes.		

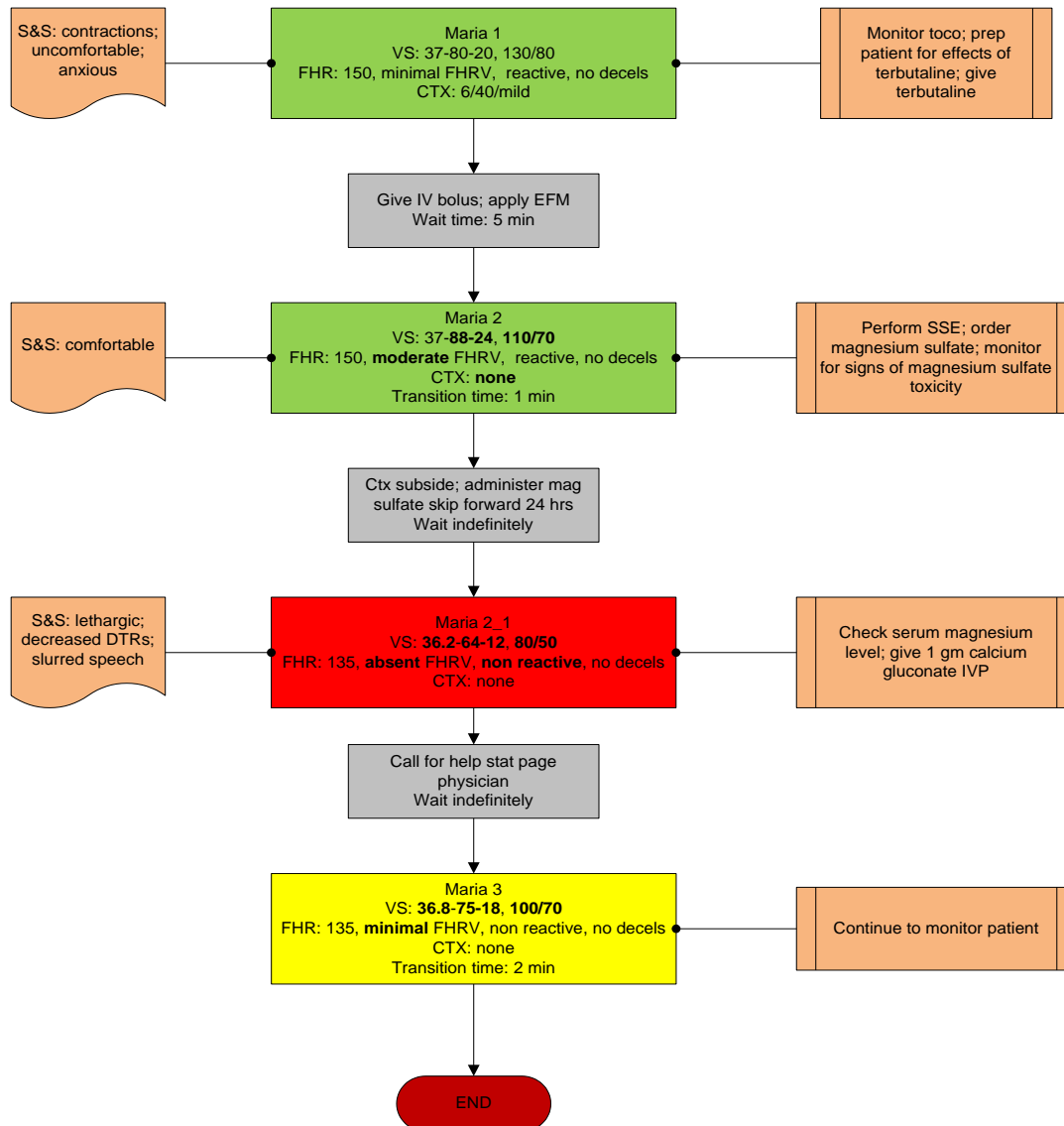

## Newborn Flowcharts

Quick Start Newborn	
1	Alice's Baby
2	Asphyxia
3	Beth's (Dona's) Baby
4	Cynthia's Baby
5	Elaine's Baby
6	Francine's Baby
7	Gloria's Baby
8	Helen's (Irene's) Baby
9	MAS
10	RDS
11	TTN

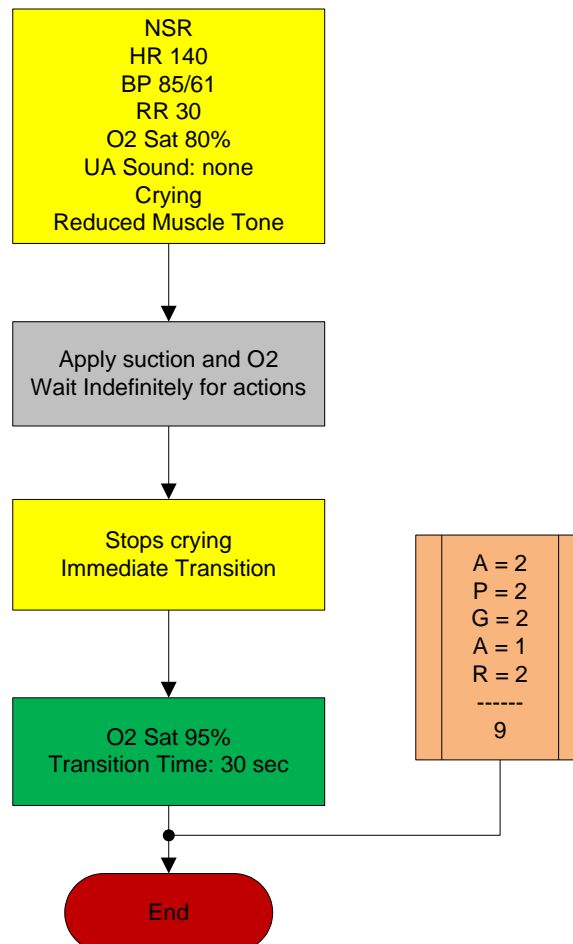 <p><b>Gaumard®</b> Simulators for Health Care Education</p>	<p>Newborn HAL® <b>Alice's Baby</b></p>
<p>Healthy baby</p>	








**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®  
**Asphyxia**

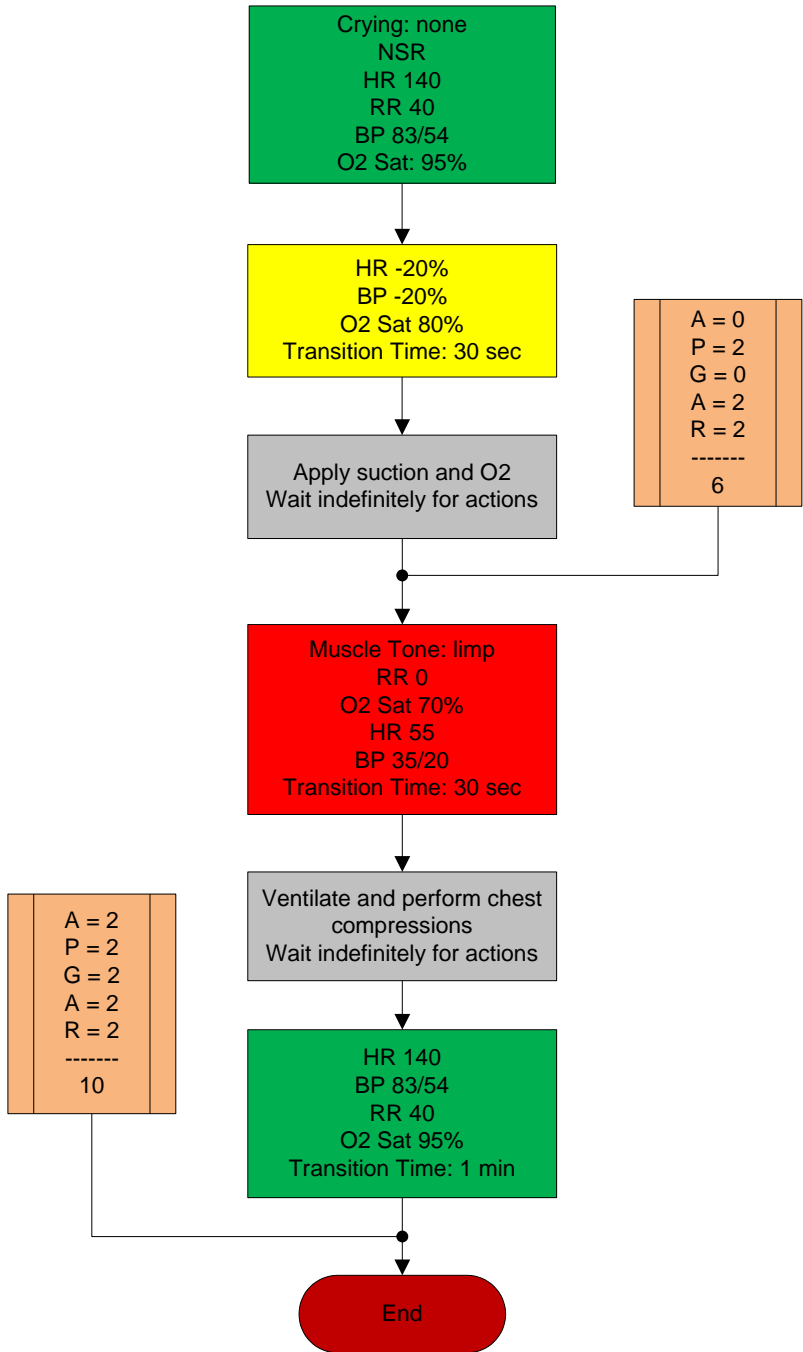Baby has an asphyxia attack and the providers need to give ventilations to help bring back the vitals to a healthy state.






**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®  
**Beth's (Donna's) Baby**

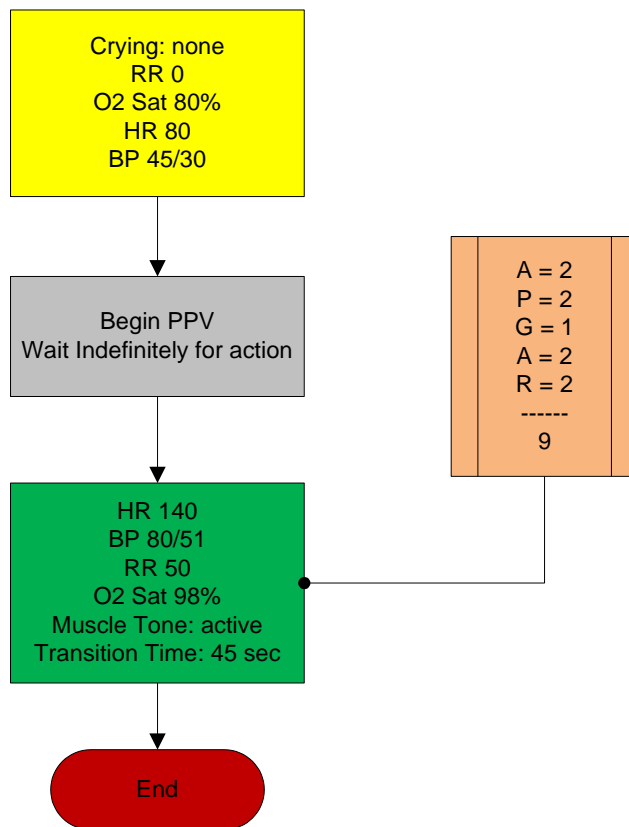Baby is born with a mild asphyxia that needs attention. Once ventilations are started, the baby's vitals go to a healthy state.






**Gaumard®**  
Simulators for Health Care Education

## Newborn HAL® Cynthia's Baby

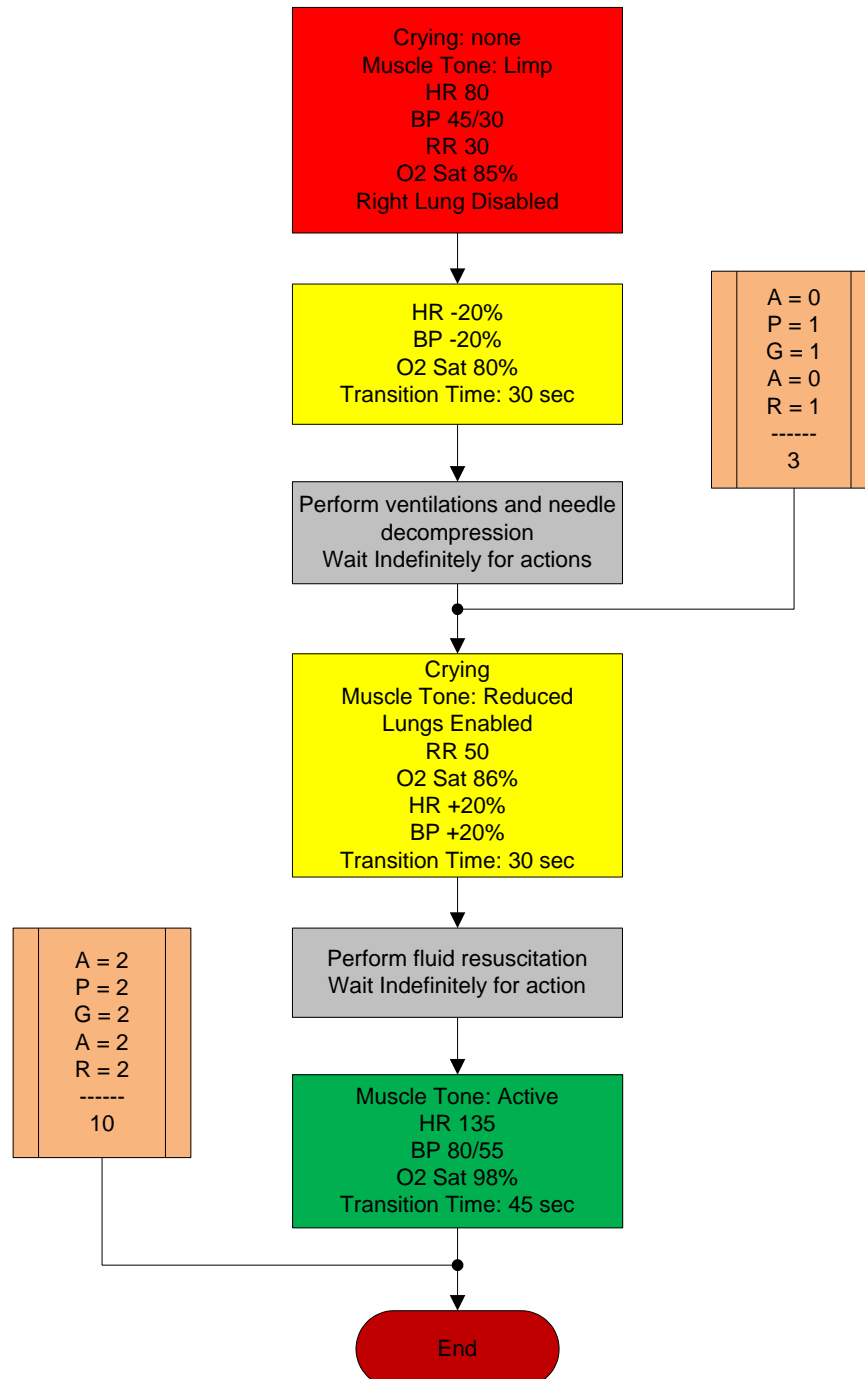Male infant with central cyanosis, limp, flaccid and requires immediate resuscitation. No spontaneous movement of right arm is noted. Stat CXR reveals a fractured right clavicle and right pneumothorax.







**Gaumard®**  
Simulators for Health Care Education

## Newborn HAL® Elaine's Baby

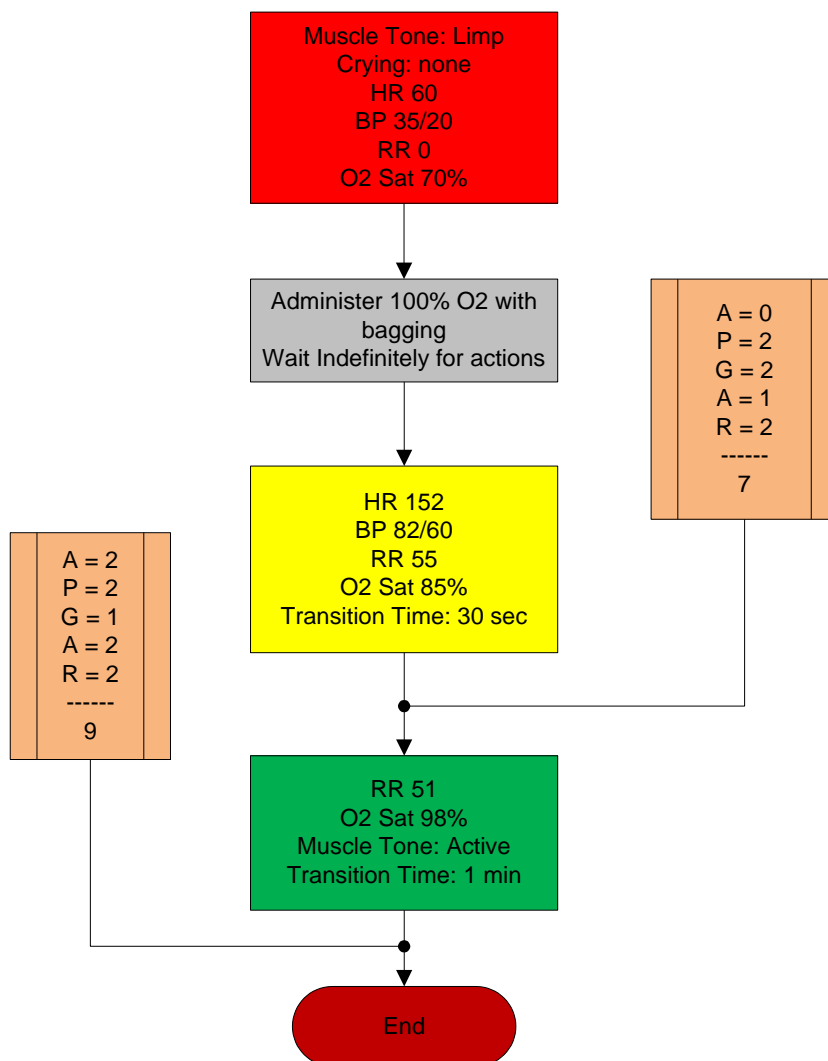This baby is born with moderate asphyxia, and will require CPR and oxygen to bring the vitals to a healthy state.






**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®

**Francines’s Baby**

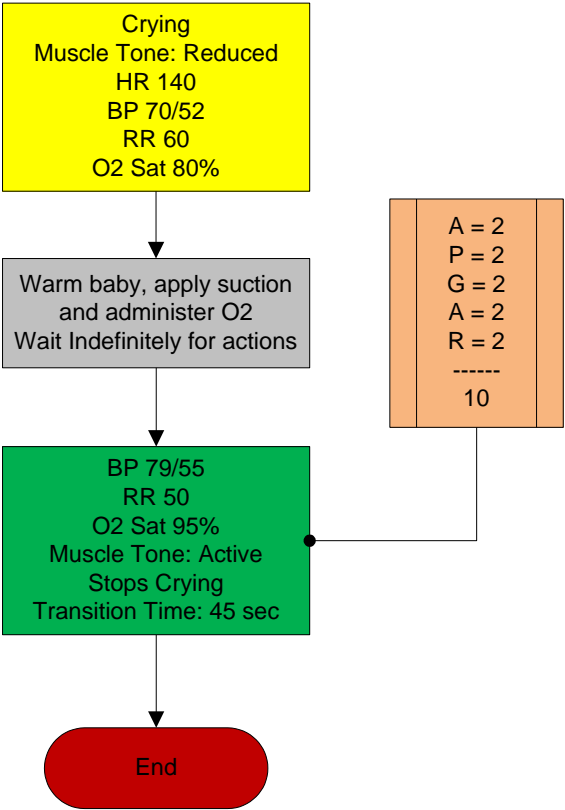This baby was born through a C-Section and is responsive but needs some attention, after a while all vitals go to a healthy state.






**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®  
**Gloria's Baby**

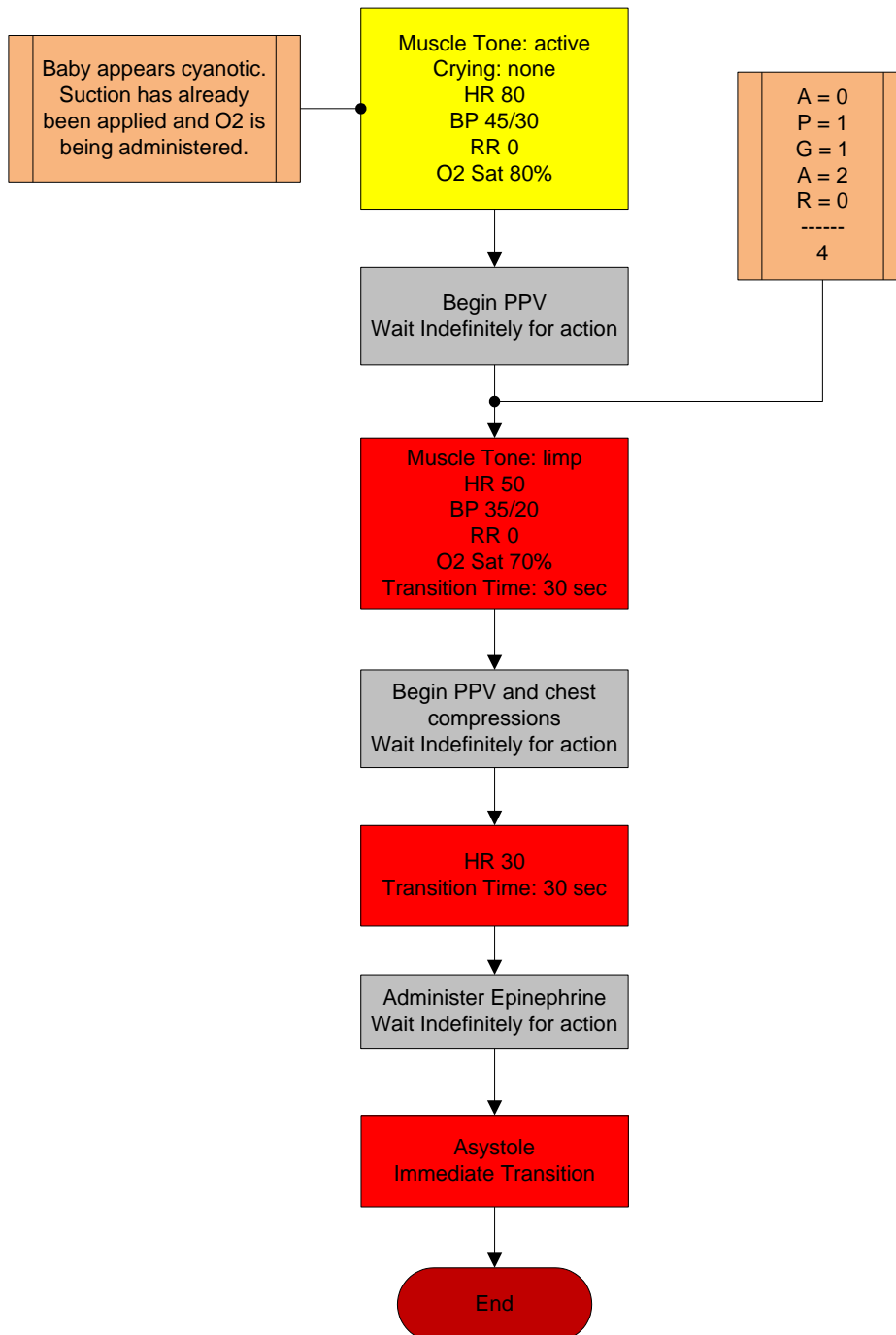This baby is born with mild asphyxia, but no matter how good the interventions are, this disastrous intrapartum complication results in neonatal death.






**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®  
**Helen's (Irene's) Baby**

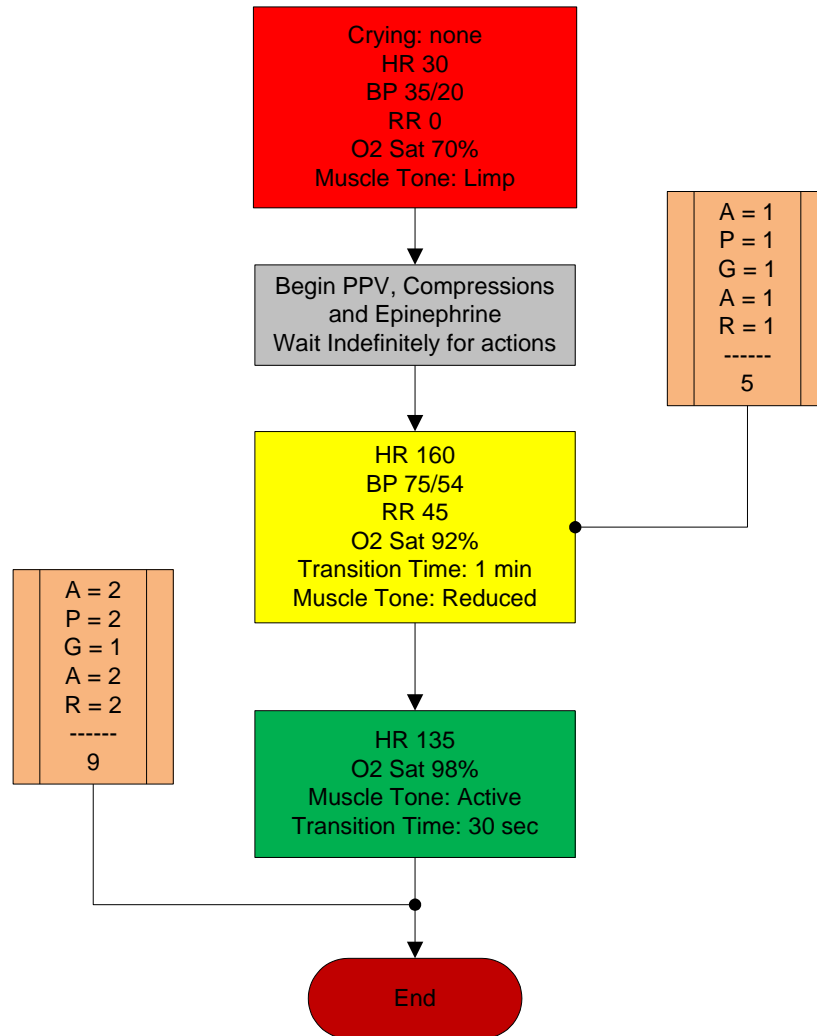This baby is born with a severe asphyxia that has to be treated immediately. After ventilations and EPI have been given, the baby's vitals go towards a good outcome.






**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®  
**MAS**

## Meconium Aspiration Syndrome

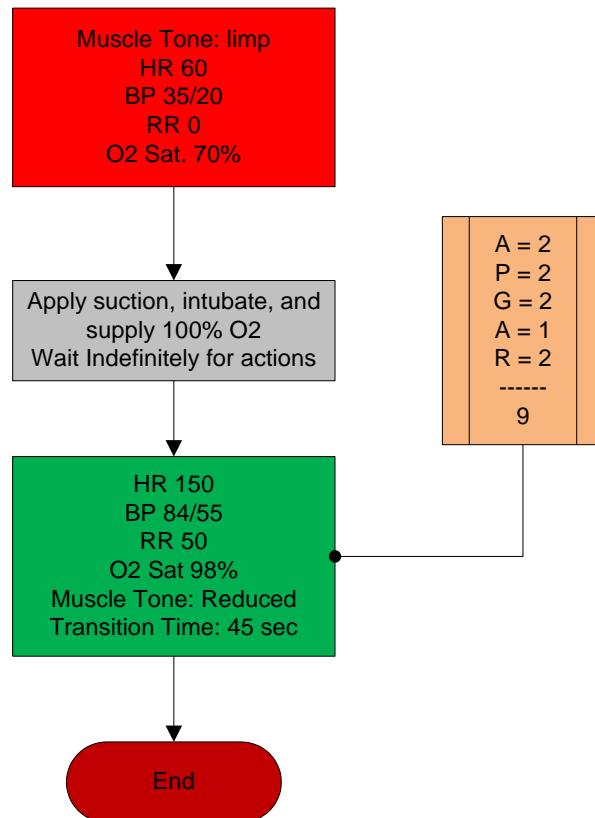






**Gaumard®**  
Simulators for Health Care Education

## Newborn HAL® RDS

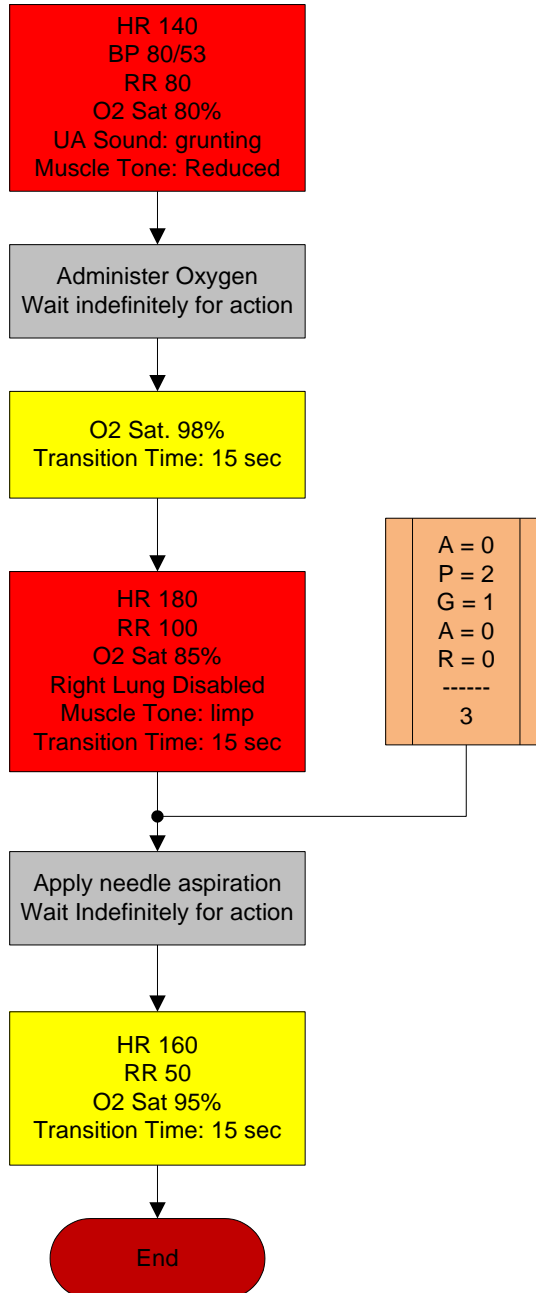Newborn with mild Respiratory Distress Syndrome gets a pneumothorax after oxygen is given.






**Gaumard®**  
Simulators for Health Care Education

Newborn HAL®  
**TTN**

## Transient Tachypnea of the Newborn

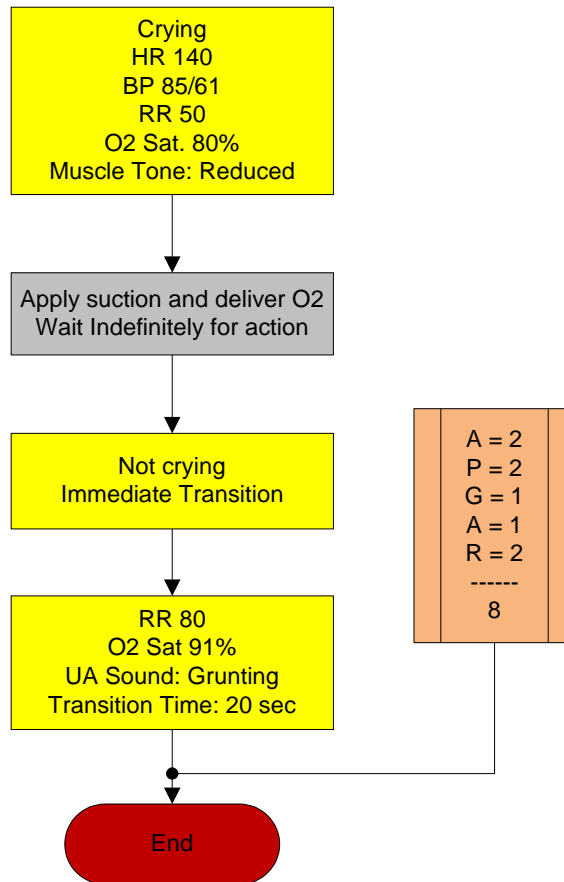

# Tips on Creating Scenarios

## Thinking in Terms of Palette Items

As described previously, Palette items represent complete or partial groups of settings that have been stored as a single item. We learned how applying partial states will hold constant all settings that are left unspecified.

Not only does it take time to customize the palette, but a very large palette becomes difficult to navigate. So, it is desirable to minimize the number of Palette Items in each Profile. To accomplish this, an experienced facilitator tries to create items that are as generally applicable as possible and can, thus, be applied to a wide range of scenarios. The key is to only include in your Palette Items the settings that are directly related to the physiological event represented by that Palette Item.

## Smart Scenarios

After reading the Details, Palette, and Scenarios sections of this guide, it should be clear how to build a scenario. You may have already tried building your own or modifying some of the factory presets. The following four guidelines will refine your ability to build the best possible scenarios.

### 1. How will the scenario begin?

The first thing to consider is the initial condition of the patient. Create a Palette Item to describe this condition. Make sure that this first step in the scenario is a complete state. That is, indicate some selection for each and every available setting on the Details page. Remember that only the settings you specify will cause a change in **NOELLE**, and all other settings will remain constant. So, by starting with a complete state, NOELLE's condition will always be the same when the scenario starts, regardless of what he was doing previously.

Likewise, the "transition duration" of the first step in the scenario should be zero, indicating that changes are applied immediately.

There is one point that can cause confusion and warrants further explanation. It is an extension of the above discussion of partial states. The issue is best illustrated through the following example:

Suppose that you are creating a Palette Item to start your scenario. In this case, you have decided that the patient will be apneic. The question is, "How should the lung sounds be set?"

Most people's first inclination is to set the lung sounds to "none." This is incorrect, despite apnea. Obviously, no lung sounds should be heard during apnea, but since you have already set respiratory rate to zero, none will be. (Sounds are synchronized to the breathing cycle.)

What you are really setting here when you choose a lung sound is the condition of the lungs, given respiratory drive. That is, if the patient's respiratory rate were changed from zero, what sound would be heard? Assuming that the lungs themselves are normal in this scenario, you would choose "normal" for the lung sound setting.

Then, as the scenario progresses, if the patient starts breathing, there will be no need to set the lung sound again. It will already be set. The same principle applies to the heart sound and other settings.

### 2. Include notes to guide the facilitator during the simulation.

It is common for scenario designers, especially those who act as facilitators, to neglect the importance of notes in the scenario. They think that they will remember the learning objectives, patient history, and other details at the time they are ready to conduct the simulation. They usually don't, especially when revisiting a scenario months after creating it.

When you add "Wait" and "Wait Indefinitely" steps to a scenario, you have an opportunity to edit the item description. Use this description field to hold notes to the facilitator. Typically, scenario designers put notes there to indicate what the provider(s) or facilitator should be doing at that point.

Further, when saving the scenario, you may edit the scenario description. This is the best place to put patient history and any other longer notes and instructions.

### 3. Assume that providers will do the right thing.

Usually, you should create a scenario with the assumption that the providers will perform correctly. As long as they do, the scenario can simply be allowed to continue.

Naturally, you must be prepared for what might happen to NOELLE when providers deviate from expectations. The consequences of such deviations can sometimes be included in the scenario, punctuated by "Wait Indefinitely" items. In other cases, the simulation will require more direct control by the facilitator via either the Palette or Details page.

### 4. Choose auto-response settings based on the scenario content and the objectives.

As you've seen, auto-responses can be used to free the facilitators attention. They also enhance realism by presenting instant reactions to the care providers. On the other hand, sometimes it is not possible or desirable to determine the responses before the simulation begins. Different environments and applications call for different settings.

Loosely structured teaching and practice is usually best done with the auto-response settings in Prompt mode. Responses must be triggered by a vigilant facilitator. Though it is slower and requires more attention, the benefit of Prompt over other modes is that the simulation can be allowed to go in any direction, and it will be possible to choose the response on a case-by-case basis.

Tightly structured teaching and assessment requires a higher degree of automation. For such applications, most facilitators choose Auto mode for the auto-response settings. The key issue is standardized timing of symptom presentation. A consistent, repeatable simulation is essential for fair assessment of that care provider in relation to others and for the broader interpretation of results in the context of training validation studies.

When in doubt, it is best to choose *Prompt* mode, in which the facilitator will be given direct control of the responses as events are detected.

# Troubleshooting

## General Troubleshooting Guide

Use the following table to find causes and solutions to a number of possible problems.

Symptom	Possible Cause	Solution
Communication never gets established or is lost ( <b>blinking communication indicator is consistently red</b> )	Battery is discharged	If NOELLE's backup battery is completely discharged, connect the charger and wait 20 minutes to power on the simulator. Leave the simulator connected while in use. <ul style="list-style-type: none"> <li>NOELLE should always be plugged in while being operated. The birth mechanism is power intensive and will drain the battery quickly.</li> <li>Newborn must be plugged in to the charger during use.</li> </ul>
	Communication Module is not detected.	Perform a full shutdown of the tablet.
	Disconnected power plug	NOELLE's internal battery is used for transportation. Always use NOELLE connected to the charger.
	Communication module RJ-45 Cable is not connected to the simulator.	Both simulators operate using a wired connection. Connect the wired communication module to the simulator and restart GaumardUI.
	All others	Close the GaumardUI software and unplug the USB communication module for at least 5 seconds, then plug it back in. Restart the software and wait for initialization
Sound quality while streaming is poor.	Sound is too low or too loud.	Sound volume at PC side is managed from PC's volume control. Simulator sound volume is managed from PC's Microphone gain control. Adjust microphone gain until simulator voice level equals user's voice intensity.  Always talk as close as possible to the microphone in order to improve quality. Using a headset is recommended.
	Respiration and other undesirable sounds are heard by instructor.	Since simulator's microphone has high sensitivity in order to capture the voice of providers, it also captures all surrounded noises on or around the simulator. This is normal and it is not a malfunction.
	Simulator is set to "Generic" "Setup/Options/Environment" menu.	Make sure to select multiple simulator environment (Setup→Options→Environment tab), and enter the Serial Number of the simulator you are using.  <b>Warning: Streaming audio is disabled in "generic mode" simulator is checked.</b>
Streaming audio does not work, tab is not displayed.	The backup battery on the simulator is depleted	Plug charger into the simulator. Verify LED light on charger indicates "charging" status. <ul style="list-style-type: none"> <li>NOELLE should always be plugged in while being operated. The birth mechanism is power intensive and will drain the backup battery quickly.</li> <li>Reconnect Newborn to the charger.</li> </ul>
GaumardUI has set the power mode to STAND-BY automatically	Wired communication module is not connected	Connect the communication module to any USB port.
"Wired communication module not found" message is displayed when GaumardUI is started	Wired communication is module not identified by the computer	Close the software and try disconnecting the communication module for at least five seconds, then plug it back in and restart the software
	Drivers not properly loaded	Perform a full shutdown procedure on the tablet.
	Is the communication indicator panel consistently yellow?	See solution above in section making reference to " <b>blinking communication indicator is consistently yellow</b> "
Chest compressions are not properly detected or not detected at all	Is the respiratory rate set to "0 / min"? Chest compressions are only detected when the respiratory rate is	Set respiration rate to zero

Symptom	Possible Cause	Solution
	set to 0 per minute (0 / min). Otherwise they are ignored	
	Simulator is not connected	Verify connection to the simulator.
	All others	See “Calibration Wizard” section inside User’s Manual
Newborn artificial ventilations are not properly detected or not detected at all	Simulator is not connected	Verify connection to the simulator.
	All others	See “Calibration Wizard” section inside User’s Manual
	Simulator not running	Simulator must be powered on and respiratory rate set at 0.
Newborn simulator’s chest does not rise with artificial ventilation (e.g. BVM)	Incorrect respiratory settings	Respiration rate and inspiration percentage must be higher than “0”. If problem persist, shut down the simulator and restart the tablet. Turn on the system then verify connection, battery and lung settings.
Newborn’s Low chest rise (or no chest rise at all) while breathing	Wrong force sensor reading	When baby is pulled with more than 35 lbs. of force, to avoid damage, the baby is released. If user is not pulling the baby it means that “Force Sensor” is offset. See the “Calibration Wizard” section inside User’s Manual and reset the force sensor to its default value (“Reset to Default Force” button)
NOELLE®’s delivery baby is disengaged prematurely	Wrong force sensor reading	Try manual “Release” from button located on the left hand site corner under “Delivery” tab. If baby releases, then calibrate “Force Sensor”. See the “Calibration Wizard” section inside User’s Manual and reset the force sensor to its default value (“Reset to Default Force” button)
NOELLE®’s delivery baby does not release at delivery end when user is pulling the baby	Dystocia is “ON”	When “Dystocia On” button is checked, the baby is not released until checking “Dystocia Off”
	Motor is “disoriented”	Reset the delivery motor by going to “Setup/Options/” menu, “NOELLE Features” tab, and then click on “Labor Motor” under “Reset...” tab. If delivery mechanism is at the very beginning in a way that it is compressing the rubber boots, the grinding noise won’t go away for a minute or so. Should that be the case either let it go (it won’t break) or simply move motor forward a few turns, and then reset.
	Initial fetus position was not specified on the software.	Lock the baby into the birthing mechanism and manually turn the baby to either ROA, LOA, LOP and ROP. Then synchronize the baby’s position in the womb with the labor position icon located on the labor tab, page 99.
NOELLE®’s delivery mechanism doesn’t come to its initial position when using the “Reset” button under the “Labor” tab, or it makes a grinding noise when reaching the end of the rails		Select “Quick Start Scenarios” when starting the software. To change profiles from within the GaumardUI. Go to “File/Profile” menu and then select “Modeled Scenarios”
Pre-built scenarios do cannot be loaded	Volume not set to user’s criterion.	Every sound has a volume control. Adjust the volume control located on the status panel to reach the desired level.
A sound is absent or is not heard at desired volume level	Cyanosis intensity not set to user’s criterion.	Set Cyanosis level to a desired level by playing with the “Set Max cyanosis level” control.

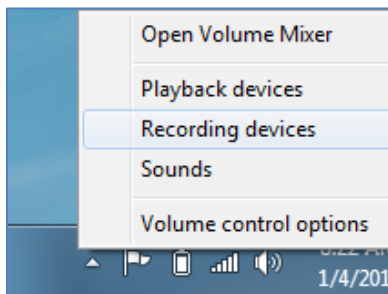
## Microphone Boost (Windows® 7)

Use the instructions below to increase the streaming audio volume. The Headset must be connected to the tablet in order to adjust the microphone volume properties. In addition, adjust the recording control on the headset's physical control to high.


1. Right click on the speaker icon located on the bottom right corner of the laptop's taskbar.




2. The speaker menu is displayed. Click on the recording devices option.



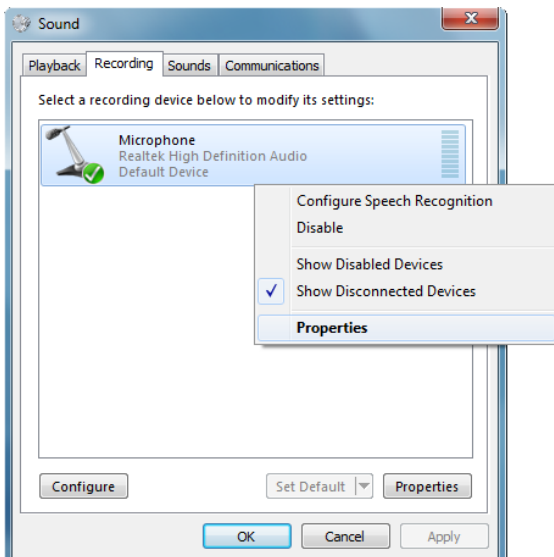3. The sound properties window and recording tab are displayed. Right-click the microphone option and select properties.




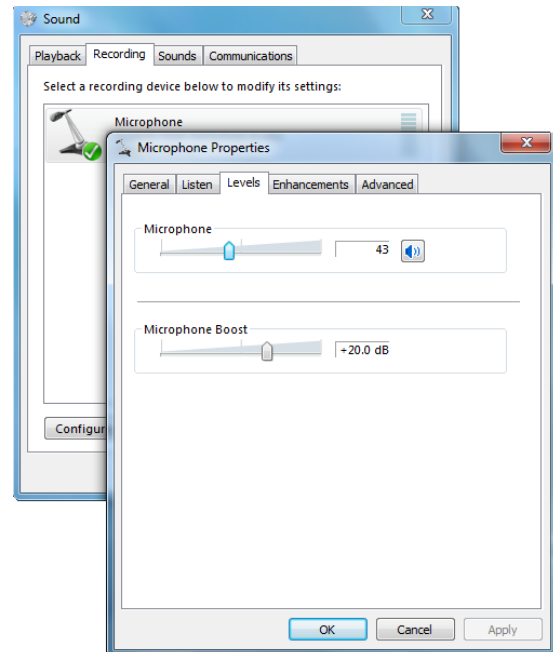4. From the microphone properties sub menu, select the Levels tab. Use the microphone control to decrease and increase the recording volume. For an additional increase in recording volume, use the microphone boost control.



Microphone boost increases volume and saturation which can decrease overall clarity. For optimal clarity, adjust the microphone volume to 100 and the microphone boost to +10.

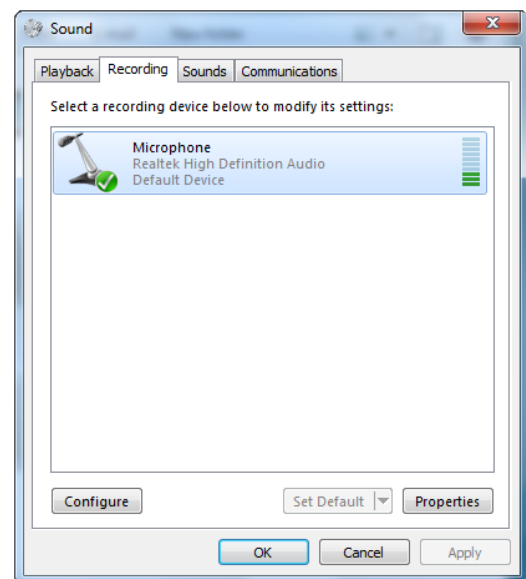5. Click OK to save the changes to the volumes on the microphone properties menu.
6. Click OK to close the Sound properties window.



- It might be necessary to re-adjust the microphone settings to accommodate environment noise.

## Connecting to the Gaumard Monitors

To establish a connection between the Gaumard Virtual Monitor software and GaumardUI, both computers must be connected to the same network **ad-hoc network**. The section below describes in detail how to configure the ad-hoc wireless settings and establish a connection between the Gaumard computer systems. Follow the guide when troubleshooting connection issues.

Vital sign information is sent and updated from the GaumardUI software via Wi-Fi ad-hoc network connection, and not from the simulator or the RF communication module.

### Procedure overview for each computer

- Configuring Static IP addresses
- Configuring the Network SSID name
- Verifying the ad-hoc connection

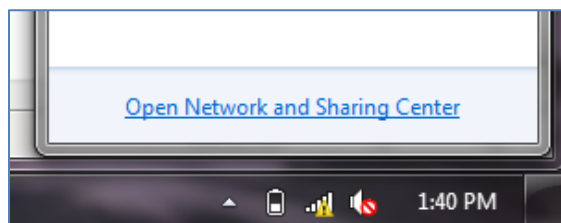
### Configuring the computer static IP

Perform the steps below to assign a static IP addresses to the Gaumard Virtual Monitor computer and the GaumardUI tablet. Please note that each computer requires a unique static IP address.


- Locate the **wireless connection icon** on the bottom right corner of the computer screen.




- Click the wireless icon and select **Open Network and Sharing Center** from the options on the menu.



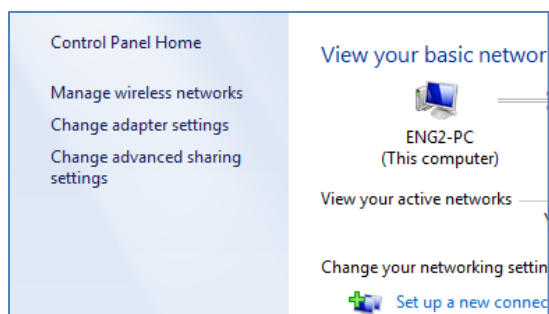
- In the **Network and Sharing Center** window, click **Change adapter settings** from the left-hand panel.



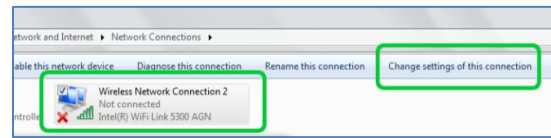- Click on the **Wireless Network Connection** icon to select it, and then click **Change settings of this connection** on the menu panel. The Properties window for the connection will be shown.



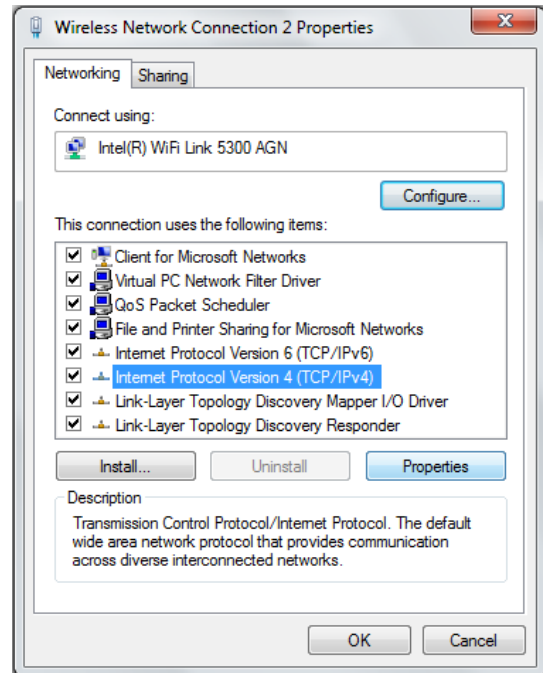
- Select **Internet Protocol Version 4 (TCP/IPv4)** to highlight it, and click Properties.





6. Configure both computers with static IP addresses and subnet masks following the IP address convention below. Each computer in the ad-hoc network must be configured with the same first three octets and only differ on the last.


Computer	IP address	Subnet Mask	SSID Network Name
Tablet 1 NOELLE	1.0.0. <u>1</u>	255.0.0.0	GaumardNet
Tablet 2 Newborn HAL	1.0.0. <u>2</u>	255.0.0.0	GaumardNet
Virtual Monitor Computer	1.0.0. <u>3</u>	255.0.0.0	GaumardNet

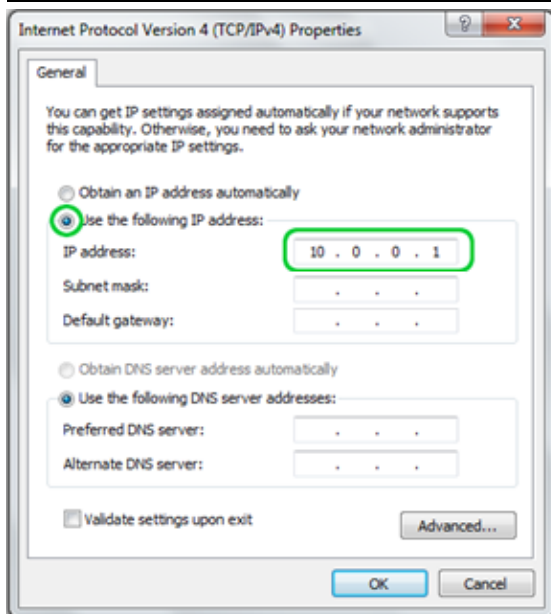

Example of IP addresses that will NOT communicate with each other:

Tablet 1 - 1.0.**10**.1      VM Computer - 1.0.**0**.3  
 Tablet 1 - **50**.0.10.1      VM Computer - **1**.0.10.3  
 Tablet 1 - 10.10.10.**1**      VM Computer - 10.10.10.**1**

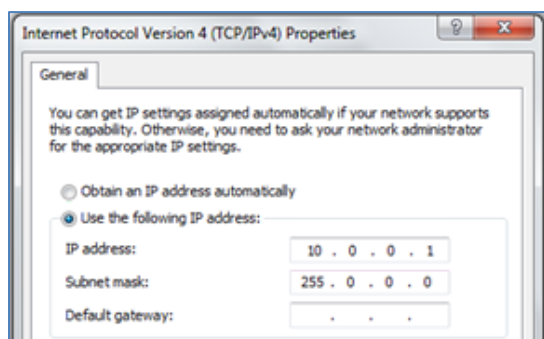
7. Click once in the **subnet mask** field to auto populate the correct address. Leave all other fields blank and click **OK** to save the changes and return to the **Network and Sharing Center** window.




## Configuring the Wireless ad-hoc Network

After each computer is configured with a unique static IP address, follow the steps below to create the **ad-hoc wireless network**. The following steps are performed identically on the GaumardUI tablet and the virtual monitor computer.

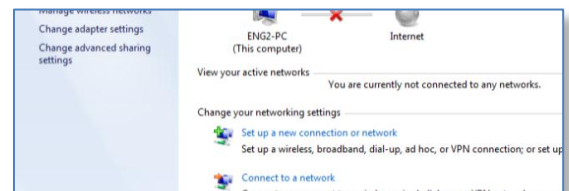1. On Network and Sharing Center and select **Set up a new connection or network**.



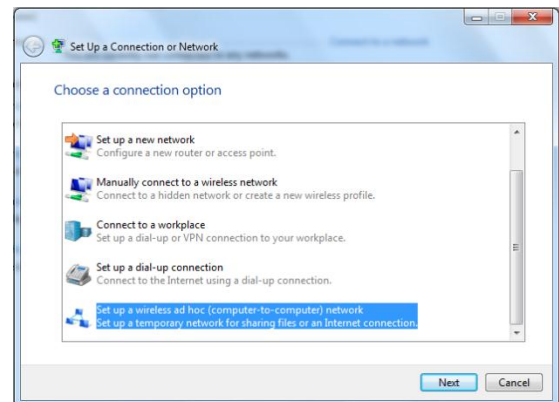
2. Select **Set up a wireless ad hoc (computer-to-computer) network** and click **Next** to continue.




3. The next screen provides some information about the ad hoc connection, click **Next**.



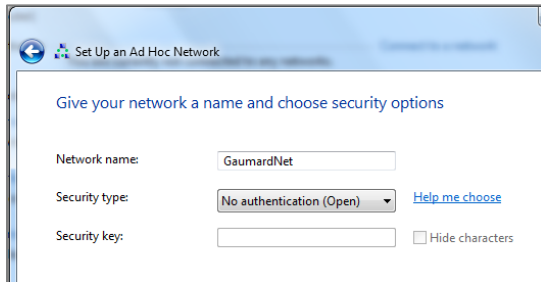1. In the **Network name (SSID)** field, type **GaumardNet**. Network names are **case-sensitive** and must be typed identically on all the computers that will belong to this ad-hoc network.




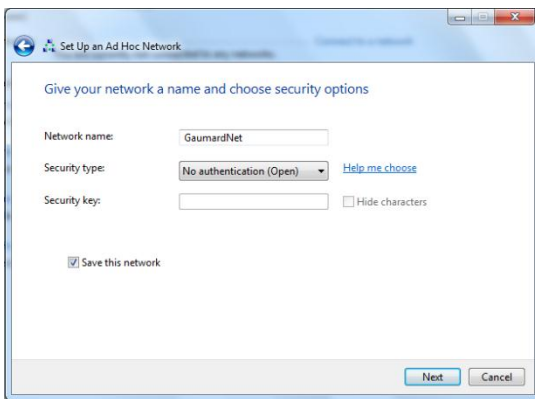2. Ensure the security type is **No authentication**. Checkmark **Save this network**, then click **Next**.




If the system warns of another network with the same name, simply power down any other computer configured with the same network name during this step and try again.

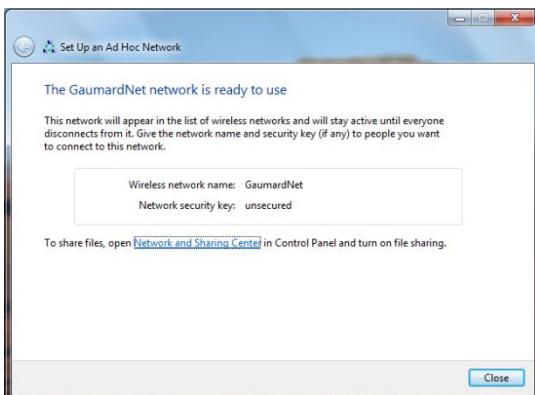4. Windows® will advise if the network is successfully completed by displaying a ready message and a summary of the network's properties. Click "Close" and restart the computer.



## Verifying the ad-hoc wireless connection

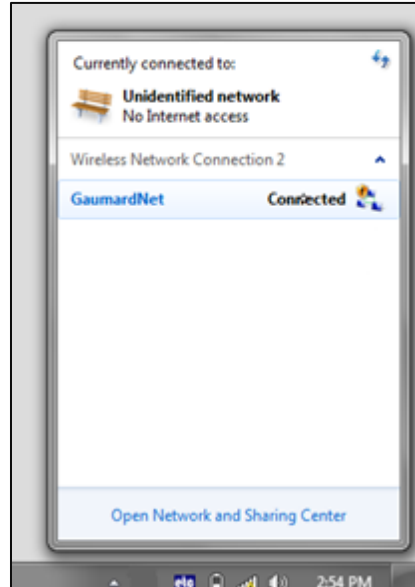Window 7 systems do not automatically connect at startup. To establish a connection to the ad-hoc wireless network, navigate to the wireless icon located at the bottom right corner of the screen, select the ad-hoc network name and click connect. Because the computers are connecting to each other, both systems must be powered on




## Connecting GaumardUI to the Gaumard Monitors software


GaumardUI sends all the vital signs information to the Gaumard Monitor software **over** the wireless ad-hoc connection configured in the previous section. Follow the guide below to establish the connection between the programs.

**CAUTION:** To avoid connection issues, always establish the ad-hoc wireless connection between the computers before activating the GaumardUI and Gaumard Monitors program.

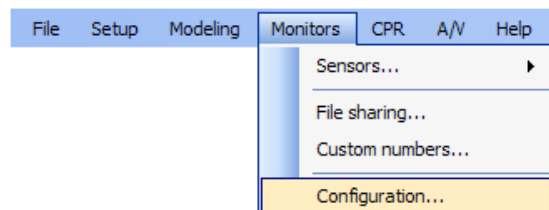
1. On the tablet computer, start the **GaumardUI** software and navigate to **Monitors>Configuration**.



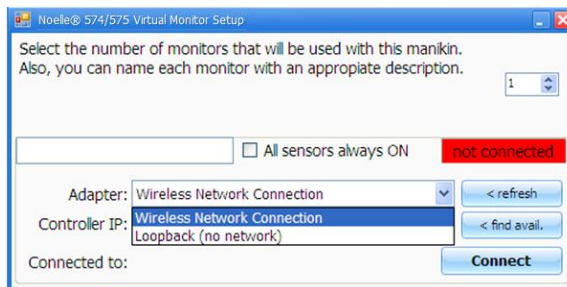- In the **NOELLE® Virtual Monitor Setup** window, access the **Adapter** dropdown menu and select **Wireless Network Connection**. If the adapter is inaccessible, click the stop button first.



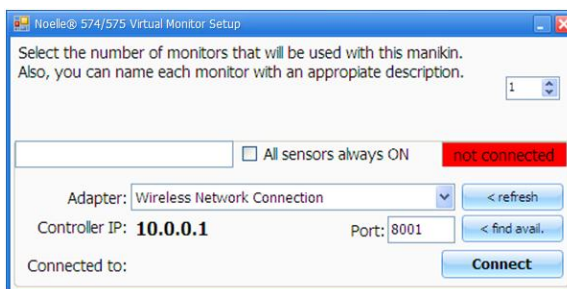
- The controller IP shown is the static IP address previously configured.




- Click **Connect** and take note of the controller IP address and port number. This information will be used on the following steps. The status icon will read **connecting** until the Gaumard Monitor software is configured.
- Activate the Gaumard Monitors software on the Virtual monitor computer.



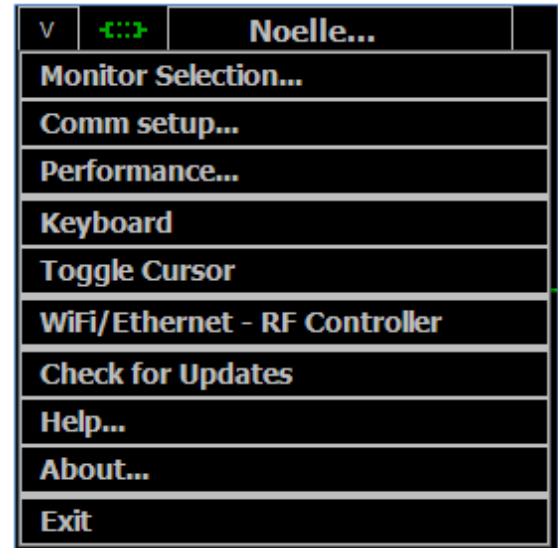
- Click on the menu labeled **V** located on the upper left corner and click **Monitor Selection**.



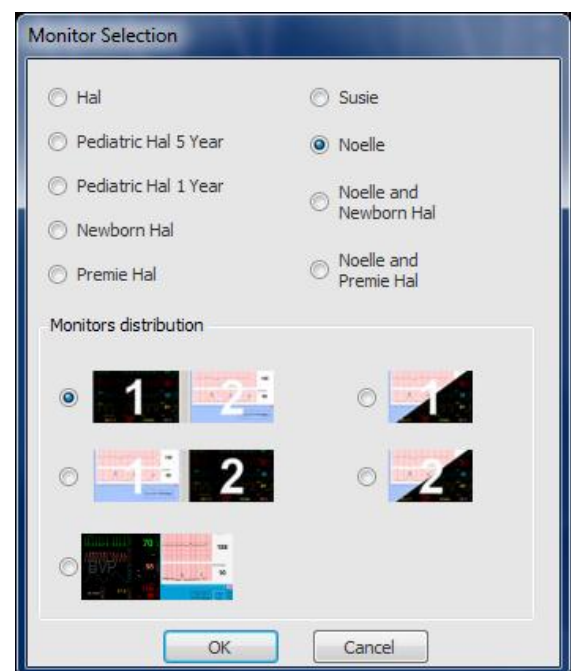- From the **Monitor Selection** menu, select the applicable simulator screen and click **OK** to save the settings.




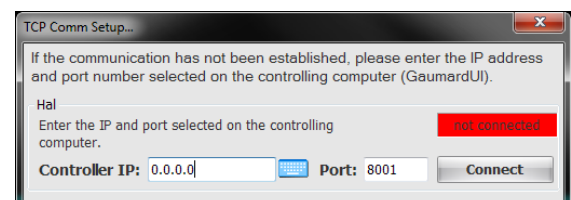- Return to the **V** menu and select **Comm... setup** option to access the **TCP Comm Setup...** window.




- Input the **Controller IP** address and **port number** as previously noted on GaumardUI's **Virtual Monitor Setup** window.


## Monitor selection - NOELLE

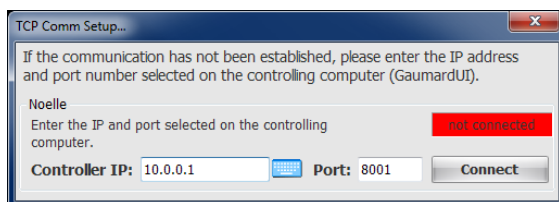


## Monitor Selection – NOELLE and Newborn HAL

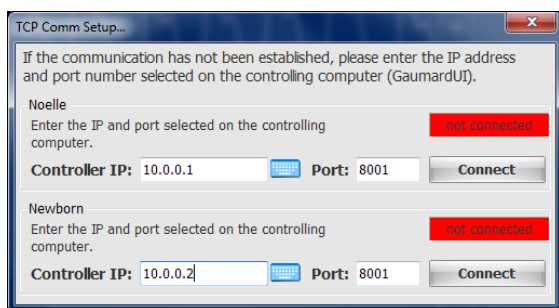

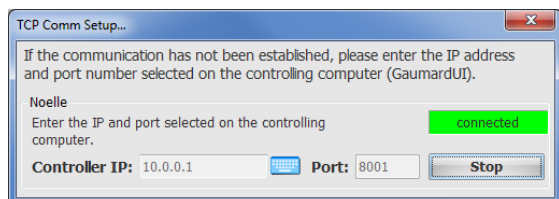10. Click **Connect** to finalize the connection.




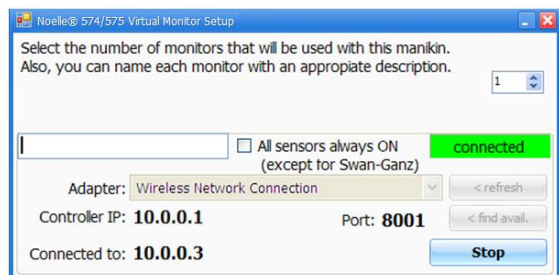11. On GaumardUI's **Virtual Monitor Setup** window the connection status will also display connected.




## Diagnostics

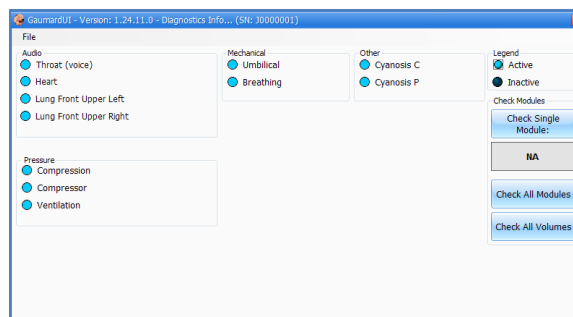The Diagnostics window can be accessed by going to the Help menu and selecting "Diagnostics". This window is very useful for troubleshooting because it gives the user feedback on all of the working modules inside the simulator. The user can click on the button that says "Check All Modules" and the software checks which modules are responding.



Also check individual modules by clicking on the module to highlight. Notice that the name of the module is displayed on the right column.

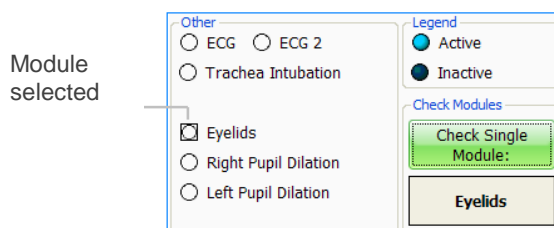Now click on the "Check Single Module" button:



Active modules report light blue, and inactive modules report black. If there is a specific module that fails to respond please contact customer support (ensure that the module that is unresponsive is not specific to an Add-On feature that is not installed on your simulator).

# Selected Consumables and Replacements Parts

## Selected Parts List

Contact Gaumard Scientific for a **complete list** of consumables and replacement parts and their prices.

Item ID	Name	Type	Description
S555.100.001.R2	A/C Virtual Monitor	R	All in one computer
S555.100.004L.D	Lower Left Arm Skin Cover	M	Lower left arm skin cover
S555.100.004L.L	Lower Left Arm Skin Cover	M	Lower left arm skin cover
S555.100.004L.M	Lower Left Arm Skin Cover	M	Lower left arm skin cover
S555.100.004R.D	Lower Right Arm Skin Cover	M	Lower right arm skin cover
S555.100.004R.L	Lower Right Arm Skin Cover	M	Lower right arm skin cover
S555.100.004R.M	Lower Right Arm Skin Cover	M	Lower right arm skin cover
S555.100.007.D	C-Section Abdominal Cover	C	NOELLE stomach cover for C-Section Exercises, dark color
S555.100.007.L	C-Section Abdominal Cover	C	NOELLE stomach cover for C-Section Exercises, light color
S555.100.007.M	C-Section Abdominal Cover	C	NOELLE stomach cover for C-Section Exercises, medium color
S555.100.008.R2.D	Abdominal Cover	R	NOELLE stomach cover with foam insert, dark color
S555.100.008.R2.L	Abdominal Cover	R	NOELLE stomach cover with foam insert, light color
S555.100.008.R2.M	Abdominal Cover	R	NOELLE stomach cover with foam insert, medium color
S555.100.010	Battery	M	Rechargeable battery
S555.100.011	Battery Charger	R	Battery charger with label
S555.100.013	Birthing Mechanism	R	Automatic Birthing Mechanism
S555.100.016.D	Birth canal	C	Dark color
S555.100.016.L	Birth canal	C	light color
S555.100.016.M	Birth canal	C	Medium color
S555.100.032	Automatic Boggy Uterus	R	Boggy Uterus for <b>automatic</b> PPH
S555.100.033	Episiotomy Trainer Set	R	Episiotomy Trainer set with vulva insert
S555.100.040.D	Articulating baby	R	Articulating birthing fetus
S555.100.040.L	Articulating baby	R	Articulating birthing fetus
S555.100.040.M	Articulating baby	R	Articulating birthing fetus
S555.100.048	Adult IV Filling Kit	A	Fluid dispensing syringe with filling tube
S555.100.060	Simulator Transport Case	R	Soft storage and transport case with wheels

Item ID	Name	Type	Description
S555.100.062	Shipping Cardboard Box	R	Box dimensions 60"x20"x30" for simulators with soft case
S555.100.080	Simulated Blood Concentrate	C	
S555.100.081	Silicone Oil	R	Oil-based Silicone lubricant
S555.100.087	Wireless Streaming Audio Headset	R	
S555.100.207	Laptop PC	R	GUI Software included
S555.100.310	Wireless keypads	R	Extra pack of five wireless keypads
S555.100.EXW	Two Year Extended Warranty	A	Extended warranty for years Two <b>AND</b> Three
S555.100.INST	In-Service Training	A	Day of in-service training and installation

C=Consumables; R=Replacements; A=Accessories; U=Upgrades; M=Replace in Miami Factory ONLY



# Replacing Common Consumable and Replacement Parts

## Birth Canal

Follow the instructions below to replace the birth canal. You may also use the steps below in case the birth canal has to be removed temporarily to adapt the episiotomy kit.

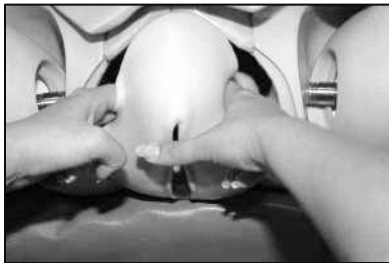
### Remove birth canal:

1. Insert your left hand at the 9 o'clock position, between the lower torso and the birth canal insert. Slowly push inward.

**WARNING:** Do not pull the birth canal to remove.



2. Repeat instructions in step 1 with right hand at the 3 o'clock position.



3. Push inward on the birth canal insert. Then place your thumbs on the top of the birth canal insert and push downward.



4. Slowly pull the birth canal insert towards you. When the urine catheter tube is visible, detach from the birth canal insert.



5. Finally, after removing the urine catheter tube, gently remove the birth canal insert.

### Attach the birth canal

1. Place birth canal insert into polythene bag and ensure the bag covers the Velcro on the outer wall of the insert as shown in picture.



2. While first inserting the plastic bag through the birth canal opening on the mannequin, gently position the birth canal insert with the approximate final alignment to the opening.



3. Push the birth canal insert into position. Place your hands on the lower half of the insert with your fingers holding the bag in place. Steady the front of the insert with your thumbs. Without removing the plastic, orient the birth canal insert so that it is properly aligned with the opening on the lower torso of the mannequin.



4. Gently begin to remove the polythene bag with one hand while holding the insert in place with the other, as pictured below. It is easier to start at the top, then work down the sides towards the bottom.



5. Place one hand towards the bottom of the birth canal insert when pulling out the bottom portion of the polythene bag.



6. With the birth canal insert in place, pull the post-partum hemorrhage tube through the hole in the bottom of the birth canal insert as shown in the picture. This is to position the hemorrhage tube in the correct location. The birth canal insert is now ready for normal use.





# Warranty

## **EXCLUSIVE ONE-YEAR LIMITED WARRANTY**

Gaumard warrants that if the accompanying Gaumard product proves to be defective in material or workmanship within one year from the date on which the product is shipped from Gaumard to the customer, Gaumard will, at Gaumard's option, repair or replace the Gaumard product.

This limited warranty covers all defects in material and workmanship in the Gaumard product, except:

1. Damage resulting from accident, misuse, abuse, neglect, or unintended use of the Gaumard product;
2. Damage resulting from failure to properly maintain the Gaumard product in accordance with Gaumard product instructions, including failure to properly clean the Gaumard product; and
3. Damage resulting from a repair or attempted repair of the Gaumard product by anyone other than Gaumard or a Gaumard representative.

**This one-year limited warranty is the sole and exclusive warranty provided by Gaumard for the accompanying Gaumard product, and Gaumard hereby explicitly disclaims the implied warranties of merchantability, satisfactory quality, and fitness for a particular purpose.** Except for the limited obligations specifically set forth in this one-year limited warranty, Gaumard will not be liable for any direct, indirect, special, incidental, or consequential damages, whether based on contract, tort, or any other legal theory regardless of whether Gaumard has been advised of the possibilities of such damages. Some jurisdictions do not allow disclaimers of implied warranties or the exclusion or limitation of consequential damages, so the above disclaimers and exclusions may not apply and the first purchaser may have other legal rights.

**This limited warranty applies only to the first purchaser of the product and is not transferable.** Any subsequent purchasers or users of the product acquire the product "as is" and this limited warranty does not apply.

**This limited warranty applies only to the products manufactured and produced by Gaumard.** This limited warranty does not apply to any products provided along with the Gaumard product that are manufactured by third-parties. For example, third-party products such as computers (desktop, laptop, tablet, or handheld) and monitors (standard or touch-screen) are not covered by this limited warranty. Gaumard does not provide any warranty, express or implied, with respect to any third-party products. Defects in third-party products are covered exclusively by the warranty, if any, provided by the third-party.

Any waiver or amendment of this warranty must be in writing and signed by an officer of Gaumard.

In the event of a perceived defect in material or workmanship of the Gaumard product, the first purchaser must:

1. Contact Gaumard and request authorization to return the Gaumard product. Do **NOT** return the Gaumard product to Gaumard without prior authorization.
2. Upon receiving authorization from Gaumard, send the Gaumard product along with copies of (1) the original bill of sale or receipt and (2) this limited warranty document to Gaumard at 14700 SW 136 Street, Miami, FL, 33196-5691 USA.
3. If the necessary repairs to the Gaumard product are covered by this limited warranty, then the first purchaser will pay only the incidental expenses associated with the repair, including any shipping, handling, and related costs for sending the product to Gaumard and for sending the product back to the first purchaser. However, if the repairs are not covered by this limited warranty, then the first purchaser will be liable for all repair costs in addition to costs of shipping and handling.

## **Extended Warranty**

In addition to the standard one year of coverage, the following support plans are available:

- Two-Year Extension (covers second and third years)
- Call for pricing (USA only)

## Contact Us

**E-mail Technical Support:** support@gaumard.com

**E-mail Sales and Customer Service:** sales@gaumard.com

**Phone:**

Toll-free in the USA: (800) 882-6655

Worldwide: 01 (305) 971-3790

**Fax:** (305) 667-6085

Before contacting Tech Support **you must:**

1. Have the simulator's Serial Number (located in the left leg under the IM site)
2. Be next to the simulator if troubleshooting is needed.

**Post:** Gaumard Scientific

14700 SW 136 Street

Miami, FL 33196-5691

USA

**Office hours:** Monday-Friday, 8:30am - 4:30pm EST (GMT -4 Summer Time)

© Gaumard Scientific Company, 2012

All rights reserved.

Gaumard®, NOELLE®, HAL® are trademarks of Gaumard Scientific Company, Inc.

Patented; Other Patents Pending

Always dispose of this product and its components in compliance with local laws and regulations.

